




Scout® *Pro* Balance Instruction Manual

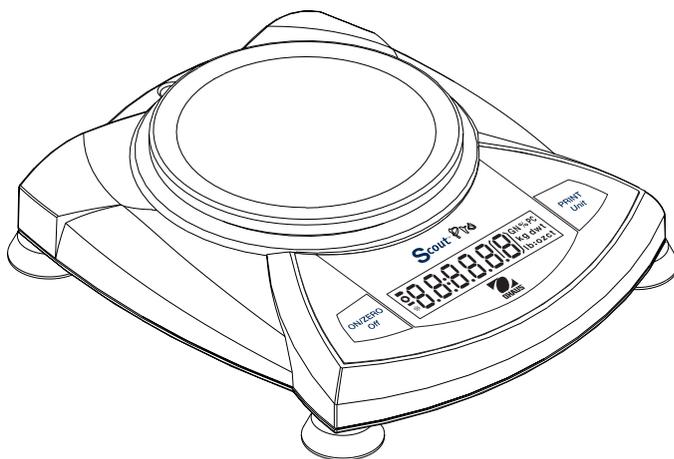

Declaration of Conformity

The undersigned, representing the following manufacturer

Ohaus Corporation
 19A Chapin Road
 P.O. Box 2033
 Pine Brook, NJ 07058
 USA

hereby declares that the following products are in conformity with the EEC directives listed below (including any and all modifications).

Balance models: SP123, SP202, SP401, SP402, SP601, SP2001, SP2001N, SP4001, SP6000, SPE123, SPE202, SPE402, SPE401, SPE601, SPE2001, SPE4001, SPE6000, SPU123, SPU202, SPU402, SPU401, SPU601, SPU2001, SPU4001 SPU6000, SPG123F, SPG202F, SPG402F, SPG401F, SPG601F, SPG2001F, SPG4001F, SPG6000F, SPS202F, SPS402F, SPS401F, SPS601F, SPS2001F, SPS4001F, SPS6000F, JS40, JS500, JS1200

Marking	Directive	Standard
	73/23/EEC Low Voltage Basse tension Baja tensión	EN60950: 1992 + A1: 1993 + A2: 1993 + A3: 1995 + A4: 1997 Safety of information technology equipment. Seguridad de los equipos de tratamiento de la información incluyendo los equipos eléctricos de oficina. Sécurité des matériels de traitement de l'information.
	Electromagnetic compatibility Compatibilité électromagnétique Compatibilidad electromagnética	Electrical equipment for measurement, control and laboratory use (Class B) Matériels électriques de mesure, de commande et de laboratoire — Prescriptions relatives à la CEM (Class B) Equipo eléctrico de medida, control y uso en laboratorio — Requisitos de compatibilidad electromagnética (Class B)

Last two digits of the year which the CE marking was affixed: 03

ISO 9001 Certificate for Ohaus Corporation – Ohaus Corporation, USA was examined and evaluated in 1994 by the Bureau Veritus Quality International (BVQI) and was awarded the ISO 9001 certificate. This certifies that Ohaus Corporation, USA, has a quality system that conforms to the international standards for quality management and quality assurance (ISO 9000 series). Repeat audits are carried out by BVQI at intervals to check that the quality system is operated in the proper manner.



Ted Xia
President

Ohaus Corporation
Pine Brook, NJ USA

Date: September 22, 2003



Urs Müller
General Manager
Ohaus Europe

Greifensee, Switzerland

Date: September 22, 2003

FCC NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CANADIAN NOTE:

This class B digital apparatus complies with Canadian ICES-003.

Cet appareil numérique de la classe B est conforme à la norme NMB-003 du Canada.



13123

AS/NZS4251.1 AS/NZS4252.1 Emission and Immunity

TABLE OF CONTENTS

1. INTRODUCTION	EN-3
Safety Precautions	EN-3
2. INSTALLATION	EN-3
Unpacking	EN-3
Installing Components	EN-4
Releasing the Shipping Lock	EN-4
Platform Installation	EN-4
Draft Shield Installation (120g Model)	EN-4
Security Bracket	EN-4
Selecting the Location	EN-5
Balances with Level Adjustment	EN-5
Connecting Power	EN-5
Battery Installation	EN-5
AC Adapter Installation	EN-5
3. OPERATION	EN-6
Overview of Controls and Display Functions	EN-6
Button Functions	EN-8
Symbols Used for Operation of the Balance	EN-9
Turning the Balance On	EN-9
Turning the Balance OFF	EN-10
Navigating the Menu	EN-10
Menu Structure	EN-11
Entering the Menu	EN-11
Accepting/Bypassing an Individual Menu Item	EN-11
Entering the .S.E.T.U.P. Menu	EN-12
Turning Display Hold, Totalize or Specific Gravity Mode On	EN-13
Exiting the .S.E.T.U.P. Menu	EN-13
Entering the .U.N.I.T. Menu	EN-13
Parts Counting	EN-14

TABLE OF CONTENTS (Cont.)

Calibration	EN-14
Span Calibration	EN-14
Linearity Calibration	EN-15
Applications	EN-16
Weighing	EN-16
Weighing with Tare	EN-16
Parts Counting	EN-17
Percent Weighing	EN-18
Establishing a New Reference Weight	EN-19
Exiting Percent Weighing	EN-19
Display Hold	EN-19
Exiting Display Hold	EN-20
Totalize	EN-20
Clear Exit Totalize	EN-21
Specific Gravity	EN-21
Clear Exit Specific Gravity	EN-22
Additional Features	EN-22
Weigh Below	EN-22
Lock Switch	EN-23
4. MAINTENANCE	EN-24
Cleaning	EN-24
Troubleshooting	EN-24
Error Codes List	EN-25
Accessories	EN-26
5. TECHNICAL DATA	EN-27
Specifications	EN-27
Warranty	EN-29

1. INTRODUCTION

The Scout *Pro* offers parts counting with auto optimization, display hold, totalize, % weighing and specific gravity mode. Models are available with ranges from 120g to 6000g.

Scout *Pro* standard features include:

- Battery or AC operation (AC adapter included)
- Density/Specific Gravity determination (certain models)
- Integral security bracket
- Programmable auto shut-off
- Span calibration masses included on certain models
- Optional USB or RS232 interface available

Safety Precautions

Please follow the safety precautions as listed:

CAUTION:



- Do not operate the balance around corrosive fumes.
- Use only the adapter provided with the balance.
- Do not try to service the Scout *Pro* balance.
- Before plugging in the balance, make sure that the voltage rating of the power adapter and the AC Adapter match.

2. INSTALLATION

Unpacking

Inform your Ohaus dealer if parts are missing.

Your Scout *Pro* package contains:


- Scout *Pro* Balance
- Warranty card
- AC Power Adapter
- Platform
- Instruction Manual
- Calibration Masses (on certain models)
- Draft shield (120g Model only)
- Specific Gravity Kit (on certain models)

Store the packaging material for future transport.

Installing Components

Releasing the Shipping Lock



On top of the balance, turn the pointer 90 degrees counter-clockwise to unlock.

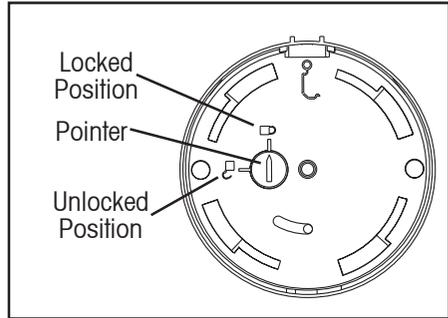

Releasing the Shipping lock.

Platform Installation



Balances with a rectangular platform are placed into the subplatform as shown and rotated counter-clockwise until it locks. Round platforms are placed straight down on subplatform.

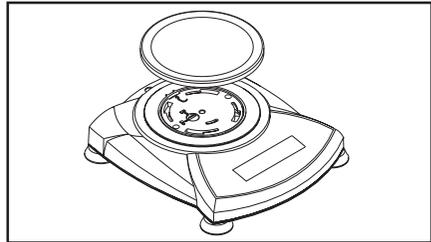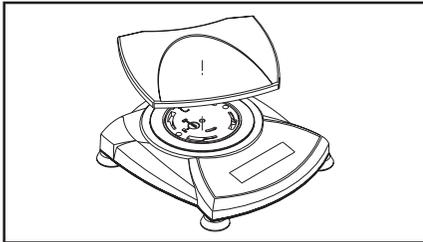

Platform Installations.

Draft Shield Installation (120g Model)

Position the keyed draft shield on top of the balance and rotate it until it faces forward.

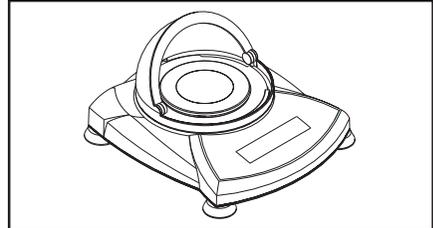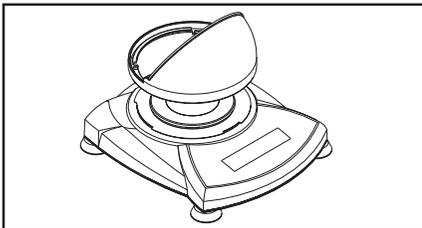

Installing Draft Shield.

Security Bracket


A security bracket is provided at the rear of the balance allowing the balance to be secured by an optional cable and lock accessory.

Selecting the Location

For best performance, the Scout *Pro* balance should be used in a clean, stable environment. Do not use the balance in environments with excessive drafts, with rapid temperature changes, near magnetic fields or near equipment that generates magnetic fields, or vibrations.

Balances with Level Adjustment

Balances containing leveling feet and a bubble level must be leveled before using. See illustration for correct leveling.

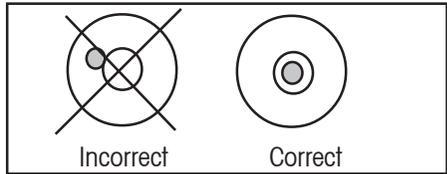

Bubble Level indicator


Connecting Power

Battery Installation

Install the Four "AA" batteries with polarity as shown in the battery compartment.

AC Adapter Installation


Plug the AC adapter into the jack at the rear of balance.

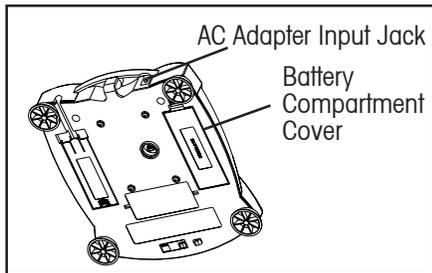

Battery and AC Power Connections

3. OPERATION

Overview of Controls and Display Functions

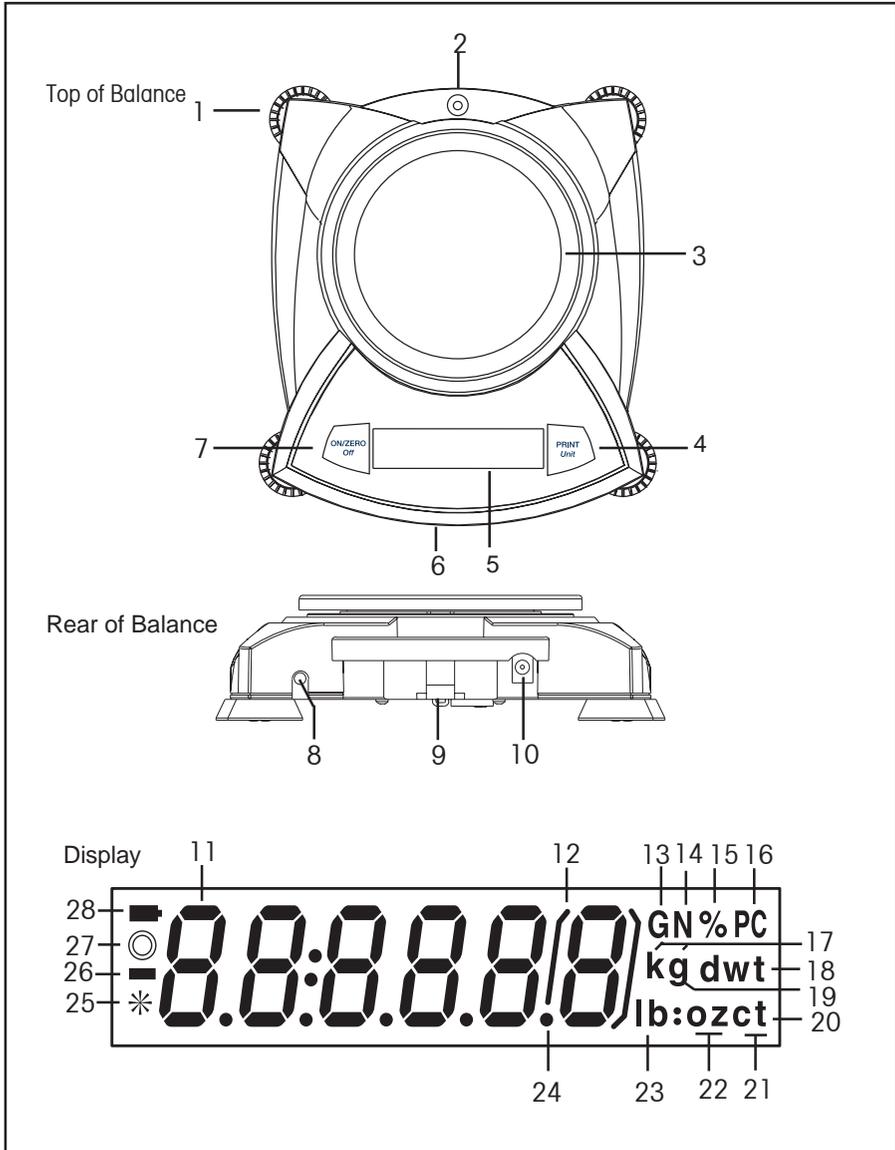

No.	Designation	Function
1.	Feet ¹	Provides leveling for certain models.
2.	Spirit Level ¹	Provides level indication.
3.	Platform	Weighing platform, either round or rectangular.
4.	PRINT <i>Unit</i> button	Prints data, scrolls through units, steps through menu options.
5.	Display	LCD display with icons.
6.	Lockswitch	Locks certain menu functions, located under balance.
7.	ON/ZERO <i>Off</i> button	On/Off, Zero, enters menu, accepts menu settings.
8.	USB or RS232 port	Optional kit for either RS232 or USB operation.
9.	Security Bracket	Part of balance for optional external cable and lock.
10.	Power Input Jack	Connector for AC adapter.
11.	7-segment LCD	Part of 6-digit LCD display.
12.	Brackets	Auxilliary indication.
13.	G	Indicates specific gravity
14.	N	Indicates weight in Newtons.
15.	%	Indicates percent weighing.
16.	PC	Indicates pieces during parts counting.
17.	kg	Indicates weight in kilograms. (certain models)
18.	dwt	(not used)
19.	g	Indicates weight in grams.
20/22.	oz †	(not used)
20.	†	Indicates weight in totalize mode (or tael in certain models).
21.	ct	(not used)
22.	oz	Indicates weight in ounces.
23/22.	lb:oz	(not used)
23.	lb	Indicates weight in pounds. (certain models)
24.	•	Decimal point.
25	*	Stability indicator, indicates stable weight and Display hold when flashing.
26.	-	Negative sign.
27.	○	(not used)
28.	■	Battery indicator flashes when battery is down to approximately 20 minutes of power remaining.

NOTES: 1. Certain models are provided with leveling feet and spirit level.

2. Unit measurement varies by model.

Button Functions

Two switches provide the necessary functions to access a given menu, select a function and to turn it on or off. Functions are listed as follows:



ON/ZERO Off Button

Primary Function (ON/ZERO)- Turns on balance. If balance is on, zeros the display.

Secondary Function (Off)- Turns balance off, **OFF** will be displayed after button is held for 3 seconds. In Display Hold or Totalize mode, a long press exits the mode without turning the balance off.

Menu Function- An extended long press (>5 seconds) during power up will cause the balance to enter the Menu mode. A short press is used to accept a setting on a display.

PRINT Unit Button

Primary Function (PRINT)- Sends print command to interface port. If Display Hold or Totalize mode is active, a short press will enter that mode.

Secondary Function (Unit)- Press and hold scrolls through units. Release on desired unit.

Menu Function- Will bypass setting shown on display.

Symbols Used for Operation of the Balance

Symbols are used to simplify the setup and operation of the balance. A description of each symbol follows:

 Press

The clock symbols adjacent to the finger symbol indicates the length of time to press a button.

 1 second momentary press.

 3 second extended press.

 5 second extended press.

ON/ZERO
Off

PRINT
Unit

Panel control buttons used to initiate actions.



Displays are shown as they actually appear on the balance. A model with 200g capacity was used for the displays shown in this manual.

●●● Indicates scrolling to a final display. The first and last displays are shown.

➡ Indicates advances to next display.

Turning the Balance On




ON/ZERO
Off

➡



●●●


0.00 g

Turning the Balance Off




ON/ZERO
Off

➡



0.00 g

➡

OFF

Navigating the Menus

Menu Structure

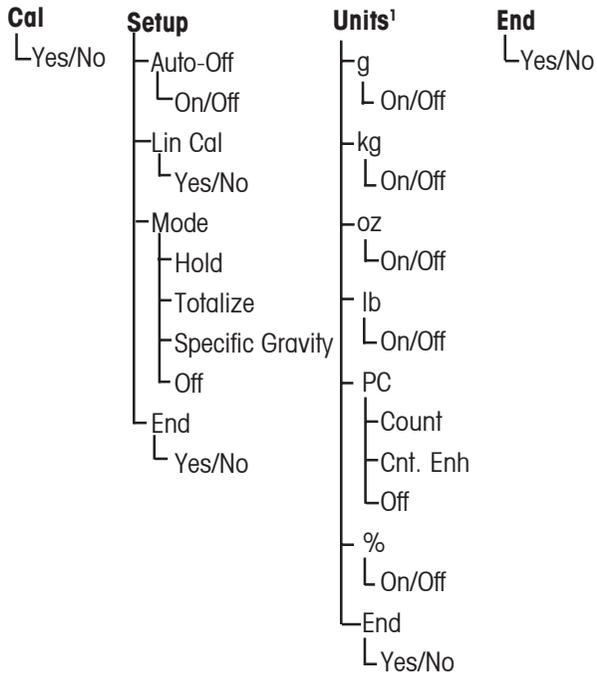

NOTES:


1. Refer to specification table for available units.

Entering the Menu

There are four main menus in the balance: **.C.A.L.**, **.S.E.T.U.P.**, **.U.N.I.T.S.** and **.E.N.D.**


Start with the balance off and the Lock Switch off (see page 23).

The sequence is shown below.

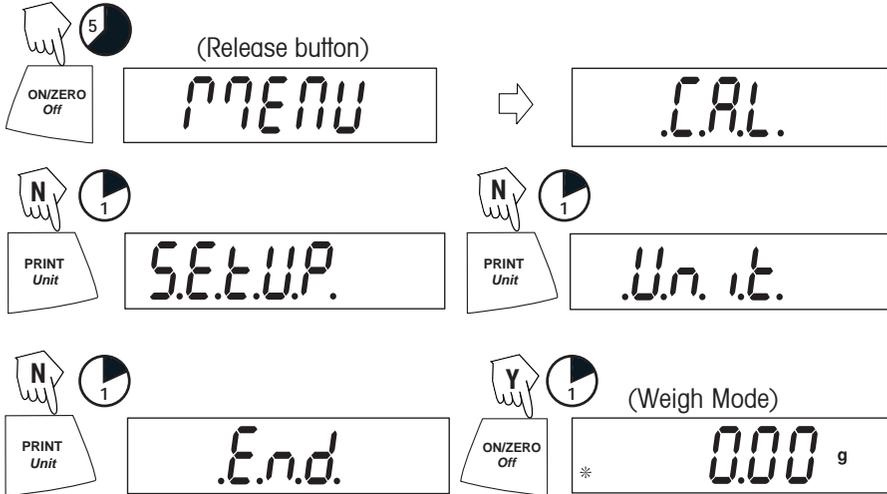

Accepting / Bypassing an Individual Menu Item


Start with menu item displayed.

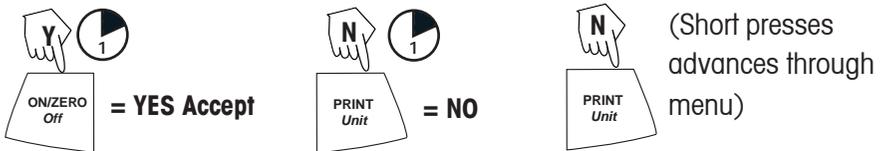

Entering the .S.E.T.U.P. Menu

The .S.E.T.U.P. menu contains Auto-Off, Linearity Calibration, Mode (Display Hold, Totalize, Specific Gravity) and END. Auto-Off can be turned on or off. Display Hold, and Totalize require entering the Mode submenu.

Starting from the .C.A.L. menu.

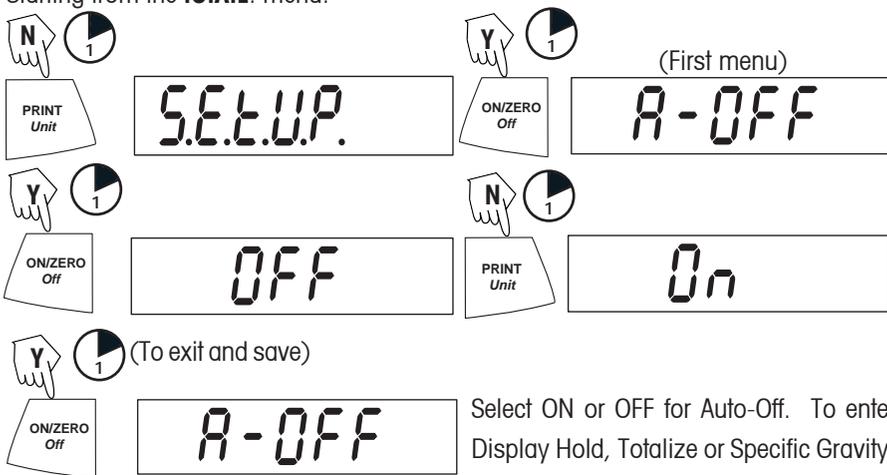


Select ON or OFF for Auto-Off. To enter Display Hold, Totalize or Specific Gravity, you must first enter the Mode sub-menu.

Turning Display Hold, Totalize or Specific Gravity Mode On

NOTE: Only one mode can be active at one time.


Start at the .S.E.T.U.P. menu.

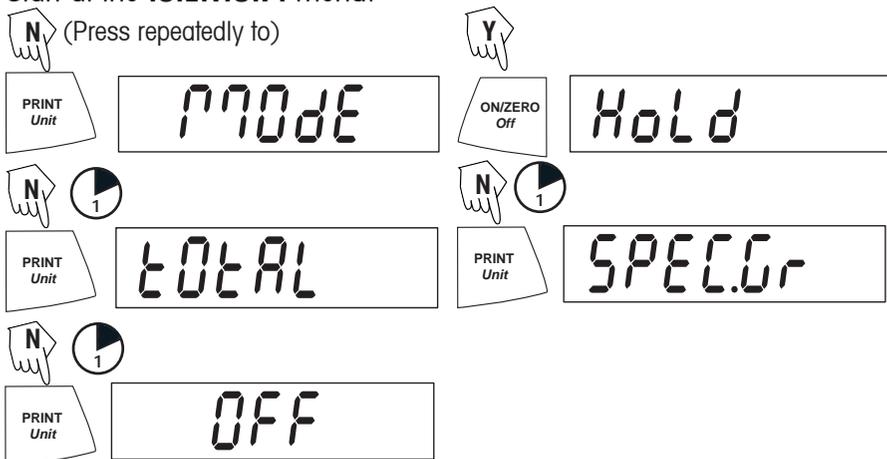

Press YES to desired mode, then continue.

Exiting the **.S.E.T.U.P.** Menu


Select ON/OFF or YES/NO to desired menu items, proceed to **.E.N.D.** menu.

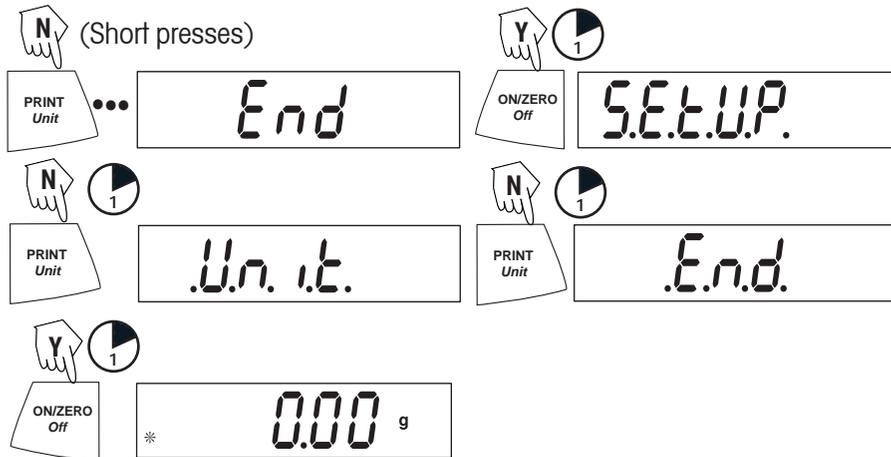

Entering the **.U.N.I.T.** Menu

The **.U.N.I.T.** menu contains units of measure, PC (parts counting), % weighing and END. Units vary with the model type. Determine which units are to be turned on or off.


Start in the **.U.N.I.T.** menu. Select either ON or OFF for each unit.

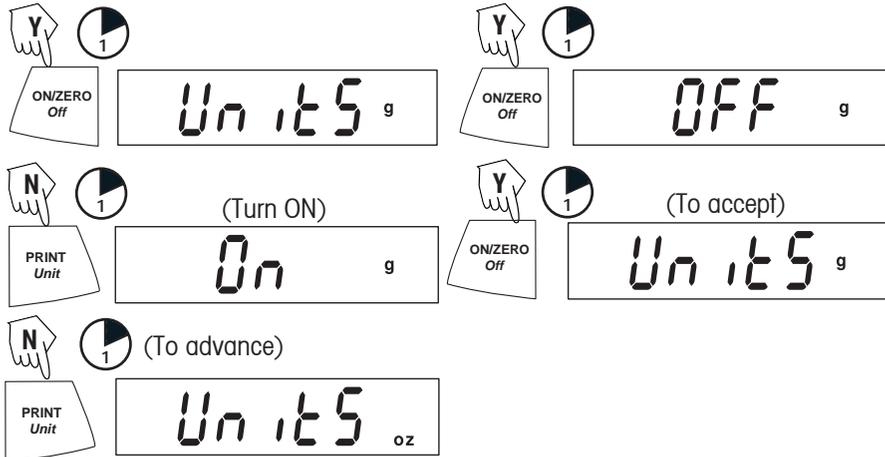

NOTE: Repeated presses of **PRINT Unit** button will go through all units, you then may select ON or OFF. Parts Counting is slightly different.

Parts Counting

Two types of counting modes are available, standard or enhanced.

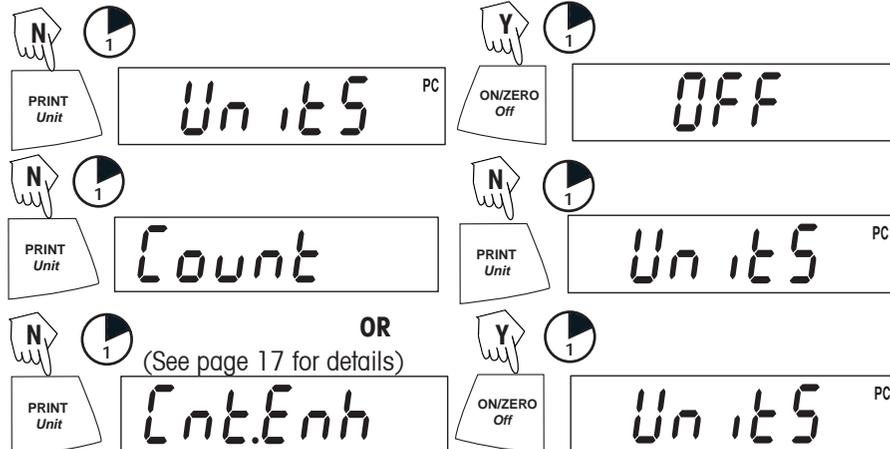

Exiting the .U.N.I.T. Menu

Use the same procedure as Exiting the .S.E.T.U.P. Menu.


Calibration

Span Calibration


Span calibration uses two calibration points, zero and a specified calibration weight. Before beginning calibration, make sure the Lock Switch is off. Clear the platform.

NOTE: Value of calibration mass depends on capacity of balance. After calibration, the balance returns to the currently selected weigh mode.

Start with the balance OFF.

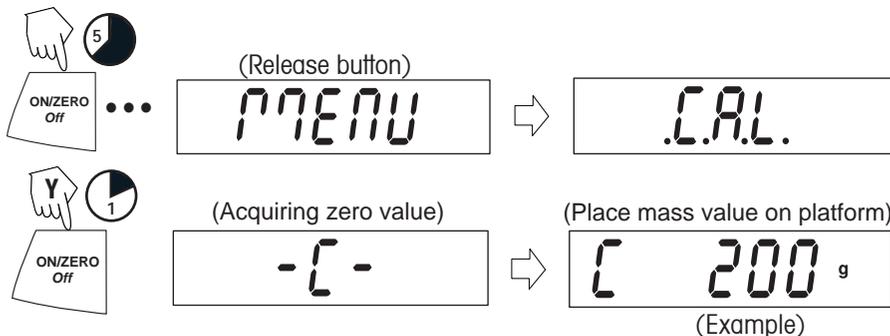


Span Calibration (Cont.)

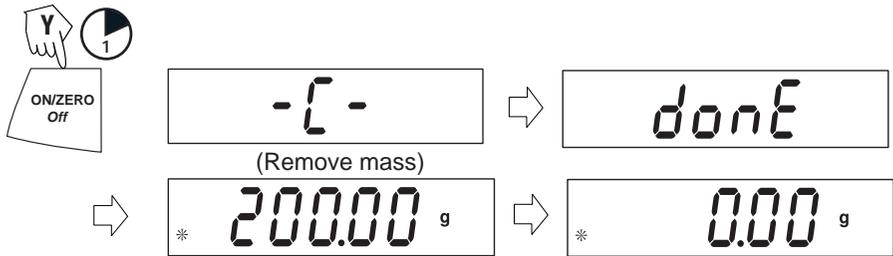

Linearity Calibration

Linearity calibration uses three calibration points; zero, mid-scale and full scale. Lin Cal must be selected and set to YES in the **.S.E.T.U.P.** Menu. Before beginning calibration, make sure the menu Lock Switch is off. Clear the platform. Start with the balance OFF.

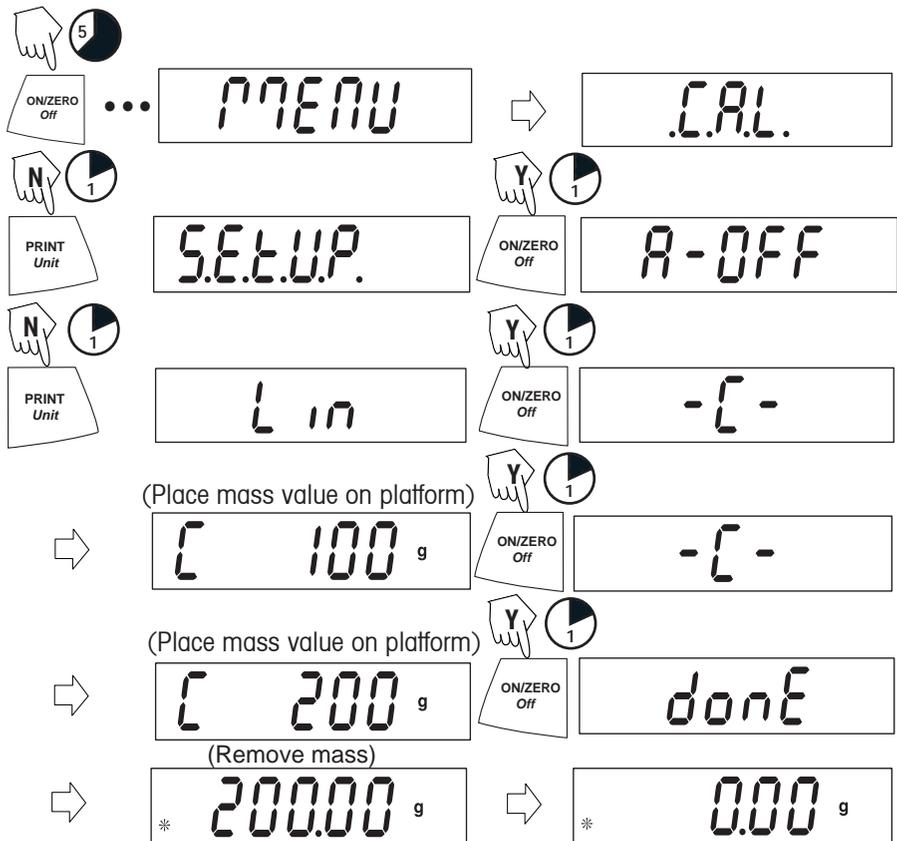


Applications

Scout Pro applications include: Weighing, Parts Counting, Percent Weighing, Display Hold, Totalize and Specific Gravity.


Weighing

Start with the balance on.

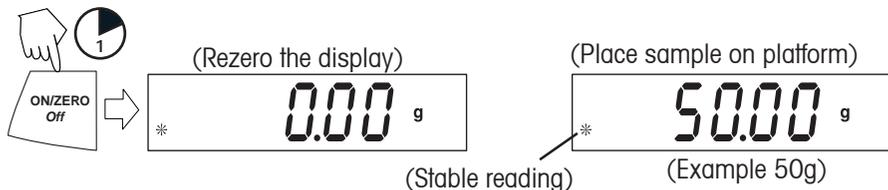

Weighing with Tare

Taring zeros the container weight. With the balance on, place an empty container on the platform. (Display example indicates a container weight of 23g.)

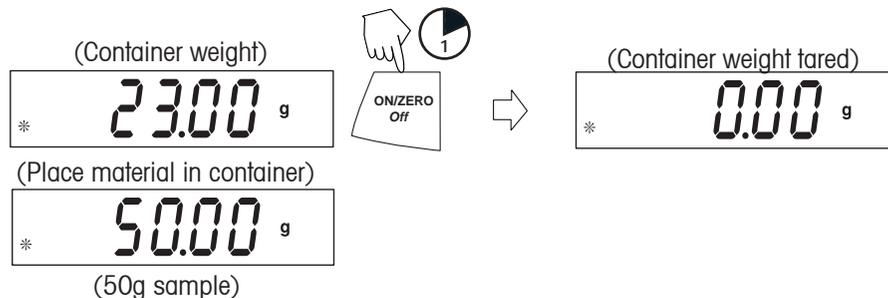

NOTE: Removing the container and material from the platform will cause the balance to display the container's weight as a negative number. The tared weight remains until **ON/ZERO Off** button is pressed again or the balance is turned off.


(Remove container with material)




Parts Counting

Parts Counting is *enabled only* when PC is turned ON in the **.U.N.I.T.** menu. In parts counting mode, there are two modes of parts counting, normal and enhanced. In normal parts counting, the balance determines the quantity based on the average weight of the parts in the original reference quantity. In the enhanced mode, additional parts can be added to the platform equal to or less than the original number. The additional reference quantity produces a more accurate average part weight.

Standard Parts Counting

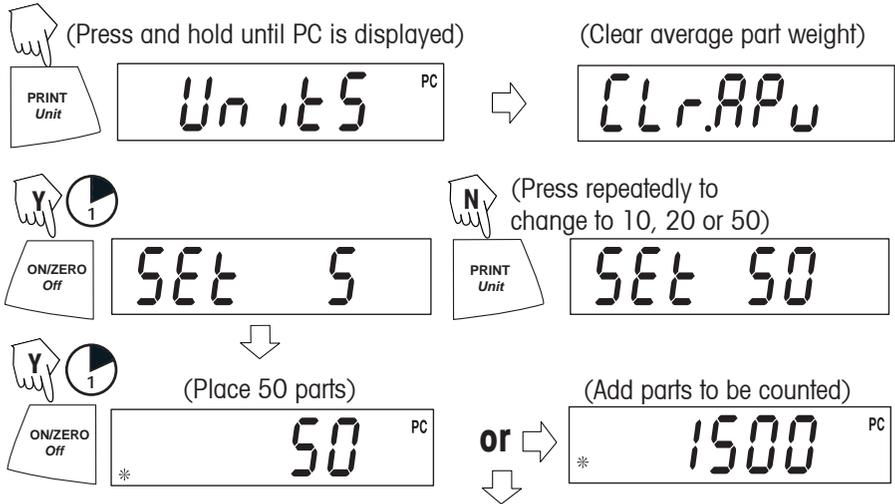


Enhanced Parts Counting

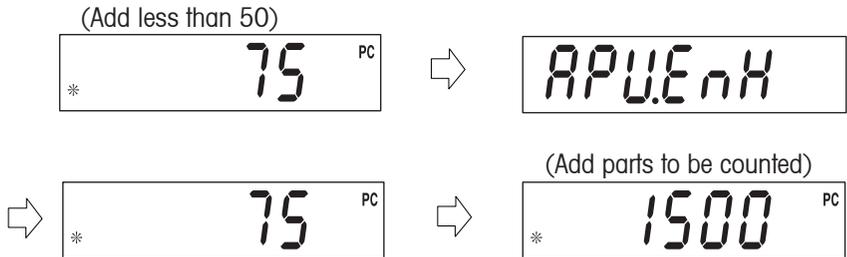

NOTE: The preceding procedure for enhanced counting can be repeated as many times as necessary providing the quantity added is less than the original entry.


Parts Counting (Cont.)

To count different parts, press and hold, until PC reappears, release Unit button.

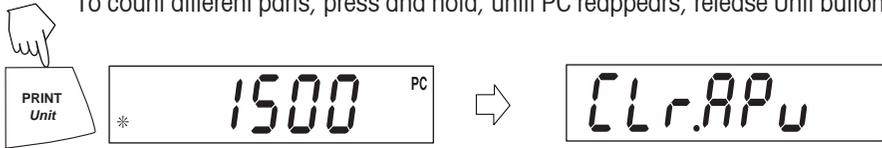

Percent Weighing

Percent Weighing is *enabled only* when Percent is turned ON in the **.U.N.I.T.** menu. Percent weighing permits placing a reference weight on the balance, then viewing other loads as a percentage of the reference. The reference weight equals 100%. Start in the weighing mode and zero the display.

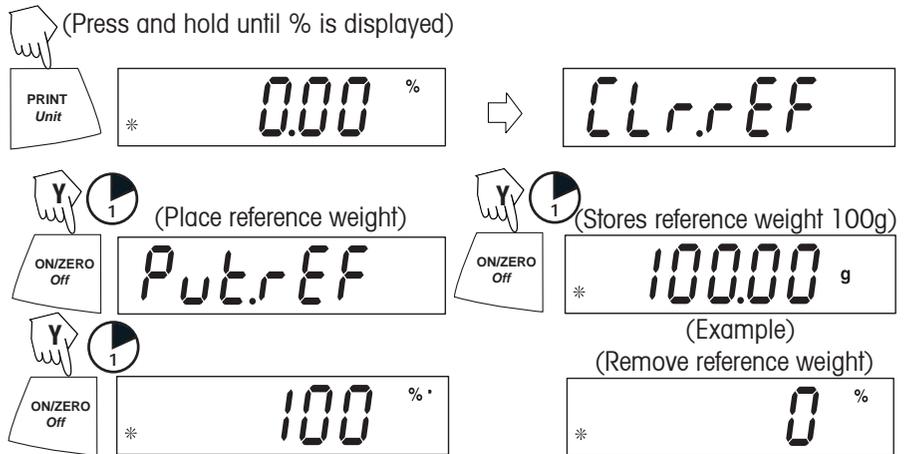


(Place load on the platform, display indicates percentage of reference weight.)



Establishing a New Reference Weight

 Press and hold until % on the display reappears, then release.

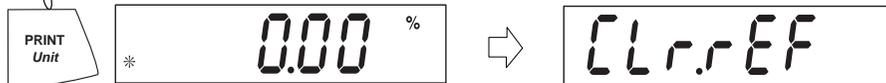

Repeat above procedure for new a reference weight.


Exiting Percent Weighing

 Press and hold until desired unit is displayed.



Display-Hold

Display-Hold is *enabled only* when Hold is turned ON in the Mode submenu in the **.S.E.T.U.P.** menu. Display-Hold mode captures and stores the highest stable value. When displayed, the stable icon will blink. **NOTE:** Units cannot be changed when in Display-Hold mode.

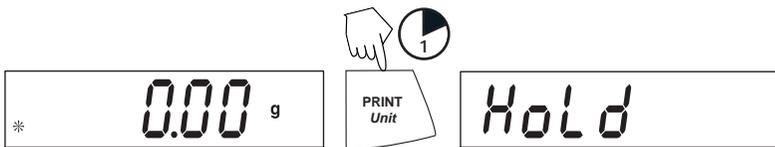

Place item(s) on platform.



Exiting Display Hold



To return to display hold, repeat above procedure.

Totalize

Totalize is *enabled only* when Total is turned ON in the Mode submenu in the **.S.E.T.U.P.** menu. Totalize allows storage of a series of weight measurements. Totalize mode has been initiated when **T** and the current unit, i.e. (g) is displayed. When totalized weight is shown, the **T** indicator will blink. **NOTE:** Units cannot be changed when in Display-Hold mode.

(Place 1st weight on platform)

(Remove weight)

(Place 2nd weight on platform)

(Total weight shown)

(Place 3rd weight on platform)


(Total weight shown)

(Remove weight)

Total weight will remain on the display until weight is removed. The total weight remains in memory. Total is limited to 999999.

Clear/Exit Totalize

Performing this next step will erase all totalized memory.

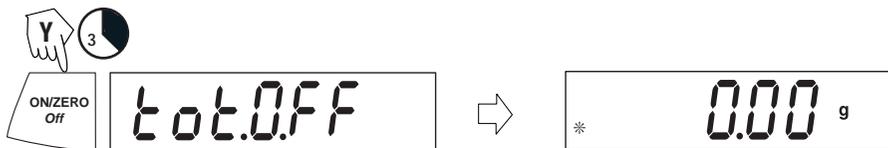


Specific Gravity (on certain models)


Specific gravity is *enabled only* when SPEC.Gr is turned ON in the Mode submenu in the **.S.E.T.U.P.** menu. The Specific Gravity mode allows calculating the specific gravity of a sample. **SPEC.GR mode must be turned ON before removing platform and installing hook.** Prepare the balance as shown. Remove the balance platform, insert the hook and install the beaker stand. Beaker is not supplied.

The sample is weighed in air suspended from the hook and then weighed in water. The formula is

$$\frac{\text{dry wt}}{\text{dry wt} - \text{wet wt}}$$

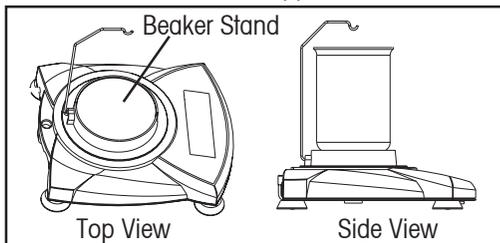

Specific Gravity Setups

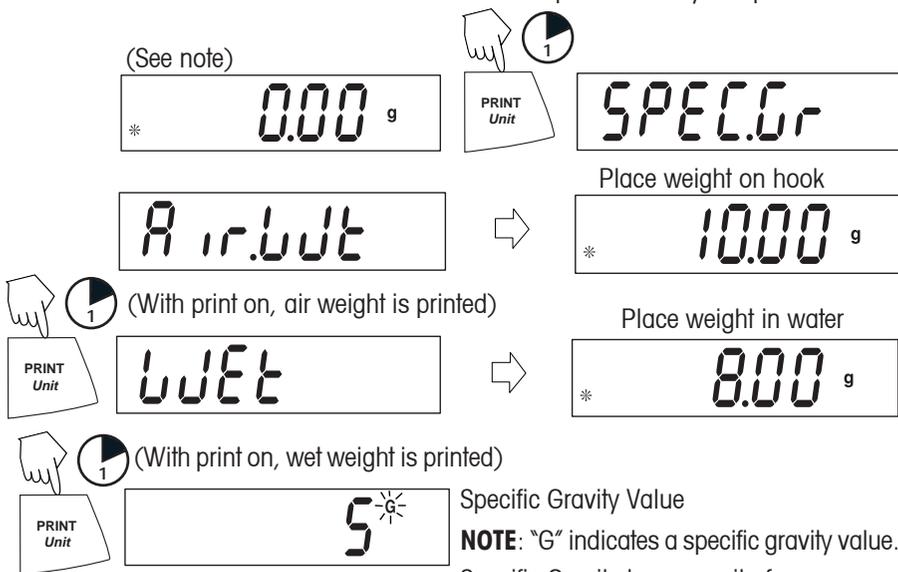


Specific Gravity Value

NOTE: "G" indicates a specific gravity value. Specific Gravity has no unit of measure.

Clear/Exit Specific Gravity


**Additional Features**Weigh Below

Scout Pro permits below balance weighing for applications such as specific gravity/density determination. The balance is normally elevated, supported on all feet and leveled. A fine wire is attached to the built-in hook at the bottom of the balance (see illustration).

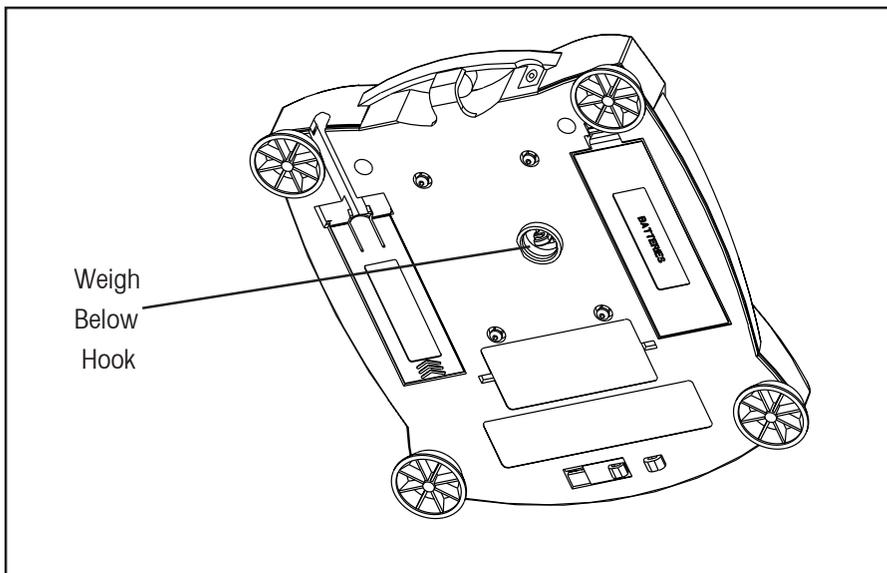

Lock Switch

The Lock Switch is used to lock out the menu mode preventing unauthorized changes. To lock the menus, position the Lock Switch next to the tab on the bottom of the balance as shown.

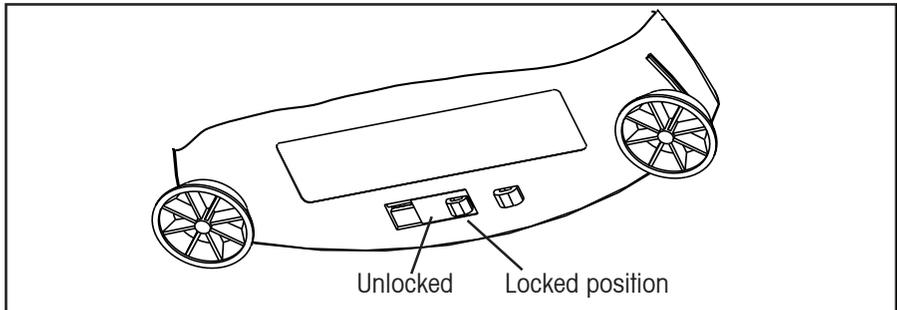

Lock Switch

4. MAINTENANCE

Cleaning

To ensure proper balance operation, keep the housing and platform clean. If necessary, a cloth dampened with a mild detergent may be used. Check under the platform for debris and remove. Keep calibration masses in a safe dry place. Unplug the AC Adapter when not in use. For long term storage, remove the batteries.

Troubleshooting

SYMPTOM	PROBABLE CAUSE	REMEDY
No display.	<ol style="list-style-type: none"> 1. Power Adapter not connected. 2. Batteries are exhausted. 	<ol style="list-style-type: none"> 1. Connect AC Adapter. 2. Replace batteries.
Battery Indicator is flashing.	Batteries are weak.	Replace batteries.
Incorrect weight reading.	<ol style="list-style-type: none"> 1. Balance out of calibration. 2. Balance was not rezeroed before weighing. 	<ol style="list-style-type: none"> 1. Calibrate the balance. 2. Press the ON/ZERO Off button with no weight on the platform, then weigh item.
Calibration procedures do not work.	Incorrect calibration masses being used.	<ol style="list-style-type: none"> 1. Use correct masses. 2. See error codes note ERR4.
Unable to display weight in a particular weighing unit.	<ol style="list-style-type: none"> 1. Weighing unit not activated in menu. 2. Mode prevents unit change. 	<ol style="list-style-type: none"> 1. Use Units menu to set desired units ON. 2. Exit Hold or Totalize mode.
Balance readings unstable.	Balance location may have drafts from air conditioning vents or vibration from other appliances nearby affecting operation.	Either move or shield the balance from external air currents or vibration.
Error code is displayed.	Various internal and external problems	Review error code table and take appropriate action.

Error Codes List

When internal and some external problems occur with the balance, the internal software will display messages as error codes. Explanations are given for possible problems. Try remedies in order indicated.

Error Codes

- Err 1 chEct** Invalid checksum data
- a) Cycle balance ON/OFF,
 - b) Return balance for servicing.
- Err 2 LoAd** Overload (>cap+90d) or Underload
- a) Remove load,
 - b) Check shipping lock,
 - c) Recalibrate.
- Err 4 WEIght** Incorrect calibration weight
- a) Use correct calibration weights,
 - b) Check environmental conditions.
- Err 5 999999** Displayed value >99999
- a) Exit application,
 - b) Rezero balance.
- Err 6 count** Parts counting error - balance exits parts counting
- a) Exit application,
 - b) Rezero balance.
- Err 8 FULL** RS232 buffer is full
- a) Cycle ON/OFF,
 - b) Check RS-232 settings in balance and computer.
- Err 9 dAtA** Internal data error. The next button press or a 5 second time-out causes the balance to turn off
- a) Cycle ON/OFF,
 - b) Return balance for servicing.

Accessories

	Part No.
RS232 Interface Kit (includes cable)	71147376
USB Interface Kit (includes cable)	71147377
Security Device	76288-01
Hard Shell Carrying Case	77256-01
Impact Printer	SF42
Impact Printer Paper 5 pack	78204-01
Cable for SF42 Printer	AS017-06
Scoops: Aluminum,	
3.62 x 4.50 x 1.0"/9.20 x 11.34 x 2.54 cm	4590-10
Black anodized, aluminum,	
3.62 x 4.50 x 1.0"/9.20 x 11.34 x 2.54 cm	4590-30
Aluminum,	
1.5 x 2.00 x 0.43"/3.81 x 5.08 x 1.11 cm	5076-00
Gold anodized aluminum,	
2.25 x 3.00 x 0.75"/5.71 x 7.62 x 1.90 cm	5077-00
Calibration Masses:	
See specification table for required masses.	
50g	53054-00
100g	51015-05
200g	51025-06
300g	51035-05
500g	51055-06
1000g	51016-06
2000g	51026-02
AC Adapters:	
120V/60Hz USA	12102320
230V/50Hz Europe	12102321
230V/60Hz Australian	12102323
230V/50Hz UK	12102322
100V/50Hz Japan	12102324
230V/60 Hz China	12104881

5. TECHNICAL DATA

Specifications

Item No.	SPx123	SPx202	SPx402	SPx401	SPx601	SPx2001	SPx4001	SPx6000
Capacity (g)	120	200	400	400	600	2000	4000	6000
Span Calibration Mass (g)	100	200	400	400	600	2000	4000	6000
Linearity Calibration Mass (g)	50	100	200	200	300	1000	2000	3000
Readability (g)	0.001	0.01		0.1				1.0
Repeatability (Std. dev.)(g)	0.003	0.01		0.1				1.0
Linearity (g)	±0.003	±0.01		±0.1				±1.0
Weighing modes	g, N, oz, %, Parts Counting				g, kg, N, oz, lb, %, Parts Counting			
Tare range	To capacity by subtraction							
Over range capacity	Capacity +90d							
Stabilization time	3 seconds							
Operating temp. range	50° - 104°F / 10° - 40°C							
Power requirements	AC Adapter (supplied) or 4 AA batteries (not included)							
Calibration	digital calibration from keypad							
Display (in/mm)	LCD (0.6 / 15 high digits)							
Pan size (in/mm)	4.7 / 120 diam.				6.5 x 5.6 / 165 x 142			
Dimensions WxHxD (in/mm)	7.5 x 2.2 x 8.3 / 192 x 54 x 210							
Net Weight (lb/kg)	2.0 / 0.9			3.5 / 1.6				

Capacity x Readability

	SPx123	SPx202	SPx402	SPx401	SPx601	SPx2001	SPx4001	SPx6000
gram g	120.00 x 0.001	200.00 x 0.01	400.00 x 0.01	400.0 x 0.1	600.0 x 0.1	2000.0 x 0.1	4000.0 x 0.1	6000 x 1
ounce oz avoirdupois	4.23290 x 0.00005	7.0550 x 0.0005	14.1095 x 0.0005	14.110 x 0.005	21.165 x 0.005	70.550 x 0.005	141.095 x 0.005	211.65 x 0.05
pound lb avoirdupois					1.3230 x 0.0005	4.4090 x 0.0005	8.8185 x 0.005	13.230 x 0.005
kilogram kg						2.0000 x 0.0001	4.0000 x 0.0001	6.000 x 0.001
newton N	1.17680 x 0.00001	1.9613 x 0.0001	3.9227 x 0.0001	3.923 x 0.001	5.884 x 0.001	19.613 x 0.001	39.227 x 0.001	58.84 x 0.01

NOTE: SPGxxxx models only contain SI units.

LIMITED WARRANTY

Ohaus products are warranted against defects in materials and workmanship from the date of delivery through the duration of the warranty period. During the warranty period Ohaus will repair, or, at its option, replace any component(s) that proves to be defective at no charge, provided that the product is returned, freight prepaid, to Ohaus.

This warranty does not apply if the product has been damaged by accident or misuse, exposed to radioactive or corrosive materials, has foreign material penetrating to the inside of the product, or as a result of service or modification by other than Ohaus. In lieu of a properly returned warranty registration card, the warranty period shall begin on the date of shipment to the authorized dealer. No other express or implied warranty is given by Ohaus Corporation. Ohaus Corporation shall not be liable for any consequential damages.

As warranty legislation differs from state to state and country to country, please contact Ohaus or your local Ohaus dealer for further details.

