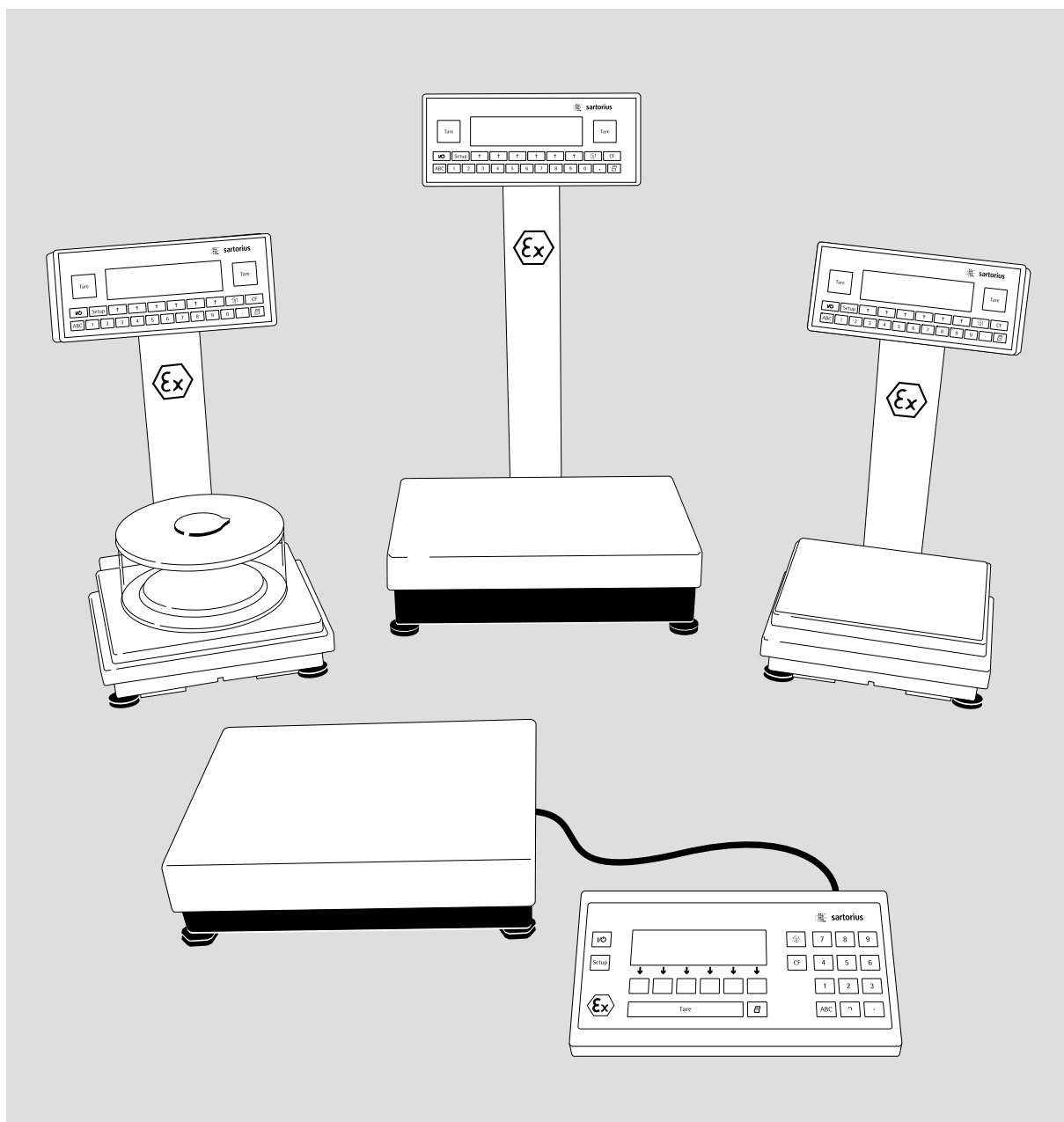


Operating Instructions

Sartorius Factory Series

FC and FCA Models

Electronic Precision Scales for Hazardous Areas/Locations



Contents

2	Contents	115	Error Codes
3	Operating Design	117	Overview
7	Configuration	117	General Views of the Scales
7	Setting the Language	121	Specifications
8	Navigating in the Setup Menu	127	Dimensions (Scale Drawings)
9	Entering the Time and Date	131	Accessories (Options)
10	Setting the Scale Functions	136	Declarations of Conformity
13	Setting the Device Parameters	138	EC Type Examination Certificates
13	Entering a Password	162	Certificates of Compliance
17	Setting the Application Parameters	168	EC Type-approval Certificate
23	Selecting the Printout Function	170	Plates and Markings
25	Printout Configuration	172	Index
28	Displaying Info		
29	MP8 Interface Emulation		
29	Factory Settings		
31	Operation		Appendix
31	Basic Weighing Function		Entering the General Password
34	Device Parameters		
36	Calibration/Adjustment		
45	ReproTEST		
46	Application Programs		
46	Toggle between Two Weight Units		
48	Counting		
51	Weighing in Percent		
54	Animal Weighing		
58	Recalculation		
63	Calculation		
65	Checkweighing		
70	Time-Controlled Functions		
72	Totalizing		
76	Formulation		
80	Statistics		
84	Extra Functions (in the Application Menu)		
84	Second Tare Memory (Preset Tare)		
86	Individual Identification Codes (ID)		
90	Manual Storage (M+)		
91	Product Data Memory		
93	“FlexPrint” Printout Function		
94	Combining Applications		
95	Examples of Application Combinations		
97	Data Output Functions		
99	Printouts		
101	Interface Description		
106	Pin Assignment Chart		
107	Cabling Diagram		

Symbols

The following symbols are used in these instructions:

- indicates steps you must perform
- indicates steps required only under certain conditions
- > describes what happens after you have performed a certain step
- indicates an item in a list
- ⚠ indicates a hazard

Important Note to Users

Make sure to carefully read and follow sections marked with this symbol – they contain important safety instructions.

Hotline

For advice on the use of these applications, just call or fax
Telephone: +49.551.308.4440
Telefax: +49.551.308.4449

Operating Design


The scales in the Factory Series consist of a weighing cell and a display and control unit. Besides an AC adapter, your scale also has an interface port for connecting a printer, computer, universal remote control switch, etc.

The display and control unit and the weighing cell can be set up separately.

Where not expressly indicated otherwise, the uses described in this manual apply to verified and verifiable scale versions* (indicated by the suffix "...CE" in the model number), as well as the standard version.

Combining Applications

You can combine the use of various application programs to meet your more complicated requirements.

Press the  key to select the desired application programs.




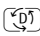

Keys

You can operate the scale either using the keys on the terminal or from a connected PC. This manual describes operation using the keys on the terminal.

Labeled Keys

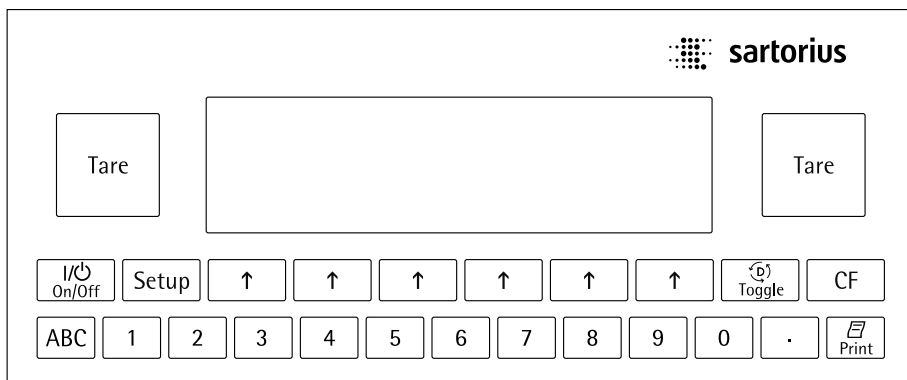
These keys always have the function indicated by their label, but are not available at all times. Availability of their functions depends on the current operating status of the scale and the menu settings.

Meaning


-  **Alphabetic keys**
Please see section on "Text Input"
-  **On/off/standby key**
Turns the scale on and off or switches it to the standby mode
-  **Menu settings**
Accesses and exits the Setup
-  **Toggles to the next application program**
-  **Clear Function**
Deletes keypad input
Interrupts a calibration and adjustment routine in progress
Quits application programs


* For the use in legal metrology in the EU and European Economic Area


FC Models:




FCA Models:

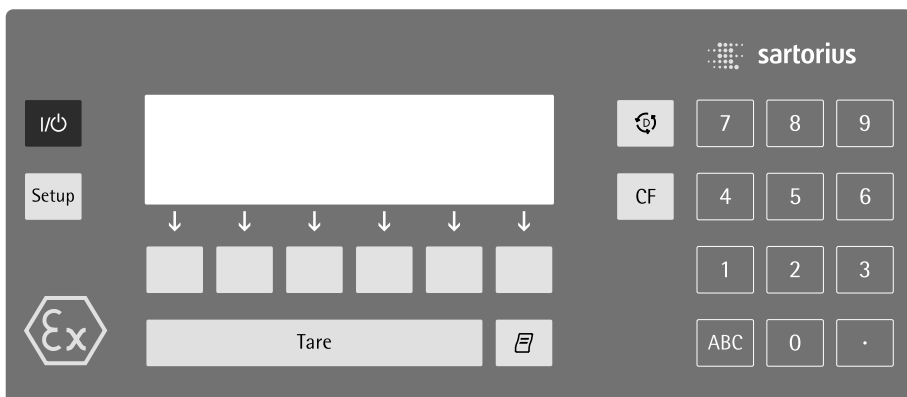

 **Print key**
Outputs displayed values or data logs to the interface port





 **Enters a decimal point**

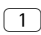 **1 ... 9 0 keys**
See the section on "Numeric Input"

 **Tares the scale**

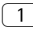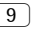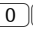
Numeric Input

To enter numbers:

press  ...   


To store numbers entered:

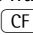press the corresponding function key directly below the soft key label

To delete an entire numeric input digit by digit: press the  key

Text Input

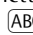● To enter numbers:
see the section on "Numeric Input"


● To enter letters or characters:
press the  key



> Letters are displayed in the bottom line for selection


● To select a different letter:
press the corresponding soft key to change the letter shown



● To select the letter/character shown:
press the corresponding function key below the soft key label
> The selected letter is shown on the display

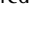○ Enter the next letter/character, if desired, as described above

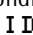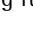○ To exit the letter input mode e.g., if the last character entered is a letter):
press the  key

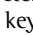● To store a word: press the corresponding function key (soft key), such as  

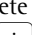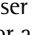● To delete an input character by character: press the  key

● To delete user data:
enter  or a space  and save

Operating Design

Keys

Your Factory scale is operated either through the keys on the display and control unit or via a connected PC. Operation through the scale keys is described in the following.

Function Keys (Soft Keys)

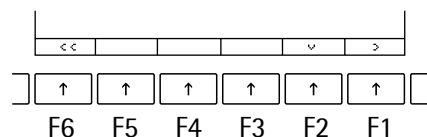
The current function of a soft key is indicated in the bottom line of the display.

The bottom line shows symbols and/or abbreviations to indicate soft key functions. The abbreviations are usually self-explanatory.

Example

i s o T S T: start calibration/adjustment
S I D: save ID

The function keys are numbered F1 through F6, from right to left.



The arrows shown in the footer indicate the following functions:

- < <** Return to Setup menu (in the Setup menu: save settings and exit the Setup program)
- <** Go back to the higher selection level
- >** Show sub-items under the active item
- ^** Move upward in the input/output window
- v** Move downward in the input/output window
- ↓** Set the selected menu parameter

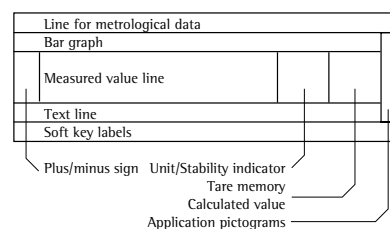
Display

There are two fundamentally different types of display:

- display of measured and calculated values
- display for menu parameter settings (setup)

Display of Measured and Calculated Values

This display is divided into nine sections.



Line for Metrological Data:

If the scale is verified for use in legal metrology, the following metrological specifications are shown here:

- Max** Maximum capacity of the scale
- Min** Minimum capacity of the scale; i.e., the minimum weight allowed when the scale is used in legal metrology
- e** Verification scale interval of the scale
- d** Readability: indicates the scale interval of the scale
- R 1** Displayed when $e = d$ through
- R 4** $e = d$

On standard scales, only **Max** and **d** are shown.

Bar Graph:

The bar graph indicates how much of the scale's capacity is "used up" by the current load; during checkweighing, **t** indicates the control limits.

The following symbols may be displayed here:

- 0%** Lower load limit
- 100%** Upper load limit
- Bar graph showing 10% intervals
- Minimum for checkweighing
- =** Target for checkweighing
- +** Maximum for checkweighing

Plus/Minus Sign, Stability Symbol:

A plus or minus sign (**+** or **-**) is shown here for a weight value (e.g., a calculated value when weighing in percent) or the **Q** symbol, indicating that the verified or verifiable scale has been zeroed or tared.

Measured Value Line:

This section shows the weighed or calculated value or alphanumeric input. Note Concerning Verified Scales Approved for Use as Legal Measuring Instruments in the EU*:

For verified scales that have a verification scale interval **e** not equal to the scale interval **d**, the last digit on the display is bordered

Unit and Stability:


When the scale reaches stability, the weight unit or calculation unit is displayed here.

The symbol may be displayed for readouts on a scale verified for legal metrology. However, these readouts can be used only for standard applications (not in legal metrology/not legal for trade).

* including the Signatories of the Agreement on the European Economic Area

Tare Memory, Calculated Values:
The symbols displayed here indicate when there is a value in one of the tare memory modules or when the value shown is a result of calculation rather than direct measurement.

These symbols are as follows:


 Calculated value


NET 1 Tare memory used by an


NET 2 Application program
(e.g., formulation, second tare)

Application Pictograms:
The pictograms displayed here indicate the application(s) selected. The pictogram is displayed inversely when the corresponding application is active.

For example, the following symbols may be displayed simultaneously:

 The counting application is active

 Checkweighing is also active

 Print

 Data record

Text Line:
Additional information is displayed here (e.g., operator guidance prompts, name of the active program, etc.).

Soft Key Labels:
The current functions of the soft keys are indicated here; during calibration/adjustment, this line shows up- and down-arrows (↕ and ↗) for selecting calibration and adjustment functions.

Display for Menu Parameter Settings (Setup)

This display is divided into three sections.

Header
Input and Output Window
Footer

Header

The header indicates the function of the current screen page. In the Setup program, the current menu path is shown here.

Example in the path "Setup/Menu":

SETUP	MENU	[]

Input and Output Window

This window contains either detailed information (e.g., on the active application) or a pick list. A selected item is displayed inversely. You can also enter information in an active field in this window using the alphanumeric keys.

Example in the path "Setup/Menu":

1	Minimum vibration
o 2	Normal vibration
3	Strong vibration
4	Extreme vibration

The o symbol in this window indicates the current menu setting.

Soft Key ID Labels:
see "Function Keys (Soft Keys)" on the previous page.

To select a parameter:

- Press ↕ or ↗ repeatedly until the desired setting is selected (displayed inversely)

- Confirm your selection by pressing the ↵ soft key

To change the alphanumeric value of a parameter:

- Press ↕ or ↗ repeatedly until the desired setting is selected (displayed inversely)
- Enter the desired character using the 0 1 ... 9 . keys or press the ABC key and enter letters or characters
- Confirm your selection by pressing the ↵ soft key
- To return to the Setup/Menu level:
- Press the << soft key

Operating Design

Input

Numeric Input

To enter numbers:

Press the **1** **2** ... **0** **.** keys

To store numbers entered:

Press the corresponding soft key (i.e., the arrow key under the appropriate abbreviation in the bottom line of the display)

To interrupt/cancel numeric input:

Press **CF**

Alphabetic Input

(see also the example given on page 33)

- To enter letters or characters:
first press the **ABC** key
> Letters are displayed in the bottom line
- To select a different letter:
press the corresponding soft key to change the letter shown (i.e., the arrow key under the letter displayed)
- To select the letter/character shown:
press the corresponding soft key
> The selected letter is shown in the display
- Enter the next letter/character, if desired, as above.
- To store a word: press the corresponding soft key (e.g., **S** **ID**)
- To delete a word: press **CF**

Parameter Settings

The parameters for configuration are in the application menu and the scale operating menu. These menus have several levels.

- To set parameters: press **Setup** and then the appropriate soft key (e.g., **A** **▢** **▢** for the application menu)
- To move within a menu level:
use the **↶** and **↷** soft keys
- To select a parameter:
 - Press **↶** or **↷** repeatedly until the desired setting is selected (displayed inversely)
 - Confirm your selection by pressing the **↵** soft key

To change the numeric value of a parameter:

- Press **↶** or **↷** repeatedly until the desired setting is selected (displayed inversely)
- Enter the desired number using the **1** **2** ... **0** **.** keys
- Confirm your selection by pressing the **↵** soft key
- To return to the Setup/Select level:
 - Press the **↶** **↶** soft key

See the chapter entitled “Configuration” for a complete description of all parameters.

To save the parameter settings and exit the Setup menu: press **Setup** or the **↶** **↶** soft key

To cancel the parameter setting operation: press **U/C**

Data Output

Your Factory scale is equipped with an interface port for connecting your choice of the following:

- Printer*
- Peripheral device (e.g., computer)*
- Universal remote control switch

Printer

You can configure the print functions to meet your individual requirements by selecting the corresponding menu code.

You can have printouts generated automatically, or by pressing **(E)**; dependent on or independent of the stability or time parameters; with or without IDs; and as standard or ISO/GMP-compliant printouts.

ISO: International Organization for Standardization

GMP: Good Manufacturing Practice

See the section on “Data Output Functions” in the chapter entitled “Operation” for a detailed description of data output options.

Interface Port

Instead of a printer, you may choose to connect a different peripheral device, e.g. a computer (PC)*. With an on-line PC you can control both the weighing cell and the display unit of the Factory scale.

Request messages are sent via the interface to initiate functions in the weighing cell and in the display unit. Some of the functions generate response messages.

See the chapter entitled “Operation” under the section on “Data Output” for a detailed description of the interface port.

Error Codes

If you press a key that has no function, or which is blocked at a certain point in an application program, this error is indicated as follows:

- a double-beep is sounded as an acoustic signal if the key has no function
- a double-beep is sounded and a message is displayed for 2 seconds in the text line if the key function is not available at that time

The response to an operator error is identical in all models of the Factory series. See the chapter entitled “Error Codes” for a detailed description.

Storing Settings

Storing Parameter Settings

The settings configured are stored in the scale’s non-volatile memory. Additionally, the factory settings can be reloaded.

Saving Parameter Settings

You can assign passwords in order to block access to:

- printing function **Confis**
- the application menu **A** **▢** **▢**
- the scale operating menu **Menu** and
- user data input functions **Input**

* using a Zener barrier (see also page 131)

Configuration

Purpose

You can configure the scale terminal to meet individual requirements by entering user data and setting parameters in the Setup program.

The Setup menu contains the following submenus:

- Balance/scale functions
- Device parameters
- Application parameters
- Printout
- Device information
- Language
- Factory settings

Setting the Language

You can choose from 5 languages for the information display:

- German
- English (factory setting)
- English with U.S. date/time format
- French
- Italian
- Spanish

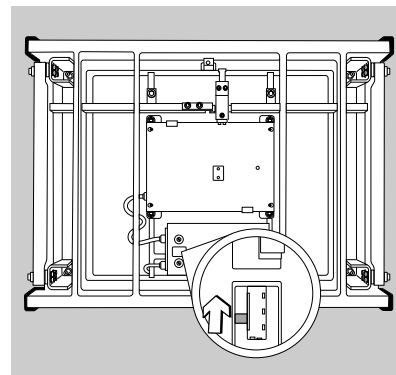
Configuring the Scale for Use in Legal Metrology

Set the menu access switch as described below to configure the following functions for use of the scale in legal metrology:

- Display: Verification scale interval **e**; lower limit of the weighing capacity **Min**
- External calibration blocked

Preparation FCA...IGG:

- Remove the load plate from the weighing platform
- Remove the cap from the back of the scale housing
- Move the switch upwards



- > When the switch is in the upper position, the Setup menu is locked and the scale can be used in legal metrology
- > When the switch is in the lower position, the menu is accessible

Preparation FC..., FCA...EDE:

- Switch located see sections entitled "Overview" and "Plates and Markings"

Example:

Setting the Language to "U.S. Mode"

Step	Press key (or follow instructions)	Display/Output
1. Select "Setup" menu	Setup	<pre> SETUP Balance/scale functions Device parameters Application parameters Printout Info << >> </pre>
2. Select "Language" and confirm	repeatedly press v soft key, then > soft key	<pre> SETUP LANGUAGE Deutsch oEnglish U.S.-Mode Français Italiano << < ^ v j </pre>
3. Select "U.S. mode"	^ soft key	<pre> SETUP LANGUAGE Deutsch oEnglish oU.S.-Mode Français Italiano << < ^ v j </pre>
4. Save language	j soft key	<pre> SETUP LANGUAGE Deutsch English oU.S.-Mode Français Italiano << < ^ v j </pre>
5. Exit the Setup menu	<< soft key	<pre> Max 6200g d= 0.01g 0% 100% 0.00 g isoTST PT1/T1 </pre>

Configuration

Navigating in the Setup Menu (Examples):
Example: Adapting the scale to “Extreme vibration”

Step	Press key(s) (or follow instructions)	Display/Output
1. Select Setup menu	Setup	<div>SETUP</div> <div>Balance/scale functions</div> <div>Device parameters</div> <div>Application parameters</div> <div>Printout</div> <div>Info</div> <div><< v ></div>
2. Confirm “Scale functions”	> soft key	<div>SETUP BAL.FUNC.</div> <div>Calibration/adjustment</div> <div>Adapt filter</div> <div>Application filter</div> <div>Stability range</div> <div>Taring</div> <div><< < v ></div>
3. Select menu item “Adapt filter” and confirm	v, then > soft key	<div>SETUP BAL.FUNC. ADAPT FILT.</div> <div>Minimum vibration</div> <div>Normal vibration</div> <div>Strong vibration</div> <div>Extreme vibration</div> <div><< < ^ v ↓</div>
4. Select menu item “Extreme vibration”	v soft key	<div>SETUP BAL.FUNC. ADAPT FILT.</div> <div>Minimum vibration</div> <div>Normal vibration</div> <div>Strong vibration</div> <div>Extreme vibration</div> <div><< < ^ ↓</div>
5. Confirm menu item “Extreme vibration”	↓ soft key	<div>SETUP BAL.FUNC. ADAPT FILT.</div> <div>Minimum vibration</div> <div>Normal vibration</div> <div>Strong vibration</div> <div>Extreme vibration</div> <div><< < ^ ↓</div>
6. If required, select further menu items	< v ^ > soft keys	
7. Save setting and exit Setup Menu	<< soft key	

Example: Entering the time and date

Step	Press key(s) (or follow instructions)	Display/Output
1. Select Setup menu; select "Device parameters"	[Setup], then ↓ soft key and → soft key	
2. Set clock	press ↓ repeatedly, then press →	
3. Enter the time	[1] [1] [.] [1] [2] [.] [3] [0]	
4. Set the time according to your local clock	↓ soft key	
5. Enter the date	[1] [3] [.] [0] [3] [.] [0] [0]	
6. Store the date	↓ soft key	
7. Enter other data, if desired	< ^ ^ > soft keys	
8. Exit Setup menu	< < soft key	

Setting the Scale Functions (BAL.FUNC.)

Purpose

This menu item enables you to configure the scale functions, i.e., to meet individual requirements by selecting predefined parameters in the Setup menu. You can block access to the menu by assigning a password.

Features

The scale functions are combined in the following groups (1st menu level):

- Calibration/adjustment
- Adapt filter
- Application filter
- Stability range
- Taring
- Auto zero
- Weight unit 1
- Zero range
- Zero range at power on
- Tare/zero at power on
- Factory settings: only wgh. param. (only the scale functions)

Factory Settings

Parameters: The factory settings are identified by the symbol "O" in the list starting on the next page.

Preparation

Show available balance/scale functions:

- Select Setup menu: press the **[Setup]** key
- > SETUP is displayed

SETUP					
Balance/scale functions					
Device parameters					
Application parameters					
Printout					
Info					
<<				v	>

- Select "Scale functions":
press the **[>]** soft key
- If you already assigned a password:
> The password prompt is displayed
- If access is blocked by a password:
enter the password using the numeric/
alphabetic keys.
- If the last character of the password is a
letter: conclude input by pressing **[ABC]**
- Confirm your password and have
the scale functions displayed:
Press the **[v]** soft key
- > Scale functions are displayed:

SETUP	BAL.FUNC.				
Calibration/adjustment					
Adapt filter					
Application filter					
Stability range					
Taring					
<<				v	>

- To select the next group:
press the **[v]** soft key (down arrow)
- To select the previous item of a group:
press the **[^]** soft key (up arrow)
- To select the next sub-item within
a group: press the **[>]** soft key
(right arrow)
- To select the previous group:
press the **[<]** soft key (left arrow)
- To confirm: press the **[v]** soft key

Extra Functions

- Exit the Setup menu:
press the **[<<]** soft key
- > Restart your application
- Print parameter settings:
 - When the balance/scale functions are displayed, press **[F]**
- > Printout (example)
Texts with more than 20 characters are cut off

SETUP

BAL.FUNC.

```
-----
Calibration/adjustm
  CAL/iso TST key fun
    Internal cal./adju
      Cal/adjustm seq
        Cal. with adjustm au
          .
            isoCAL-function
              Off
            Start autom. adjus
              isoCAL
            Print GLP/GMP adju
              Automatic if GLP is
                selected
              Parameter for exte
                Wt. ID (W ID):

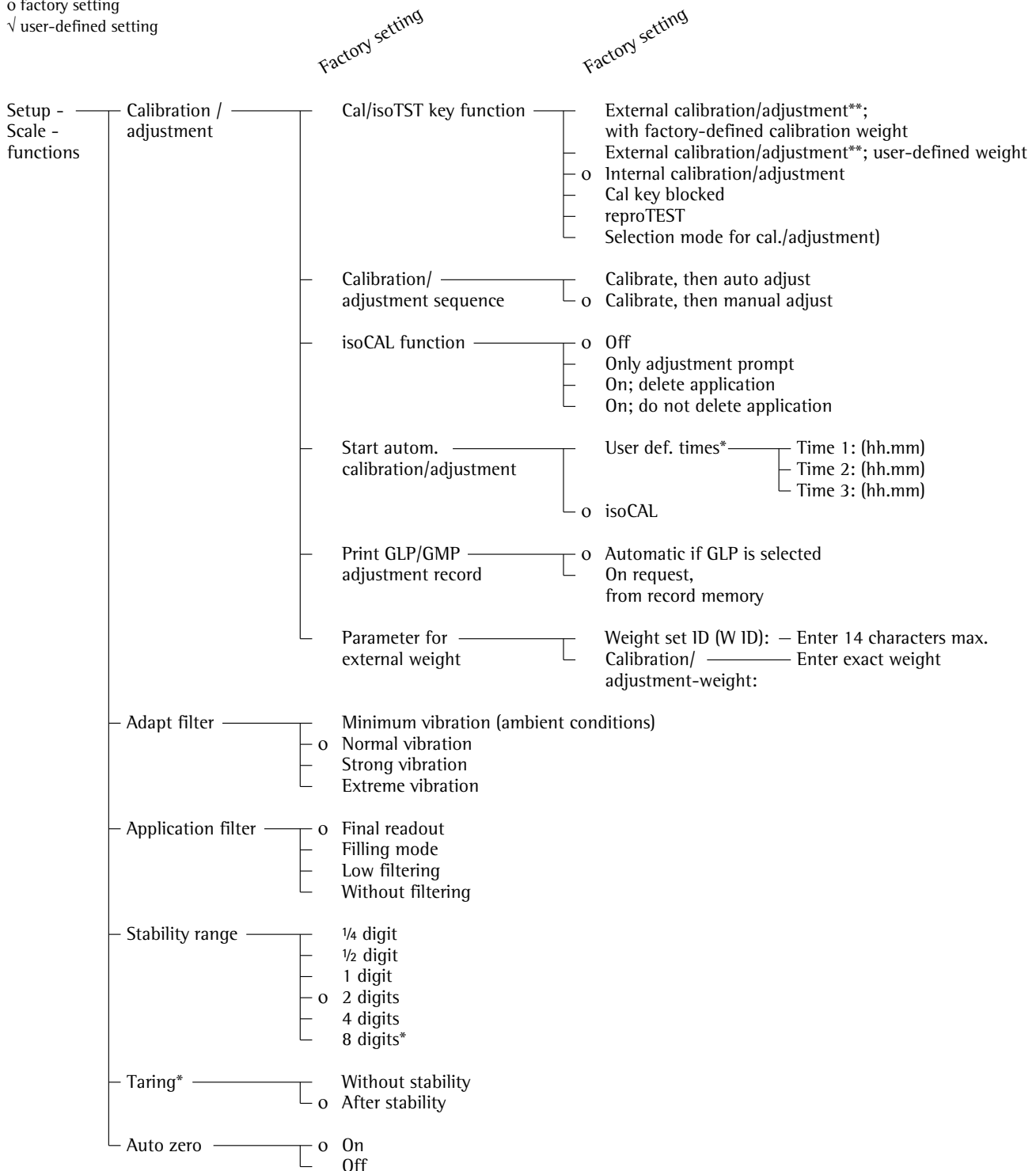
                Cal./adjust.-wt:
                  5000.00 g
            Adapt filter
              Normal vibration
            Application filter
              Filling mode
            Stability range
              2 digits
            Taring
              After stability
            Auto zero
              Off

            Weight unit 1
              Grams /g
etc.
```

Setting the Scale Functions (Overview)

o factory setting


√ user-defined setting



* = not applicable to verified scales used in legal metrology in the European Economic Area

** = verified scales can only be externally calibrated, not adjusted (you can only check the external weight)

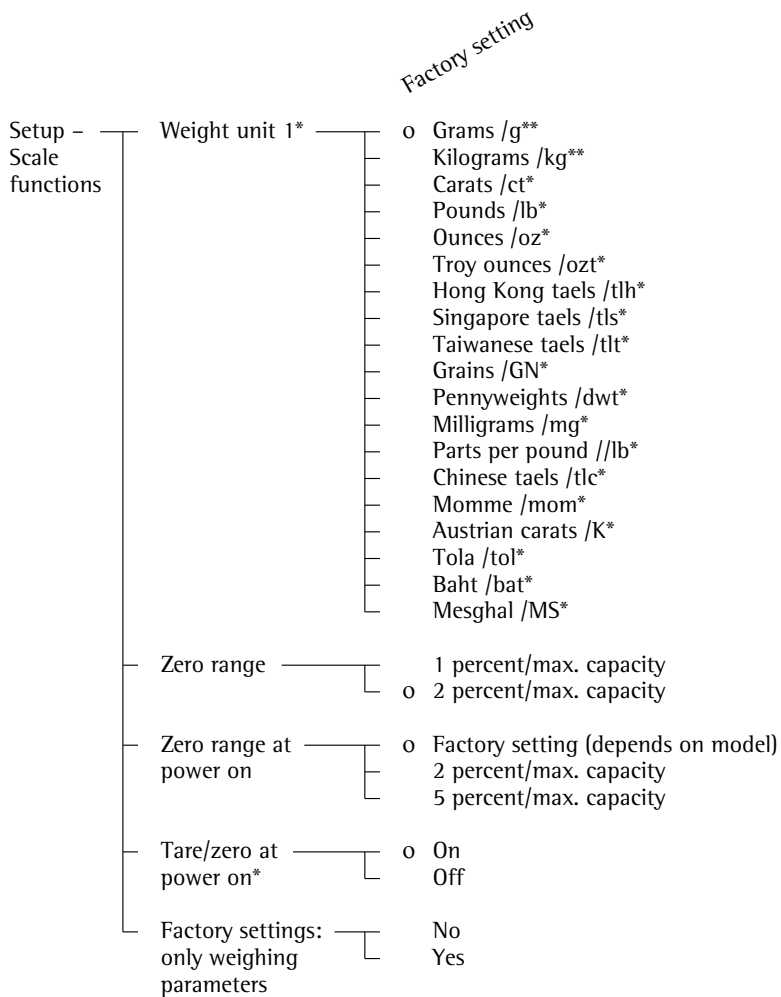Setting the Scale Functions (Bal Func.)



* = not applicable to verified scales used in legal metrology in the European Economic Area

** = factory setting depends on weighing range: – ≤ 33 kg: grams
– ≥ 34 kg: kilograms

Setting the Device Parameters (DEVICE)

Purpose

This menu item enables you to configure the scale to meet individual requirements by selecting predefined menu parameters in the Setup menu. You can block access to the menu by assigning a password.

Features

The device parameters are combined in the following groups (1st menu level):

- Password
- User ID
- Clock
- Interfaces
- Display
- Keys
- Extra functions
- Factory settings: only device parameters

Factory Settings

Parameters: The factory settings are identified by the symbol “o” in the list starting on the page after next.

Preparation

Display available device parameters

- Select the Setup menu: press **Setup**

> SETUP is displayed:

SETUP					
Balance/scale functions					
Device parameters					
Application parameters					
Printout					
Info					
<<					>>

- Select “Device parameters”:
use the **↓** and **↑** soft keys

If no password has been assigned, anyone can access the Setup menu device parameters.

If a password has already been assigned:

- > The password prompt is displayed
- If access is blocked by a password:
enter the password using the numeric and/or alphabetic keys
- If the last character of the password is a letter: conclude input by pressing the **ABC** key
- Press **↓** to confirm the password

> Device parameters are now displayed:

SETUP					
DEVICE					
Password					
User ID					
Clock					
Interfaces					
Display					
<<					>>

- To select the next group:
press the **↓** soft key (down arrow)
- To select the previous menu item of a group: press **↑** soft key (up arrow)
- To select the next sub-item within a group: press the **→** soft key (right arrow)
- To select the previous group:
press the **←** soft key (left arrow)
- Press **↓** soft key to confirm the selected menu item

Entering or Changing a Password

- Let's assume that a password with 8 characters max. has already been assigned to access the Setup device parameters

- Select the Setup menu: press **Setup**

> SETUP is displayed

- Select parameters:
Use the **↓** and **→** soft keys

> The password prompt is displayed:

SETUP					
PASSW. CHECK					
Enter password:					
<<					>>

- Enter the password

- Press the **↓** soft key to confirm your password and view the device parameters

- Write down your password here for easy reference:
Password =.....

If you assign a password and then forget what the word is:

- Enter the General Password (see Appendix)

- Press the **↓** soft key to confirm and display the password

> The parameters are displayed

- Select the device parameter “Password”:
If necessary, repeatedly press **↓** or **↑**, until you see

> **Password:** and any existing password

SETUP					
DEVICE					
PASSWORD					
Password:					
ABC123					
ESC					↓

Setting the Device Parameters (DEVICE)

- New password: Enter the numbers and/or letters for the new password (8 characters max.)
If "none" is displayed, this means no password has been assigned to delete the user password: Press and confirm

- To confirm: press the ↵ soft key

- Exit the Setup menu:
press the ⏮ soft key
> Restart the application

Extra Functions

- Exit the Setup menu:
press the ⏮ soft key
> Restart the application
- Print the parameter settings:
 - If the device parameters are displayed:
press
- > Printout (example)

SETUP

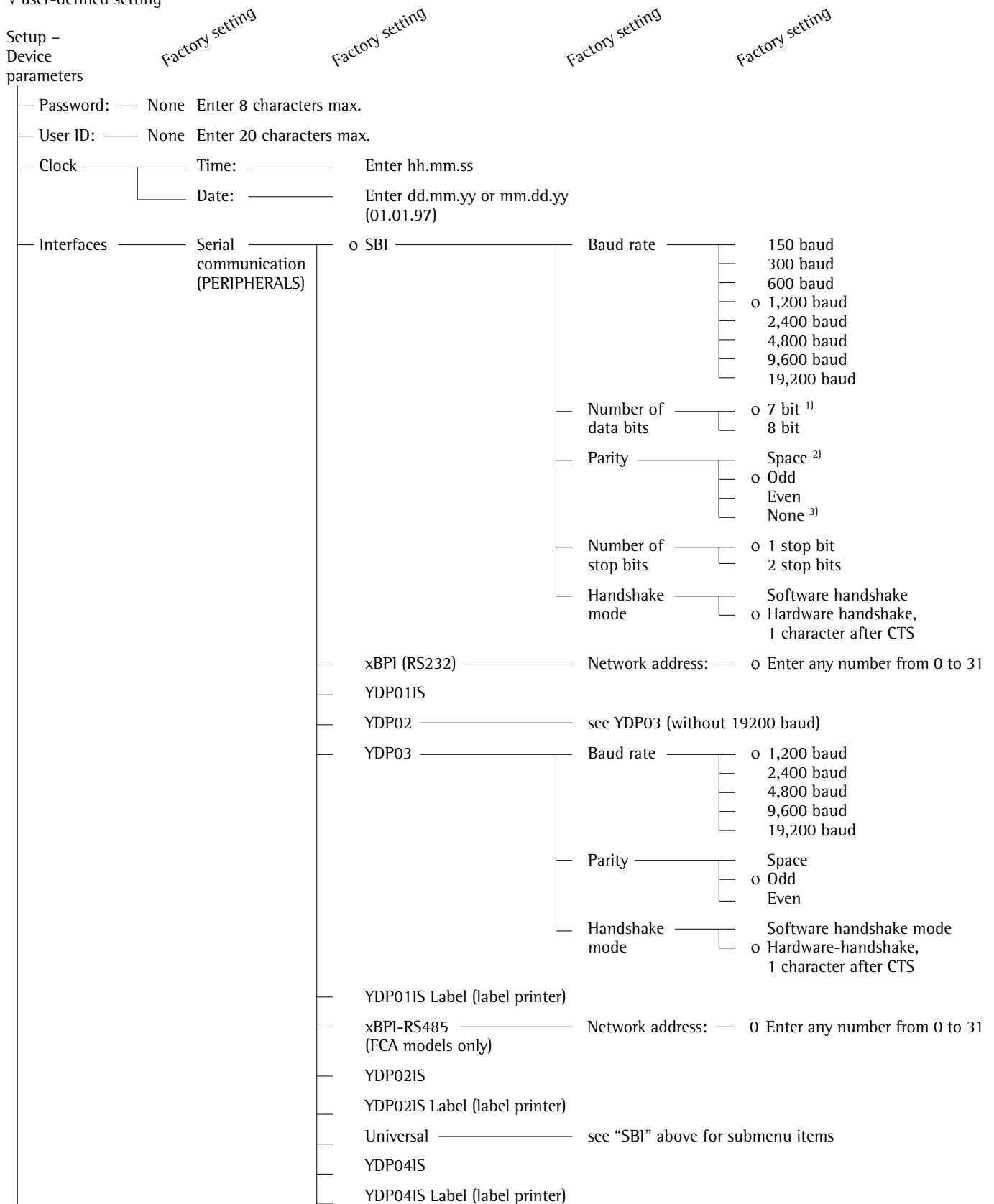
DEVICE

```
-----  
User ID  
User ID:  
  
Interfaces  
Serial communicati  
SBI  
Baud rate  
1200 baud  
Number of data b  
7 data bits  
Parity  
Odd  
Number of stop b  
1 stop bit  
Handshake-mode  
Hardware handshake  
after 1 char  
Function external  
Print key  
Function control  
Output  
Display  
Contrast  
2  
Background  
White  
Digit size  
10mm + bar graph  
+text display  
Application symbo  
On  
Keys  
CF function in ap  
Clear all applicati  
CF function for i  
Delete last charact  
Block key functio  
All keys unblocke  
  
etc.
```

Device Parameters (Overview)

- o factory setting
- √ user-defined setting

Setup – Device parameters



— For the display, keys and extra functions, see next pages

¹⁾ not if “None” parity is selected

²⁾ only if 7 data bits selected

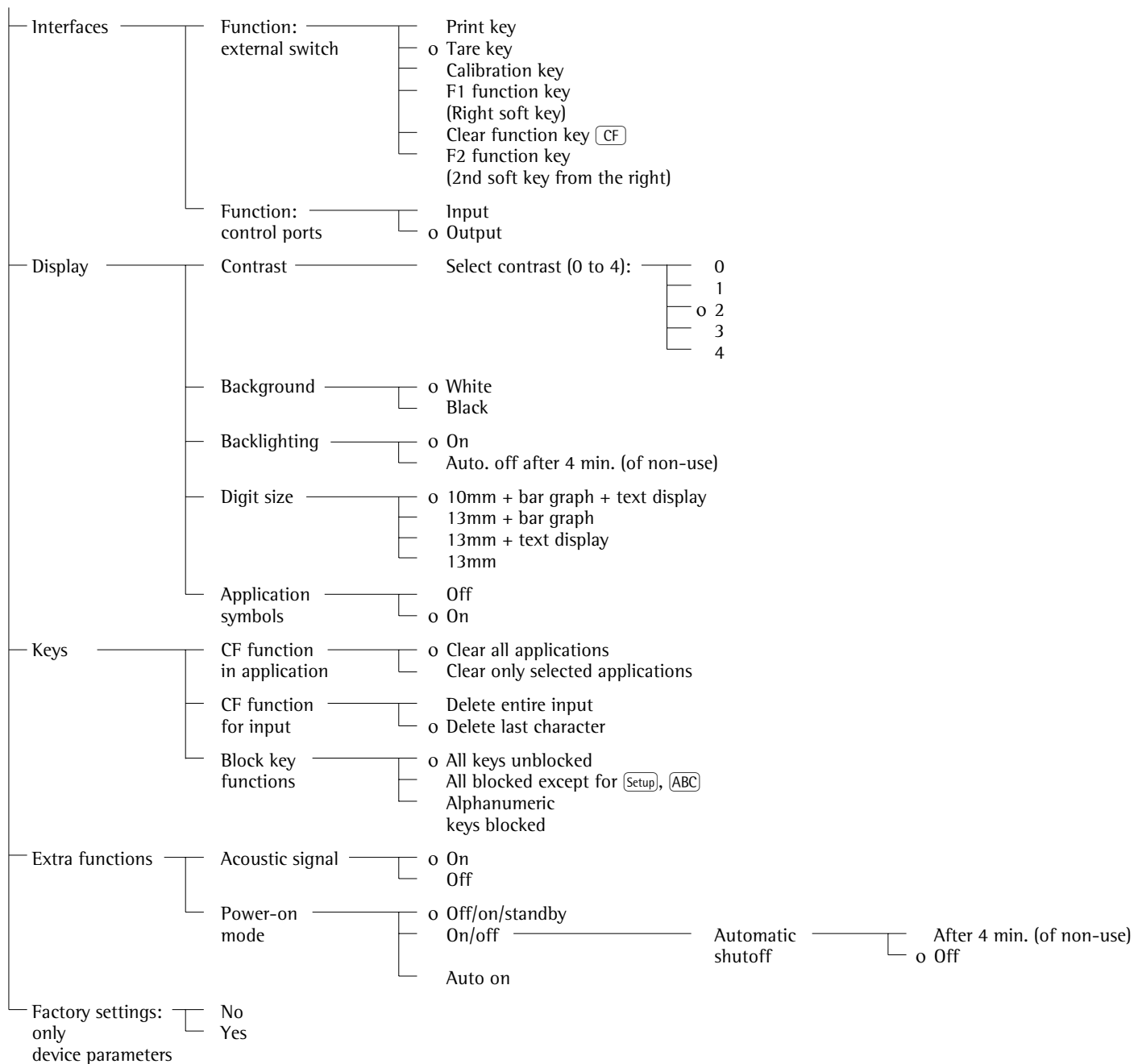
³⁾ only if 8 data bits selected

Setting the Device Parameters (DEVICE)

Setup –
Device
parameters

Factory setting

Factory setting



Setting the Application Parameters (Application)

Purpose

This menu item enables you to configure the scale, i.e., adapt the scale to your individual requirements by selecting from a list of parameter options in a menu. You can block access to this menu by assigning a password.

Features

The simple weighing function is available at all times. You can select one from each of the following application groups. This means a number of combinations are possible.

Application 1 (basic settings)

- Toggle weight units
- Counting
- Weighing in percent
- Animal weighing (averaging)
- Recalculation
- Calculation
- Density determination
- Differential weighing

Application 2 (control functions)

- Checkweighing
- Time-controlled functions

Application 3 (data records)

- Totalizing
- Formulation
- Statistics

In addition, you can assign 2 extra functions to each of the soft keys, in some cases (depending on the Setup configuration):

- Second tare memory
- Identification codes
- Manual storage M+ key
- Product data memory

Auto-start application the scale is switched on

Factory settings:
only application parameters

Factory Settings

The factory settings are identified by the symbol "o" in the list starting on the next page.

Preparation

Display available application parameters:

- Select the Setup menu:
press the **[Setup]** key
- > SETUP is displayed

SETUP					
Balance/scale functions					
Device parameters					
Application parameters					
Printout					
Info					
<<				v	>

- Select parameters: repeatedly press the **v** and **>** soft keys

If you have already assigned a password:

- > The password prompt is displayed:
- If access is blocked by a password:
enter the password using the numeric/alphabetic keys
- If the last character of the password is a letter: conclude input by pressing **[ABC]**
- Confirm your password and have the application parameters displayed:
press the **↓** soft key
- > The application menu is displayed:

SETUP APPLICATION					
Application 1 (basic settings)					
Application 2 (control functions)					
Application 3 (data records)					
Extra function (F4)					
Extra function (F5)					
<<		<		v	>

- To select the next group:
press the **v** soft key (down arrow)
- To select the previous item of a group:
press the **^** soft key (up arrow)
- To select the next sub-item within a group: press the **>** soft key (right arrow)
- To select the previous group:
press the **<** soft key (left arrow)
- To confirm: press the **↓** soft key

Extra Functions

- Exit the Setup menu:
press the **<<** soft key
- > Restart your application
- Print parameter settings:
 - When the scale functions are displayed, press **[E]**
 - > Printout (example)
Texts with more than 20 characters are truncated

SETUP


APPLICATION

Application 1 (basic settings)	Off
Application 2 (control functions)	Off
Application 3 (data records)	Off
Extra function (F4)	Off
2nd tare memory	No
Container tare weight	No
Automatic printout	Off
Extra function (F5)	Off
Auto-start app. when scale is switched on	Off

Setting the Application Parameters (Application)


Application Parameters (Overview)

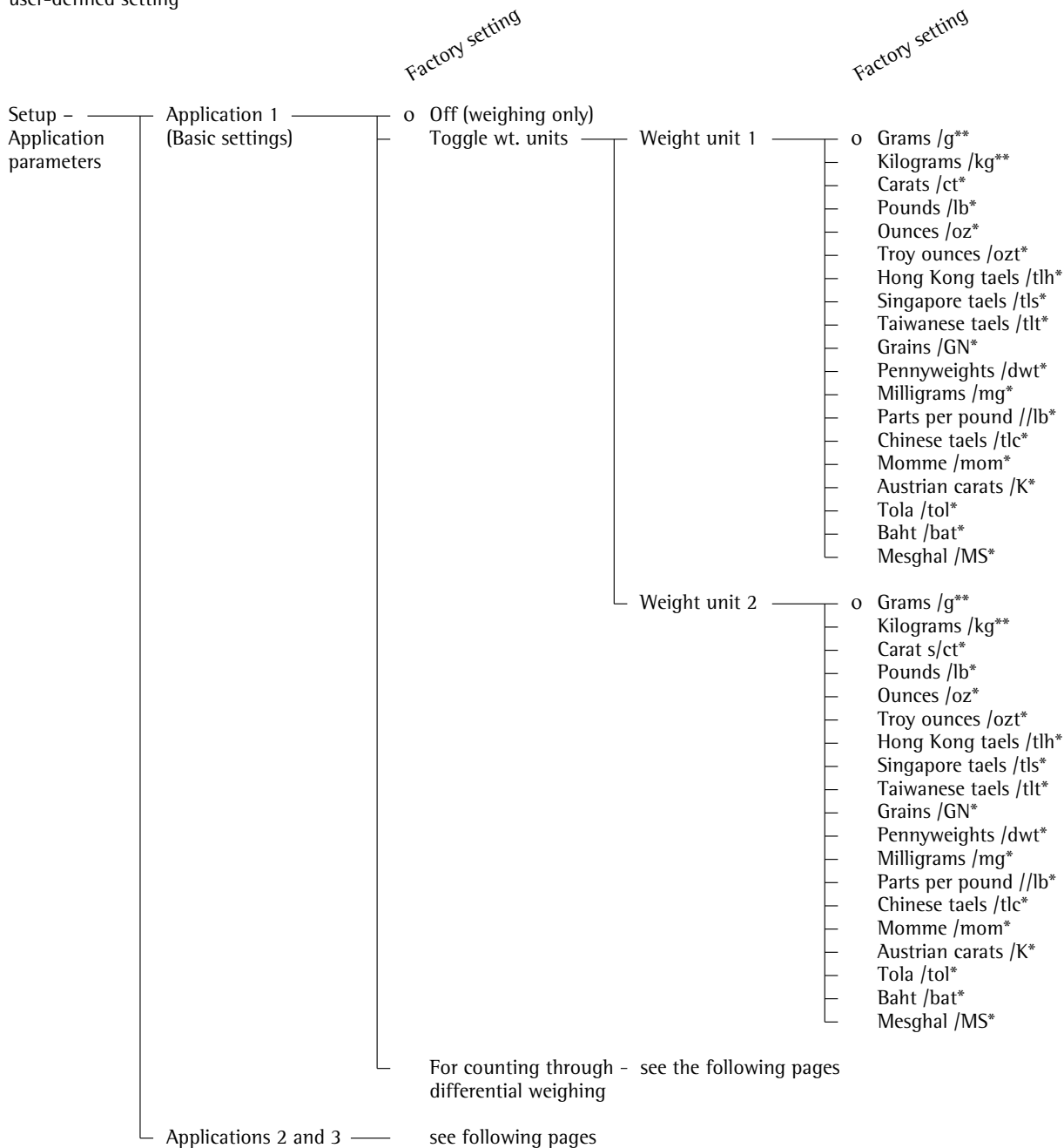- o factory settings
- √ user-defined setting




* = not applicable to verified scales used in the European Economic Area
 ** = factory setting depends on weighing range: – ≤ 33 kg: grams
 – ≥ 34 kg: kilograms

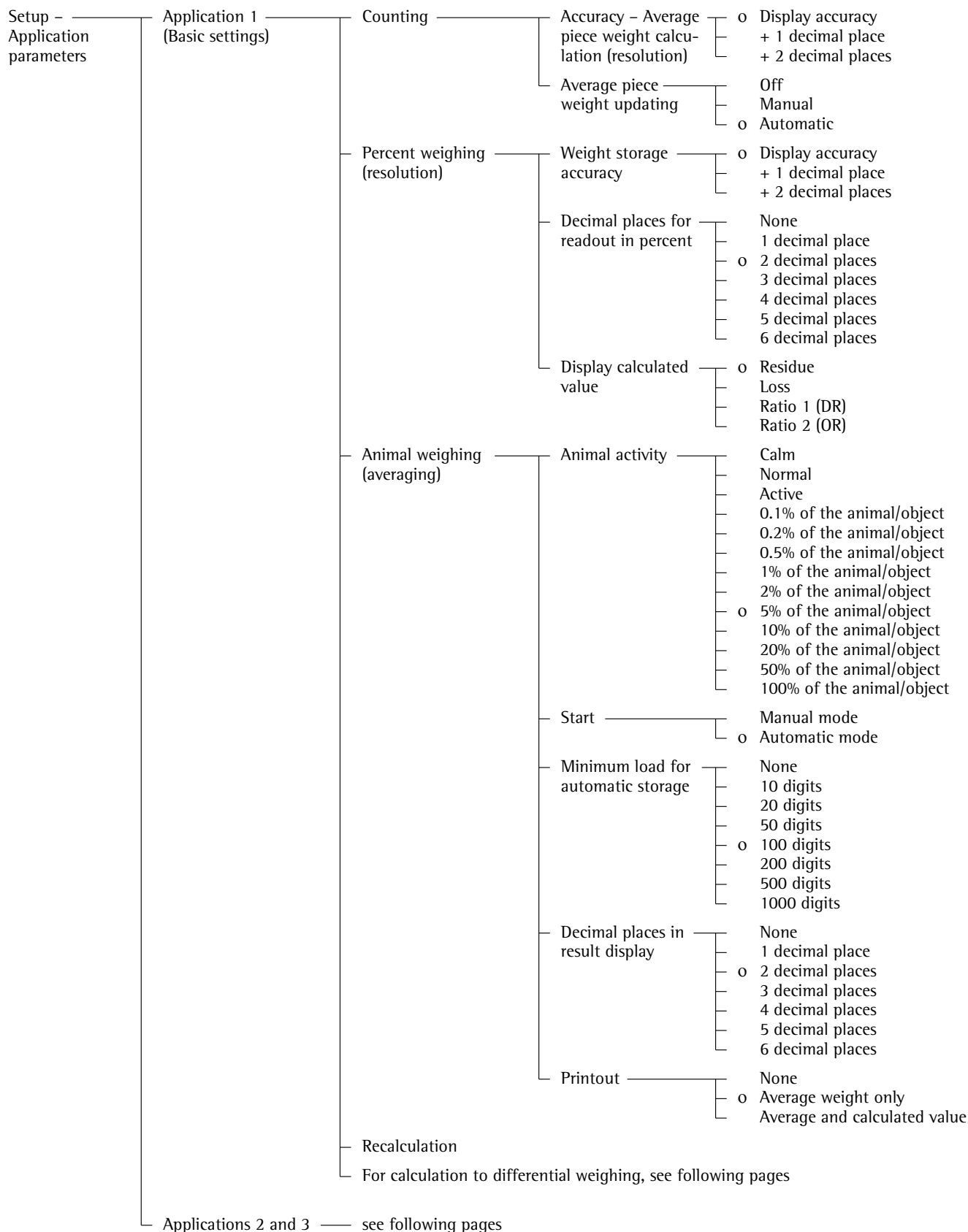
Factory setting



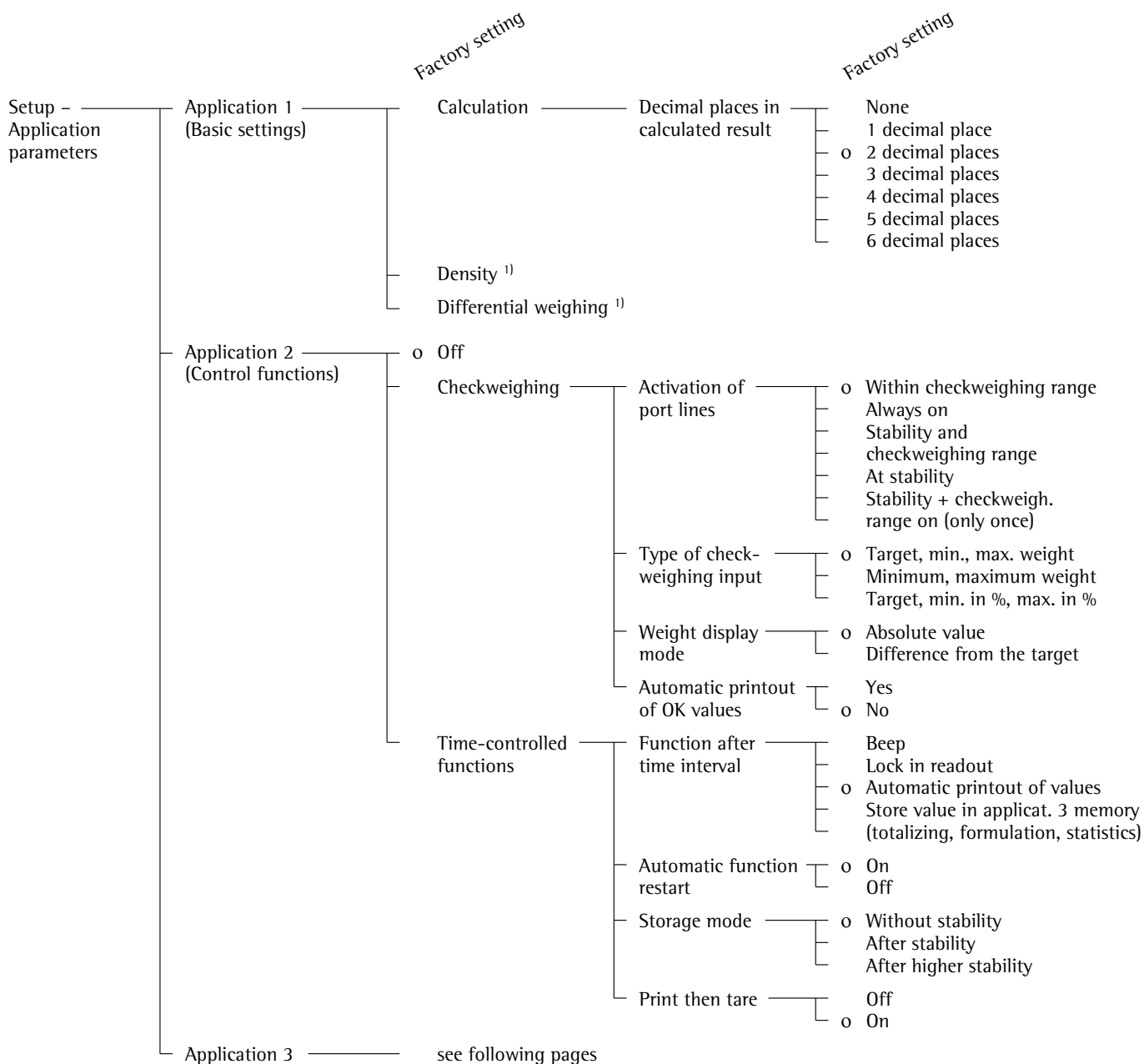
Setting the Application Parameters (Application)



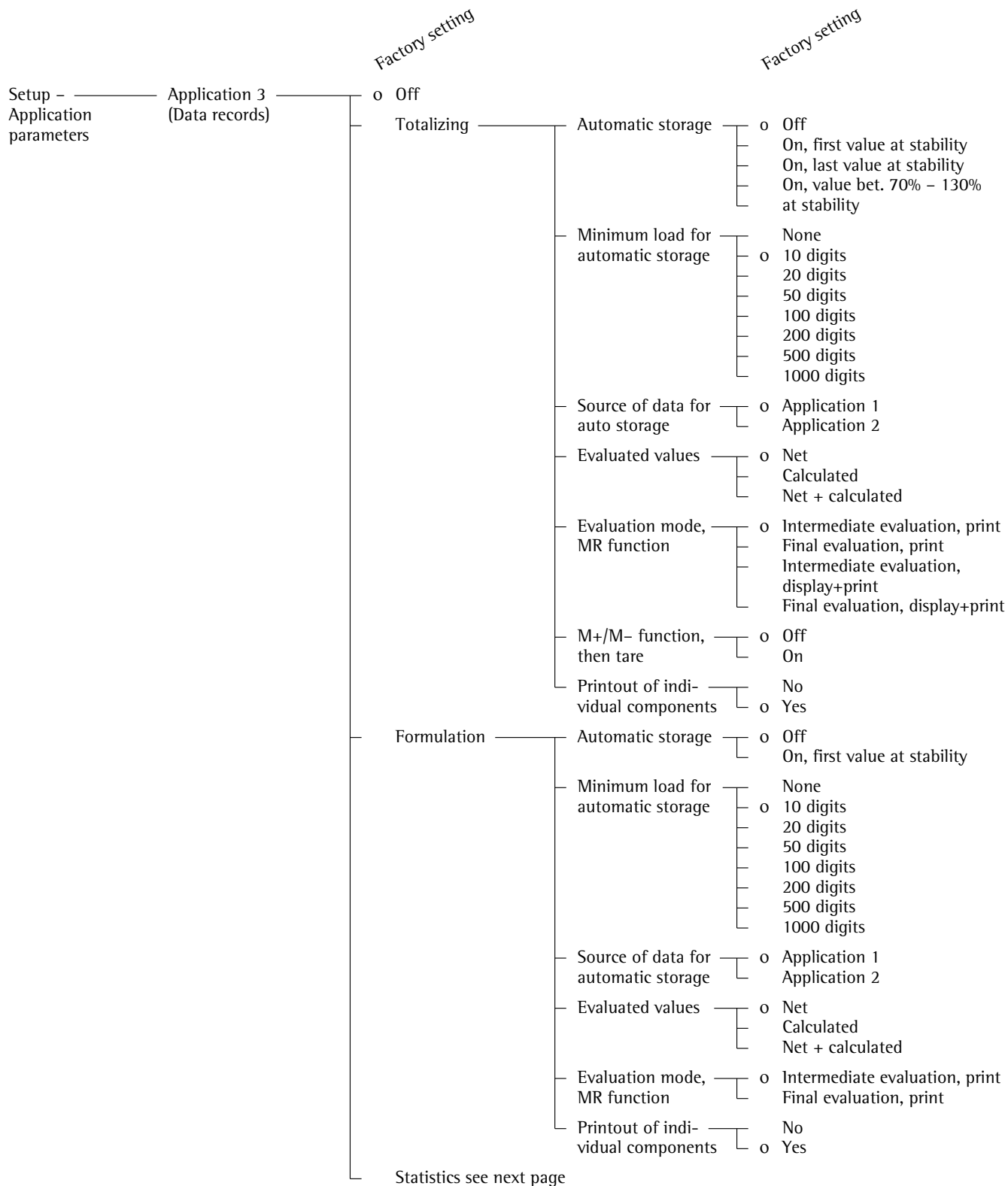¹⁾=For details on using the “Differential Weighing” application, please refer to the operating instructions for Master^{pro} LA balances, available on request from Sartorius or by Internet download (www.sartorius.com → Laboratory Mechatronics → Downloads; file name: “LA.pdf”).



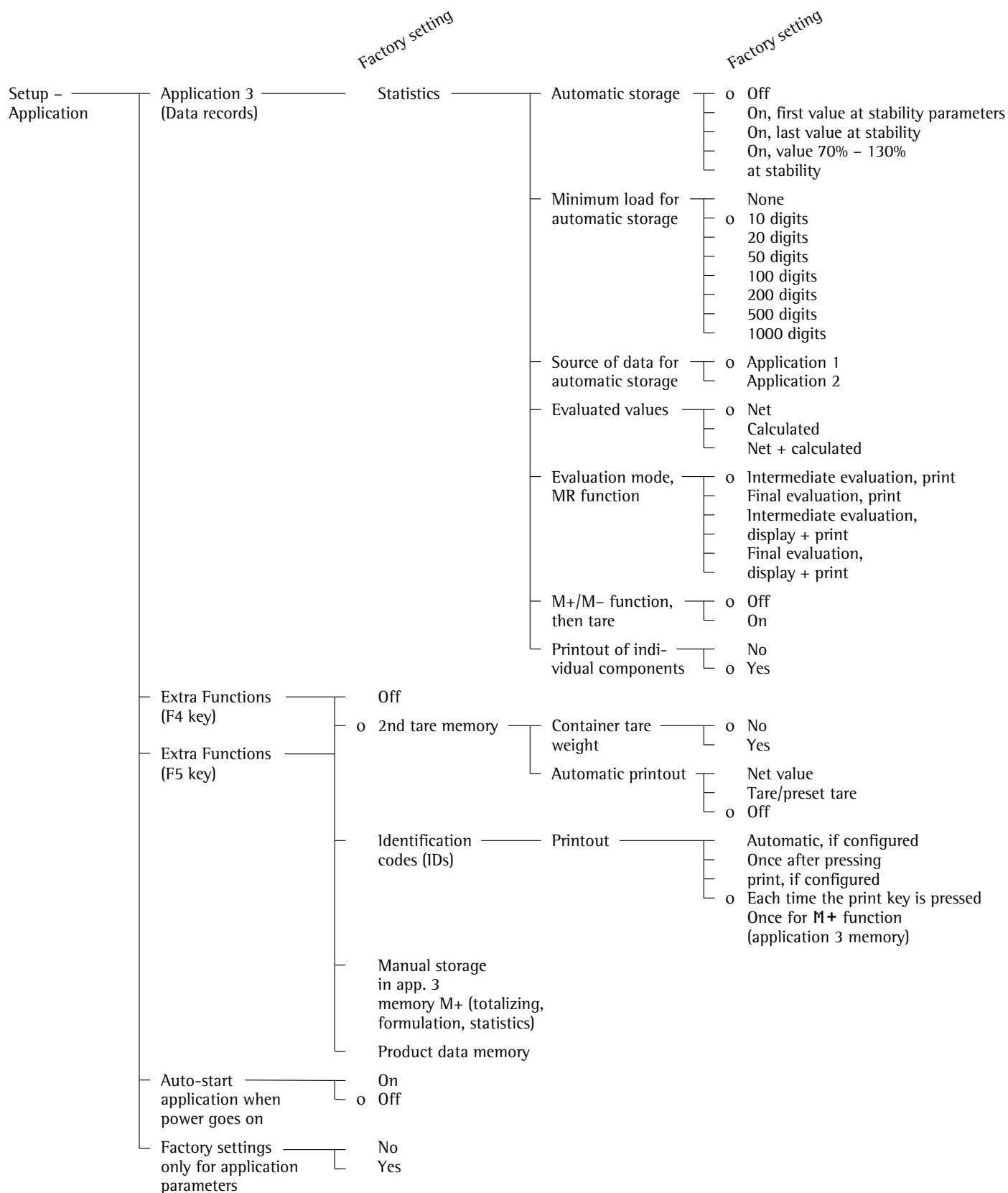Setting the Application Parameters (Application)



Selecting the Printout Function (PRINTOUT)

Purpose

This menu item enables you to configure the printout to meet your individual requirements by selecting predefined menu parameters in the Setup menu. Printouts of weights and other measured or calculated values and IDs enable you to document your data. You can select the particular data you wish to print. To prevent changes to your settings, you can block access to the menu by assigning a password.

Features

The device parameters are combined in the following groups (1st menu level):

- Application-defined output
 - Configured Printout
 - FlexPrint
- Automatic output of displayed values
- Output to interface port
- Line format
- ISO/GLP/GMP printout
- Identification (identifier)
- Factory settings – printout only

Factory Settings

Parameters: The factory settings are identified by the symbol “o” in the list on the next page.

Preparation

Display available printout parameters

- Select the Setup menu: press **[Setup]**
- > SETUP is displayed:

SETUP
Balance/scale functions
Device parameters
Application parameters
Printout
Info
<< < > >>

- Select “Printout”: use the **↓** and **→** soft keys

If no password has been assigned, anyone can access the printout parameters in the Setup menu

If a password has already been assigned:

- > The password prompt is displayed
- If access is blocked by a password: enter the password using the numeric and/or alphabetic keys
- If the last character of the password is a letter: conclude input by pressing the **[ABC]** key

- Press **↓** to confirm the password

- > Printout parameters are now displayed:

SETUP	PRINTOUT
Application-defined output	
Automatic output of displayed value	
Output to interface ports	
Line format	
ISO/GLP/GMP printout	
<< < > >>	

- To select the next group: press the **↓** soft key (down arrow)
- To select the previous item of a group: press the **↑** soft key (up arrow)
- To select the next sub-item within a group: press the **→** soft key (right arrow)
- To select the previous group: press the **←** soft key (left arrow)
- To confirm: press the **↓** soft key

Extra Functions

- Exit the Setup menu: press the **←** soft key
- > Restart your application
- Print parameter settings:
 - When the printout parameters are displayed, press **[F7]**
- > Printout (Example)

SETUP

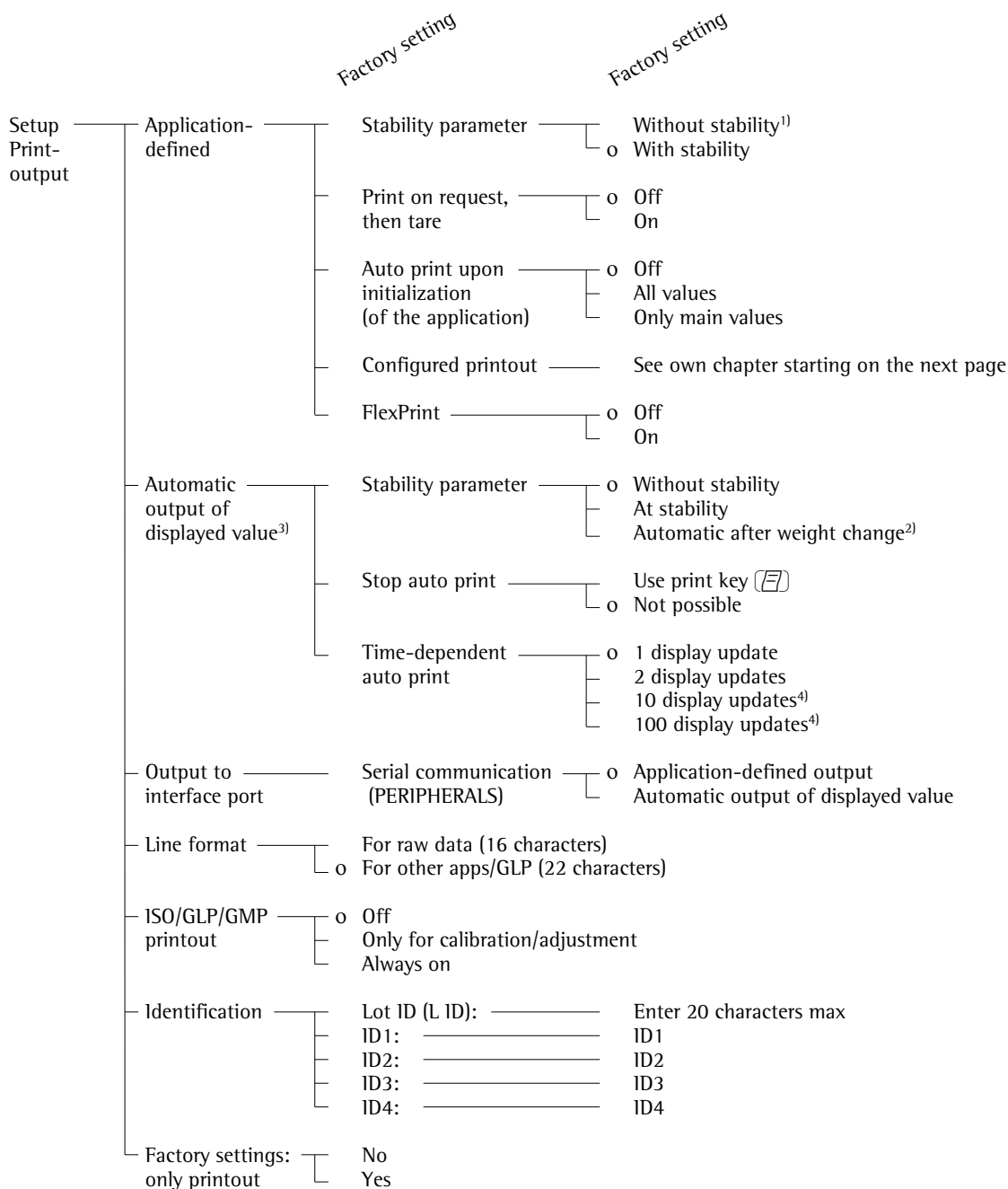
PRINTOUT

```
-----
Application defined
Stability paramete
With Stability
Print on request t
Off
Auto print upon in
All values
Configured printou
Indiv.: Printout
Comp.: Printout
Total: Printout
FlexPrint
Off
Automatic Output of
Stability paramet
Without stabili
Stop auto print
Not possible
Time-dependent aut
1 display update
Output to interface
Serial communicat:
Application-defined
output
Line format
For other apps/GLP
(22 characters)
ISO/GLP/GMP printou
Off
Identification
Lot (L ID):
ID1:
ID1
etc
```

Selecting the Printout Function (PRINTOUT)

Printout Parameters (Overview)

- o factory setting
- √ user-defined setting



¹⁾ = When the scale is used in legal metrology, this setting may be configured for control purposes only; printout of values is not permitted

²⁾ = Autoprint when load change is >10 d and stability is reached: no printout until residual difference in load value >5 d

³⁾ = To enable automatic output of the displayed value, the corresponding option ("Automatic output of displayed value") must be activated under the next menu item at this level, "Output to interface port"

⁴⁾ = These settings may not be changed on verified scales

Printout Configuration

Purpose

This menu item enables you to configure individual printout formats. With the formulation, totalizing and statistics applications, you can also define the values to be included on the total printout when the MR key is pressed.

Under “Setup: Printout: Application-defined output: Configured printout”, you can configure individual, component or total data records that contain the items available for printouts in each application. Configure these printouts after you have configured the applications, because some entries in the data record depend on the particular application.

Features

- Maximum items in a data record: 60
- Separate configuration of printout formats for individual weights, components, total, backweighing and statistics

- Individual printout generation: press the \boxed{E} key

Automatic printout of application data: results from animal weighing or density application (Setup menu: Application 1: Density: Printout: All data) OK values from checkweighing application, time-controlled printouts, 2nd tare memory

- Component printout:
For results from totalizing, formulation or statistics applications, press **M+** or **M-** (Setup: Application 3: ..., Printout of individual components: On)
- Total printout:
For totalizing, formulation or statistics applications, press **MR**
- Backweighing printouts or records: automatically generated after backweighing or manually by pressing the \boxed{E} key when the result is displayed at the end of backweighing
- Statistics printout or output:
To generate, press the \boxed{E} key when the statistics are displayed

Printouts for Differential Weighing:
These printouts can be generated as standard or configured (user-defined) reports.

You can configure the following printouts:

- Individual printout
- Backweighing printout
- Statistics printouts

Printouts are generated in one of two ways:

- at the request of the user by pressing the \boxed{E} key (print on request)
- automatically, if configured in the Setup menu [Application parameters: Application 1: Differential weighing: Generate printout: Auto]

You can turn off automatic printout generation in the Setup menu [Application parameters: Application 1: Differential weighing: Generate printout: None]

Data records are deleted after you have switched to a different application or activated or de-activated an extra function in the application parameters of the Setup menu

- A new pick list for a data record is created based on the currently active application programs and extra functions

- Printout items can be deleted individually

- No printout is generated when the following setting is configured:

Setup: Printout: Line format:
For raw data (16 characters)

- Print item “Form feed” for footer:
Advance to beginning of next label in the “YDP01IS-Label” and “YDP02IS-Label” interface mode

Extra Functions

- Exit printout configuration:
press $\leftarrow \leftarrow$ soft key
- > Restart application

Printing “Select” and “List” Settings

- **LIST**: print the currently selected list
- **Select**: printout items that can still be selected

- When the select bar is on **LIST** or **Select**: press the \boxed{E} key
- > Printout (Example)

```
BACKW. PRINT.LIST
=====
Sample date
Net initial wt.
Backweighed res
Loss in %
=====
etc.
```

Printout Configuration

Practical Example:

Configure an Individual Printout for Counting Application to Include Dotted Line, Date/Time, Piece Count and Net Weight

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Application 1: Counting

Exit the Setup menu: press the << soft key

Then call Setup again and select: Printout: Application-defined output: Configured printout

Step	Press key(s) (or follow instructions)	Display/Output
1. Select Setup menu, "Printout"	[Setup], then v repeatedly and > soft key	<pre> SETUP PRINTOUT Application-defined output Automatic output of displayed value Output to interface ports Line format ISO/GLP/GMP printout << < v > </pre>
2. Confirm "Application-defined output"	> soft key	<pre> SETUP PRINTOUT APPLICATION Stability parameter Print on request then tare Auto print upon initialization Configured printout FlexPrint << < v > </pre>
3. Select and confirm "Configured printout"	v soft key 3 x > soft key	<pre> PRINTOUT APPLICATION CONFIG Indiv.: Printout f. app./weighing << < v > </pre>
4. Confirm "Indiv. printout"	> soft key	<pre> LIST INDIV.PRT SELECTION Blank line Form feed Date/time Time << < v > </pre>
5. Select "Dotted line"	>, v, ↓ soft keys	<pre> LIST INDIV.PRT SELECTION Blank line Form feed Date/time Time << < v ↓ </pre>
6. Select "Date/time"	v soft key twice, then ↓ soft key	<pre> LIST INDIV.PRT SELECTION Date/time Form feed Time GLP header << < ^ v ↓ </pre>
7. Select "Piece count"	v repeatedly, then ↓ soft key	<pre> LIST INDIV.PRT SELECTION Date/time Gross (G#) Ref. quantity Ref. weight Target << < ^ v ↓ </pre>
8. Select "Net weight"	^ soft key repeatedly, then ↓ soft key	<pre> LIST INDIV.PRT SELECTION ID1 ID2 ID3 ID4 Gross (G#) << < ^ v ↓ </pre>
9. Exit "Printout configuration"	<< soft key	
10. Perform weighing operations, then print	[E]	<pre> ----- 14.01.2000 09:19 Qnt + 598 pcs N + 2003.13 g </pre>

Data Items for the Printout:

Parameter	Display text	Indiv.	Comp.	Total
Blank line**	Blank line	x	x	x
Dotted line**	-----	x	x	x
Form feed*	Form feed	x	x	x
Date/Time*	Date/Time	x	x	x
Time with seconds*	Time	x	x	x
GLP header*	GLP header	x	x	x
GLP footer*	GLP footer	x	x	x
Sample ID*	S ID	x	x	x
ID 1*	ID1	x	x	x
ID 2*	ID2	x	x	x
ID 3*	ID3	x	x	x
ID 4*	ID4	x	x	x
Net weight*	Net (N)	x		
Gross weight*	Gross (G#)	x	x	x
Preset tare/ Tare 1 weight*	Tare1 (T1/PT1)	x	x	x
With "Counting" application:				
Reference quantity	Ref. quantity	x	x	x
Reference weight	Ref. weight	x	x	x
Piece count	Piece Count	x		
With "Weighing in percent" application:				
Reference percentage	Ref. percent	x	x	x
Reference weight	Ref. weight	x	x	x
Percentage	Percent	x		
With "Animal Weighing" application:				
Number of weighing operations	No. of weights	x	x	x
Calculation factor	Factor	x	x	x
Average animal weight	Mean value	x		
Average animal weight calculated	Mean factor	x		
With "Calculation" application:				
Equation	Equation	x	x	x
Calculated result	Calc. result		x	

* = Items are available independently of the applications selected

** = Items are available independently of the applications selected and can be selected more than once

Data Items for the Printout:

Parameter	Display text	Indiv.	Comp.	Total
With the "Check-weighing" application:				
Target value	Target	x	x	x
Minimum value	Minimum	x	x	x
Maximum value	Maximum	x	x	x
With "Time-controlled Functions" application:				
Time/interval	Time/interval	x		
With the "Totalizing" application:				
No. of weights	No. of wts.		x	x
Weight of transaction	Trans. wt.		x	
Weight total	Wt. total			x
Number of calculated values	No. of calc. val.		x	x
Calculated value transactions	Calc. val. trans.		x	
Total of calculated values	Total calc.			x
Nominal no. of weighing operations	Nom. no. wghs.		x	x
With "Formulation" application:				
Number of components	Number		x	x
Net component	Net component		x	
Components calculated	Net transact.		x	
Total net components	Net/comp. calc.			x
Total calc. components	Tot. comp. calc.			x
Preset tare/ Tare 2 weight	Tare2	x	x	x
Nominal no. of weighing operations	Nom. no. wghs.		x	x
With the "Statistics" application:				
No. of weights	No. of wts.		x	x
Weight of trans.	Trans. wt.		x	
Mean weight	Average wt.			x
Standard deviation – weight	Std. dev. wt.			x
Variation coefficient – weight	Var. coeff. wt.			x
Weight total	Wt. total			x
Minimum weight	Min. wt.			x
Maximum weight	Max. wt.			x
Difference – weight	Diff. wt.			x
No. of calc. values	No. of calc. val.		x	x
Calc. value – transactions	Calc. val. trans.		x	
Mean calc. value	Mean calc. val			x
Standard derivation calculated values	Std. dev. calc.			x
Variation coefficient – calc. values	Var. coeff. calc.			x
Total – calc. values	Total calc.			x
Minimum – calc. values	Min. calc.			x
Maximum – calc. values	Max. calc.			x
Difference – calc. values	Diff. calc.			x
Nominal no. of weighing operations	Nom. no. wghs		x	x

Displaying Info

Purpose

This menu item enables you to have information displayed about the specific scale ("device"), as well as "FlexPrint" information.

Displaying Device Information

- Select the Setup menu:
press the **[Setup]** key
- > "SETUP" is displayed:

SETUP
Balance/scale functions
Device parameters
Application parameters
Printout
Info
<< < > >>

- Select "Device information":
Repeatedly press the **↓** soft key, then press the **→** soft key
- Press **→** soft key to confirm "Device information"
- > Device information is displayed:

SETUP	INFO	DEVICE
Version no.:	01-45-01	
Wgh. sys. ver. #:	00-20-11	
Model:	FCA64EDE-HX	
Serial no.:	91205355	
<< < > >>		

- Print device information:
Press the **[F]** key

- > Printout (Example)

```

-----
23.02.2000      13:02
Model      FCA64EDE-HX
Ser. no.    91205355
Vers. no.   01-45-01
(Version of the operating program)
ID          BECKER123
(User-ID)
-----

```

```

-----
SETUP
INFO
DEVICE
-----
Version-no.:
              01-45-01
(Version of the operating program)
Wgh. sys. vers:
              00-20-11
(Version no. of the weighing cell)
Model:
              FCA64EDE-HX
Serial no.:
              91205355
-----

```

- Return to SETUP overview:
press the **<** soft key
- Exit Setup menu:
press the **<<** soft key
- > Original settings are restored

Display Flexprint Information

- Select the Setup menu:
press the **[Setup]** key
- > "SETUP" is displayed:

SETUP
Balance/scale functions
Device parameters
Application parameters
Printout
Info
<< < > >>

- Select "Info": press the **↓** soft key repeatedly and then the **→** soft key

SETUP		INFO	
Device information			
FlexInfo			
<<	<	>	>>

- Select "FlexInfo": press the **↓** soft key and then the **→** soft key
- > The FlexPrint information is displayed, with print instruction file name, software ID and version number:

SETUP	INFO	FLEXINFO
PDIRECT	ID---	U.---
PGMPF00T	ID403	U.000001
PGMPHEAD	ID403	U.000001
<< < > >>		

- To select a particular print file name with software ID (for example, ID403), if desired: press key **↓** or **↑** as required

- > If the display shows **ID---** :
The weight block for legal metrology is not printed by this print file.

- > Display of version number:
U. xx. xx. xx
Created by Sartorius:
U. S. xx. xx. xx

- Return to SETUP overview:
press the **<** soft key

- Exit Setup menu:
press the **<<** soft key
- > The device returns to the previous mode

Date of Manufacture

The month and year of manufacture are encoded in the serial number as follows:

Y M M x x x x x

Y Year

1	2000-2006
2	2007-2013
3	2014-2020
4	2021-2027
5	2028-2034
6	2035-2041
7	2042-2048
8	2049-2055
9	2056-2062

The first digit represents a 7-year period as indicated in the table above. The next 2 digits represent the month. The months are numbered consecutively, starting with 13, over the entire 7-year period. Thus the number representing the month also indicates the specific year of manufacture.

2000 13-24
2001 25-36
... etc.

Example:
113xxxxx → January 2000

The individual devices are numbered consecutively in the last 5 digits, starting from 00000 again at the beginning of each month.

MP8 Interface Emulation

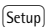
Purpose

With the MP8 interface emulation function, you can connect peripheral devices of the MP8 generation that have separate AC power supplies (such as the 73822... Data Control terminal, a YFC..., or a YDI 50 Z Data Input dedicated keyboard, for example) to your scale.

Features

- The scale can be used only to determine weights.
- The interface communicates exclusively in the MP8 binary protocol.
- Select the application program and the program index for MP8, as well as individual application parameters, in the Setup menu.
- See next page for menu setting options

Preparation

- Activate the MP8 emulation mode*:
 - Press 
 - Select the **factory settings** and confirm: press the \vee soft key repeatedly and then the \triangleright soft key
 - Select **Reset to MP8:** press \vee soft key and \triangleright soft key
 - Select **Yes** and press \downarrow to confirm
- > The terminal is restarted

* follow the same procedure to return from MP8 emulation back to factory setting

Factory Settings

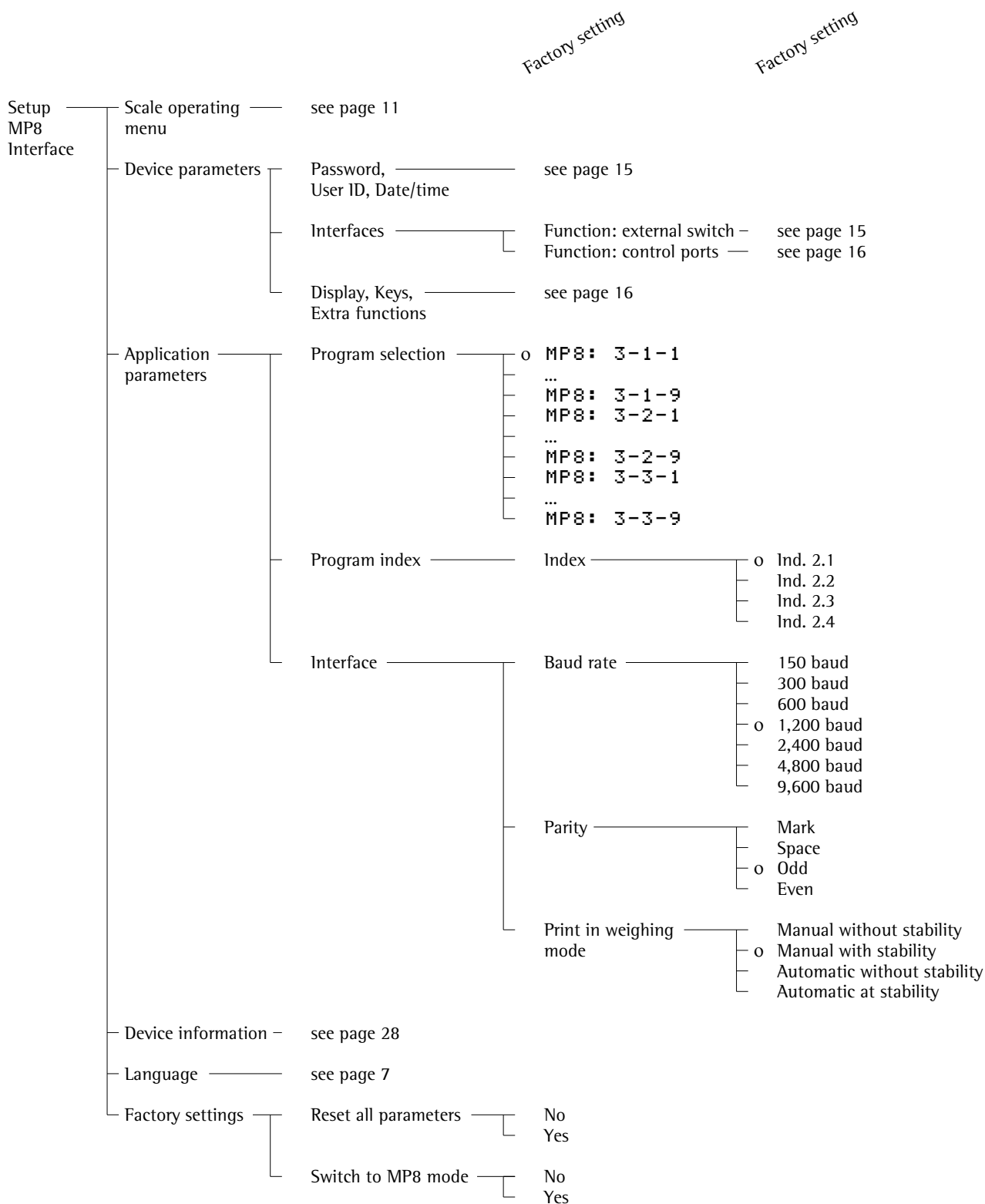
Each parameter category has a factory setting. To restore the factory settings, select this item in the Setup menu and select **YES** to confirm.

The following parameters are not restored to factory settings when you activate this function:

- Language
- Password
- Display contrast
- Time (clock)

MP8 Interface Emulation Parameters (Overview)

- o factory setting
- √ user-defined setting



Operation

Basic Weighing Function

Purpose

The basic weighing function is always accessible and can be used alone or in combination with an application program (Toggle between Weight Units, Counting, Weighing in Percent, etc.).

Features

- Taring the scale
- Assigning IDs to weights
- Printing weights
- Printing ID codes for weights

Factory Settings

Tare: **After stability**

Manual/auto print mode:

Manual with stability

Line format:

**For other apps/GLP
(22 characters)**

Soft Key Functions


isoTST	Initiate calibration/ adjustment routine
isoCAL	Press to start isoCAL routine
S ID	Store ID entered

Under-Scale Weighing


A port for an under-scale weighing hanger is located on the bottom of the scale.

FC... BBE, FC... CCE:

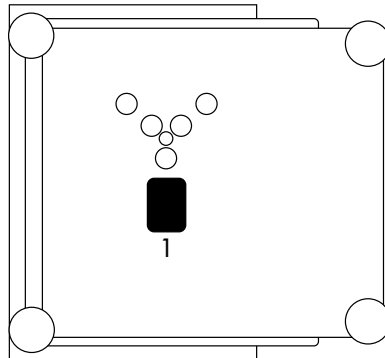- Open cover plate **(1)** on the bottom of the scale




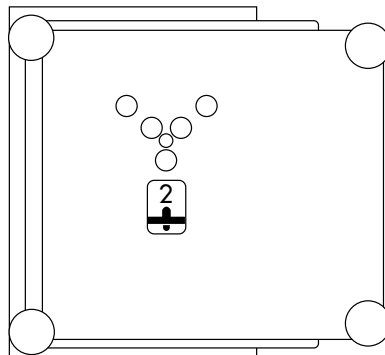- Attach the sample (e.g., using a suspension wire) to the hook **(2)**.



FCG... EDE:

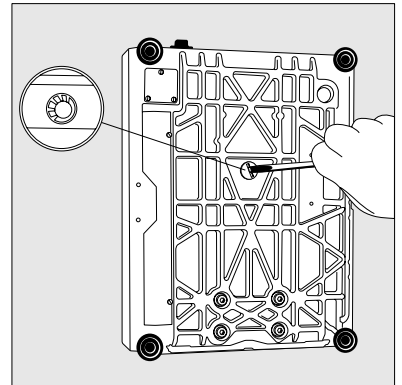- Use a suitable screwdriver to remove the cover plate from the bottom of the scale.

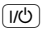

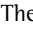


- Fasten the hook to the cross-bar (see "Accessories").

- If necessary, install a shield for protection against drafts

⚠ Important Note Concerning Verified Scales Approved for Use as Legal Measuring Instruments in the EU*: The under-scale weighing port may not be opened or used when an approved scale is being operated as a legal measuring instrument

Preparation

- Turn on the scale: press 
> The Sartorius logo is displayed
- To tare the scale, if desired: press 
- > The  symbol is displayed when a verified scale is tared or zeroed (± 0.25 digits).

⚠ Important Note Concerning Verified Scales Approved for Use as Legal Measuring Instruments in the EU*: The type-approval certificate for verification applies only to non-automatic weighing instruments; for automatic operation with or without auxiliary measuring devices, you must comply with the regulations of your country applicable to the place of installation of your scale.

- You must calibrate the scale at the place of installation before using it as a legal measuring instrument (see the section entitled "Calibration and Adjustment" in this chapter)

- The temperature range indicated on the verification ID label must not be exceeded during operation

* including the Signatories of the Agreement on the European Economic Area

Additional Functions

In addition to the functions:

- alphanumeric input
- taring (not during alphanumeric input)
- printing

you can also access the following functions from the weighing application:

- Calibration
 - Press **isoTST**
 - > See the section on “Calibration/Adjustment” for further instructions.
- Setup
 - Press **Setup**
 - > See the chapter entitled “Configuration” for further instructions.
- Turning Off the Scale
 - Press **1/0**
 - > The scale shuts off
 - > The display goes blank

Practical Examples

Example W1: Simple weighing

Step	Key (or instruction)	Display/Output
1. If necessary, tare the scale	<div>Tare</div>	<div>Max6200 g d= 0.01g 0% 100% 0.00 g COUNTING: nRef = 10 pcs isoTST PT1/T1 Start</div>
2. Enter sample ID	see Example W2	
3. Determine sample weight (Example)	Place load on scale	<div>Max6200 g d= 0.01g 0% 100% 223.00 g COUNTING: nRef = 10 pcs isoTST PT1/T1 Start</div>
4. Print weight value	<div>⏏</div>	<div>S-ID ABC123 N + 2231.56 g</div>

Example W2

Enter "ABC123" as sample ID

△ Note:

- The sample ID generally applies to one weighing operation only
- The ID is deleted after data output

Step	Key (or instruction)	Display/Output
Initial status (scale unloaded) (ID can also be entered while scale is loaded)		
1. Select alphabetic input	(ABC)	
2. Select the required letter group	ABCDEF soft key	
3. Enter the letter "A" (To delete a letter: (CF))	A soft key (CF)	
4. Select the letter group and enter "B"	ABCDEF soft key B soft key	
5. Select the letter group and enter "C" (If only letters are entered, conclude input: (ABC))	ABCDEF soft key C soft key (ABC)	
6. Enter the numbers 1, 2 and 3	(1) (2) (3)	
7. Store the ID (max. 20 characters) - The next printout will include the sample ID	S ID soft key	

Device Parameters

Password

You can enter a password to block access to the operating menu and to the functions for ID code input and exact calibration weight input.

For details, see “Setting the Device Parameters” in the chapter entitled “Configuration”.

User ID

You can enter your own personal password (20 characters max.)

Clock

ISO/GMP printouts in particular must be generated with the date and time stamp of the specific measurement. This date and time stamp is optional on other printouts.

For details, see “Setting the Device Parameters” in the chapter entitled “Configuration”.

Interface

Purpose

This item enables you to set the parameters for the following interfaces:

- Serial interface port
- External switch function

Serial communications port

You can set the serial communications port to use for the following modes:

- SBI
- XBPI
- YDP01IS
- YDP02
- YDP03
- YDP01IS-Label
- XBPI-RS485*
- YDP02IS
- YDP02IS-Label
- Universal
- YDP04IS
- YDP04IS-Label

Universal Remote Switch

You can connect an external universal remote switch (foot switch) to one of the two serial ports. Then you can assign one of the following functions to be performed when the switch is activated:

- Print key
- Tare key
- Cal key
- F1 function key
- CF key
- F2 function key





For further information on the pin assignment chart, see “Pin Assignment Charts” in the chapter entitled “Data Output Functions”.

Control Port Function

You can connect either a checkweighing display or an external universal switch to the serial communications port on the balance (factory setting).

To do so, you need to configure the interface for **input** or **output** in the operating menu under “Device parameters: Interfaces.”

Pin Assignment Chart of the Female Interface Connector

Pin	Function: Input
0	 key; see “Universal switch”
F	 key
E	Soft key 6 (Cal)
G	Soft key (F1)
D	 key
F	 key
Pin	Function: Output
0	“External switch” (see above)
F	Control port 1: lighter
F	Control port 2: equal
F	Control port 3: heavier
F	Control port 4: “set”

For further information on the pin assignment chart, see the section on “Pin Assignment Charts” in the chapter entitled “Overview”.

* FCA models only

Display

You can configure the display for your individual needs.

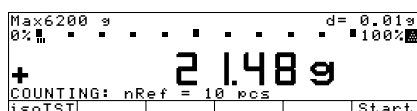
The contrast can be adjusted in 5 levels:

Contrast


Characters can be displayed in black on white or vice versa: **Background**




You can blank out either the bar graph or the text line or both **Digit size**




10mm + bar graph + text display

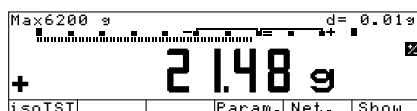

13mm + bar graph

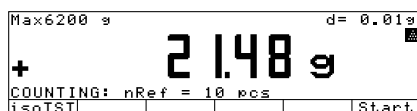

13mm + text display

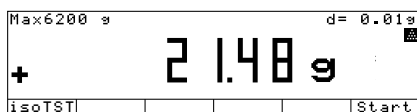

13mm

You can blank out the display of application symbols **Application symbols**

Keys

You can assign different functions to the **[CF]** key for deleting input and applications.

When you delete applications, you can delete either the data stored for all applications or just selected data. **CF function in application**

When you delete input, you can either delete all the data input in a field, or only the last character entered. **CF function for inputs**

You can block key functions; you can choose whether to block all keys (except **[I/O]** and **[Setup]**) or just the alphanumeric keys.

Block key functions

Extra Functions

Acoustic Signal

An acoustic signal is emitted when you press a key. When the key pressed is allowed, the signal is a single beep-tone; when it is not allowed, this is signaled by a double-beep (key does not initiate a function). In the Setup menu, you can configure whether

- the acoustic signal should sound (**On**)
- the acoustic signal should not sound (**Off**)

Power-On Mode

You can configure the scale so that when a power supply is connected,

- the scale is off (**Off/on/standby**), or
- the scale switches on automatically (**Auto on**)

You can also configure the scale to go into the standby mode (**Off/on/standby**) when it is turned off.

After you turn on the scale, a self-test of the functions is run (**TEST** is displayed in the text line and the bar graph is shown)

Calibration/Adjustment "isoTEST"*

Purpose

Calibration is the determination of the difference between the weight readout and the true weight (mass) of a sample. Calibration does not entail making any changes within the scale.

Adjustment is the correction of this difference between the measured value displayed and the true weight (mass) of the sample, or the reduction of the difference to an allowable level within maximum permissible error limits.

Using Verified Scales as Legal Measuring Instruments in the EU**:

Before using your scale as a legal measuring instrument, you must perform "internal calibration" at the place of installation after the warmup period.

Available Features

You can start the isoTEST function at the press of a key to check a scale used as measuring, inspection and test equipment at any time. The scale is calibrated and any deviation is displayed. Press the **Start** soft key to start adjustment. If you do not wish to have the scale adjusted, press the **End** soft key to cancel the isoTEST.

Your scale can be calibrated externally (Scale menu: CAL/isoTST key function; menu item **Ext. cal./adj.;** **factory-def. wt.** or **Ext. cal./adj.;** **user-defined wt.**) or internally (**Internal cal./adjustment**).

External calibration can be performed

- using a pre-set weight value
Ext. cal./adj.;
factory-def. wt., or
- with a user-defined weight
Ext. cal./adj.;
user-defined wt.

The adjustment can be performed

- automatically following calibration:
Cal., then auto
adjust. or
- if desired, the adjustment operation can be started manually after calibration:
Cal., then manual
adjust

You can also configure whether the calibration mode

- will be activated according to the specific setting (external/internal) or
- can be selected by the user after pressing the **isoTST** soft key:
Selection mode.

You can have the scale automatically display an adjustment prompt after a certain time interval has elapsed since the last calibration/adjustment or when the ambient temperature changes by a defined amount.

You can configure the scale to perform calibration and adjustment automatically (isoCAL) when the pre-set time(s) and/or temperature limit is reached:

On and reset applica-
tion and On without
resetting applica-

You can have the calibration/adjustment results documented on an ISO/GMP-compliant printout or on a block printout with up to 50 adjustment sequences; see page 100.

Factory Settings

Calibration/adjustment mode:
Selection mode

Calibration/adjustment sequence:
Calibrate, then auto
adjust

Automatic calibration/adjustment function: **On without**
resetting app.

Start automatic adjustment: **isoCAL**

Generate GLP/GMP-compliant record:
Automatic if GLP is
selected

External Calibration in Verified Scales of Accuracy Class (II)

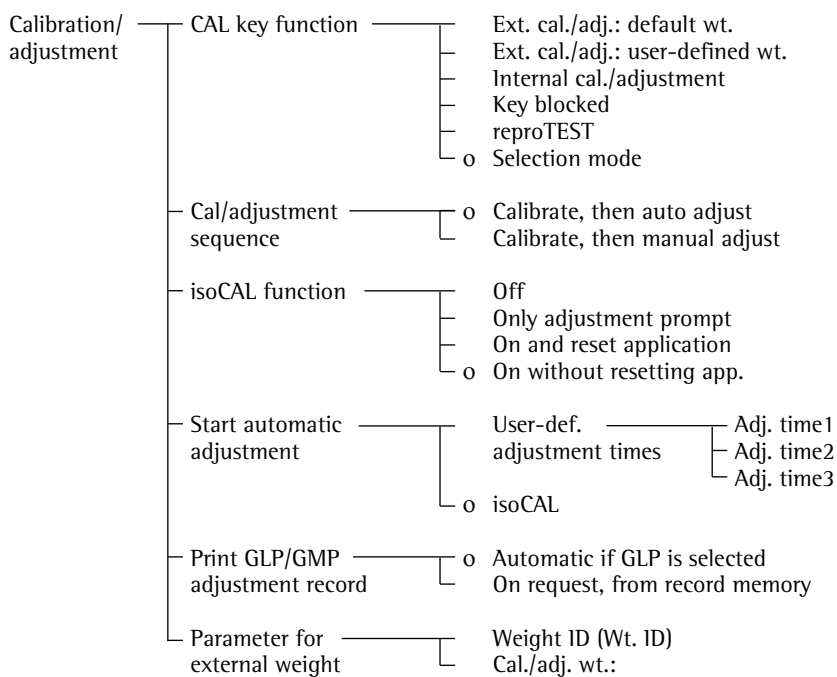
- External calibration is blocked when the scale is used in legal metrology
- > External calibration can only be released after removing the verification control seal, in which case the validity of the verification becomes void and the scale must be re-verified
- External calibration can now be performed

* isoTEST = TEST in the U.S. and Canada

** including the Signatories of the Agreement on the European Economic Area

Preparation

- Select the balance function for “calibration/adjustment”: press **Setup**
- To select the **Balance/scale functions**: press the **>** soft key
- Select **Calibration/adjustment**: press the **>** soft key





o = factory setting

- Save settings and exit Setup menu: press the **<<** soft key

Calibration/Adjustment

Preparation

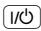Set the parameters for calibration and adjustment; e.g., with manual calibration/adjustment, isoCAL off

Step	Press key(s) (or follow instructions)	Display/Output
1. Switch on the scale, if not already on		Sartorius logo and self-test <div> Max 6200g d= 0.01g 0% 100% 0.00 g isoTST PT1/T1 </div>
2. Select the Setup menu		<div> SETUP Balance/scale functions Device parameters Application parameters Printout Device information << >> </div>
3. Select "Balance/scale functions"	> soft key	<div> SETUP BAL.FUNC. Calibration/adjustment Adapt filter Application filter Stability range Taring << < > >> </div>
4. Select "Calibration/adjustment"	> soft key	<div> SETUP BAL.FUNC. CAL./ADJ. CAL/isoTST key function Cal/adjustment sequence isoCAL function Print GLP/GMP adjustment record Parameter for external weight << < > >> </div>
5. Select CAL/isoTST key function	> soft key	<div> BAL.FUNC. CAL./ADJ. CAL KEY Parameter for external weight Ext. cal./adj.: user-defined wt. Key blocked Selection mode << < > >> </div> <p>o = last setting selected</p>
6. Select desired function and confirm (e.g., "Ext. cal./adj.; factory-def. wt.")	^ soft key, repeatedly, if necessary ↓ soft key	<div> BAL.FUNC. CAL./ADJ. CAL KEY Ext. cal./adj.: factory-def. wt. Ext. cal./adj.: user-defined wt. Key blocked Selection mode << < > >> </div>
7. Exit CAL/isoTST key function	< soft key	<div> SETUP BAL.FUNC. CAL./ADJ. CAL/isoTST key function Cal/adjustment sequence isoCAL function Start automatic adjustment Print GLP/GMP adjustment record << < > >> </div>
8. Select "Cal./adjustment sequence"	v soft key	<div> SETUP BAL.FUNC. CAL./ADJ. CAL/isoTST key function Cal/adjustment sequence isoCAL function Start automatic adjustment Print GLP/GMP adjustment record << < > >> </div>

Step	Press key(s) (or follow instructions)	Display/Output
9. Confirm calibration and adjustment sequence	> soft key	<div>BAL.FUNC. CAL./ADJ. CAL/ADJ SEQ</div> <div>oCalibrate, then auto adjust</div> <div>oCalibrate, then manual adjust</div> <div><< < v ↓</div>
		o = last setting selected
10. Select other settings, if desired and confirm (e.g., Calibration with manual adjustment)	v and ↓ soft keys	<div>BAL.FUNC. CAL./ADJ. CAL/ADJ SEQ</div> <div>oCalibrate, then auto adjust</div> <div>oCalibrate, then manual adjust</div> <div><< < ^ ↓</div>
11. Exit Cal./adjustment sequence	< soft key	<div>SETUP BAL.FUNC. CAL./ADJ.</div> <div>CAL/isoTST key function</div> <div>Cal/adjustment sequence</div> <div>isoCAL function</div> <div>Print GLP/GMP adjustment record</div> <div>Parameter for external weight</div> <div><< < ^ v ></div>
12. Save settings and exit the Setup menu	< < soft key	<div>Max 6200g d= 0.01g</div> <div>0% 100%</div> <div>0.00 g</div> <div>isoTST PT1/T1</div>

Calibration/Adjustment

Internal Calibration/Adjustment

First set either **Internal cal./adjustment** or **Selection mode** (factory setting) in the Setup: Scale menu.

Inside the scale housing is a built-in motorized calibration weight.

The internal calibration/adjustment sequence is as follows:

- Select the calibration function:
Press the **isoTST** soft key twice
- > The internal calibration weight is applied automatically
- > The scale is calibrated
- > If the setting **Calibrate, then auto adjust** is selected in the Scale menu, the scale is now automatically adjusted
- > If the setting **Calibrate, then manual adjust** is selected in the Scale menu, the internal calibration routine is now ended without adjusting the scale (see “Calibration and Adjustment Sequence”, next column)
- > The internal calibration weight is removed
- > (ISO/GMP printout: see page 100)

Calibration and Adjustment Sequence

In the Setup menu, you can configure the scale so that:

- calibration is always followed automatically by adjustment
Calibrate, then auto adjust (factory setting) or
- you have the choice of ending the sequence or starting adjustment after calibration
Calibrate, then manual adjust

If no deviation is determined in calibration, or the deviation is within the tolerance limits dictated by the degree of accuracy you require, it is not necessary to adjust the scale. In this case, you can end the calibration/adjustment sequence after calibration.

There are 2 soft keys active at this point:

- **Start** to start adjustment
- **End** to end the sequence

Selecting the Calibration/Adjustment Parameter

The setting **Selection mode** must be selected in the Setup menu (factory setting).

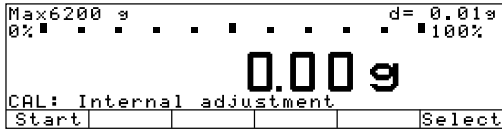

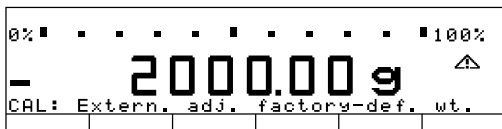
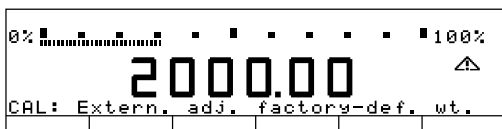
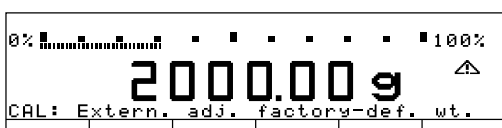
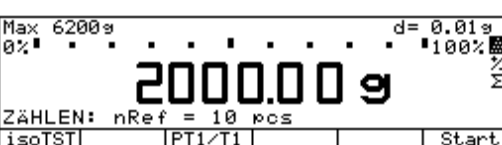
After pressing the **isoTST** soft key, you can choose from among the following settings by pressing the **Select** soft key:

- External calibration/adjustment with the preset calibration weight:
**Ext. cal./adj.;
factory-def. wt.**
- External calibration/adjustment with a calibration weight determined by the user: **Ext. cal./adj.;
user-defined wt.**
- Internal calibration/adjustment
**Internal cal./
adjustment**
- Reproducibility test **reproTEST**

- Start the desired routine:
Press the **isoTST** soft key again

In the selection mode: Perform external calibration followed by automatic adjustment with the factory-set weight

Configuration:
factory settings

Step	Key(s) (or instruction)	Display/Output
1. Select Calibration	isoTST soft key	
2. Select external calibration/adjustment with factory-defined weight (for scales of accuracy class II , only "external adjustment" is possible)	Select soft key repeatedly	
3. Start external calibration/adjustment	Start soft key	
4. Place the weight on the scale (e.g., 2,000.00 g) Minus sign -: Weight too low Plus sign + Weight too high no plus/minus sign: Weight o.k.	Place weight on scale	
This is displayed after calibration, for approx. 10 seconds: (on verified scales, the difference between the displayed weight and the true weight (mass) is displayed)		
5. Unload the scale (ISO/GMP printout: see page 100)		

Calibration/Adjustment

External Calibration/Adjustment* with a User-Defined Calibration Weight

First set either **Ext. cal./adj.:** user-defined wt. or **Selection mode** (factory setting) in the Setup: Scale menu.

You can define a weight for calibration/adjustment. External calibration/adjustment must be performed with weights that are traceable to a national standard and that have error limits which are at least 1/3 of the required tolerance of the display accuracy. The defined weight must equal at least 10% of the maximum scale capacity.

See page 41 for the external calibration/adjustment sequence. For this example, select external calibration/adjustment with a user-defined weight.

The scale has a factory-set weight value (see “Specifications”).

To reset a user-defined calibration weight to the original factory setting:

- Enter the factory-defined value manually (see “Specifications”)

Define the Calibration Weight

Step	Press key(s) (or follow instructions)	Display/Output
1. Select Setup menu	Setup	
2. Select “Balance/scale functions”	> soft key	
3. Select “Calibration/adjustment”	> soft key	
4. Select parameter for external weight	∨ soft key 5 x > soft key	
5. Select “Cal./adj. wt.”	∨ soft key	<p>3000 = last setting selected</p>
6. Enter calibration weight (e.g., 5000 g) and save	5 0 0 0	
7. Save the calibration weight	↓ soft key	
8. Exit the Setup menu	<< soft key	

isoCAL:

Automatic Calibration and Adjustment after a Change in Temperature

First set either **On and reset the application** or **On without resetting the app.** (factory setting) in the Setup: Scale menu.

The “isoCAL” display automatically begins flashing if the ambient temperature changes in relation to the temperature at the time of the last calibration/adjustment, or after a defined time interval has elapsed. The scale is telling you that it wants to adjust itself.

This automatic calibration prompt is triggered when:

- The change in temperature is greater than 10 Kelvin
- The scale status does not correspond to Setup configurations
- No number or letter input is active
- The load has not been changed within the last 2 minutes
- The scale has not been operated within the last 2 minutes
- The load on the scale does not exceed 2% of the maximum capacity

When these requirements are met, **C** is displayed in the measured value line.

If the scale is not operated and the load is not changed, internal calibration and adjustment starts after 15 seconds have elapsed.

Automatic Calibration and Adjustment at Specific Time Settings

Select either **On and reset application** or **On without resetting app.** (factory setting) in the Setup menu.

In the Setup: Input menu, you can now enter up to three different times of day for automatic calibration/adjustment. The scale will display the flashing calibration prompt (“isoCAL”). Calibration/adjustment is not performed if the scale is off or in the Setup mode at the time set for calibration.

If the scale is being operated at the time set for automatic calibration/adjustment, the calibration/adjustment sequence is prompted afterward.

If time settings are selected for automatic calibration/adjustment, the time and temperature criteria for the isoCAL function are switched off.

Automatic calibration/adjustment is prompted at fixed times when:

- The user-defined time is reached
- The scale status does not correspond to Setup configurations
- No alphanumeric input is active (e.g., equation for calculation)
- The load has not been changed within the last 2 minutes
- The scale has not been operated within the last 2 minutes
- The load on the scale does not exceed 2% of the maximum capacity

When these requirements are met, **C** is displayed in the measured value line.

If the scale is not operated and the load is not changed, internal calibration and adjustment starts after 15 seconds have elapsed.

In the Setup menu, you can configure the scale so that after calibration and adjustment

- the application program is restarted **On and reset the application**
- the application program remains at its previous status **On without resetting the app.**

Also in Setup, you can configure the scale so that it displays a calibration prompt, but does not perform the calibration/adjustment functions automatically

Only adjustment prompt

* including the Signatories of the Agreement on the European Economic Area

Calibration/Adjustment

Printout of Calibration/Adjustment Data

You can have the results of a calibration/ adjustment procedure printed out. You can configure whether the printout is generated as soon as the procedure is completed, or whether a number of calibration/ adjustment procedures (up to 50) are collected for a block printout.

Block Printout of Calibration/ Adjustment Results

With the following menu configuration under Setup: Balance/scale functions: Calibration/adjustment, you can store the data from up to 50 calibration/ adjustment procedures and have them printed on request:

- isoCAL printout
On request, from record memory

When the memory contains 50 data records:

- additional records are output automatically

If at least one block printout data record has been configured, the following soft keys are available after you press the **Cal** soft key:

Info The number of records is displayed in the text line
PrtPro Print accumulated records
DelPro Delete accumulated records; records can only be deleted after a printout has been generated.
If a password has been assigned in the Setup: Input menu, you must enter either the configured password or the General Password before you can delete the records.

For internal calibration/adjustment, the initialization mode of the procedure is displayed in the **Start** line.

```
-----
13.05.1997      09:17
      SARTORIUS
Mod.      FC6CCE-HX
Ser. no.   60419914
Ver. no.   01-35-18
ID
-----

24.04.1997      12:03
Start:      manual
Diff. +     0.01 g
External calibration
      completed

25.04.1997      12:10
Start:      isoCAL/temp
Diff. +     0.01 g
Internal adjustment
      completed
Diff. +     0.00 g

25.04.1997      18:30
Start:      Adj.time
Diff. +     0.01 g
Internal adjustment
      completed
Diff. +     0.00 g

26.04.1997      9:37
Start:      manual
Diff. +     0.01 g
Internal adjustment
      completed
Diff. +     0.00 g

27.04.1997      11:53
Start:      Ext.cal.
W ID
Nom. +     2000.00 g
Diff. +     0.01 g
External calibration
      completed
Diff. +     0.00 g

-----
13.05.1997      09:17
Name:
-----
```

GLP header

List of Calibration/Adjustment Procedures:

Example 1:
External calibration

Example 2:
isoCAL triggered by difference
in temperature

Example 3:
isoCAL at defined time

Example 4:
Internal calibration/adjustment
triggered manually

Example 5:
External calibration/adjustment

GLP footer

Determination of the Repeatability (reproTEST)

Definition

Repeatability is the ability of the scale to display identical readouts when it is loaded several times with the same weight under constant ambient conditions (also called “reproducibility”).

The standard deviation for a given number of measurements is used to quantify the repeatability.

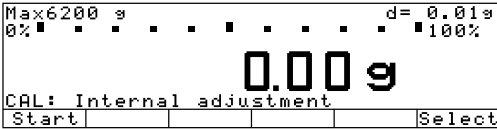
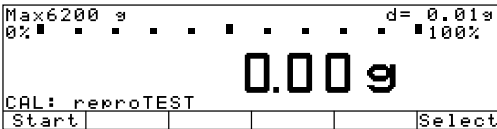
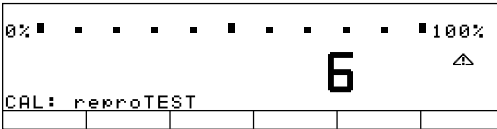
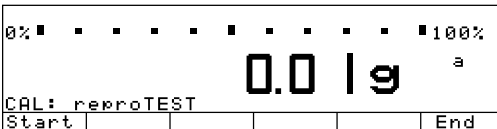
Purpose

The “reproTEST” function automatically determines the repeatability of results (based on 6 individual measurements). In this way, the scale determines one of the most important quantities in relation to the place of installation. The results are displayed with the scale’s accuracy.

Preparation

- Turn on the scale: press **1/0**
 - > The Sartorius logo is displayed
 - > The scale performs a self-test
- Select reproTEST in the Setup menu: Press **Setup**
- Select Calibration/adjustment: CAL-key function: press the **> v** keys
- Select either reproTEST or Selection mode (factory setting): see “Configuration”
- Exit the Setup menu: Press the **<<** soft key

Check the Reproducibility of the Scale

Step	Key(s) (or instruction)	Display/Output
1. If reproTEST is set: and proceed with step 4. If Selection mode is set:	isoTST soft key isoTST soft key	
2. Select reproTEST	Select soft key	
3. Start reproTEST	Start soft key	
4. Number of measurements is displayed; 6 measurements will now be performed The standard deviation is displayed		 
5. End reproTEST or restart reproTEST	End soft key Start soft key	

Application Programs

Soft Key Functions

Start Start application program

Weight Toggle to basic weighing functions

Auto-Start Application When the Power Goes On

In the Setup menu, you can select whether the application that is active before you turn off the power will automatically start when the power is turned on again (Setup: Application parameters: Auto-start app. when power goes on: On).

Using Verified Scales as Legal Measuring Instruments in the EU*: All application programs can be selected on scales verified for use in legal metrology. Non-metric values can be indicated as follows:

- Percent = %
- Piece counting (Counting) = pcs
- Computed value = o, Δ

Toggle between Two Weight Units

Purpose

With this application program you can switch the display of a weight value back and forth between two weight units by pressing a soft key.

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formula-tion, statistics).

Available Features

- Toggling the displayed weight
- Other features as for the basic weighing function

Factory Settings

Weight unit 1: **Grams** /g

Weight unit 2: **Kilograms** /Kg

Preparation

Scales used as legal measuring instruments: grams and kilograms are the only weight units available.

Standard scales: The following weight units are available in both ranges:

Unit	Conversion factor	Display/ Printout	Line for metrological data
Grams	1.00000000000	g	g
Kilograms	0.00100000000	kg	kg
Carats	5.00000000000	ct	ct
Pounds	0.00220462260	lb	lb
Ounces	0.03527396200	oz	oz
Troy ounces	0.03215074700	ozt	ozt
Hong Kong tael	0.02671725000	tlh	tlh
Singapore tael	0.02645544638	tls	tls
Taiwanese tael	0.02666666000	tlt	tlt
Grains	15.43235835000	GN	GN
Pennyweights	0.64301493100	dwt	dwt
Milligrams	1000.00000000000	mg	mg
Parts per pound	1.12876677120	/lb	lb
Chinese tael	0.02645547175	tlc	tlc
Mommes	0.26670000000	mom	M
Austrian carats	5.00000000000	K	K
Tola	0.08573333810	tol	tol
Baht	0.06578947437	bat	bat
Mesghal	0.21700000000	MS	MS

- Turn on the scale: Press \square
> Sartorius logo is displayed
- Select the “Toggle weight units” program in the Setup menu: press \square
- Select **Application parameters**: press the ∇ soft key 2x, then the \triangleright soft key once
- **Application 1 (basic settings)**: press the \triangleright soft key
- Select **Toggle wt. units**: Δ or ∇ soft key (repeatedly)
- Confirm **Toggle wt. units**: \triangleright soft key
- Select and confirm:
 - **Weight unit 1**: see above
 - **Weight unit 2**: see abovesee also the “Application Menu (Overview)” in the chapter entitled “Configuration”
- Save settings and exit the Setup menu: Press the \ll soft key

* including the Signatories of the Agreement on the European Economic Area

Additional Functions

In addition to the functions for:

- alphanumeric input,
- taring (not during alphanumeric input),
- printing (NUM print; S ID),

you can also access the following functions from this application:

Calibration/Adjustment

- Press the **i s o T S T** soft key
- > See “Calibration/Adjustment” for further instructions

Toggling to the Next Application

- Press **(Q)**
- > See the section on the corresponding application program for further instructions

Setup (setting parameters)

- Press **Setup**
- > See “Configuration” for further instructions

Turning Off the Scale

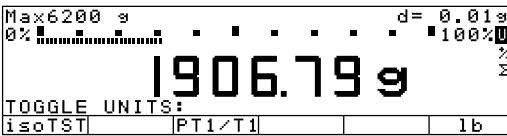
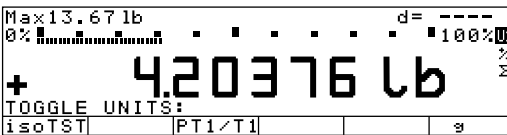
- Press **(U)**
- > The scale shuts off

Practical Example

Toggle the Display From Grams [g] (1st Unit) to Pounds [lb] (2nd Unit)

Settings (changes in the factory settings required for this example):

Setup: App: Application 1: Toggle wt. units: Weight unit 2: Pounds/lb


Step	Key (or instruction)	Display/Output
1. Delete previous setting if necessary (U1: Weight unit 1)	(CF)	
2. Change weight unit to Pounds [lb] (U2: Weight unit 2)	1 b soft key	
3. Change weight unit to Grams [g]	g soft key	

Purpose

With the Counting application, you can determine the number of pieces of approximately equal weight.

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics).

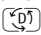Features

- Optional scale configuration in Setup for automatically initializing this application and loading the most recent reference sample quantity “nRef” and average piece weight “wRef” when you switch on the scale (this is the automatic setting when the scale is initialized; Setup menu: Application parameters: Auto-start app. when power goes on: On).
- Reference sample quantity “nRef” entered manually
- Average piece weight “wRef” entered manually
- Storage of the current weight value for the preset reference sample quantity “nRef”, to be loaded when the Counting program is initialized
- Setting the accuracy when the reference sample weight “wRef” is stored for calculating a piece count
- Automatic output of the quantity and sample weight via the data interface port after initialization or reference sample updating while running the Counting program (Printout: Application-defined printout: Auto print upon initialization: All values)
- Toggling between piece count and weight value by pressing the **Count.** or **Weigh.** soft key
- Toggling between counting and other applications using the  key (for example, checkweighing)

Factory Settings

Accuracy when calculating piece weight:
Display accuracy

Reference sample updating:
Automatic

Soft Key Functions

- | | |
|---------------|---|
| nRef | Store value input as reference sample quantity |
| wRef | Store input value as reference sample weight |
| Update | Reference updating criteria met; reference updating can be performed to optimize the accuracy |
| Count. | Toggle to the Counting application |
| Weigh. | Toggle to the Weighing mode |
| Start | Storage of the current weight value for the preset reference sample quantity |

Preparation

To calculate a piece count, the average weight of one piece must be known. This average piece weight can be entered into the Counting program in one of three ways:

- Enter the average piece weight using the numeric keys and store it;
- The last reference sample quantity entered is loaded and displayed when you turn on the scale. Place the same number of parts on the scale and initialize the Counting program;
- When the automatic initialization parameter (see previous page) is on (Setup: Printout: Application-defined output: Auto print upon initialization: All values), the scale goes into the "Counting" mode when you turn it on and loads the last average piece weight and corresponding reference sample quantity that were entered or calculated.

Reference Sample Updating

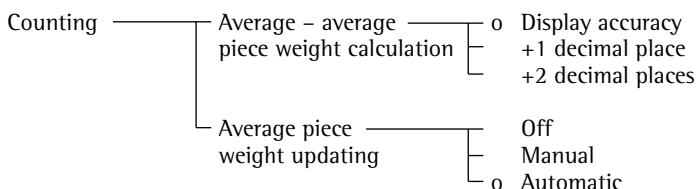
You can have the average piece weight updated during counting (with the piece count displayed) if "AWP update" is set to "manual" or "automatic" in the Setup menu. Manual updating can only be performed when the **Update** soft key is displayed. Reference sample updating must be completed before using an application program from Application 3.

The **Update** soft key is displayed when:

- the scale has reached stability (stability symbol displayed)
- the current piece count is less than double the original piece count
- the current piece count is less than 100
- the internally calculated piece count (e.g., 17.24) differs from the nearest whole number (here: 17) by less than 0.3

Reference updating can be repeated several times with an approximately doubled piece count.

- To perform reference updating: Press the **Update** soft key
- Turn on the scale: Press **1/0**
 - > Sartorius logo is displayed, self-test is performed
- Select the "Counting" program in the Setup menu: Press **Setup**
- Select **Application parameters**: press the **v** soft key 2x, then the **>** soft key once
- Select **Application 1 (basic settings)**: press the **>** soft key
- Select **Counting**: **^** or **v** soft key, repeatedly
- Confirm **Counting**: **>** soft key



o = factory setting

see also the "Application Menu (Overview)" in the chapter entitled "Configuration"

- Save settings and exit the Setup menu: Press the **<<** soft key

Additional Functions

In addition to the functions for:

- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:

Calibration/Adjustment

- Press the **isOTST** soft key
- > See "Calibration/Adjustment" for further instructions

Toggling to the Next Application

- Press **↔**
- > See the section on the corresponding application program for further instructions

Setup (setting parameters)

- Press **Setup**
- > See "Configuration" for further instructions

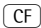
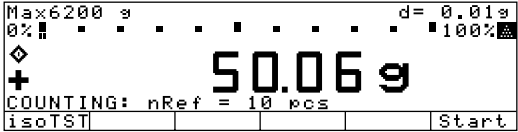

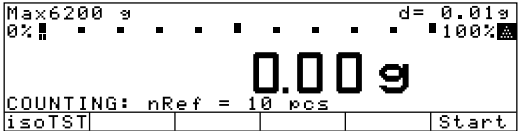
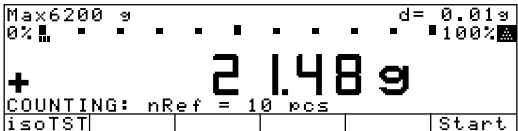
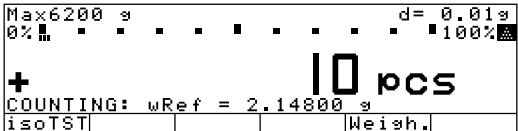
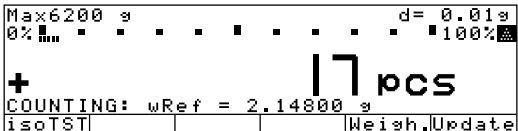
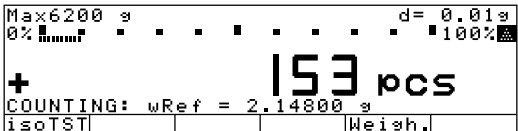

Turning Off the Scale

- Press **1/0**
- > The scale shuts off

Counting

Practical Example
Determining an Unknown Piece Count; Weighing in the Preset Reference Sample Quantity

Settings (changes in the factory settings required for this example):
Setup: Application parameters: Application 1: Counting: Average piece weight updating: Manual
Setup: Printout: Application-defined output: Autoprint upon initialization: All values

Step	Key (or instruction)	Display/Output
1. Delete previous setting if necessary		
2. Prepare a container for the parts to be counted	Place the empty container on the scale	
3. Tare the scale		
4. Place the reference sample quantity on the scale (example: 10 pcs)	Place the displayed number of parts in the container	
5. Determine the average piece weight (number of decimal places displayed depends on the scale model)	Start soft key	 <div><div>nRef + 10 pcs</div><div>wRef + 2.14800 g</div></div>
6. If necessary, increase the number of parts and perform reference sample updating (example: 7 additional pieces)	Place additional parts in the container Update soft key	 <div><div>nRef + 17 pcs</div><div>wRef + 2.14800 g</div></div>
7. Weigh uncounted parts	Place parts to be counted in the container	
8. If desired, print total piece count (here: 153 pcs)		<div>Qnt + 153 pcs</div>

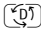
Weighing in Percent %

Purpose

This application program allows you to obtain weight readouts in percent which are in proportion to a reference weight. Alternatively, you can have the value displayed as a difference in percent between the weight on the scale and the reference weight, or as a special ratio 1 or ratio 2.

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formula-tion, statistics).

Available Features

- Reference percentage “pRef” loaded from long-term memory when you turn on the scale
- Optional scale configuration in Setup for automatically initializing this application and loading the most recent reference percentage “pRef” entered with reference weight “Wxx%” when you turn on the scale
- Value displayed as:
 - Residual quantity (portion)
 - Difference (deviation)
 - Ratio1
 - Ratio2depending on the Setup menu code selected.
- Reference percentage “nRef” entered manually
- Store the current weight as the reference percentage weight “Wxx%” for initializing the weighing-in-percent application program
- Reference weight “Wxx%” entered manually
- Storage parameter (rounding-off factor) for storing the reference weight “W100%” in percentage calculation can be configured
- Configuration of decimal places displayed with a percentage
- Optional configuration for having the reference weight “Wxx%” and reference percentage automatically output via the data interface port after initialization of the weighing-in-percent program (print application parameters)
- Toggle the display between percentage and weight readout by pressing the **W e i g h .** soft key
- Toggle between the weighing-in-percent program and other applications (e.g., checkweighing) by pressing 

Factory Settings

Storage parameter:

Display accuracy

Digits displayed with percentage:

2 digits

Display calculated value:

Residue

Soft Key Functions



- | | |
|---------------|---|
| pRef | Store value input as reference percentage |
| Wxx% | Store input value as reference sample weight |
| Perc. | Toggle to the Weighing-in-percent application |
| New | Store next value |
| Weigh. | Toggle to the Weighing mode |
| Start | Store the current weight value for the preset reference sample quantity |


Preparation

To calculate a value in percent, the reference percentage must be known.

This value can be entered into the weighing-in-percent program in one of three ways:

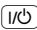- The last reference percentage entered is loaded and displayed when you turn on the scale.
Place the corresponding weight on the scale and initialize the weighing-in-percent program;
- With automatic initialization switched on (see previous page), the scale goes into the “weighing in percent” mode when you turn on the power and loads the last reference percentage entered as well as the corresponding reference weight (Setup: Printout: Application-defined output: Auto print upon initialization: All values);
- Enter the reference weight using the numeric keys and store it (W××% soft key).

- Turn on the scale: press 
 - > Sartorius logo is displayed, self-test is performed
- Select the “Weighing in percent” application in the Setup menu: press 
- Select **Application parameters**: press the ∇ soft key 2x, then the \triangleright soft key once
- Select **Application 1 (basic settings)**: press the \triangleright soft key
- Select **Percent weigh.**: Δ or ∇ soft key repeatedly
- Confirm **Percent weigh.**: \triangleright soft key
- Select and confirm:

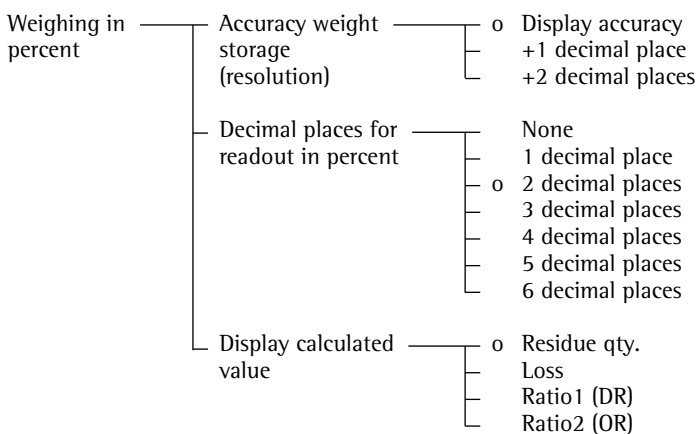

o = factory setting

see also the “Application Menu (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: Press soft key $\triangleleft \triangleleft$

Equations

The following equations are used for the corresponding calculations:

$$\text{Residue} = \text{Current weight} \div 100\% \text{ weight} \times 100\%$$

$$\text{Loss} = (\text{Current weight} - 100\% \text{ weight}) \div 100\% \text{ weight} \times 100\%$$

$$\text{Ratio1} = (100\% \text{ weight} - \text{current weight}) \div \text{current weight} \times 100\%$$

$$\text{Ratio2} = 100\% \text{ weight} \div \text{current weight} \times 100\%$$

Additional Functions

In addition to the functions for:

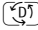
- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:


Calibration/Adjustment

- Press the **ISO TEST** soft key
 - > See “Calibration/Adjustment” for further instructions


Toggle to the Next Application

- Press 
 - > See the section on the corresponding application program for further instructions

Setup (setting parameters)

- Press 
 - > See “Configuration” for further instructions

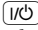Turning Off the Scale

- Press 
 - > The scale shuts off

Examples

Practical Example: Weighing in Percent with Reference Weight Taken from Weight on Scale

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Application 1: Weighing in percent

Setup: Printout: Application-defined output: Autoprint upon initialization: All values

Step	Key (or instruction)	Display/Output
1. Delete previous setting if necessary		
2. Prepare a container for the parts	Place the empty container on the scale	
3. Tare the scale		
4. Place the reference weight on the scale (here: 1821.48 g = 100%)	Place weight equal to reference weight in the container	
5. Initialize the scale	Start soft key	<div> <p>pRef + 100 %</p> <p>Wxx% + 1821.48 g</p> </div>
6. Unload the scale	Remove reference sample from the container	
7. Determine the percentage of an unknown weight	Place sample to be measured in the container	
8. If desired, print percentage (here: 98.37%)		<div> <p>Prc + 98.37 %</p> </div>

Animal Weighing


Purpose

Use this program to determine the weights of unstable samples (e.g., live animals) or to determine weights under unstable ambient conditions. In this program, the scale calculates the weight as the average of a defined number of individual weighing operations. These weighing operations are also known as “subweighing operations.”

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics).

Available Features

- Animal weighing started manually or automatically
- Automatic start:
 - when a defined threshold has been exceeded (Minimum load threshold: None; 10; 20; ...; 500; 1,000 display increments)
 - when three successive subweights lie within a user-defined tolerance range (calm; normal; active; 0.1%; 0.2%; ...; 50%; 100% of the animal/object)
- Manual start:
 - also possible when the load is under the minimum load threshold
 - when three successive subweights lie within a user-defined tolerance range (calm; normal; active; 0.1%; 0.2%; ...; 50%; 100% of the animal/object)
 - Optional scale configuration in the Setup menu for automatically initializing this application when you turn on the scale (Setup: Printout: Auto-start application when power goes on: On)

- Number of weighing operations for calculation of an average (**mDef**) can be set before the beginning of each animal weighing operation
- The factor for calculation of the result can be set before the beginning of each animal weighing operation
- The number of subweighs remaining to be performed is indicated in the text display during weighing
- Arithmetic average displayed as a result in the pre-set weight unit (identified by the  symbol).
- Optional multiplication of the arithmetic average by a user-defined factor **Mu1**.
A circle “o” is displayed as weight unit and **Mu1 = xxx** is shown in the text line
- Toggling between the weight and the calculated value by pressing the **xNet** soft key and the **xRes** soft key
- Automatic output of results via the interface port:
 - Number of weighing operations **mDef**
 - Multiplication factor **Mu1**
- Automatic output of results (printout) via the interface port:
 - Weighing result **xNet**
 - Calculated result **xRes**
 The following options have to be set:
Printout: Application-defined output:
Auto print upon initialization: All values
- The unload threshold is equal to one-half the minimum scale capacity
- Return to weighing mode by unloading the scale; i.e., when the load is below the unload threshold

Factory Settings

Animal activity:
5% of the animal/object

Start: **Automatic**

Minimum load for automatic storage:
100 display increments

Decimal places in result display:
2 decimal places

Printout:
Average weight only

Soft Key Functions

New	Automatic start: <ul style="list-style-type: none"> – Unload scale and weigh next animal, if desired – Press key to start next subweigh Manual start: Start next subweigh
mDef	Store user-defined number of subweighs for averaging
Mu1	Store user-defined factor as multiplication factor for calculated the arithmetic mean
xNet	Toggle to the animal weighing application
xRes	Toggle to the calculated animal weighing result
Start	Activate Animal weighing

Printout for Animal Weighing

Upon completion of the averaging process, you can have the results printed out automatically. You can also have both the weight and the calculated result printed.

```

mDef          10
Mul           0.00347
xNet  +      153.00 g
xRes  +        5.30 o
  
```

mDef: Number of subweighing operations for averaging

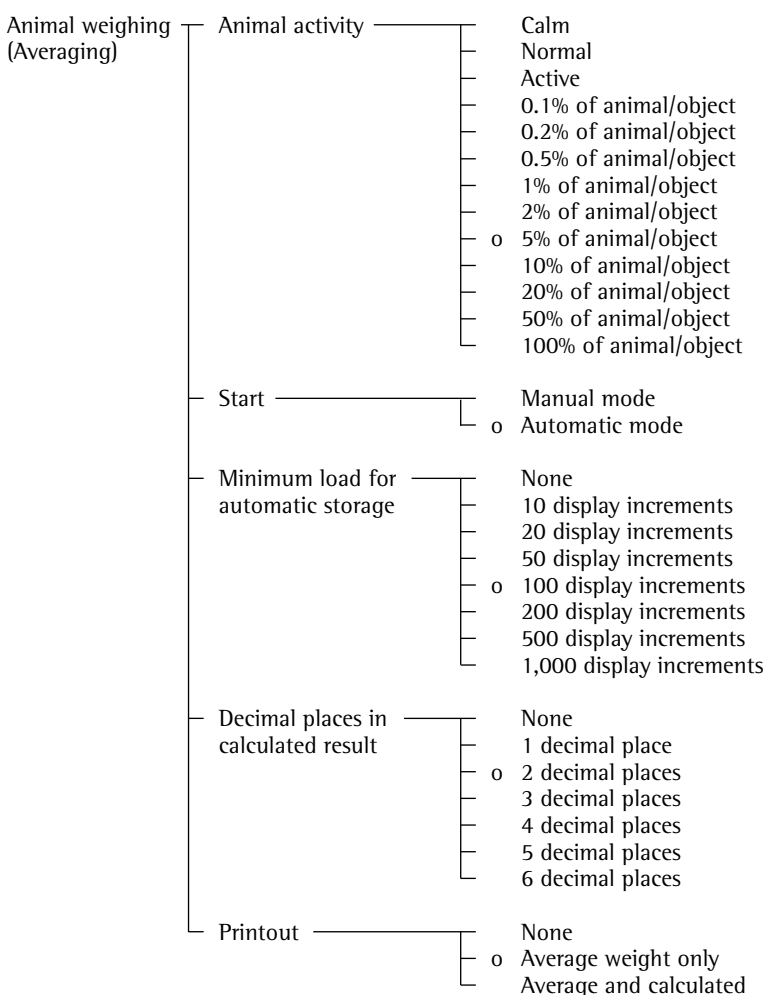
Mul: Multiplication factor

xNet: Result of averaging

xRes: Calculated result

Preparation

- Turn on the scale: press **[I/O]**
> Sartorius logo is displayed
- Select the “Animal weighing” application in the Setup menu: press **[Setup]**
- Select **Application parameters**: press the **▼** soft key 2x, then the **>** soft key once
- Select **Application 1 (basic settings)**: press the **>** soft key
- Select **Animal weigh.**: **▲** or **▼** soft key repeatedly
- Confirm **Animal weigh.**: **>** soft key



o = factory setting

see also the “Application Menu (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: press the **<<** soft key

Additional Functions

In addition to the functions for:

- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:

Calibration/Adjustment

- Press the **isoTST** soft key
> See “Calibration/Adjustment” for further instructions

Toggling to the Next Application

- Press **[↔]**
> See the section on the corresponding application program for further instructions

Setup (setting parameters)

- Press **[Setup]**
> See “Configuration” for further instructions

Turning Off the Scale

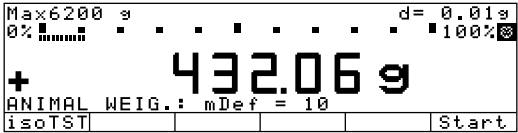
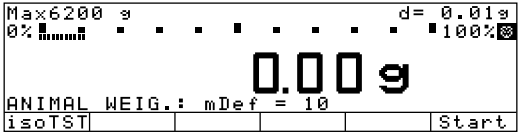
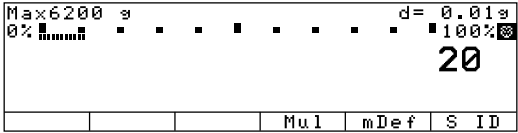
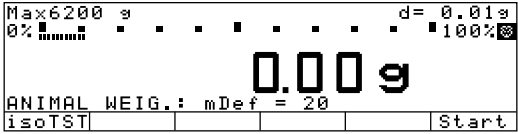
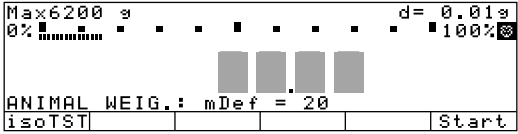
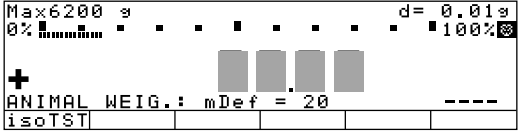
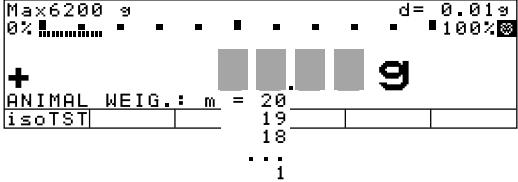
- Press **[I/O]**
> The scale shuts off

Animal Weighing

Practical Example

Determining Animal Weight With Automatic Start of 20 Subweighing Operations for Averaging; Automatic Printout of the Number of Subweighing Operations and of the Animal Weight

Settings (changes in the factory settings required for this example):
Setup: Application parameters: Application 1: Animal weighing: Animal activity: Active
Setup: Application parameters: Application 1: Animal weighing: Decimal places in calculated result display: 2 decimal places
Setup: Application parameters: Application 1: Animal weighing: Printout: Average and calculated values
Setup: Printout: Application-defined output: Auto print upon initialization: All values

Step	Key (or instruction)	Display/Output
1. Delete previous setting if necessary	CF	
2. Prepare a container (cage)	Place empty cage on the scale	
3. Tare the scale	Tare	
4. Enter number of subweighing operations for averaging	2 0	
5. Save number	mDef soft key	
6. Weigh the first animal	Place 1st animal in cage	weight value fluctuates due to animal activity 
7. Start automatic animal weighing	Start soft key	
The scale delays starting the subweighing operation until three successive subweights lie within the range defined for an "active" animal	When this criterion is met, the subweighing series begins	

Step	Key (or instruction)	Display/Output
After 20 subweighing operations the arithmetic average (xNet) is displayed		<div>Max6200 g d= 0.01g</div> <div>0% ██████ ■ ■ ■ ■ ■ ■ ■ ■ 100%</div> <div>+ 69.72 g</div> <div>ANIMAL WEIG.: xNet</div> <div>isoTST xRes New</div>
(mDef: no. of subweighs Mul: calculation factor xNet: arithm. average, net value)		<div>mDef 20</div> <div>Mul 1</div> <div>xNet + 69.72 g</div> <div>xRes + 69.72 o</div>
8. Unload the scale	Remove animal from cage	<div>Max6200 g d= 0.01g</div> <div>0% ██████ ■ ■ ■ ■ ■ ■ ■ ■ 100%</div> <div>0.00 g</div> <div>ANIMAL WEIG.: mDef = 20 Auto</div> <div>isoTST</div>
9. If desired, weigh next animal	Place animal in cage	<div>Max6200 g d= 0.01g</div> <div>0% ██████ ■ ■ ■ ■ ■ ■ ■ ■ 100%</div> <div>+ ■ ■ ■ ■</div> <div>ANIMAL WEIG.: mDef = 20 ----</div> <div>isoTST</div>
Next weighing series begins automatically		<div>Max6200 g d= 0.01g</div> <div>0% ██████ ■ ■ ■ ■ ■ ■ ■ ■ 100%</div> <div>+ ■ ■ ■ ■ g</div> <div>ANIMAL WEIG.: m = 20</div> <div>isoTST 19</div> <div> 18</div> <div> ...</div> <div> 1</div>

Recalculation

Purpose

With this application program you can compensate for over-poured components in formulation.

If a component is over-poured when weighing in the individual formulation components, the mixture already poured cannot be used in its current composition. To avoid having to discard the materials weighed, you can adjust the proportions of the formulation to compensate for the over-pour.

When you use this application, the recalculation procedure is mainly performed by the scale.

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics).

Available Features

- Individual components (up to 99) weighed in with a readout showing from “0” to the desired component weight
- Transaction counter shows the next component expected
- Weighed components are stored, followed by automatic printout and taring
- Additive weighing of components with printout
- Toggle the display between component weight and total formulation weight (additive mode) after first component is stored
- Stored component weight displayed as true net weight for 2 seconds
- Enter a divisor before or during component weighing. For example, if the formulation has a total weight of 100 g, enter the divisor 10 to weigh in a total formulation of 1,000 g
- If a component is over-poured, you can use the recalculation function to change the amount of this component indicated in the formulation by using plus or minus keys or numeric input. A factor is then calculated by which all components amounts will then adjusted
- Recalculation factor displayed in the text line, with a warning symbol if the factor is not equal to 1
- All components displayed with number and the amount (by weight) to be added in follow-on filling. Components displayed in sequence by the scale
- Display of actual net weight during follow-on filling
- After the amounts of the components already weighed have been corrected, weighing continues according to the adjusted formulation amount. The readout is recalculated (updated) according to the divisor
- You can repeat the over-pour correction procedure as often as necessary, in case other components are over-poured
- After follow-on (corrective) filling, the total amount differs from that given for the formulation, but the proportion of components in relation to each other is the same
- You can have the weight printed after each measurement
- Choose whether the current component weight or tare value is printed after each measurement
- Individual component weights are printed as “Compxx.”
- Press **[CF]** to exit the application program. The component memory is cleared and the sum of components printed as “S-Comp.”
- Toggle between the recalculation program and other applications (e.g., checkweighing) by pressing **[<D>]**.

Factory Settings of the Parameters

Print application parameters (automatic output of application parameters):

All values



Line format:

for other apps/GLP
(22 characters)

Soft Key Functions

Comp.xx	Store component
Add.xx	Store component in additive weighing mode
Div.	Store divisor before or during component weighing
Recalc	Start correction procedure for recalculation
→Add./→Comp.	Toggle display between component weight and total weight (additive mode)
Comp.	Store numeric input for recalculation
Minus	Set value given for the formulation
Plus	Set value given for the formulation

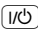Preparation

- Turn on the scale: press 
- > The Sartorius logo is displayed
- Select the “Recalculation” application in the Setup menu: press 
- Select **Application parameters**: press the \vee soft key 2x, then the \triangleright soft key once
- Select **Application 1 (basic settings)**: press the \triangleright soft key
- Select **Recalculation**: \wedge or \vee soft key repeatedly
- Confirm **Recalculation**: \leftarrow soft key
- Save settings and exit the Setup menu: press the $\leftarrow\leftarrow$ soft key

Additional Functions

In addition to the functions for:

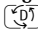
- alphanumeric input,
- taring (not during alphanumeric input), and
- printing,

you can also access the following functions from this application:


Calibration/Adjustment

- Press the **ISO TST** soft key
- > See “Calibration/Adjustment” for further instructions


Tagging to the Next Application

- Press 
- > See the section on the corresponding application program for further instructions

Setup (setting parameters)

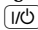- Press 
- > See “Configuration” for further instructions

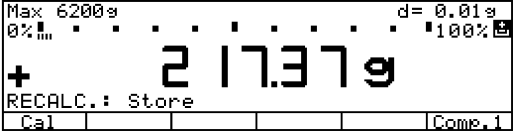
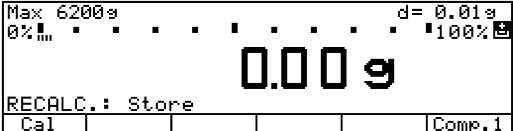
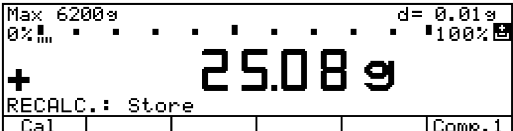
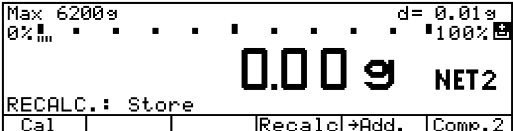

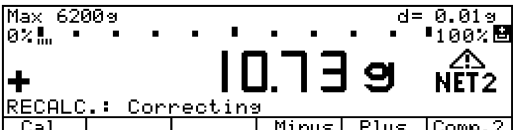
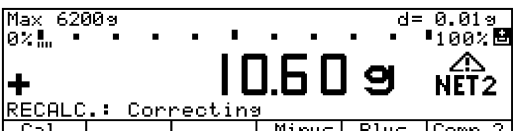
Turning Off the Scale

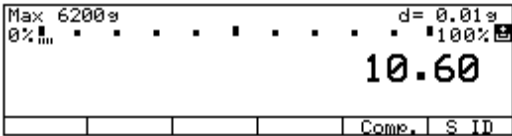
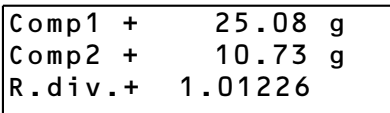
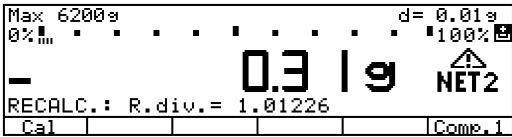
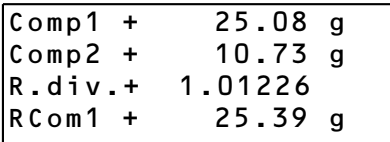
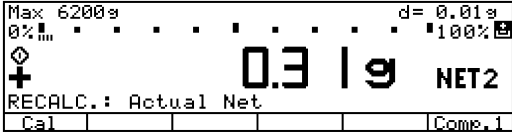

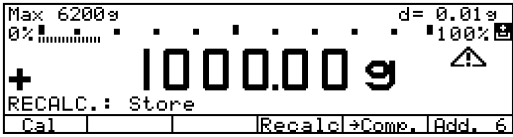
- Press 
- > The scale shuts off

Recalculation

Practical Example
When weighing in formulation components, the second component is over-poured.

Settings (changes in the factory settings required for this example):
Setup: App(lication parameters): Application 1: Recalculation

Step	Key (or instruction)	Display/Output
1. Delete previously stored values, if necessary	CF	
2. Place container for filling components on the scale	Place empty container on the scale	
3. Tare	Tare	
4. Add the first component	Weigh the first component into the container	
5. Store component	Press the Comp. 1	<div>Comp1 + 25.08 g</div> 
6. Add the second component	Weigh the second component into the container	
7. Start recalculation, because 10.73 g were poured rather than 10.60 g	Recalc soft key	
8. Either press the minus key to correct the value ...	Minus soft key repeatedly	

Step	Key (or instruction)	Display/Output
... or enter the desired value	1 0 . 6 0	
9. Confirm the new value	Comp. soft key	
Follow-on filling amount for first component is displayed		
10. Follow-on filling of 1st component	Weigh the first component up to 0	
and store	Comp. 1 soft key	
The true net value is displayed for 2 seconds		
11. Weigh in further components, if called for in the formulation	Repeat steps 4 and 5 as needed	
12. Toggle to the additive mode, if required	→Add. soft key	
13. Add further components, as required ... (here, e.g., up to the total weight of the formulation: 1,000 g)	Add components to container	

Recalculation

Step	Key (or instruction)	Display/Output
14. ... and store (here, e.g., the 6th component)	Add . 6 soft key	<div>Comp1 + 25.08 g</div> <div>Comp2 + 10.73 g</div> <div>R.div.+ 1.01226</div> <div>RCom1 + 25.39 g</div> <div>Comp3 + 22.03 g</div> <div>Comp4 + 31.49 g</div> <div>Comp5 + 107.50 g</div> <div>Comp6 + 812.61 g</div>
The true net value (of the 6th component) is displayed for 2 seconds		<div>Max 6200g d= 0.01g</div> <div>0% 100%</div> <div>+ 812.61g NET2</div> <div>RECALC.: Actual Net</div> <div>Cal Comp. 6</div>
Then the total weight is displayed		<div>Max 6200g d= 0.01g</div> <div>0% 100%</div> <div>+ 1000.00g NET2</div> <div>RECALC.: Store</div> <div>Cal Recalc+Comp. Add. 7</div>
15. End the weighing procedure Total weight is printed	CF	<div>Comp1 + 25.08 g</div> <div>Comp2 + 10.73 g</div> <div>R.div.+ 1.01226</div> <div>RCom1 + 25.39 g</div> <div>Comp3 + 22.03 g</div> <div>Comp4 + 31.49 g</div> <div>Comp5 + 107.50 g</div> <div>Comp6 + 812.61 g</div> <div>Tot.cp+ 1009.75 g</div>
Total weight is displayed Component memory is cleared		<div>Max 6200g d= 0.01g</div> <div>0% 100%</div> <div>+ 1009.75g</div> <div>RECALC.: Store</div> <div>Cal Comp. 1</div>


Calculation

Purpose

With this application program you can calculate a weight value using an algebraic equation. This can be used, for example, to determine the gsm weight (grams per square meter) of paper.

You can use this application program in combination with a program chosen from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics).

Available Features

- You can store an equation and configure the Setup menu to initialize this program automatically with the stored equation (Setup: ... Auto-start application when power goes on: On)
- The  symbol is displayed to indicate a calculated value. The equation used is displayed in the text line
- If no equation was entered, the weight value is displayed
- Toggle between the weight readout, equation input and display of the calculated result by pressing the **Start** or **Weigh** soft key (or press **CF** to toggle between weight and calculated value)
- There are four operators (+, -, *, /) and one factor (weight value) available when you enter an equation
- Max. equation length: 28 characters
- Pressing **CF** will delete either the equation or the last character entered, depending on the configuration in the Setup menu (Setup: ... Keyboard: CF function for input: Delete last character; see also "Configuration")
- The calculated result is displayed with the number of decimal places configured in the Setup menu. Not all decimal places are displayed if the result is longer than the display allows. If there are more digits before the decimal point than the display can show, an error message is displayed.
- The equation is stored in non-volatile memory

Factory Settings

Decimal places in calculated result:
2 decimal places

Soft Key Functions

Equat. Toggle to equation

- +** Enter an addition operator in the equation
- Enter a subtraction operator in the equation
- *** Enter a multiplication operator in the equation
- /** Enter a division operator in the equation

Start Start calculation

Weigh. Toggle to the weighing mode

Weight Enter a weight operand in the equation


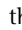




Printout for Calculation

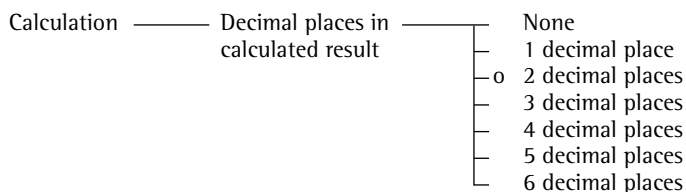
The calculation result is printed.

Res + 693.88 o

Res: Result of calculation with equation

Preparation

- Turn on the scale: press **ON**
- > Sartorius logo is displayed
- Select the "Calculation" application program in the Setup menu: press **Setup**
- Select **Application parameters**: press the  soft key 2x, then the  soft key once
- Select **Application 1 (basic settings)**: press the  soft key
- Select **Calculation**:  or  soft key repeatedly
- Confirm **Calculation**:  soft key



o = factory setting

see also the "Application Menu (Overview)" in the chapter entitled "Configuration"

- Save settings and exit the Setup menu: press  soft key

Additional Functions

In addition to the functions for:

- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:

Calibration/Adjustment

- Press the **ISO TEST** soft key
- > See "Calibration/Adjustment" for further instructions

Toggling to the Next Application

- Press **DP**
- > See the section on the corresponding application program for further instructions

Setup (setting parameters)

- Press **Setup**
- > See "Configuring the Scale" for further instructions

Turning Off the Scale

- Press **OFF**
- > The scale shuts off

Calculation

Practical Example

Calculate the gsm weight of paper: determine the gsm of a sheet of A4 paper with the dimensions $0.210\text{ m} \times 0.297\text{ m} = 0.06237\text{ m}^2$. The gsm weight is a product of the division of the weight by the surface area.

Settings (changes in the factory settings required for this example):
Setup: App: Application 1: Calculation

Step	Key (or instruction)	Display/Output
1. Turn on the scale and configure the settings as indicated above		
2. Delete previous setting if necessary		
3. Tare the scale		
4. Select equation input	Equat. soft key	
5. Enter weight value Enter division sign Enter the surface area of a sheet of A4 paper	Weight soft key / soft key 	
6. Turn on the calculated result display	Start soft key	
7. Determine the gsm weight	Place A4 sheet on the scale	

Checkweighing ✂

Purpose

This program is used to check whether a sample corresponds to a pre-set target value or is within a specific tolerance range. In addition to the display in the measured value line, the results are shown on the bar graph and can be routed through the interface port via control lines for further electronic processing.

You can use this application in combination with any program chosen from Application 1 (such as counting, weighing in percent) and one from Application 3 (totalizing, formulation, statistics).

Available Features

- Optional configuration in the Setup menu for long-term storage of target value and tolerance limits
- Optional scale configuration in the Setup menu for automatically initializing this application and loading the values stored in long-term memory for the target value and the upper and lower tolerance limits when you turn on the scale
- You can perform checkweighing
 - without entering a target value, but only upper and lower tolerance limits;
 - as differential checkweighing;
 - with symmetric or asymmetric limits which can be entered as percentages
- Enter the target value and limits by placing a load on the scale or using the numeric keys
- Control parameter in entering target and tolerance values, so that the upper limit \geq the target \geq the lower limit ≥ 1 display increment
- Accuracy of a weight readout or keyboard input as target/tolerance values corresponds to the display accuracy

- Optional scale configuration in the Setup menu for automatic output to the interface port (print application parameters) of the target value and tolerance limits when initialization is completed (... Auto print upon initialization: All values).
- Control range for the scale's data output port lines is 30% to 170% of the target value
- Optional configuration in the Setup menu for activation of control lines dependent on weight value (weight value within checkweighing range, stability reached)
- Toggling the display between weight readout and checkweighing display by pressing the corresponding soft key. If the weight value exceeds tolerances, the line for measured values shows the weight while the checkweighing display shows »LL« for »too low« or »HH« for »too high.«
- Press the **Show** soft key to display target value and tolerance limits in the text line after initializing the application.
- Weight value in bar graph displayed in relation to upper and lower limits and target value
- »OK« transaction counter displayed in the text line (e.g., **n = 4**), if selected (Checkweighing: Automatic printout of OK values: Yes). This counter shows the number of measured values that lie within the tolerance range.
- Optional automatic printout of the weight value when it is within the control range at stability

After an automatic printout, the printing of OK values is blocked. Before you can generate the next printout, you must unblock the scale by unloading it (weight must be under 30% of the target) or by placing a load on the scale (bringing the weight up to at least 170% of the target).
- Initialization parameters are overwritten after the **Param.** soft key has been pressed
- Press **[CF]** to delete the initialization parameters and end the Checkweighing program

Factory Settings

Activation of port lines:
Within checkweighing range

Type of checkweighing input:
Target, minimum, maximum weight

Weight display mode:
Absolute value

Automatic printout of OK values: **No**

Soft Key Functions

Param.	Begin input of target and tolerance values
Show	Display target and tolerance values in succession during checkweighing
LLHH	Toggle to checkweighing display (»LL« for too light and »HH« for too heavy)
Diff.	Display difference between current value and target
Net	Display net weight
Start	Start checkweighing

Auto Print Checkweighing

With the over/under checkweighing application, you can have the result printed automatically as soon as the weight lies within a defined range.

N	+	153.00 g
Setp	+	180.03 g
Min	+	160.05 g
Max	+	200.06 g
N	+	165.14 g

N:	Net weight
Setp:	Target weight
Min:	Lower limit
Max:	Upper limit
N:	Printout of "OK" values

Checkweighing

Preparation

The checkweighing program often requires a target value for comparison with the current value. This target has a tolerance range, which is defined by absolute weight values. The tolerance range is defined as either an absolute value or a percentage with upper and lower limits. The values defining the limits can be symmetric or asymmetric to the target value. These values can be entered either by storing weights on the scale or by key input.

There are four control lines, called data output port lines, which are activated as follows: (see also the diagram at the right):

- lighter
- equal
- heavier
- set

The control range spans 30% to 170% of the target value. You can configure this parameter in the Setup menu (... Application 2: Checkweighing: Activation of port lines:) to select whether the control lines are:

- activated within the checkweighing range
- always on
- activated at stability within the checkweighing range
- activated at stability
- activated at stability within the checkweighing range
-> once

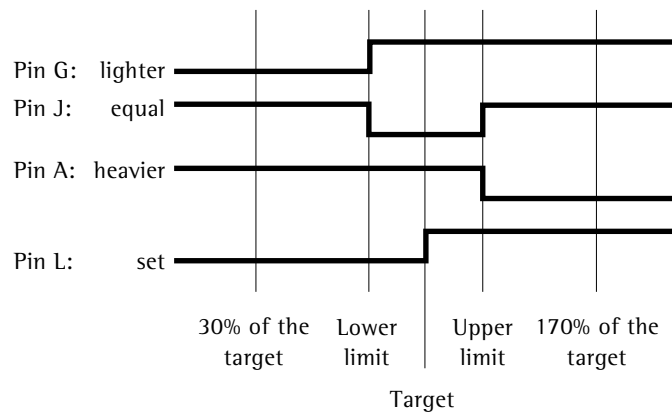
This makes it possible, for example, to connect a simple indicator for the weighing results (e.g., three different colors, one each for the weighing results: too light, O.K., too heavy).

Response of Control Lines During Checkweighing

Configurations:

- always on
- activated at stability


12-pin via Zener barrier



Configurations:

- activated within checkweighing range
- activated at stability within checkweighing range
- activated at stability within checkweighing range

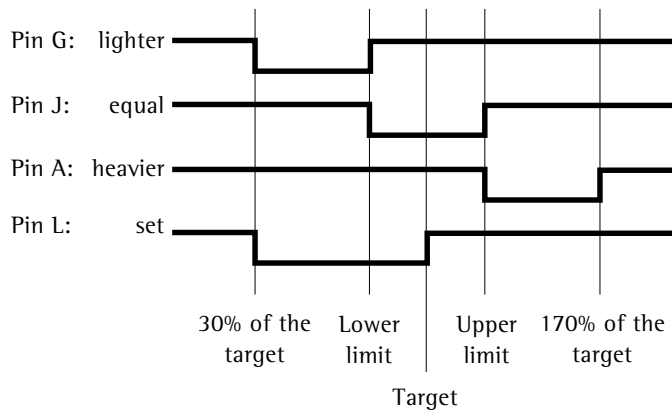12-pin via Zener barrier





For further information about the pin assignment, see "Pin Assignment Charts."


Output port specifications

- When not in use, the voltage level is high: $>2.4\text{ V}/+2\text{ mA}$.
- When activated, the voltage level is low: $<0.4\text{ V}/-2\text{ mA}$.

⚠ The output ports are not protected against short circuits!

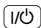Preparation

- Turn on the scale: press 
- > The Sartorius logo is displayed; a self-test is performed
- Select the “Checkweighing” application in the Setup menu: press 
- Select the **Application parameters**: press the \vee soft key 2x, then the \rightarrow soft key
- Select **Application 2 (control functions)**: Press the \vee soft key, then the \rightarrow soft key
- Select **Checkweighing**: press the \wedge or \vee soft key, repeatedly, if necessary
- Confirm **Checkweighing**: press the \rightarrow soft key

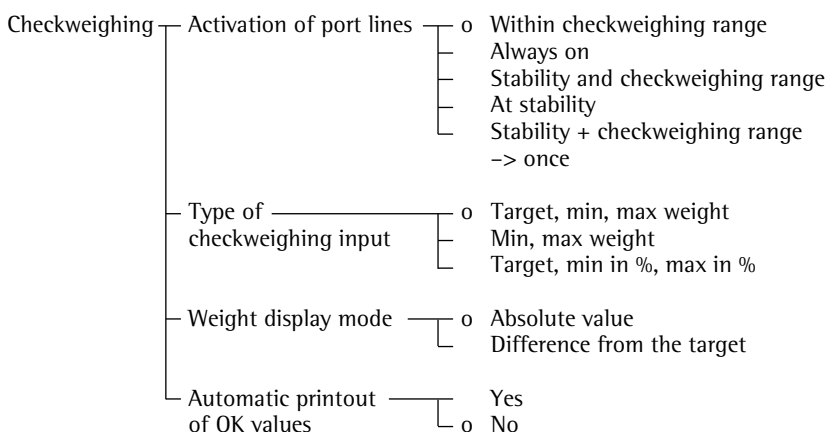

o = factory setting

see also the “Application Parameters (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: press the $\leftarrow \leftarrow$ soft key

Additional Functions

In addition to the functions for:

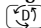
- alphanumeric input, (not during initialization),
- taring (not during alphanumeric input)
- printing,

you can also access the following functions from this application:


Calibration/Adjustment

- Press the **isoTST** soft key
- > See “Calibration/Adjustment” for further instructions


Toggling to Another Application

- Press 
- > See the section on the corresponding application program for further instructions

Setup (Setting Parameters)

- Press 
- > See “Configuration” for further instructions

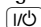Turning Off the Scale

- Press 
- > The scale shuts off





Checkweighing ⚠

Practical Example

Checkweighing samples of 170 g, with an allowable tolerance of -5 g and +10 g. Printout of upper and lower tolerance limits. Weighed values are printed out automatically when stability is reached and the weight value is within the checkweighing range.

Settings (changes in the factory settings required for this example):
Setup: Application parameters: Application 2: Checkweighing: Automatic printout of OK values: On

Step	Press key(s) (or follow instructions)	Display/Output
1. Turn on the scale and configure the settings as indicated above	[I/O]	
2. Delete previous values, if necessary	[CF]	
3. Prepare a container for the samples	Place empty container on the scale	<div>Max6200 g d= 0.01g 100% 0% + 210.00 g CHECKWEIGH: Initialize isoTST Param. Start</div>
4. Tare the scale	[Tare]	<div>Max6200 g d= 0.01g 100% 0% 0.00 g CHECKWEIGH: Initialize isoTST Param. Start</div>
5. Enter initialization values	Param. soft key	<div>CHECKWEIGH: 0.00 g Target: Setp= + 0.00 g Minimum: Min = + 0.00 g Maximum: Max = + 0.00 g << v J</div>
6. Accept target value (in this example: 170 g)	Place ideal sample in container	<div>CHECKWEIGH: 170.00 g Target: Setp= + 0.00 g Minimum: Min = + 0.00 g Maximum: Max = + 0.00 g << v J</div>
7. Save target and unload the scale	J soft key Remove ideal sample from the scale	<div>CHECKWEIGH: 0.00 g Target: Setp= + 170.00 g Minimum: Min = + 0.00 g Maximum: Max = + 0.00 g << ^ v J</div>
8. Enter value for lower limit (170 g - 5 g) and save	[1] [6] [5] J soft key	<div>CHECKWEIGH: 0.00 g Target: Setp= + 170.00 g Minimum: Min = + 165.00 g Maximum: Max = + 0.00 g << ^ J</div>

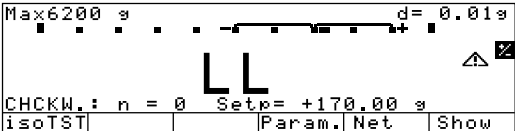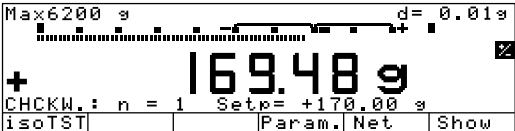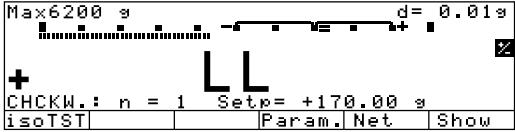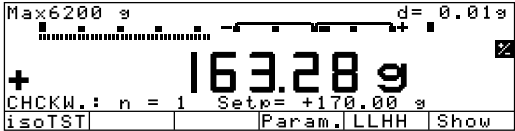Step	Press key(s) (or follow instructions)	Display/Output									
9. Enter value for upper limit (170 g + 10 g) and save	1 8 0	 <table border="1"> <tr> <td>Setp</td> <td>+</td> <td>170.00 g</td> </tr> <tr> <td>Min</td> <td>+</td> <td>165.00 g</td> </tr> <tr> <td>Max</td> <td>+</td> <td>180.00 g</td> </tr> </table>	Setp	+	170.00 g	Min	+	165.00 g	Max	+	180.00 g
Setp	+	170.00 g									
Min	+	165.00 g									
Max	+	180.00 g									
10. Weigh sample (in this case: 169.48 g)	Place sample in container	 <table border="1"> <tr> <td>N</td> <td>+</td> <td>169.48 g</td> </tr> </table>	N	+	169.48 g						
N	+	169.48 g									
If the weight value is too low:											
11. In this case, switch to net value (for ex., a weight of 163.28 g)	Net soft key										
12. Weigh next sample (if any)	Place sample in container										

Time-Controlled Functions

Purpose

With this application program, you can configure the scale to perform certain functions (such as automatic printout of values, store value in totalization memory) at a given time or after a set interval.

You can use this application in combination with any program chosen from Application 1 (such as counting, weighing in percent) and one from Application 3 (such as totalizing, formulation).

Features

- Time-controlled activation of scale functions:
 - one time only, at a given time (**Setting=** is displayed in the text line)
 - repeatedly, at given intervals (**Interval=** is displayed in the text line before the function is started, and **Repeat =** is displayed after the function has been started)
- Functions that can be time-controlled include:
 - Acoustic signal (beep)
 - Lock in readout
 - Automatic printout of values
 - Store values for totalizing, formulation or statistics
- Print time in addition to weight value
- Store value depending on the stability parameter
- Tare the scale after printout of weight values
- Press the corresponding soft key to cancel time-controlled functions

Factory Settings

Function after time interval:
Automatic printout of values

Automatic function restart: **On**

Storage mode:
Without stability

Print then tare: **On**

Soft Key Functions

Stop	Stop the application
Quit	Confirm performed function (e.g., »Lock in readout« or »Beep«)
Interval	Store input interval for time-controlled functions
Set.	Store input time for one-time performance of function
Start	Start timer function


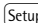
Printout for Time-Controlled Functions


If the "Automatic printout of values" parameter is set, the time and weight (or other value) are printed out.

Time: 10:15:00
N + 150.00 g

Time: Time that the values were stored
N: Net weight

Preparation

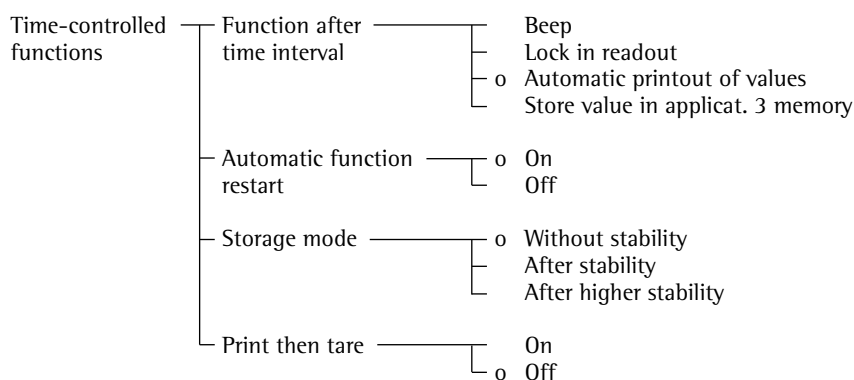- Turn on the scale: press 
 - > The Sartorius logo is displayed
- Select the "Time-controlled functions" application in the Setup menu: press 
- Select the **Application parameters**: press the ∇ key 2x, then the \triangleright soft key
- Select **Application 2 (control functions)**: press the ∇ soft key, then the \triangleright soft key
- Select **Time-controlled functions**: press the \wedge or ∇ soft key
- Confirm **Time-controlled functions**: press the \triangleright soft key



o = factory setting

see also the "Application Parameters (Overview)" in the chapter entitled "Configuration"

- Save settings and exit the Setup menu: press the \ll soft key

Print Net Values without Printout of Time

Select the Setup menu:

Setup: Printout: Application-defined output: Auto print upon initialization: Off

Practical Example

Document the evaporated amount of a sample with defined surface, temperature and air pressure at preset intervals of 1 minute, 30 seconds.

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Application 2: Time-controlled functions

Setup: Balance/scale functions: Taring: Without stability

Setup: Printout: Application-defined output: Stability parameter: Without stability

Step	Press key(s) (or follow instructions)	Display/Output
1. Turn on the scale and configure the settings as indicated above		
2. Delete stored values, if necessary		
3. Place container with sample on the scale and tare		
4. Enter time interval: 1 minute, 30 seconds		
5. Store time interval	Interval soft key	
6. Begin documentation (Time remaining until the next printout is displayed in the text line)	Start soft key	
Printout of evaporated amount every 1/2 minutes		<div> Time: 15:19:50 N - 0.37 g Time: 15:21:20 N - 0.33 g Time: 15:22:50 N - 0.30 g Time: 15:24:20 N - 0.40 g </div>
7. Stop the documentation procedure	Stop soft key	

Totalizing Σ

Purpose

This application program runs as a cumulative memory function.

You can use this application in combination with any program chosen from Application 1 (such as counting, weighing in percent) and one from Application 2 (checkweighing, time-controlled functions).

Features

- Totalization of weight values and calculated values
- Optional configuration in the Setup menu for simultaneous storage of net and calculated values
- Optional configuration in the Setup menu for loading weight values and calculated values either from Application 1 (e.g., counting, weighing in percent) or from Application 2 (checkweighing, time-controlled functions)
- Totalization memory for up to 65,535 values
- Simultaneous display in the text line of transaction counter and, e.g., the current total
- Optional configuration in the Setup menu for having the scale tare automatically after a value is stored in the totalization memory, if no preset tare has been entered
- Manual input of the number of individual weighing operations and confirmation using the **nDef** soft key (target no. of operation nDef). Result printed and memory cleared after printout of nDef.
- Optional configuration in the Setup menu to add the current weight, with display accuracy, to the current total by pressing the **M+** soft key and generate a printout of the result
- Optional configuration in the Setup menu for stability-dependent storage of the measured value: **Balance/scale functions, Stability range**
- Optional automatic storage of measured values
Storage of measured value is indicated by $\rightarrow \pm$; $\pm \pm$ indicates that you can place a load on the scale.
- Minimum load threshold for automatic storage
- Press the **M-** soft key to delete the last value added to the totalizing memory. The transaction counter value is reduced by one and a printout is generated.
- Press the **MR** soft key for information about number of transactions and the current total. In the Setup menu, you can define whether the information is displayed and printed, or only printed, and whether the information comprises an intermediate or final evaluation
- In the Info window you can choose which value is displayed in the text line during weighing
- Printout of the end result independent of which program is configured for Application 1 or Application 2. Configure the Setup menu to define which values are included on the printout (printout of individual components)
- Press the key identified by **MR** (soft key label) for a printout of an intermediate evaluation after each addition or a final evaluation
- If you end the totalization process by pressing **[CF]** without having first pressed the **MR** soft key for a printout, a final evaluation is printed when you press **[CF]**
- Optional configuration in the Setup menu to clear the totalizing memory and reset the transaction counter by pressing **[CF]** or after an evaluation is printed out
- Totalization data and transaction counter data are stored in non-volatile memory
- Continue totalization after turning the scale off and back on

Factory Settings of Parameters

Automatic storage: **Off**

Minimum load for automatic storage: **10 digits**

Source of data for auto storage: **Application 1**

Evaluated values: **Net**

Evaluation mode, MR key function: **Intermediate evaluation, print**

M+/M- function, then tare: **Off**

Printout of individual components: **Yes**

Stability range: **2 digits**

Application-defined output: **Print then tare: Off**

Soft Key Functions

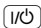
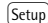
- M+** Add weight values or application values to the total in the totalization memory. The component or transaction counter value increases by one each time you press this key.
- M-** Delete the last value added to memory. The transaction counter value decreases by 1. You cannot delete previous values by repeatedly pressing this key.
- MR** Print or display an intermediate or final evaluation
- nDef** Store the input number of components


Printout for Totalizing

The transaction or component counter is printed in front of each measured value (weight). When an intermediate or final evaluation is printed out, all results up to this point are included.

n 5
Total +151.67 g

Preparation

- Turn on the scale: press 
- > The Sartorius logo is displayed; a self-test is performed
- Select the “Totalizing” application program in the Setup menu: press 
- Select the **Application parameters**: press the ∇ soft key 2x, then the \rightarrow soft key
- Select **Application 3 (data records)**: press the ∇ soft key 2x, then the \rightarrow soft key once
- Select **Totalizing**: press the \wedge or ∇ soft key
- Confirm **Totalizing**: press the \rightarrow soft key

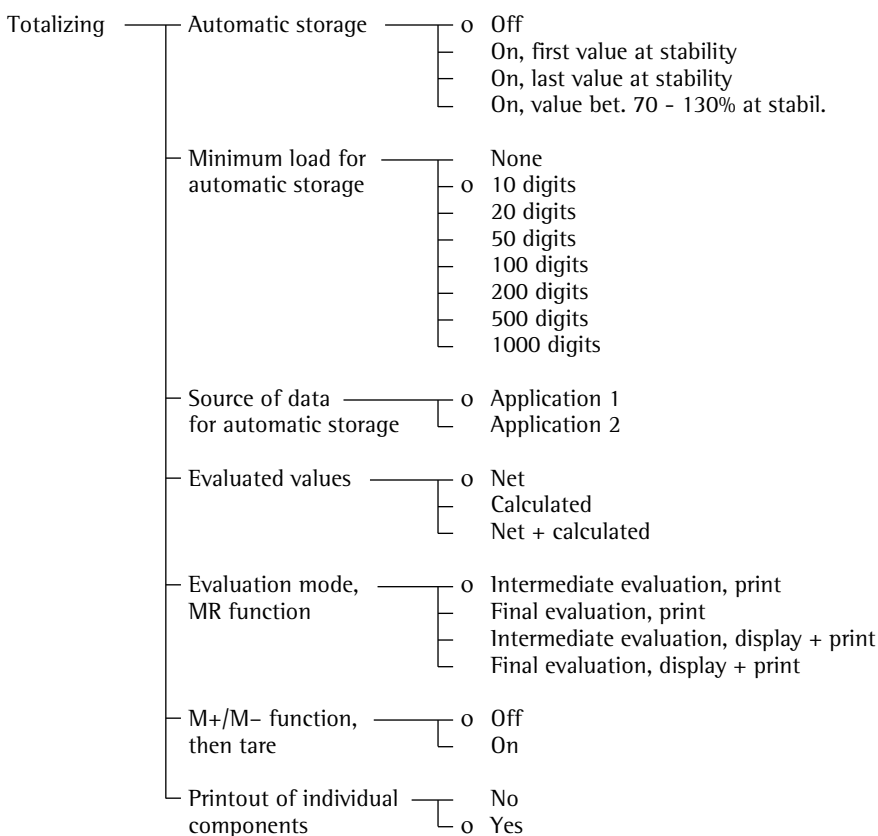

o = factory setting

see also the “Application Parameters (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: press the $\Leftarrow \Leftarrow$ soft key

Additional Functions

In addition to the functions for:

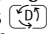
- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:


Calibration/Adjustment

- Press the **isoTST** soft key
- > See “Calibration/Adjustment” for further instructions


Toggleing to Another Application

- Press 
- > See the section on the corresponding application program for further instructions

Setup (Setting Parameters)

- Press 
- > See “Configuration” for further instructions

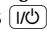Turning Off the Scale

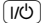
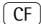

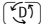
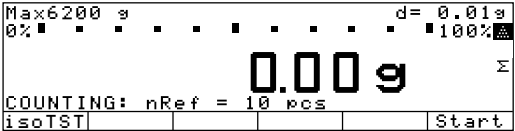
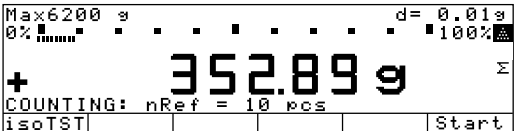
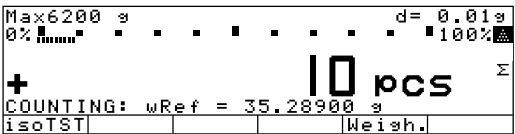
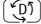
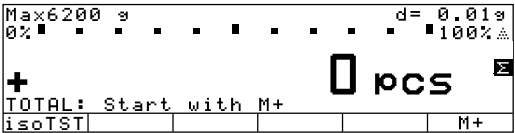
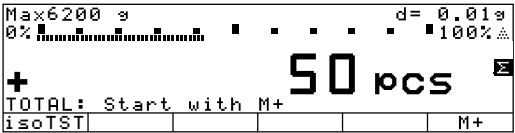
- Press 
- > The scale shuts off



Totalizing Σ

Practical Example

Totalize counted pieces

Settings (changes in the factory settings required for this example):
Setup: Application parameters: Application 1: Counting
Setup: Application parameters: Application 2: Off
Setup: Application parameters: Application 3: Totalizing: Evaluated values: Net + calculated
Setup: Application parameters: Application 3: Totalizing: Evaluation mode, MR function: Final evaluation, display + print

Step	Press key(s) (or follow instructions)	Display/Output
1. Turn on the scale and configure the settings as indicated above		
2. Delete old totalization data, if necessary		
3. Tare the scale		
4. Toggle to Application 1: Counting		
5. Place the displayed number of parts on the scale (here: 10 pcs)	Place parts to be counted on the scale	
6. Initialize the Counting application	Start soft key	 <div><div>nRef10 pcs</div><div>wRef35.28900 g</div></div>
7. Remove the reference sample quantity and toggle to Totalizing	Unload the scale 	
8. Place a number of parts on the scale (here: 50 pcs)	Place parts on the scale	

Step	Press key(s) (or follow instructions)	Display/Output
9. Store piece count	M+ soft key	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 50 pcs TOTAL: n=1 Σ tot= + 100 pcs isoTST MR M- M+ </div> </div> <div> ----- 16.01.1997 11:06 n + 1 N + 1764.45 g Qnt + 50 pcs ----- </div>
10. Unload the scale	Remove parts from the scale	
11. Place another load of parts on the scale (e.g., 60 pcs)	Place parts on the scale	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 60 pcs TOTAL: n=1 Σ tot= + 110 pcs isoTST MR M- M+ </div> </div>
12. Add piece count to stored total	M+ soft key	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 60 pcs TOTAL: n=2 Σ tot= + 170 pcs isoTST MR M- M+ </div> </div> <div> n 2 N + 2117.34 g Qnt + 60 pcs </div>
13. Repeat steps 10 and 11 as required		
14. Display final evaluation ("Info" window) (here: 5 weighing operations; total weight: 8751.67 g; total quantity: 248) The  indicates which value is displayed in the text line; you can change this selection	MR soft key	<div> TOTAL: Net: n = 5 oNet: Σ = + 8751.67 g Calculated: n = 5 Calculated: Σ = + 248 pcs << v ↓ </div> <div> ----- n 5 Total + 8751.67 g Total + 248 pcs 16.01.1997 11:16 ----- </div>
15. Print final evaluation		

Formulation ↕

Purpose

With this application program you can add weight values and calculated values as components of a formula to a totalizing memory.

You can use this application in combination with any program chosen from Application 1* (such as counting, weighing in percent) and one from Application 2 (checkweighing, time-controlled functions) as well as with the extra functions.

*= not with recalculation or 2nd tare memory; cannot be activated during formulation

Available Features

- Totalization of weight values and calculated values
- Weigh in different components toward zero to a total amount defined by pressing the **nom** soft key and entering the value through the numeric keys
- Simultaneous storage of net and calculated values, if available
- Optional configuration in the Setup menu for loading weight values and calculated values either from Application 1 (e.g., counting, weighing in percent) or from Application 2 (checkweighing, time-controlled functions)
- Totalizing memory for up to 65,535 values
- Transaction counter and current total displayed in the text line
- Scale tared after a value has been stored
- Manual input of the number of individual weighing operations and confirmation using the **nDef** soft key (target no. of operations nDef). Result printed and memory cleared after printout of nDef.
- Optional configuration in the Setup menu to add the current weight, with display accuracy, to the current total by pressing the **M+** soft key and generating a printout of the result

- Optional configuration in the Setup menu for stability-dependent storage of the measured value: Balance/scale functions, Stability range
- Optional automatic storage of measured values
Storage of measured value is indicated by **→←**; **↕↕** indicates that you can place a load on the scale.
- Minimum load threshold for automatic storage
- Press the **M-** soft key to delete the last value added to the totalizing memory. The transaction counter value is reduced by one and a printout is generated.
- Printout of an evaluation of results, depending on the Application 1 or Application 2 parameters. Configure the Setup menu to define the information included on this printout.
- Press the key identified by **MR** (soft key) for a printout of an intermediate evaluation after each addition or a final evaluation
- A final evaluation is printed when the formulation routine is ended by pressing **[CF]**, if no final evaluation was generated by pressing **MR**
- Optional configuration in the Setup menu to clear the totalizing memory and reset the transaction counter by pressing **[CF]** or after an evaluation is printed out
- Totalization data and transaction counter data are stored in the non-volatile memory
- Continue formulation after turning the scale off and back on

Factory Settings of the Parameters

Automatic storage: **Off**

Minimum load for automatic storage:
10 digits

Source of data for auto storage:
Application 1

Evaluated values: **Net**

Evaluation mode, **MR** key function:
Intermediate evaluation, print

Printout of individual components: **Yes**

Stability range: **2 digits**
Printout: Application-defined output:
Print on request then tare: **Off**

Soft Key Functions

M+	Add weight values or application values to the total in the totalizing memory. The component or transaction counter value increases by one each time you press this key.
M-	Delete the last value added to memory. The transaction counter value decreases by 1. You cannot delete previous values by repeatedly pressing this key.
MR	Print or display an intermediate or final evaluation
nDef	Store the input number of components
nom	Press to enter target component weight using the numeric keys



Printout of Formulation Report


When an intermediate or final evaluation is printed out, all results up to this point are included.

Comp2 + 42.38 g
Tot.cp+184.89 g

Comp2: Weight of the 2nd component
Tot.cp: Total of all components

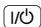Preparation

- Turn on the scale: press 
- > The Sartorius logo is displayed; a self-test is performed
- Select the "Formulation" application in the Setup menu: press 
- Select **Application parameters**: press the ∇ soft key 2x, then the \rightarrow soft key once
- Select **Application 3 (data records)**: press the ∇ soft key 2x, then the \rightarrow soft key once
- Select **Formulation**: press the \wedge or ∇ soft key
- Confirm **Formulation**: press the \rightarrow soft key

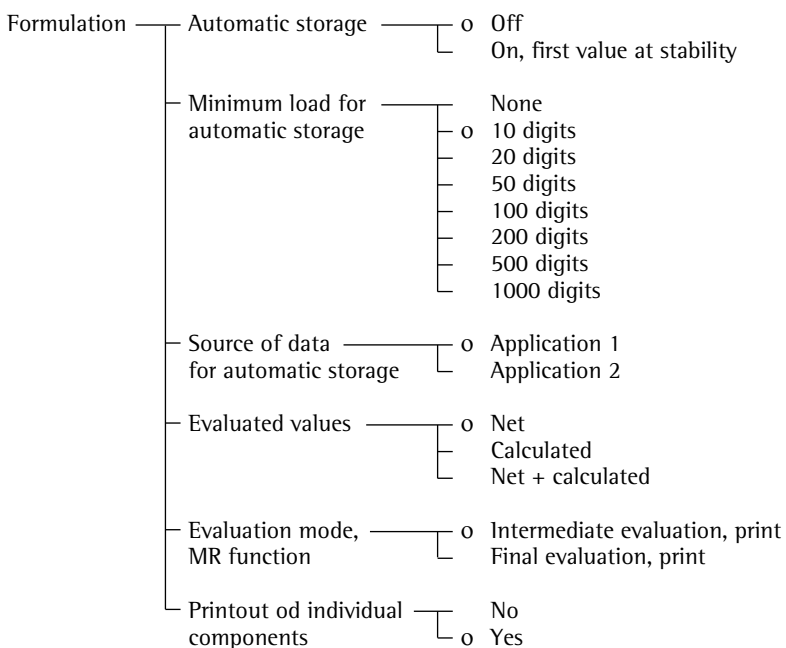

o = factory setting

see also the "Application Parameters (Overview)" in the chapter entitled "Configuration"

- Save settings and exit the Setup menu: press the $\leftarrow \leftarrow$ soft key

Additional Functions

In addition to the functions for:

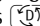
- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:


Calibration/Adjustment

- Press the **isoTST** soft key
- > See "Calibration/Adjustment" for further instructions


Toggleing to Another Application

- Press 
- > See the section on the corresponding application program for further instructions

Setup (Setting Parameters)

- Press 
- > See "Configuration" for further instructions

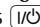Turning Off the Scale

- Press 
- > The scale shuts off

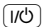
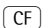

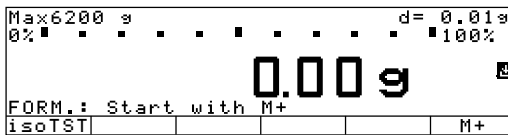






Practical Example


Weighing in Components

Setup: Application parameters: Application 3: Formulation: Automatic storage: On, first value at stability

Setup: Application parameters: Application 3: Formulation: Evaluation mode, MR function: Final evaluation,

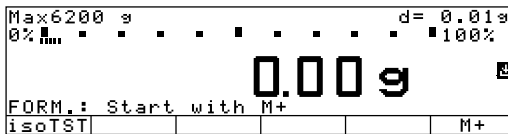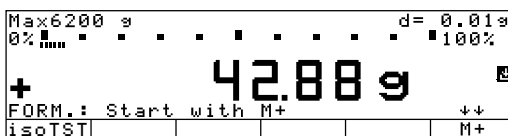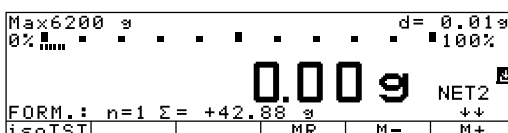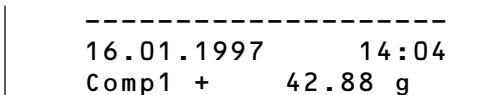Setup: Application parameters: Application 3: Formulation: Evaluation mode, MR function: Final evaluation, print

Step	Press key(s) (or follow instructions)	Display/Output
1. Turn on the scale and configure the settings as indicated above		
2. Delete old formulation data		
3. Tare the scale		
4. Place the empty container on the scale (here: 180.59 g)	Place load on the scale	
5. Tare the scale		
6. Weigh in the first component (here: 42.88 g)	Place components in container	
7. Store components in the formulation memory Scale is tared automatically	M+ soft key	
Component are printed out automatically		

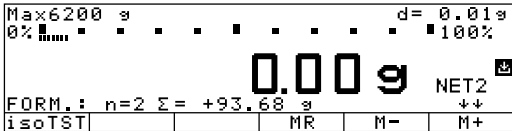Step	Press key(s) (or follow instructions)	Display/Output
<p>8. Weigh in next component (here: 50.80 g)</p> <p>Components are stored in the Totalization memory at stability and printed out</p> <p>Scale is tared automatically</p>	Place components in container	<div>Comp2 + 50.80 g</div> 
9. Repeat step 7 as required		
10. Print final evaluation (here: with total weight of all components: 212.43 g)	MR soft key	<div> <pre> ----- n 2 Tot.cp+ 212.43 g 16.01.1997 14:10 ----- </pre> </div>
11. Delete old formulation data, if necessary	CF	

Purpose

With this application program, you can have weight values and calculated values totalized and statistically evaluated.

The values determined for the evaluation are:

- average (mean value)
- standard deviation
- variation coefficient
- sum of all values
- lowest value (minimum)
- highest value (maximum)
- difference between the minimum and the maximum

You can use the statistics application in combination with any program chosen from Application 1 (such as counting, weighing in percent) and one from Application 2 (checkweighing, time-controlled functions) as well as with the extra functions.

Features

- Storage of weight values and calculated values
- Simultaneous storage of net and calculated values
- Optional configuration in the Setup menu for loading weight values and calculated values either from Application 1 (such as counting, weighing in percent) or from Application 2 (checkweighing, time-controlled functions)
- Totalizing memory for up to 65,535 values
- Simultaneous display in the text line of the transaction counter and, e.g., the current total
- Optional configuration in the Setup menu for having the scale tare automatically after a value has been stored in the totalizing memory
- Manual input of the number of individual weighing operations and confirmation using the **nDef** soft key (target no. of operations nDef). Result printed and memory cleared after printout of nDef.

- Optional configuration in the Setup menu to add the current weight, with display accuracy, to the current total by pressing the **M+** soft key and generate a printout of the result
- Optional configuration in the Setup menu for stability-dependent storage of the measured value: Balance/scale functions, Stability range
- Optional automatic storage of measured values
Storage of measured value is indicated by **→←**; **→→** indicates that you can place a load on the scale.
- Minimum load threshold for automatic storage
- Press the **M–** soft key to delete the last value added to the totalizing memory. The transaction counter value is reduced by one and a printout is generated.
- Press the **MR** soft key for information about number of transactions and the current total. By configuring the Setup menu, you can define whether the information is displayed and printed, or only printed, and whether the information comprises an intermediate or final evaluation
- In the Info window you can use the **↵**, **↩** (□) soft keys to choose which value will be displayed in the text line during weighing
- Printout of the final result depending on the Application 1 or Application 2 parameters. Configure the Setup menu to define which values are included on the printout (printout of individual components)
- Press **MR** for a printout of an intermediate evaluation after each addition or a final evaluation
- A final evaluation is printed when the statistics routine is ended by pressing **(CF)**, if no final evaluation was generated by pressing **MR**
- Optional configuration in the Setup menu to clear the totalizing memory and reset the transaction counter by pressing **(CF)** or after an evaluation is printed out
- Totalization data and transaction counter data is stored in the non-volatile memory
- Continue totalization after turning the scale off and back on

Factory Settings

Automatic storage: **Off**

Minimum load for automatic storage:
10 digits

Source of data for auto storage:
Application 1

Evaluated values: **Net**

Evaluation mode, MR key function:
Intermediate evaluation, print

M+/M– function, then tare: **Off**

Printout of individual components: **Yes**

Stability range: **2 digits**

Application-defined output:
Print on request then tare: **Off**

Soft Key Functions

M+ Add weight values or application values to the total in the totalizing memory. The component or transaction counter value increases by one each time you press this key.

M– Delete the last value added to memory. The transaction counter value decreases by 1. You cannot delete previous values by repeatedly pressing this key.

MR Print or display an intermediate or final evaluation

nDef Store the input number of components

Printout of Statistics


The transaction or component counter is printed in front of each measured value (weight). When an intermediate or final evaluation is printed out, all results up to this point are included.

n	5
Total	+ 151.67 g
Avg.	+ 33.0 pcs
s	+ 3.2 pcs
srel	+ 9.70 %
Total	+ 165 pcs
Min	+ 29 pcs
Max	+ 37 pcs
Diff	+ 8 pcs

n:	Transaction counter
Total:	Sum of all values
Mean:	Average
s:	Standard deviation
srel:	Variation coefficient
Total:	Sum of all values
Min:	Minimum
Max:	Maximum
Diff:	Difference between minimum and maximum

Preparation

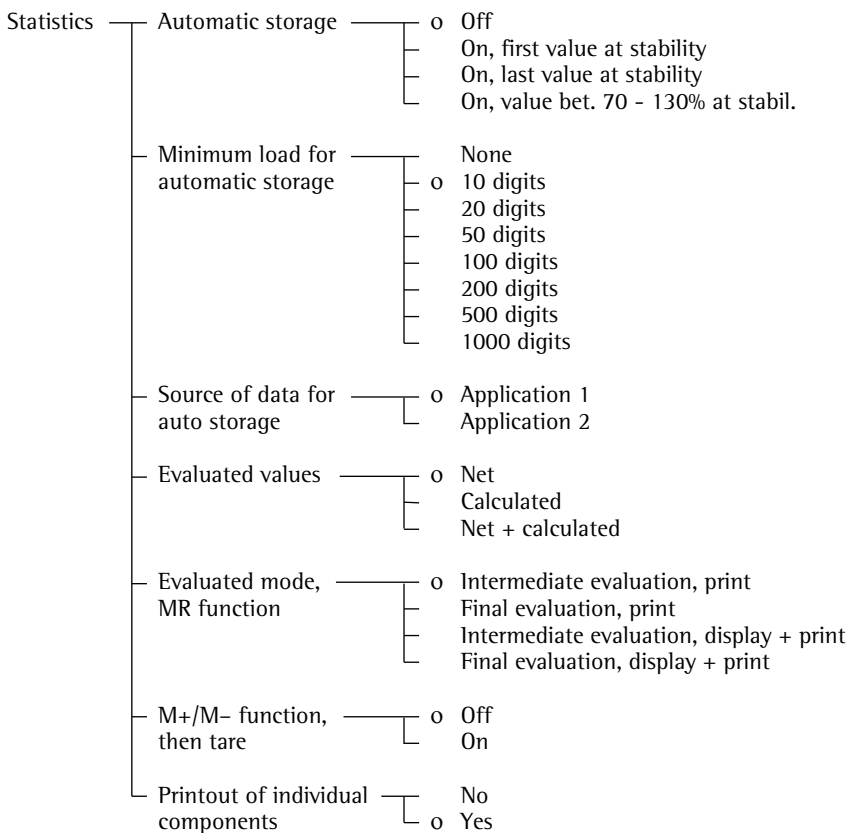- Turn on the scale: press **[I/O]**
- > The Sartorius logo is displayed
- Select the “Statistics” application in the Setup menu: press **[Setup]**
- Select **Application parameters**: press the **▼** soft key 2x, then the **→** soft key once
- Select **Application 3 (data records)**: press the **▼** soft key 2x, then the **→** soft key once
- Select **Statistics**: press the **▲** or the **▼** soft key
- Select **Statistics**: press the **→** soft key



o = factory setting

see also the “Application Parameters (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: press the **◀◀** soft key

Additional Functions

In addition to the functions for:

- alphanumeric input,
- taring (not during alphanumeric input),
- printing,

you can also access the following functions from this application:

Calibration/Adjustment

- Press the **ISO TEST** soft key
- > See “Calibration/Adjustment” for further instructions

Toggleing to Another Application

- Press **[↔]**
- > See the section on the corresponding application program for further instructions

Setup (Setting Parameters)

- Press **[Setup]**
- > See “Configuration” for further instructions

Turning Off the Scale

- Press **[I/O]**
- > The scale shuts off

Statistics

Practical Example




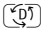
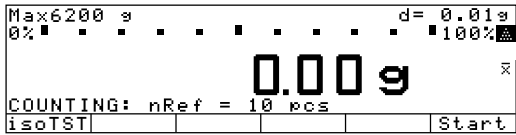
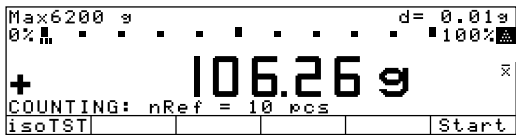
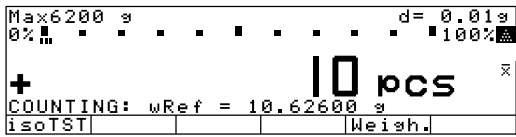
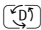
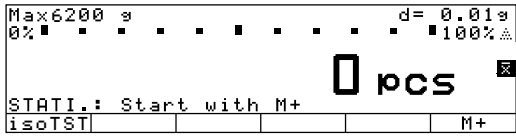
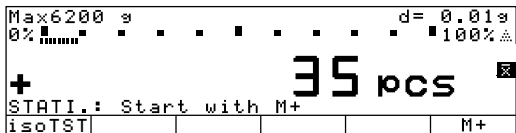
Totalize counted pieces and print out statistics



Settings (changes in the factory settings required for this example):

Setup: Application parameters: Application 1: Counting: Average piece weight updating: Manual

Setup: Application parameters: Application 3: Statistics: Evaluated values: Calculated

Setup: Application parameters: Application 3: Statistics: Evaluation mode, MR function: Final evaluation, display + print

Step	Press key(s) (or follow instructions)	Display/Output
1. Turn on the scale and configure the settings as indicated above		
2. Delete old statistics data, if necessary		
3. Tare the scale		
4. Toggle to Application 1: Counting		
5. Place the displayed number of parts on the scale (here: 10 pcs)	Place parts to be counted on the scale	
6. Initialize the Counting application	Start soft key	 <div> nRef 10 pcs wRef 10.62600 g </div>
7. Remove the reference sample quantity and toggle to Statistics	Unload the scale 	
8. Place a number of parts on the scale (here: 35 pcs)	Place parts on the scale	

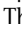Step	Press key(s) (or follow instructions)	Display/Output
9. Store piece count	M+ soft key	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 35 pcs STATI.:n=1 Qnt + 35 pcs isotST MR M- M+ </div> </div> <div> ----- 16.01.1997 11:06 n + 1 Qnt + 35 pcs </div>
10. Unload the scale	Remove parts from the scale	
11. Place another load of parts on the scale (e.g., 29 pcs)	Place parts on the scale	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 29 pcs STATI.:n=1 Qnt + 35 pcs isotST MR M- M+ </div> </div>
12. Add piece count to stored total	M+ soft key	<div> <div> Max6200 g 0% d= 0.01g 100% </div> <div> + 29 pcs STATI.:n=2 Qnt + 29 pcs isotST MR M- M+ </div> </div> <div> n + 2 Qnt + 29 pcs </div>
13. Repeat steps 11 and 12 as required		
14. Display final evaluation (»Info« window) (here: 5 weighing operations; total quantity: 165) The  indicates which value is displayed in the text line; you can change this selection	MR soft key	<div> STATI.: Calculated: n 5 Calculated: X = + 33.0 pcs Calculated: s = + 3.2 pcs Calculated: srel = + 9.70 % oCalculated: Σ = + 165 pcs << v ↓ </div>
15. Print final evaluation		<div> ----- n 5 Avg. + 33.0 pcs s + 3.2 pcs srel + 9.70 % Total + 165 pcs Min + 29 pcs Max + 37 pcs Diff + 8 pcs 16.01.1997 11:16 ----- </div>

Extra Functions (in the Application Menu)

Second Tare Memory (Preset Tare)

Purpose

With this function, you can store the weight currently on the scale as a tare weight, or use the numeric keys to enter a number for a preset tare weight.

You can use this function in combination with a program from Application 1 (such as counting, weighing in percent), one from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics) as well as with the extra functions.

Features

- Store a weight on the scale in the second tare memory (without numeric input)
- Store a numeric value in the second tare memory (input using the numeric keys)
- Identify a net value as **NET** when there is a value stored in the second tare memory
- You can assign this function to the fourth or fifth soft key (from the right); i.e., F4 or F5.
The soft key designation for this function is: **PT1/T1**
- Optional configuration in the Setup menu for storing the current weight readout as the container tare weight. Any load subsequently placed on the scale that is more than 70% of the tare weight is automatically recognized as a container and the scale is tared automatically.
- Automatic printout when a value is stored or input (see “Configuration”)
- Press **[CF]** to delete the (preset) tare value

Factory Settings

Container tare weight: **No**

Automatic printout: **Off**

Soft Key Functions

PT1/T1 Store weight as tare value

PT1 Store input value

Printout of the Data in the 2nd Tare Memory

The printout shows either:

- Net value **N1**,
- Tare weight **T1**, or
- Manually entered tare value **PT1**

N1	63.48 g
T1	138.73 g
PT1	150.00 g


N1: Net weight (value) when a weight is stored in the tare memory

T1: Tare weight

PT1: Preset tare value entered using the numeric keys

Preparation

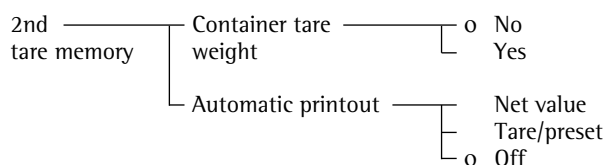- Turn on the scale: press **[ON]**
> The Sartorius logo is displayed; a self-test is performed
- Select Extra function (F4) or Extra function (F5) in the Setup menu: press **[Setup]**
- Select **Application parameters**: press the **▼** soft key 2x, then the **➤** soft key once
- Select **Extra function (F4) or Extra function (F5)**: press the **▼** soft key 3x (or 4x), then press the **➤** soft key once
- Select **2nd tare memory**
- Confirm **2nd tare memory**



O = factory setting

see also the “Application Parameters (Overview)” in the chapter entitled “Configuration”

- Save settings and exit the Setup menu: press the **◀◀** soft key

Second Tare Memory in Legal Metrology

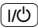
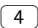
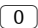
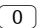
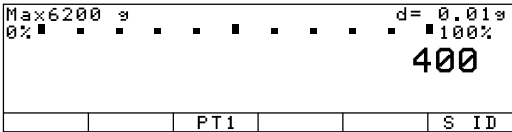
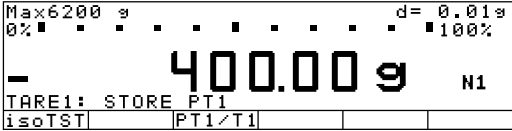
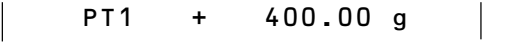

- Press the **◀ i ▶ PT1** soft key to enter information about the tare value using the number keys.
- The PT1 tare value is printed out with the net value.

Practical Example

Determine the Contents of Bottles: Bottle weight = 400 g.

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Extra function (F4): 2nd tare memory: Automatic printout: Tare/preset tare

Step	Press key(s) (or follow instructions)	Display/Output
1. If necessary: turn on the scale and enter the settings given above		
2. Enter bottle weight (example 400 g)	  	
3. Store tare value	PT 1 soft key	 
4. Determine net weight of bottles (in this case: net contents = 650 g)	Place filled bottles on the scale	

Individual Identification Codes (ID)

Purpose

With this function, you can assign IDs to values for documentation and printouts.

You can use this function in combination with any program from Application 1 (such as counting, weighing in percent), one from Application 2 (checkweighing, time-controlled functions) and one from Application 3 (totalizing, formulation, statistics) as well as with the other extra functions.

Features

- Store up to 4 IDs; these can be stored, changed or deleted individually.
- Each ID consists of a name and a value; both can be defined by the user.
- ID designations are configured as follows: Setup: Printout: Identification codes
- Each ID code can have up to 20 characters; when you enter the value later, however, no more than 15 characters of this ID are displayed.
- The ID values are entered while the application program is active; press the **ID** soft key to toggle to the ID input mode.
- Each ID value can have up to 20 characters.
- Access 1 of the 4 IDs directly using the numeric keys. The other three can only be accessed by pressing the **ID** soft key to toggle to the ID input mode.
- You can assign this function to the fourth or fifth soft key (from the right); i.e., F4 or F5.
- You can configure when the ID will be included on the printout (see “Preparation” on the next page).
- You can configure the position of IDs on the individual or total printout.
- The ID code is printed flush left; the value flush right. If the name and value together are too long for one line, the data is printed on two lines.
- Optional configuration in the Setup menu to delete a single character when entering an identification code by pressing **CF**. Setup: Device parameters: Keys: CF function for input: Delete last character
- Press the **Delete** soft key to delete an ID

Factory Settings of the ID Names

ID1: **ID 1**
ID2: **ID 2**
ID3: **ID 3**
ID4: **ID 4**

Factory Settings for ID Codes

No values set

Factory Settings

Printout:

Each time the print key is pressed

Soft Key Functions

ID Toggle to “Identification codes” menu

Delete Delete input of selected ID







Printout of ID Codes

Up to 4 (stored) identification codes are printed out.

ID1	Lot no. 1234
ID2	Daimler/Chrysler
ID3	Screws M4x6
ID4	Jack Smith

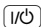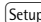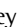ID1: Identification 1 (ID 1)
ID2: Identification 2 (ID 2)
ID3: Identification 3 (ID 3)
ID4: Identification 4 (ID 4)

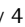
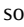

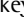

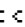
Preparation

- Turn on the scale: press 
> The Sartorius logo is displayed
- Select Extra function (F4) or Extra function (F5) in the Setup menu: press 
- Select **Application parameters**: press the  soft key 2x,
then the  soft key once
- Select **Extra function(F4)** or **Extra function(F5)**:
press the  soft key 3x (or 4x), then the  soft key once
- Select **Identification codes**
- Confirm **Identification codes**

Identification	—	Printout	—		Automatic, if configured
			—		Once after pressing print, if configured
			o		Each time the print key is pressed
			—		Once for M+ function (app. 3 memory)

o = factory setting

see also the “Application Parameters (Overview)” in the chapter entitled “Configuration”

- Save settings for the printout: press the  soft key 4x
- Enter ID name: Select “Printout”: press the  soft key, then the  soft key
- Select “Identification #”: press the  soft key 5x, then the  soft key once
- Select **ID 1**
- Enter name for **ID 1** and confirm: use the numeric keys for numbers and/or the soft keys to enter letters
- Enter names for **ID 2**, **ID 3** and **ID 4**, if desired
- Save settings and exit the Setup menu: press the  soft key

Example

See next page

Practical Example

Include company address and sample lot number on the printout. Each ID line begins with the name.
Print this ID for each net value.

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Extra function (F4): Identification codes

Setup: Input: ID1: Company

Setup: Input: ID2: Location

Setup: Input: ID3: Street

Setup: Input: ID4: Lot

Step	Press key(s) (or follow instructions)	Display/Output
1. If necessary, turn on the scale		
2. Select "Extra Function (F4)" in the Setup menu	 ↓ soft key 2x, then → soft key once ↓ soft key 3x, then → soft key once	<div> <div>SETUP</div> <div>APPLICATION EXT.FCT.F4</div> <div>o0ff</div> <div>2nd tare memory</div> <div>Identification codes</div> <div>Man. store in app.3 memory (M+)</div> <div>Product data memory</div> <div><< < > >></div> </div>
3. Select "Identification codes"	↓ or ↑ soft key; repeatedly, if necessary	<div> <div>SETUP</div> <div>APPLICATION EXT.FCT.F4</div> <div>o0ff</div> <div>2nd tare memory</div> <div>Identification codes</div> <div>Man. store in app.3 memory (M+)</div> <div>Product data memory</div> <div><< < ^ v ></div> </div>
4. Confirm "Identification codes" and exit this menu item	→ soft key; then < soft key 3 times	<div> <div>APPLICATION EXT.FCT.F4 IDENTIFIER</div> <div>Printout</div> <div><< < ></div> </div>
5. Select ID1 (Printout: Identifier)	↓ or → soft key ↓ soft key 5x, then → soft key, then ↓ soft key	<div> <div>SETUP PRINTOUT IDENTIFIER</div> <div>Lot (L ID):</div> <div>ID1: ID1</div> <div>ID2: ID2</div> <div>ID3: ID3</div> <div>ID4: ID4</div> <div><< < ^ v ></div> </div>
6. Enter name for ID 1 (in this case: COMPANY) and confirm	... see also page 33 , ↓ soft key	<div> <div>SETUP PRINTOUT IDENTIFIER</div> <div>Lot (L ID):</div> <div>ID1: COMPANY</div> <div>ID2: ID2</div> <div>ID3: ID3</div> <div>ID4: ID4</div> <div>ABCDEF GHIJKL MNOPQR STUVWX YZ/-?:#*"%</div> </div>

Step	Press key(s) (or follow instructions)	Display/Output																					
7. Repeat steps 6 and 7 for: ID2: LOCATION ID3: STREET ID4: LOT		<table><tr><th>SETUP</th><th>PRINTOUT</th><th>IDENTIFIER</th></tr><tr><td colspan="3">Lot (L ID):</td></tr><tr><td>ID1:</td><td></td><td>COMPANY</td></tr><tr><td>ID2:</td><td></td><td>LOCATION</td></tr><tr><td>ID3:</td><td></td><td>STREET</td></tr><tr><td>ID4:</td><td></td><td>LOT</td></tr><tr><td><<</td><td><</td><td>^</td></tr></table>	SETUP	PRINTOUT	IDENTIFIER	Lot (L ID):			ID1:		COMPANY	ID2:		LOCATION	ID3:		STREET	ID4:		LOT	<<	<	^
SETUP	PRINTOUT	IDENTIFIER																					
Lot (L ID):																							
ID1:		COMPANY																					
ID2:		LOCATION																					
ID3:		STREET																					
ID4:		LOT																					
<<	<	^																					
8. Save settings, exit the Setup menu and select input mode for IDs	<< soft key ID soft key	<table><tr><td>ID:</td><td></td></tr><tr><td>COMPANY</td><td></td></tr><tr><td>LOCATION</td><td></td></tr><tr><td>STREET</td><td></td></tr><tr><td>LOT</td><td></td></tr><tr><td><<</td><td>Delete</td></tr><tr><td></td><td>v</td></tr><tr><td></td><td>↓</td></tr></table>	ID:		COMPANY		LOCATION		STREET		LOT		<<	Delete		v		↓					
ID:																							
COMPANY																							
LOCATION																							
STREET																							
LOT																							
<<	Delete																						
	v																						
	↓																						
9. Enter name of company (such as Sartorius)	ABC ... see also page 33	<table><tr><td>ID:</td><td></td></tr><tr><td>COMPANY</td><td>SARTORIUS</td></tr><tr><td>LOCATION</td><td></td></tr><tr><td>STREET</td><td></td></tr><tr><td>LOT</td><td></td></tr><tr><td colspan="2">ABCDEF GHIJKL MNOPQR STUVWX YZ/--?:;#"%&</td></tr></table>	ID:		COMPANY	SARTORIUS	LOCATION		STREET		LOT		ABCDEF GHIJKL MNOPQR STUVWX YZ/--?:;#"%&										
ID:																							
COMPANY	SARTORIUS																						
LOCATION																							
STREET																							
LOT																							
ABCDEF GHIJKL MNOPQR STUVWX YZ/--?:;#"%&																							
10. Confirm input	↓ soft key	<table><tr><td>ID:</td><td></td></tr><tr><td>COMPANY</td><td>SARTORIUS</td></tr><tr><td>LOCATION</td><td></td></tr><tr><td>STREET</td><td></td></tr><tr><td>LOT</td><td></td></tr><tr><td><<</td><td>Delete</td></tr><tr><td></td><td>^</td></tr><tr><td></td><td>v</td></tr><tr><td></td><td>↓</td></tr></table>	ID:		COMPANY	SARTORIUS	LOCATION		STREET		LOT		<<	Delete		^		v		↓			
ID:																							
COMPANY	SARTORIUS																						
LOCATION																							
STREET																							
LOT																							
<<	Delete																						
	^																						
	v																						
	↓																						
11. Repeat steps 10 and 11 for LOCATION: GOETTINGEN STREET: WEENDER LANDSTRASSE LOT: 15		<table><tr><td>ID:</td><td></td></tr><tr><td>COMPANY</td><td>SARTORIUS</td></tr><tr><td>LOCATION</td><td>GOETTINGEN</td></tr><tr><td>STREET</td><td>WEENDER LANDSTRASSE</td></tr><tr><td>LOT</td><td>15</td></tr><tr><td><<</td><td>Delete</td></tr><tr><td></td><td>^</td></tr><tr><td></td><td>v</td></tr><tr><td></td><td>↓</td></tr></table>	ID:		COMPANY	SARTORIUS	LOCATION	GOETTINGEN	STREET	WEENDER LANDSTRASSE	LOT	15	<<	Delete		^		v		↓			
ID:																							
COMPANY	SARTORIUS																						
LOCATION	GOETTINGEN																						
STREET	WEENDER LANDSTRASSE																						
LOT	15																						
<<	Delete																						
	^																						
	v																						
	↓																						
12. Place the first sample on the scale (ex.: weight of 110.53214 g)	Place load on scale	<table><tr><td colspan="2">Max 6200 g</td><td colspan="2">d = 0.01 g</td></tr><tr><td>0%</td><td>100%</td><td colspan="2"></td></tr><tr><td colspan="4">+ 210.53 g</td></tr><tr><td>isoTST</td><td>ID</td><td></td><td></td></tr></table>	Max 6200 g		d = 0.01 g		0%	100%			+ 210.53 g				isoTST	ID							
Max 6200 g		d = 0.01 g																					
0%	100%																						
+ 210.53 g																							
isoTST	ID																						
13. Print weight (if desired, perform further weighing operations and print results)		<table><tr><td>COMPANY</td><td>SARTORIUS</td></tr><tr><td>LOCATION</td><td>GOETTINGEN</td></tr><tr><td>STREET</td><td></td></tr><tr><td>WEENDER LANDSTRASSE</td><td></td></tr><tr><td>LOT</td><td>15</td></tr><tr><td>N</td><td>+110.53214 g</td></tr></table>	COMPANY	SARTORIUS	LOCATION	GOETTINGEN	STREET		WEENDER LANDSTRASSE		LOT	15	N	+110.53214 g									
COMPANY	SARTORIUS																						
LOCATION	GOETTINGEN																						
STREET																							
WEENDER LANDSTRASSE																							
LOT	15																						
N	+110.53214 g																						
14. When the weighing is completed delete each ID individually	ID soft key Delete soft key 4 times	<table><tr><td>ID:</td><td></td></tr><tr><td>COMPANY</td><td>SARTORIUS</td></tr><tr><td>LOCATION</td><td>GOETTINGEN</td></tr><tr><td>STREET</td><td>WEENDER LANDSTRASSE</td></tr><tr><td>LOT</td><td>15</td></tr><tr><td><<</td><td>Delete</td></tr><tr><td></td><td>v</td></tr><tr><td></td><td>↓</td></tr></table>	ID:		COMPANY	SARTORIUS	LOCATION	GOETTINGEN	STREET	WEENDER LANDSTRASSE	LOT	15	<<	Delete		v		↓					
ID:																							
COMPANY	SARTORIUS																						
LOCATION	GOETTINGEN																						
STREET	WEENDER LANDSTRASSE																						
LOT	15																						
<<	Delete																						
	v																						
	↓																						

Manual Storage (M+)

Purpose

With this function you can load weight values and calculation results directly from Application 1 (e.g., counting, weighing in percent) or Application 2 (checkweighing, time-controlled functions) into Application 3 (totalizing, formulation, statistics).



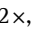




Available Features

- You can assign this function to the fourth or fifth soft key (from the right), i.e. F4 or F5. The soft key designation for this function is: **M+**
- An Application 3 program (totalizing, formulation or statistics) must be running so you can display and print the result

Factory Settings

There are no optional parameters

Preparation

- Turn on the scale: press 
> Sartorius logo is displayed
- Select Extra function (F4) or Extra function (F5) in the Setup menu: press 
- Select **Application parameters**: press the  soft key 2x, then the  soft key
- Select **Extra func. (F4)** or **Extra func. (F5)**:
press the  soft key 3x (or 4x), then the  soft key once
- Select **Man. store in app.3 memory (M+)**
- Confirm **Man. store in app.3 memory (M+)** see also the “Application Menu (Overview)” in the chapter entitled “Configuration”
- Save settings and exit the Setup menu: press the  soft key

Product Data Memory

Purpose

With this function you can enter, store and load data records for initialization of applications, including user-defined data.

You can use this function in combination with a program from Application 1 (e.g., counting, weighing in percent), one from Application 2 (checkweighing, time-controlled functions) and extra functions F4 and F5 (identifiers, second tare memory).

Available Features

- Store up to 300 data records.
- Data records can be created, stored or deleted individually.
- Press the **ProDat** soft key to display data records.
- Define a name for each data record of up to 15 alphanumeric characters; the desired location is displayed in the product data memory.
- Optional configuration in the Setup menu to delete a single character when entering a data record name by pressing **(CF)**. Setup: ... Keypad: CF function for input: Delete last character.
- Data records are displayed in alphabetical order.
- Initialization data set for an application (e.g., wRef, nRef) is saved when you select the Store option. This data is loaded from the product data memory when you access this memory from the corresponding application.
- Use alphanumeric input to search for and display individual data records.
- You can assign this function to the fourth or fifth soft key (from the right), i.e. F4 or F5.
- Error messages are displayed in the text line in plain English.
- Press the **Delete** soft key to delete a data record.

Loading Stored Data:

Data for the block printout is stored in battery-backed memory. The first time you put the scale into operation, it must remain connected to power for a full day. This data remains in memory for approx. 3 months after the equipment is disconnected from AC power.

Factory Settings

No user-definable parameters.

Soft Key Functions

ProDat Toggle to data record display

Delete Delete selected data record

Load Overwrite the initialization data with the selected data record

Change Change the data in the stored data record

New Create a new data record (after entering a data record name).

Store Store the current application data under the selected data record name. If data already exist for this data record, a prompt asks whether this data should be overwritten.

No Answer no to cancel a “delete” or “overwrite” operation

Yes Answer yes to perform the “delete” or “overwrite” operation

Preparation

- Turn on the scale: press **(I/O)**
> Sartorius logo is displayed
- Select Extra function (F4) or Extra function (F5) in the Setup menu: press **(Setup)**
- Select **Application parameters**: press the **↵** soft key 2x, then the **➤** soft key once
- Select **Extra function (F4) or Extra function (F5)**: press the **↵** soft key 3x (or 4x), then the **➤** soft key once
- Select **Product data memory**
- Confirm **Product data memory** see also the “Application Menu (Overview)” in the chapter entitled “Configuration”
- Save settings and exit the Setup menu: press the **◀◀** soft key

Product Data Memory

Practical Example

Create a New Data Record for Initializing the Checkweighing Program, Including: Target Value, Minimum, Maximum

Settings (changes in the factory settings required for this example):
Setup: Application parameters: Extra function (F4): Product data memory
Setup: Application parameters: Application 2: Checkweighing

Step	Key (or instruction)	Display/Output																												
1. If necessary: turn on the scale and enter the settings given above																														
2. In the Checkweighing application, toggle to the input mode for target, minimum and maximum values	Param. soft key	<table><tr><td colspan="2">CHECKWEIGH:</td><td>0.00 g</td><td>▲</td></tr><tr><td>Target:</td><td>Setp= +</td><td>0.00 g</td><td></td></tr><tr><td>Minimum:</td><td>Min = +</td><td>0.00 g</td><td></td></tr><tr><td>Maximum:</td><td>Max = +</td><td>0.00 g</td><td></td></tr><tr><td><<</td><td></td><td></td><td>v J</td></tr></table>	CHECKWEIGH:		0.00 g	▲	Target:	Setp= +	0.00 g		Minimum:	Min = +	0.00 g		Maximum:	Max = +	0.00 g		<<			v J								
CHECKWEIGH:		0.00 g	▲																											
Target:	Setp= +	0.00 g																												
Minimum:	Min = +	0.00 g																												
Maximum:	Max = +	0.00 g																												
<<			v J																											
3. Enter target: 170 g; minimum: 165 g; maximum: 180 g	see the Practical Example for Checkweighing, steps 5 through 9	<table><tr><td colspan="2">CHECKWEIGH:</td><td>170.00 g</td><td>▲</td></tr><tr><td>Target:</td><td>Setp= +</td><td>170.00 g</td><td></td></tr><tr><td>Minimum:</td><td>Min = +</td><td>165.00 g</td><td></td></tr><tr><td>Maximum:</td><td>Max = +</td><td>180.00 g</td><td></td></tr><tr><td><<</td><td></td><td></td><td>v J</td></tr></table>	CHECKWEIGH:		170.00 g	▲	Target:	Setp= +	170.00 g		Minimum:	Min = +	165.00 g		Maximum:	Max = +	180.00 g		<<			v J								
CHECKWEIGH:		170.00 g	▲																											
Target:	Setp= +	170.00 g																												
Minimum:	Min = +	165.00 g																												
Maximum:	Max = +	180.00 g																												
<<			v J																											
4. Toggle to display of product data (existing data records are displayed; in this example, 3 data records have been stored)	ProDat soft key	<table><tr><td colspan="2">PROD. DATA:</td><td colspan="2">PERCENT WGH</td></tr><tr><td>PERCENT WGH40</td><td>Wxx%</td><td>68.75 g</td><td></td></tr><tr><td>CALCULATION8</td><td>Ref</td><td>100 %</td><td></td></tr><tr><td>COUNTING13</td><td></td><td></td><td></td></tr><tr><td><<</td><td>Delete Load</td><td></td><td>v Store</td></tr></table>	PROD. DATA:		PERCENT WGH		PERCENT WGH40	Wxx%	68.75 g		CALCULATION8	Ref	100 %		COUNTING13				<<	Delete Load		v Store								
PROD. DATA:		PERCENT WGH																												
PERCENT WGH40	Wxx%	68.75 g																												
CALCULATION8	Ref	100 %																												
COUNTING13																														
<<	Delete Load		v Store																											
5. Enter a name for the new data record (here: CHW01)	 ABCDEF soft key, C soft key GHIJKL soft key, H soft key STUVWX soft key, W soft key 	<table><tr><td colspan="2">PROD. DATA:</td><td colspan="2"></td></tr><tr><td>CHW01</td><td></td><td></td><td></td></tr><tr><td><<</td><td><</td><td></td><td>New</td></tr></table>	PROD. DATA:				CHW01				<<	<		New																
PROD. DATA:																														
CHW01																														
<<	<		New																											
6. Store current Checkweighing parameters as a data record	New soft key	<table><tr><td colspan="2">PROD. DATA:</td><td colspan="2">NEW: KW01</td></tr><tr><td>CHECKWEIGH</td><td>Setp= +</td><td>170.00 g</td><td></td></tr><tr><td></td><td>Min = +</td><td>165.00 g</td><td></td></tr><tr><td></td><td>Max = +</td><td>180.00 g</td><td></td></tr><tr><td></td><td>Lim=</td><td>0 %</td><td></td></tr><tr><td></td><td>Lim+=</td><td>0 %</td><td></td></tr><tr><td><<</td><td><</td><td></td><td>Store</td></tr></table>	PROD. DATA:		NEW: KW01		CHECKWEIGH	Setp= +	170.00 g			Min = +	165.00 g			Max = +	180.00 g			Lim=	0 %			Lim+=	0 %		<<	<		Store
PROD. DATA:		NEW: KW01																												
CHECKWEIGH	Setp= +	170.00 g																												
	Min = +	165.00 g																												
	Max = +	180.00 g																												
	Lim=	0 %																												
	Lim+=	0 %																												
<<	<		Store																											
7. Confirm	Store soft key	<table><tr><td colspan="2">PROD. DATA:</td><td colspan="2">Data stored</td></tr><tr><td>CHW01</td><td>Setp= +</td><td>170.00 g</td><td></td></tr><tr><td>PERCENT WGH40</td><td>Min = +</td><td>165.00 g</td><td></td></tr><tr><td>CALCULATION8</td><td>Max = +</td><td>180.00 g</td><td></td></tr><tr><td>COUNTING13</td><td>Lim=</td><td>0 %</td><td></td></tr><tr><td></td><td>Lim+=</td><td>0 %</td><td></td></tr><tr><td><<</td><td>Delete Load</td><td>v</td><td>Change</td></tr></table>	PROD. DATA:		Data stored		CHW01	Setp= +	170.00 g		PERCENT WGH40	Min = +	165.00 g		CALCULATION8	Max = +	180.00 g		COUNTING13	Lim=	0 %			Lim+=	0 %		<<	Delete Load	v	Change
PROD. DATA:		Data stored																												
CHW01	Setp= +	170.00 g																												
PERCENT WGH40	Min = +	165.00 g																												
CALCULATION8	Max = +	180.00 g																												
COUNTING13	Lim=	0 %																												
	Lim+=	0 %																												
<<	Delete Load	v	Change																											
8. Exit data record display	<< soft key	<table><tr><td colspan="2">Max 6200 g</td><td colspan="2">d= 0.01g</td></tr><tr><td colspan="2">169.48g</td><td colspan="2"></td></tr><tr><td colspan="2">CHCKW.: n = 0</td><td colspan="2">Setp= +170.00 g</td></tr><tr><td>Cal</td><td>ProDat</td><td>Param.</td><td>Net Show</td></tr></table>	Max 6200 g		d= 0.01g		169.48g				CHCKW.: n = 0		Setp= +170.00 g		Cal	ProDat	Param.	Net Show												
Max 6200 g		d= 0.01g																												
169.48g																														
CHCKW.: n = 0		Setp= +170.00 g																												
Cal	ProDat	Param.	Net Show																											

"FlexPrint" Printout Function

Purpose


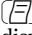
The YAD02IS "Nice Label Express" software from Sartorius lets you load user-defined label printing formats and the corresponding print instruction files in your scale. This software lets you connect any of a number of printers, equipped with a variety of printer fonts, to the RS-232 interface on your scale.


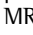
Features

With the "FlexPrint" option activated:

- Print command generates configured printout (if print instruction file exists; see table, next page). Print command generates default printout (if print instruction file does not exist).
- The function that generates an automatic printout upon initialization of an application cannot be used. Initialization data can be output only to a print instruction file.

The following items are output only as standard printouts:

- Calibration/adjustment
- SETUP printouts
-  key in "Identifier:" display page
-  key on "Product data memory:" display page

-  on "Parameter" display page for "Checkweighing" application: print function carried out
- MR function carried out when  is pressed during evaluation (info window) in Totalizing and Statistics applications.
- To recall the file names, software ID and version numbers, for FlexPrint, see the section on "Basic Settings," Info Display."

- Printout for legal metrology: Weight blocks (special weight value formats that are acceptable in legal metrology) are designated by the following line, which is printed in both header and footer and cannot be edited:
"----- CE [M] -----"

Examples of Weight Block Printouts

Without tare:

```
----- CE [M] -----
N   +           348.65 kg
----- CE [M] -----
```


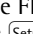
With tare:

```
----- CE [M] -----
G   +           459.70 kg
N   +           348.65 kg
T1  +           111.05 kg
----- CE [M] -----
```

Scale tare (2nd tare memory):



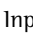
```
----- CE [M] -----
G   +           124.45 kg
N   +           100.00 kg
T1  +           24.00 kg
T2  +             0.45 kg
----- CE [M] -----
```

Preparation

- Turn on the scale: press the  key
- Configure FlexPrint in Setup: press the  key
- Select **Printout**: press ∇ soft key 3 times, and then \rightarrow soft key
- Confirm **Application-defined output**: press \rightarrow soft key
- Select **FlexPrint**: press ∇ soft key 3 times and then \rightarrow soft key
- Select **On**: press ∇ soft key and then \rightarrow soft key
- Save settings and exit Setup: press $\leftarrow \leftarrow$ soft key

Print events

Printouts generated using the "Nice Label Express" software are divided into the following groups:

Event	Explanation	File name for event group:
1.  key with individual values	Print key	PPRINT
2.  key with text input	Input and  key	PDIRECT
3. GLP /GMP header	GLP header	PGMPHEAD
4. GLP /GMP footer	GLP footer	PGMPFOOT
5. Results, Application 1	Animal weighing, MR-CF	PA1RES
6. Results, Application 2	OK values, time-controlled print	PA2RES
7. Results, Application 3	MR, MR-CF	PA3RES
8. Components, Application 1	M+ printout	PA1COMP
9. Components, Application 3	M+/M- printout	PA3COMP

Combining Applications

The following table summarizes the possibilities for combination of the application programs described here. Each line stands for one combination. The weighing function is generally available, and does not have to be combined with a calculating function.

Application 1 (basic function)	Application 2 (control function)	Application 3 (documenting function)
Counting	–	Totalizing
Counting	–	Formulation
Counting	–	Statistics
Weighing in percent	–	Totalizing
Weighing in percent	–	Formulation
Weighing in percent	–	Statistics
Animal weighing	–	Totalizing
Animal weighing	–	Statistics
Recalculation	–	Totalizing
Recalculation	–	Statistics
Calculating	–	Totalizing
Calculating	–	Formulation
Calculating	–	Statistics
Density determination	–	Statistics
Density determination	Time-controlled functions	Statistics
Differential Weighing	–	–
–	Checkweighing	Totalizing
–	Checkweighing	Formulation
–	Checkweighing	Statistics
Counting	Checkweighing	Totalizing
Counting	Checkweighing	Formulation
Counting	Checkweighing	Statistics
Weighing in percent	Checkweighing	Totalizing
Weighing in percent	Checkweighing	Formulation
Weighing in percent	Checkweighing	Statistics
Recalculation	Checkweighing	Totalizing
Recalculation	Checkweighing	Statistics
Calculating	Checkweighing	Totalizing
Calculating	Checkweighing	Formulation
Calculating	Checkweighing	Statistics
–	Time-controlled functions	Totalizing
–	Time-controlled functions	Formulation
–	Time-controlled functions	Statistics
Counting	Time-controlled functions	Totalizing
Counting	Time-controlled functions	Formulation
Counting	Time-controlled functions	Statistics
Weighing in percent	Time-controlled functions	Totalizing
Weighing in percent	Time-controlled functions	Formulation
Weighing in percent	Time-controlled functions	Statistics
Animal weighing	Time-controlled functions	Totalizing
Animal weighing	Time-controlled functions	Statistics
Recalculation	Checkweighing	Totalizing
Recalculation	Checkweighing	Statistics
Calculating	Time-controlled functions	Totalizing
Calculating	Time-controlled functions	Formulation
Calculating	Time-controlled functions	Statistics

Examples of Application Combinations

Example: Counting and checkweighing with statistical evaluation

You want to check a piece count, and have the results that lie within the tolerance range statistically evaluated and printed as a ISO/GMP-compliant record.

Settings (changes in the factory settings required for this example):

Setup: Application parameters: Application 1: Counting

Setup: Application parameters: Application 2: Checkweighing

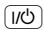
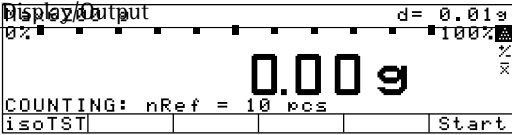
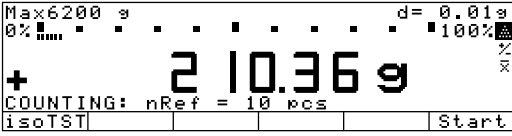
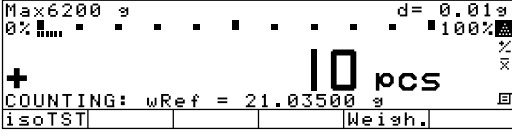
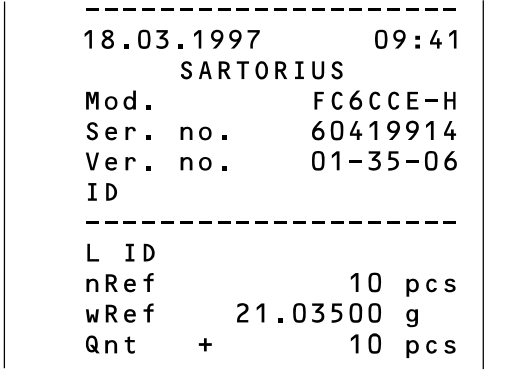
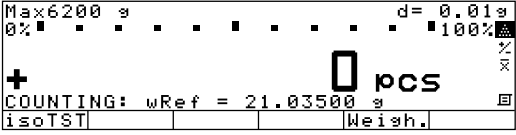
Setup: Application parameters: Application 3: Statistics: Automatic storage: On, first value at stability

Setup: Application parameters: Application 3: Statistics: Source of data for auto storage: Application 2

Setup: Application parameters: Application 3: Statistics: Evaluated value: Calculated

Setup: Application parameters: Application 3: Statistics: Evaluation mode, MR function: Intermediate evaluation, display+print

Setup: Application parameters: Basic application: Printout configuration: ISO/GLP/GMP printout: Always

Step	Key (or instruction)	Display/Output
1. If necessary: turn on the scale and enter the settings given above		
2. Place reference sample quantity	Place parts on the scale on the scale	
3. Initialize the scale	Start soft key	 
4. Remove reference sample quantity	Unload the scale	

Step	Key (or instruction)	Display/Output																												
5. Initialize Checkweighing Toggle to Checkweighing																														
6. Enter target, minimum and maximum values (here: target: 10 pcs; minimum: 7 pcs; maximum: 12 pcs)	Param. soft key , soft key , soft key 																													
7. Store input	soft key	 <table><tr><td>Setp</td><td>+</td><td>10 pcs</td></tr><tr><td>Min</td><td>+</td><td>7 pcs</td></tr><tr><td>Max</td><td>+</td><td>12 pcs</td></tr></table>	Setp	+	10 pcs	Min	+	7 pcs	Max	+	12 pcs																			
Setp	+	10 pcs																												
Min	+	7 pcs																												
Max	+	12 pcs																												
8. Determine first unknown quantity	Place uncounted parts on the scale																													
9. Toggle to Statistics																														
10. Initialize automatic storage	M+ soft key	 <table><tr><td>n</td><td>1</td></tr><tr><td>Qnt</td><td>+ 9 pcs</td></tr></table>	n	1	Qnt	+ 9 pcs																								
n	1																													
Qnt	+ 9 pcs																													
11. Determine further unknown quantities Printout is generated automatically	Place parts to be counted on the scale	<table><tr><td>n</td><td>4</td></tr><tr><td>Qnt</td><td>+ 9 pcs</td></tr><tr><td colspan="2">-----</td></tr><tr><td>Avg.</td><td>+ 10.0 pcs</td></tr><tr><td>s</td><td>+ 0.8 pcs</td></tr><tr><td>srel</td><td>+ 8.00 %</td></tr><tr><td>Total</td><td>40 pcs</td></tr><tr><td>Min</td><td>+ 9 pcs</td></tr><tr><td>Max</td><td>+ 11 pcs</td></tr><tr><td>Diff</td><td>+ 2 pcs</td></tr><tr><td colspan="2">-----</td></tr><tr><td>18.03.1997</td><td>10:26</td></tr><tr><td colspan="2">Name:</td></tr><tr><td colspan="2">-----</td></tr></table>	n	4	Qnt	+ 9 pcs	-----		Avg.	+ 10.0 pcs	s	+ 0.8 pcs	srel	+ 8.00 %	Total	40 pcs	Min	+ 9 pcs	Max	+ 11 pcs	Diff	+ 2 pcs	-----		18.03.1997	10:26	Name:		-----	
n	4																													
Qnt	+ 9 pcs																													

Avg.	+ 10.0 pcs																													
s	+ 0.8 pcs																													
srel	+ 8.00 %																													
Total	40 pcs																													
Min	+ 9 pcs																													
Max	+ 11 pcs																													
Diff	+ 2 pcs																													

18.03.1997	10:26																													
Name:																														

12. End weighing series Statistics are evaluated Final GMP printout is generated																														
13. Delete initialization of the last application																														

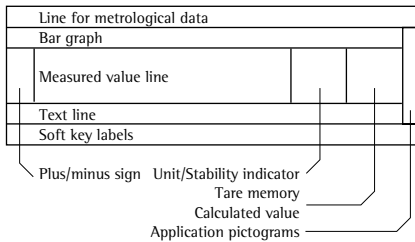
Data Output Functions

There are 3 options for data output:

- Output to the display and control unit
- Output to a printer (generate a printout)
- Output to a peripheral device (e.g., computer) via the interface port

Output to the Display and Control Unit

The display is divided into 9 sections. Information about the scale, the application being used and the sample weighed is output in the following sections:



- Line for metrological data
- Bar graph
- Plus/minus sign
- Measured value line
- Weight unit display, stability symbol display
- Data in tare memory; calculated value
- Application symbol display
- Text line
- Soft key labels

Line for Metrological Data (on scales verified for legal metrology)

This line shows:

- Max 6200 g** - Maximum scale capacity (e.g., 6,200 g)
- Min 0.5 g** - Minimum scale capacity; the weight must not go below this limit when the scale is used in legal metrology
- e = 0.1 g** - Verification interval of the scale; irrelevant if the scale is not used in legal metrology (e.g., 0.1 g)
- d = 0.01 g** - Readability: Indicates the actual scale interval (display increment of the scale) (e.g., 0.01 g)

Bar Graph (overview display)

In the bar graph, weighing results are displayed either



- as a percentage of the maximum scale capacity, or
- in relation to a target value, with tolerance limits indicated.

You can turn off (blank) the bar graph display
(Setup: App: Basic settings: Display: Digit size)

Plus/Minus Sign

This section shows:

- "Busy" symbol
- + -** - Plus or minus sign
- Zero symbol (indicating the scale has been zeroed)

125.03
35
=W* 18.3*0.9

kg
pcs

▲
NET1 NET2 PT

U1 ▲ % ⊗ ⊕

⊗ ⊕

Σ ⊕ ⊗

⊗

⊕

COUNTING: nRef = 10 pcs

Ref.wt. too light

isoTST PT1/T1 S-ID M+

<< < ^ v > ↓

Measured Value Line

This line shows:

- The current weight value (bordered values are invalid in legal metrology)
- Calculated values (e.g., piece count)
- User input (e.g., lot number, equation)

Weight Unit Display, Stability Symbol

This section shows:

- The current weight unit (e.g., kg)
- Designation of other values (e.g., "pcs")

Tare Memory, Calculated Value

This section shows:

- Indication that value is calculated (not valid in legal metrology)
- Indication that the tare memory contains application data

Application Symbols

This column shows:

- Symbol for Application 1 (toggling between weight units, counting, weighing in percent, animal weighing, calculation)
- Symbol for Application 2 (checkweighing, time-controlled functions)
- Symbol for Application 3 (totalizing, formulation, statistics)
- Symbol for current print job
- Symbol for ISO/GMP printout

Text Line

This line contains:

- Explanatory text about the application program (e.g., about "Counting")
- Explanation of error codes

Soft Key Labels

This line shows:

- Texts (abbreviations) to indicate the function assigned to each key
- Symbol for selecting and confirming parameter settings (see also "Operating Design")

Scale Information

In the Setup menu, you can select **Setup: ... Info** for a display of scale information. The display includes:

- Software version number
- Scale version number
- Scale model
- Scale serial number

SETUP	INFO
Version no.:	01-45-01
Bal. ver. no.:	00-20-13
Model:	FC6CCE-HX
Serial no.:	70604025
<<	

Printouts

Purpose

This function enables you to print out weights, other measured values and IDs. You can format the printout to meet different requirements.

Features

Line format: You can configure a data ID code with up to 6 characters at the beginning of each of the values to be printed

Weight ID: You can configure an extra line for identification of each weighed or calculated value using the code **S ID**

Print application parameters: You can generate a printout of the values configured for initialization of an application before printing the measured results

ISO/GMP-compliant printout: You can print out parameters relating to the weighing conditions

Print animal weights: Application-defined, automatic printout of animal weights or of animal weights plus calculated weights after averaging

Configuring Printout Formats

For a number of application programs, you need to set initialization values. All values upon initialization or only the main values can be automatically printed as soon as you have configured this in the Setup menu: **Auto print upon initialization**

Weights and calculated values can be printed as numeric values either with a preceding data ID code (numeric value with 22 characters) or without one (numeric value only 16 characters). See also the section on **Line format** in the chapter entitled "Data Output Functions."

You can generate an ISO/GMP print-out always or only for calibration/adjustment or turn off this option. See also page 100.

Generating an ISO/GLP/GMP Printout
In the Setup menu, you have a choice of three settings:

- No ISO/GLP/GMP printout generated (**Off**)
- ISO/GLP/GMP printout generated only for calibration/adjustment (**Only for calibration/adjustment**)
- Every printout is an ISO/GLP/GMP-compliant report (**Always on**)

Auto print checkweighing results: automatic printout of a weight when it lies within the preset limits at stability

Auto print with time-controlled functions: automatic printout of weights after a preset time has elapsed or at a defined time
Printout of intermediate or final evaluation from the application 3 memory (totalizing, formulation and statistics); generate by pressing the **MR** soft key

Generating Printouts Acceptable for Legal Metrology:
You can configure the Setup menu of the scale to generate data records that are acceptable for legal metrology on a Sartorius printer:

- YDP02
- YDP03
- YDP01IS
- YDP01IS Label
- YDP02IS
- YDP02IS Label
- YDP04IS
- YDP04IS Label

ISO/GMP-compliant Printout/Record

You can have the parameters pertaining to weighing conditions printed before (GMP header) and after (GMP footer) the values from the weighing series. These parameters include:

- Date
- Time at the beginning of a weighing series
- Scale manufacturer
- Scale model
- Model serial number
- Software version
- Lot number (weighing series no.)
- Time at the conclusion of the weighing series
- Field for operator signature

Operating the Scale with an ISO/GMP-capable Documentation Device (Printer)

ISO/GMP-compliant documentation requires a computer with special software. Contact Sartorius for a detailed description for creating this software.

Setting:
Setup: Printout: ISO/GMP printout:
Always on

The record is output to a Sartorius YDP03-OCE Data Printer or a computer.

End GMP printout:

- Press **CF**
End GMP printout while application is active:
This requires the following settings:
Setup: Device: Keys: CF function in application: Clear only selected applications

- Press **CF**
> Text line: CF selected: clear application

- Press the **GLP** soft key

```
-----  
17.01.1997      16:12  
      SARTORIUS  
Mod.      FC6CCE-HX  
Ser. no.    60419914  
Ver. no.    01-35-18  
ID      12345678901234  
-----  
L ID  12345678901234  
nRef      10 pcs  
wRef      1.35274 g  
Qnt      +      235 pcs  
Qnt      +      4721 pcs  
S ID  12345678901234  
Qnt      +      567 pcs  
-----  
17.01.1997      16:13  
Name :  
-----
```

Dotted line
Date/time
Scale manufacturer
Scale model
Scale serial number
Software vers. (display and control unit)
Scale ID no.
Dotted line
Weighing series no.
Application initialization value
Application initialization value
Counting result
Counting result
ID for counting result
Counting result
Dotted line
Date/time
Field for operator signature
Blank line
Dotted line

Record of Internal
Calibration/Adjustment:

```
-----  
17.01.1997      16:24  
      SARTORIUS  
Mod.      FC6CCE-HX  
Ser. no.    60419914  
Ver. no.    01-35-18  
ID  
-----  
L ID  12345678901234  
Internal calibration  
Start:      manual  
Diff. +      0.006 g  
Internal calibration  
completed  
Diff. +      0.000 g  
-----  
17.01.1997      16:25  
Name :  
-----
```

Dotted line
Date/time
Scale manufacturer
Scale model
Scale serial number
Software vers. (display and control unit)
Scale ID no.
Dotted line
Weighing series no.
Calibration adjustment mode
Beginning mode for calibration/adjustment
Difference after calibration/adjustment
Confirmation of completed
calibration/adjustment routine
Difference between current and target
values after calibration
Dotted line
Date/time
Field for operator signature
Blank line
Dotted line

Interface Description

Purpose

Your Factory scale comes equipped with an interface port for connection to a computer or other peripheral device.

You can use an on-line computer to change, start and/or monitor the functions of the scale and the application programs. The interface port also has four data output port lines for the over/under checkweighing program.

⚠ Warning When Using Prewired RS-232 Connecting Cables!

RS-232 cables purchased from other manufacturers often have incorrect pin assignments for use with Sartorius scales. Be sure to check the pin assignment against the chart on page 106 before connecting the cable, and disconnect any lines marked "Internally Connected". Failure to do so may damage or even completely ruin your scale and/or peripheral device.

Available Features

Type of interface:	Serial interface
Operating mode:	Full duplex
Standard:	RS-232
Transmission rates:	150; 300; 600; 1,200; 2,400; 4,800; 9,600; 19,200 baud
Parity:	Space, odd, even
Character format:	1 start bit, 7-bit ASCII, parity, 1 or 2 stop bits
Handshake:	2-wire interface: via software (XON/XOFF); 4-wire interface: via hardware handshake lines (CTS/DTR)
Operating mode:	SBI, xBPI*
Network address**:	0, 1, 2, ..., 30, 31
Data output format of the scale:	16 or 22 characters

* xBPI operating mode: 9,600 baud, 8 bits, odd parity, 1 stop bit

** Network address is only valid in the XBPI mode

Factory Settings:

Transmission rate:	1,200 baud
Parity:	Odd
Stop bits:	1 stop bit
Handshake:	Hardware 1 character after CTS
Operating mode:	SBI
Network address:	0
Print manually/automatically:	Manual after stability
Stop automatic printing:	Not possible
Automatic printout, time-dependent:	After 1 display update
Tare after ind. printout:	Off
Application initialization values:	Off
Line format:	For other applications/GLP (22 characters)

Preparation

- See page 106 for the pin assignment chart

Interface Description

Line Format (Data Output Format)

You can output the values displayed in the measured value line and the weight unit with or without adata ID code

Example: Without data ID code
+ 253 pcs

Example: With data ID code
Qnt + 253 pcs

Configure this parameter in the Setup menu (Setup: Basic settings: Printout configuration: Line format).

The output with data ID code has 16 characters; without data ID code, 22 characters.

Output Format With 16 Characters

Display segments that are not activated are output as spaces. Characters without a decimal point are output without a decimal point.

The following characters can be output, depending on the characters displayed on the scale:

Normal Operation

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	+	*	D	D	D	D	D	D	D	D	*	U	U	U	CR	LF
or	—		*	*	*		
or	*		*	*	*	*	*	*	*	*						

- *: Space
- D: Digit or letter
- U: Unit symbol
- CR: Carriage return
- LF: Line feed

Special Codes

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	*	*	*	*	*	*	—	—	*	*	*	*	*	*	CR	LF
or							H	H								
or							L	L								
or							C									

- *: Space
- : Weight
- H: Overload
- H H: Overload in checkweighing
- L: Underload
- L L: Underload in checkweighing
- C: Calibration/adjustment

Error Codes

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	*	*	*	E	r	r	*	*/#	#	#	*	*	*	*	CR	LF

- *: Space
- # # #: Error code number

Data output example: + 1255.7 g

Position	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16
	+	*	*	*	1	2	5	5	.	7	*	g	*	*	CR	LF

Position 1: Plus or minus sign or space
 Position 2: Space
 Position 3–10: Weight with a decimal point; leading zeros = space
 Position 11: Space
 Position 12–14: Unit symbol or space
 Position 15: Carriage return
 Position 16: Line feed

Data Output With ID Code

When data with an ID code is output, the ID code consisting of 6 characters precedes the data with the 16-character format. These 6 characters identify the subsequent value.

	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22
KI	I	I	I	I	I	+	*	D	D	D	D	D	D	D	D	*	U	U	U	U	CR	LF
	*	*	*	*	*	–	*	*	*	*		
						*	*	*	*	*	*	*	*	*	*	*						

I: ID code character¹⁾
 *: Space
 D: Digit or letter
 U: Unit symbol¹⁾ see “Toggle between Weight Units”
 CR: Carriage return
 LF: Line feed

¹⁾ depends on scale type; e.g., not all units and characters are available on scales verified for use in legal metrology

Special Codes

penn. 1000																					
--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--

*: Space
 –: Weight
 H: Overload
 H H: Overload in checkweighing
 L: Underload
 L L: Underload in checkweighing
 C: Calibration/adjustment

Error Codes

1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22
S	t	a	t	*	*	*	*	*	E	r	r	*	#	#	#	*	*	*	*	CR	LF

*: Space
 # # #: Error code number

ID code characters I¹⁾

S t a t	Status
I D	Identifier
L I D	Weighing series no.
W I D	Weight set number
N o m .	Exact calibration weight
S I D	Sample ID
N U M	Numeric input
T 1	Application tare memory 1
N	Net weight (T1 = 0)
N 1	Net weight (T1 ≠ 0)
Q n t	Quantity
P r c	Percentage
n R e f	Reference sample quantity
p R e f	Reference percentage
w R e f	Average piece weight
W x x %	Reference percentage weight
m D e f	Target value for animal weighing
M u l	Multiplication factor for animal weighing
x – N e t	Result in animal weighing
x – R e s	Calculated result in animal weighing
R e s	Result using equation (Calculation)
S e t p	Target value for checkweighing
M i n	Lower limit for checkweighing
M a x	Upper limit for checkweighing
T i m e	Time that a value was stored
C o m p x x	No. of components in formulation
T o t . c p	Total weight in formulation
n	Transaction counter
T o t a l	Sum of all values
A v g	Average in statistics
s	Standard deviation
s r e l	Variation coefficient
D i f f	Difference between maximum and minimum

Data Input Format

You can connect a computer to your scale to send commands via the scale interface port to control scale functions and applications. The commands sent are control commands and may have different formats; e.g., control commands can have up to 26 characters. Each character must be transmitted according to the settings configured in the Setup menu for data transmission.

Format for Control Commands


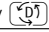
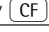
Format 1:	Esc	!	CR	LF						
Format 2:	Esc	!	#	_	CR	LF				
Format 3:	Esc	!	#	<t		(max. 20 <t) <t	_	CR	LF	
Format 4:	Esc	!	#	<t		(max. 8 <t) <t	_	CR	LF	

Esc:	Escape	_:	Underline (ASCII: 95)
!:	Command character	CR:	Carriage RETURN (optional)
#:	Number	LF:	Line FEED (optional)
<t:	Number or letter	max:	depends on command character: i.e. parameter: once the max. length is reached, input received is cut off, rather than discarded as with keyboard input

Format 1

!	Meaning
I	Weighing mode 1
L	Weighing mode 2
M	Weighing mode 3
N	Weighing mode 4
O	Block keys
P	Print
R	Unblock keys
S	Restart
T	Tare and zero
Z	Internal calibration/adjustment
Q	Acoustic signal

Format 2

!#	Meaning
f3	Zero
f4	Tare (without zeroing)
kF1	Soft key 1* Function depends on setting in application
...	
kF6	Soft key 6* program
kF7	Function key 
kF8	Function key 
s3	Function key 
x0	Perform internal calibration**
x1	Print scale model
x2	Print weighing platform serial number
x3	Print weighing platform software version
x4	Print display and control unit software version
x5	Print (GMP) scale ID number
x6	Print weight set ("inventory") number
x7	Print weighing series number

Format 3 (not allowed in the Setup menu)

!#	Meaning
z5	Input (GMP) scale ID number
z6	Input weight set ("inventory") number
z7	Input weighing series number

Format 4

!	Meaning
t	Text input in display

* numbered from right to left

** built-in calibration weight required

Synchronization

During data communication between the scale and an on-line device (computer), messages consisting of ASCII characters are transmitted via the interface. For error-free data communication, the parameters for baud rate, parity, handshake mode and character format must be the same for both units.

You can set these parameters in the Setup menu so that they match those of the on-line device. You can also define parameters in the scale to make data output dependent on various conditions. The conditions that can be configured are described under each of the application program descriptions.

If you do not plug a peripheral device into the scale interface port, no error message will be generated.

Handshake

The scale interface (Sartorius Scale Interface = SBI) has transmit and receive buffers. You can define the handshake parameter in the Setup menu:

- Hardware handshake (CTS/DTR)
- Software handshake (XON, XOFF)

Hardware Handshake

With a 4-wire interface, 1 more character can be transmitted after CTS (Clear to Send).

Software Handshake

The software handshake is controlled via XON and XOFF. When a device is switched on, XON must be transmitted to enable any connected device to communicate.

When the software handshake is configured in the Setup menu, the hardware handshake becomes active after the software handshake.

The data transmission sequence is as follows:

Scale	---	byte	---	Computer
(transmitting	---	byte	---	(receiving
device)	---	byte	---	device)
	---	byte	---	
	<---	XOFF	---	
	---	byte	---	
	---	byte	---	
	...			
	(Pause)			
	...			
	<---	XON	---	
	---	byte	---	
	---	byte	---	
	---	byte	---	
	---	byte	---	

Transmitting Device:

Once XOFF has been received, it prevents further transmission of characters. When XON is received, it re-enables the transmitting device to send data.

Receiving Device:

To prevent too many control commands from being received at one time, XON is not transmitted until the buffer is almost empty.

Activating Data Output

You can define the data output parameter so that output is activated either when a print command is received or automatically and synchronous with the scale display or at defined intervals (see application program descriptions and auto-print setting).

Data Output by Print Command

The print command can be transmitted by pressing \overline{P} or by a software command (Esc P).

Automatic Data Output

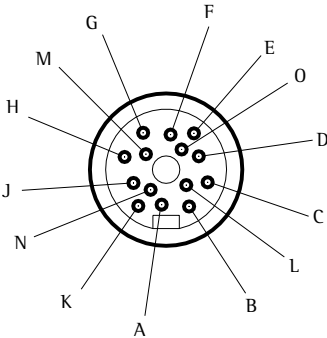
In the "auto print" operating mode, data are output to the interface port without a print command. You can choose to have data output automatically at defined print intervals with or without the stability parameter. Whichever parameter you select, the data will be output as the readouts appear on the scale display. The display update frequency depends on both the model of the scale and the current operating status.

If you select the auto print setting, data will be transmitted immediately the moment you turn on the scale. In the Setup menu you can configure whether this automatic output can be stopped and started by pressing \overline{P} .


Pin Assignment Chart

Female Interface Connector:
14-contact round connector, with screw-lock hardware for cable gland

Pin Assignment Chart
14-contact:
Scale/YAC01FC-X display and control unit



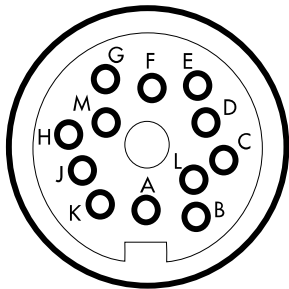12-contact:
Zener Barrier
IDI01-Z, YDI02-Z, YDI03-Z



14-contact Round connector	12-contact Round connector	RS-232 signal (SBI and xBPI)	RS-485 signal ¹⁾ (xBPI)
G	A ³⁾	Control output “heavier”	Control output “heavier”
K	B	Data output (TxD)	RxD – TxD – N
J	C	Data input (RxD)	RxD – TxD – P
N	D	Data Terminal Ready (DTR) –	
M	E	Signal GND	Signal GND
F	G ³⁾	Control output “lighter”	Control output “lighter”
A	H	Clear to Send (CTS)	–
E	J ³⁾	Control output “equal”	Control output “equal”
O	–	Universal switch ²⁾	Universal switch ²⁾
D	L ³⁾	Control output “set”	Control output “set”

Connect low-ohmic shield to the connector case

¹⁾ RS-485 interface available on request
²⁾ See “Universal Remote Switch” in the section “Device Parameters” for more information on the switch functions
³⁾ Control output available only for YDI03-Z

⚠ Important Note:
Only electrical equipment with a maximum voltage rating V_m of 250 V is permitted to be connected to the Zener barrier. The voltage rating V_Z of this Zener barrier is 12 V.


Cabling Diagram

Cabling Diagram (Adapter Cable for PC)

(Adapter cable YCC01-03ISM5 – round – DB25-PC)


Diagram for interfacing a computer via a Zener barrier to the scale using the RS-232C/V24 standard and cables up to 15 m (50 ft) long

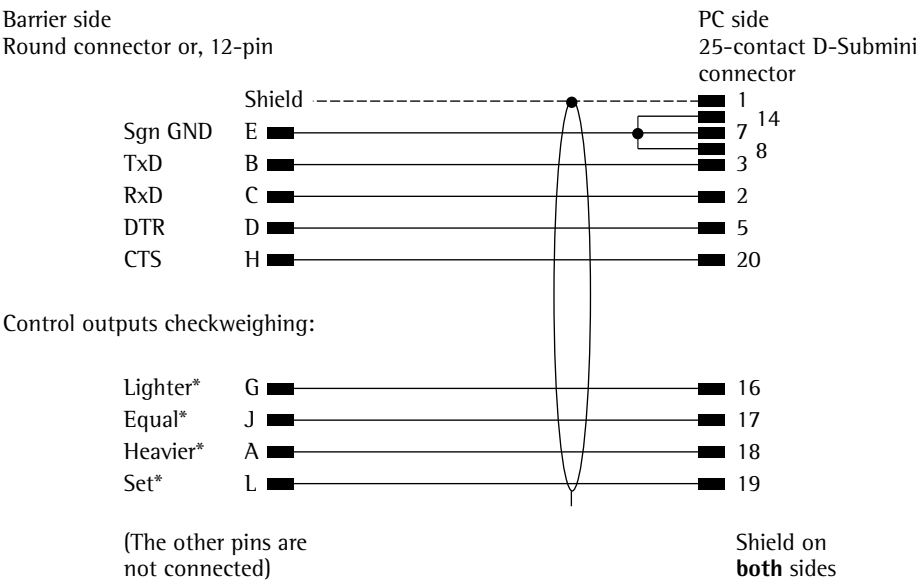Cabling diagram:
Connection assignments for the cable from the Zener barrier to an RS-232 PC interface



* only with YDI03-Z Zener barrier

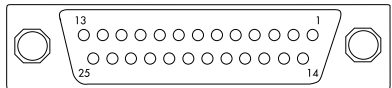Female interface connector:
25-contact, D-Submini DB25S with screw lock hardware
Pin labeling of the 25-contact D-SUB connector




Connector – front view



Male connector used:
(please use connectors with the same specifications):
25-pin D-Submini DB25S, with integrated shielded cable clamp assembly (Amp type 826 985-1C) and fastening screws (Amp type 164 868-1)

Warning: When Using Pre-wired RS-232 Connecting Cables!
RS-232 cables purchased from other manufacturers often have incorrect pin assignments for use with Sartorius scales. Be sure to check the pin assignment against the chart below before connecting the cable, and disconnect any lines marked “Internally Connected”. Failure to do so may damage or even completely ruin your scale and/or peripheral device.

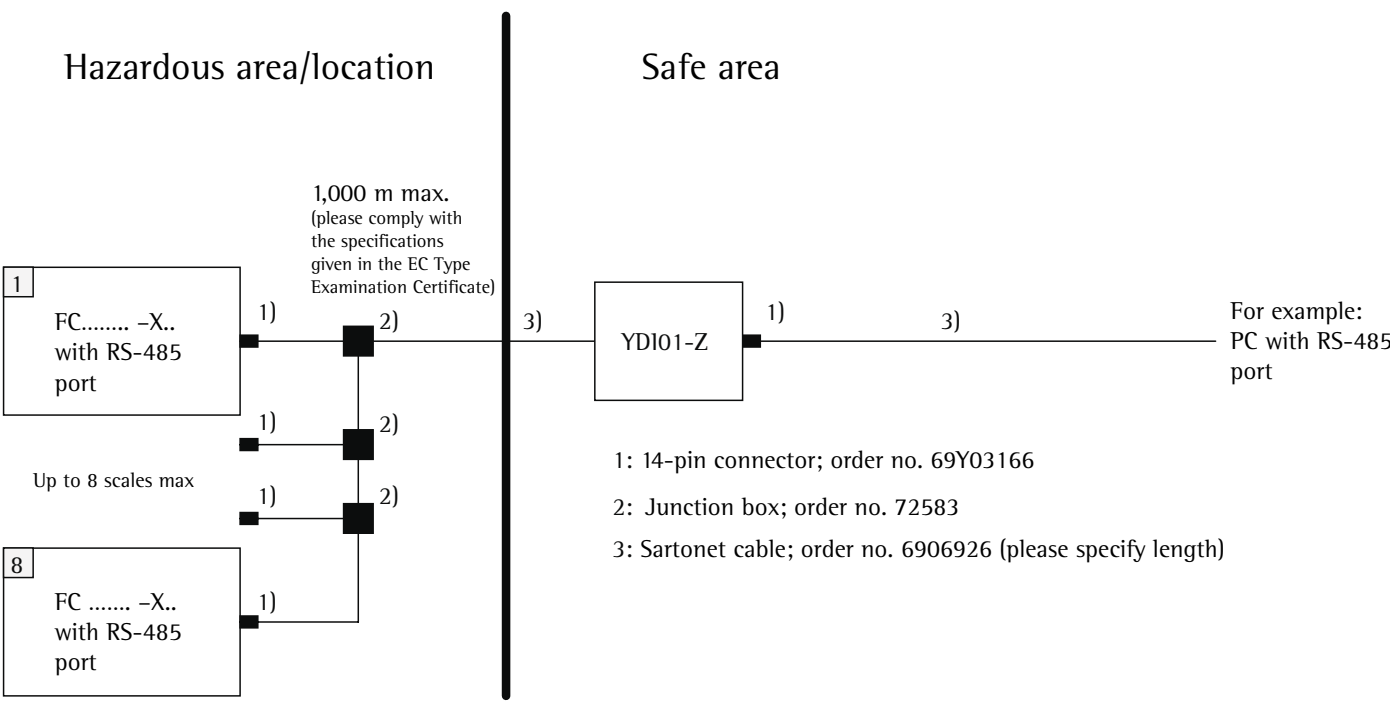Cabling Diagram: RS-485



Cabling Diagram: RS232

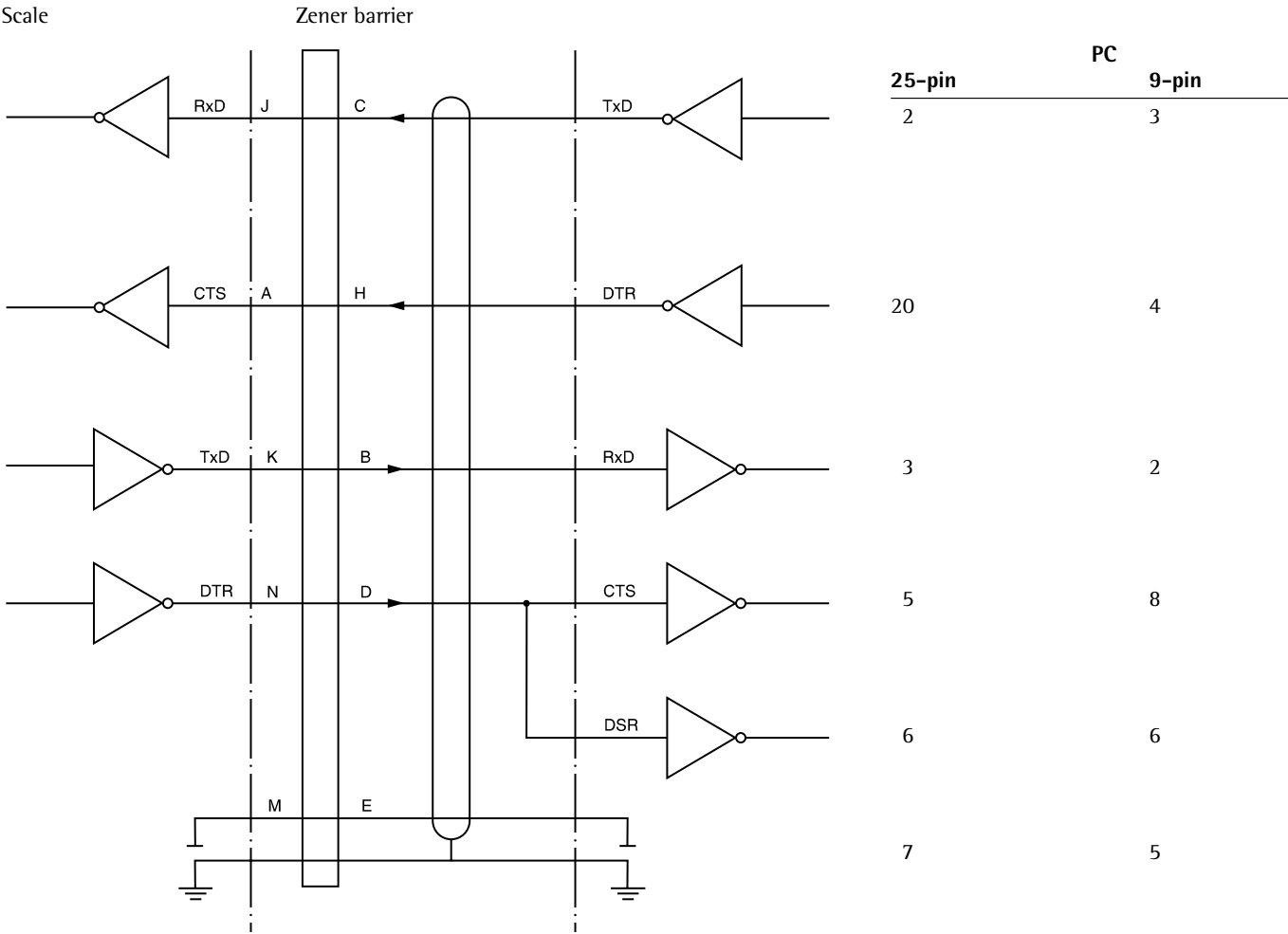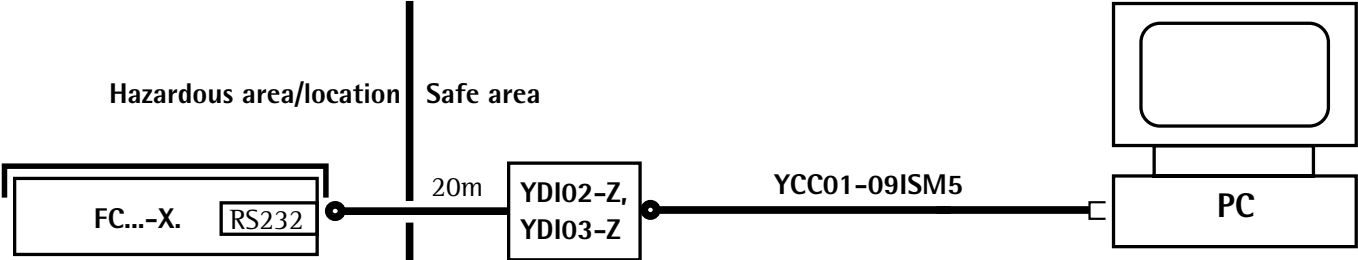


⚠ Important Note
Only electrical equipment with a maximum voltage rating V_m of 250 V is permitted to be connected to the Zener barrier. The voltage rating V_z of this Zener barrier is 12 V.

Hazardous area/location




Cabling Diagram

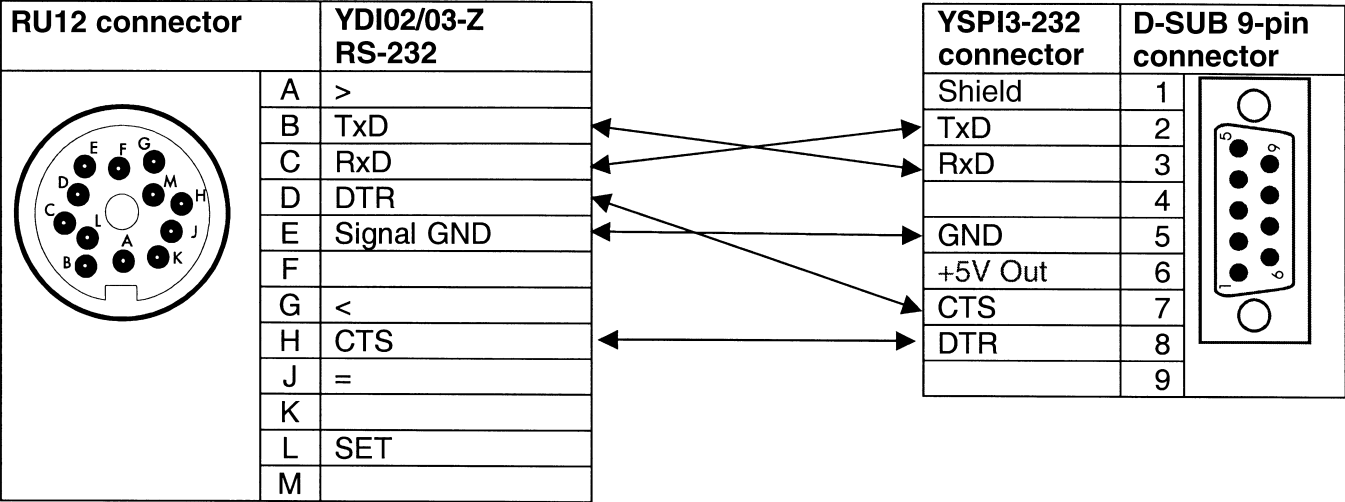Connecting Cable (between YDI02/03-Z and YSPI3-232): **VF2464**



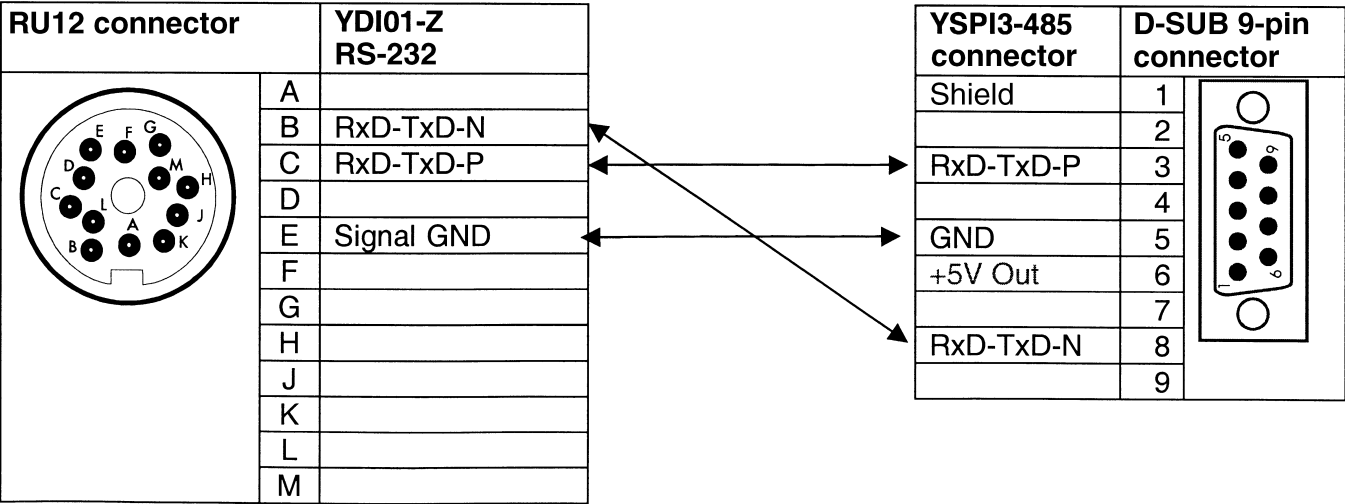Connecting Cable (between YDI01-Z and YSPI3-485): **VF2465**



Circuit:	V1	V2	V3	V4		
Uo	12.6 V	12.6 V	8.6 V	12.6V		
Io	133 mA	133 mA	187 mA	150 mA		
Po	1.46 W	1.46 W	1.51 W	1.68 W		
Co	1µF	1µF	4 µF	1 µF		
Lo	300µH	300 µH	300 µH	300 µH		


Power supply
YPS02-Z.R


II (2) G
[EEx ib] IIC
II (1) D

Non hazardous area


--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--	--

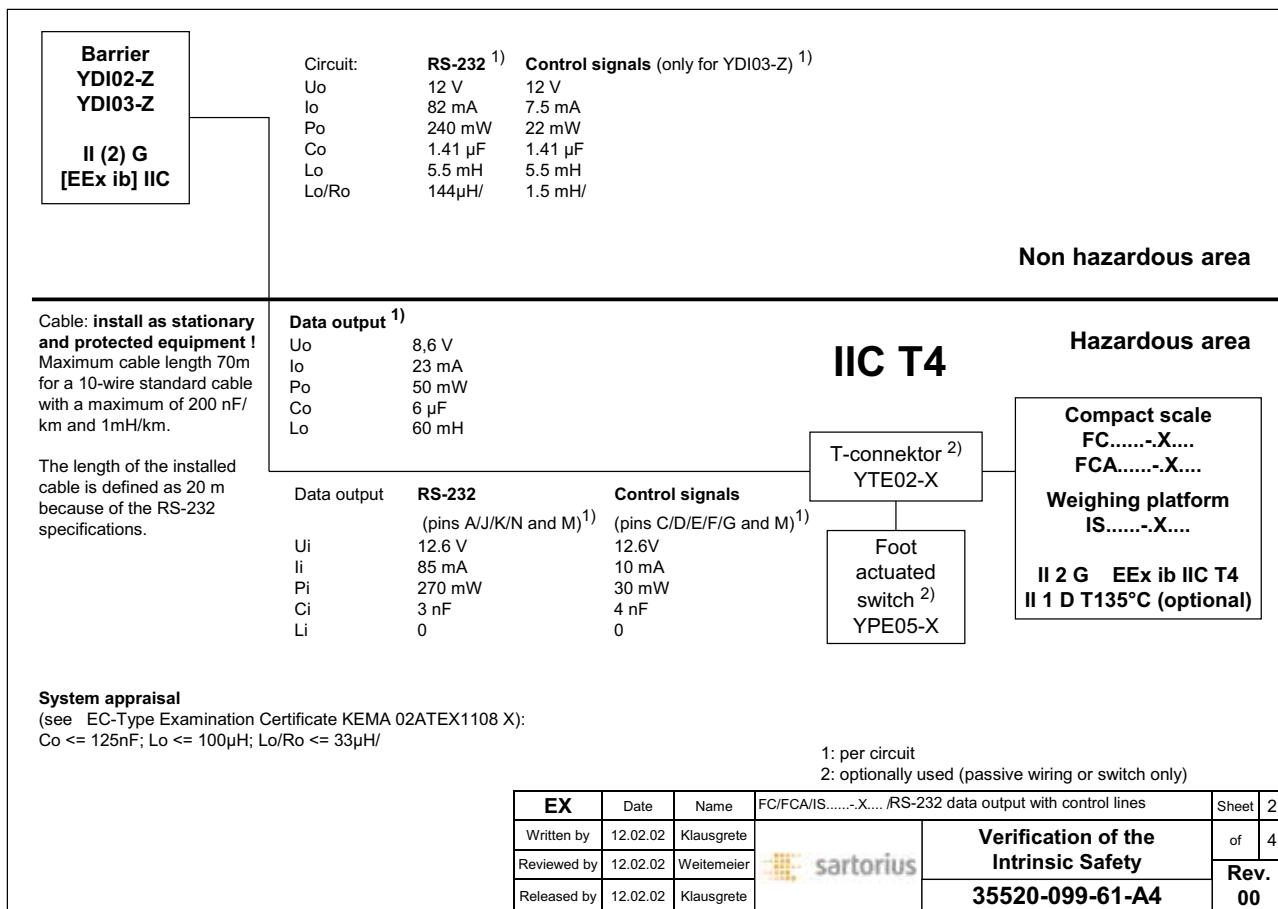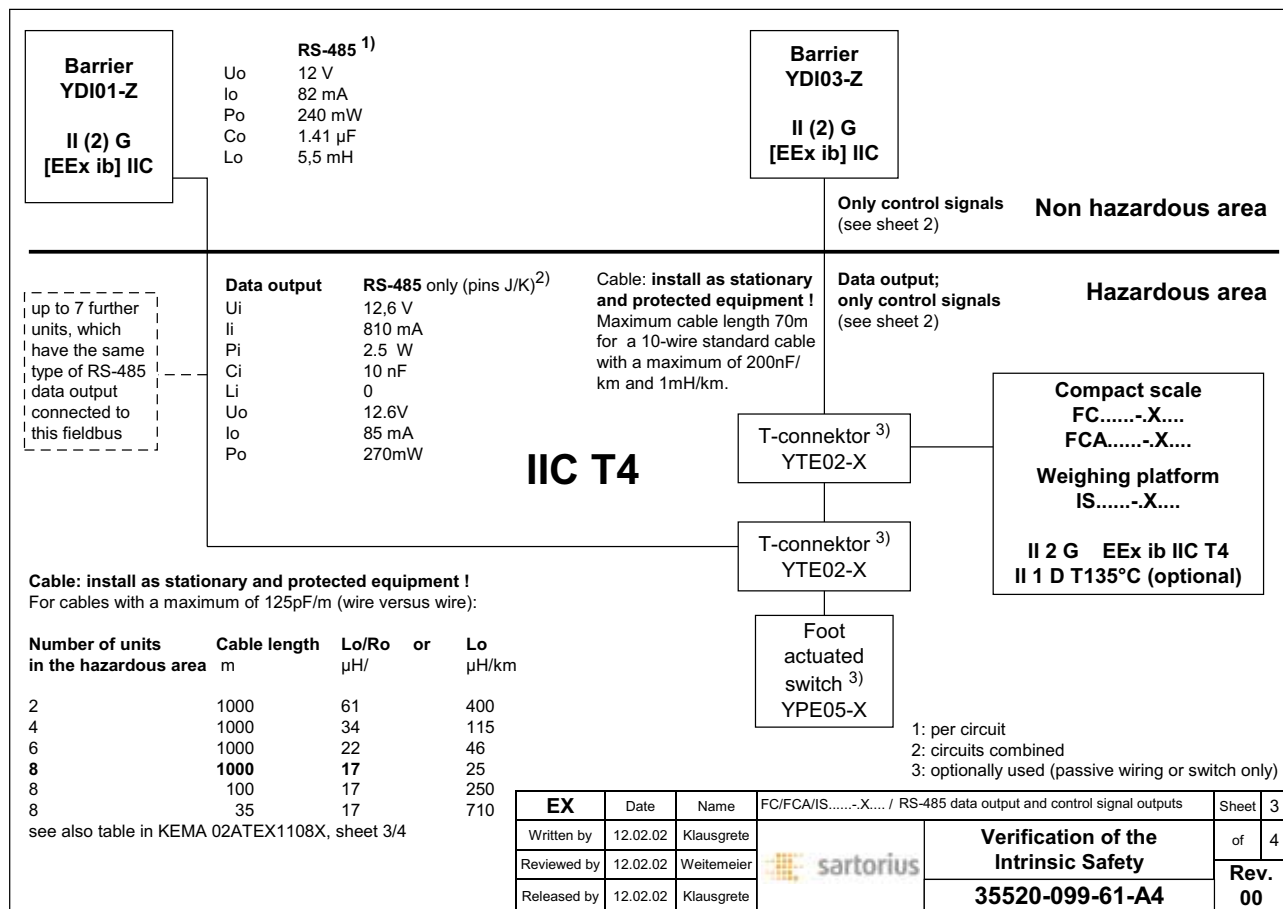Cabling Diagram





Cabling Diagram

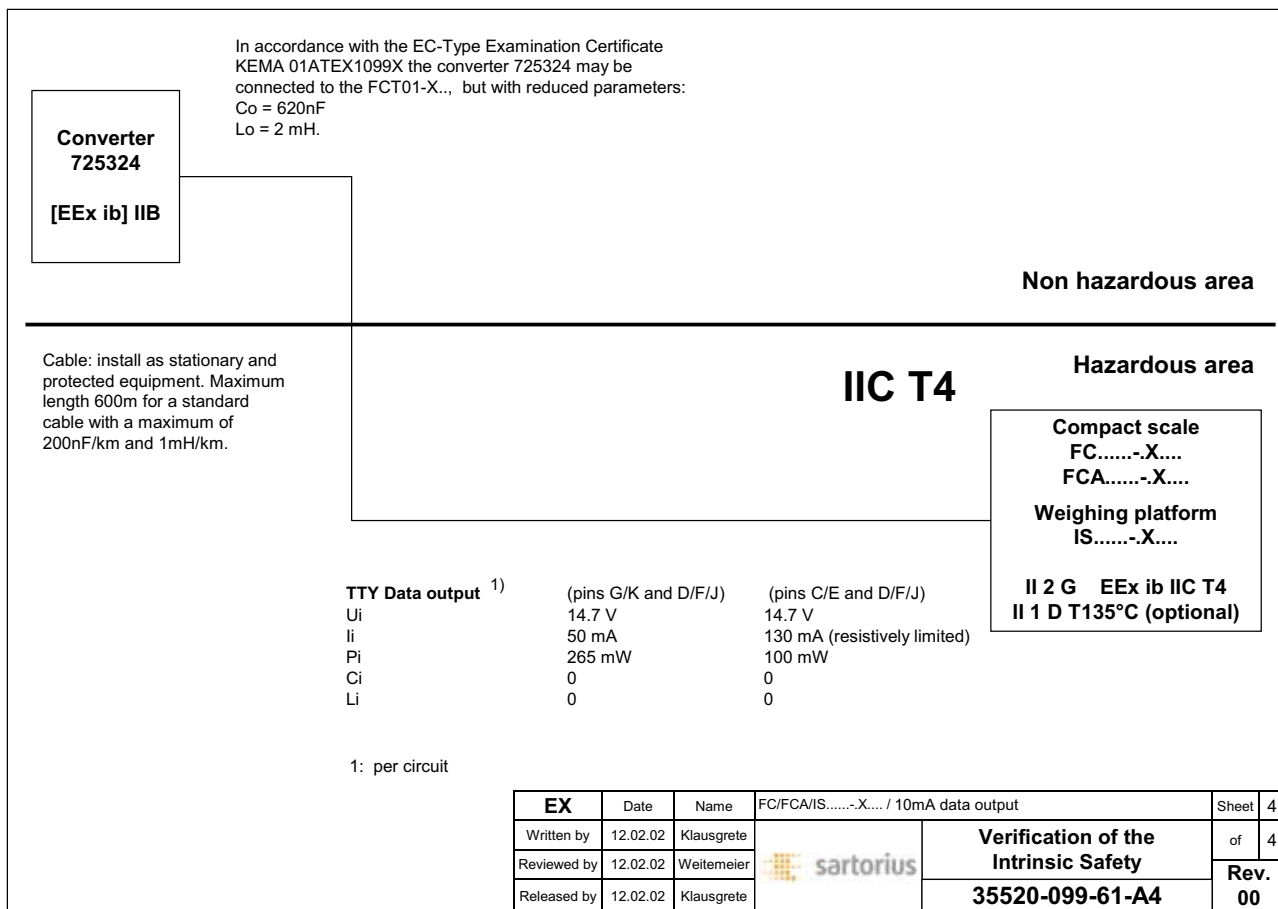

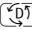
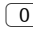
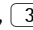
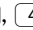
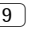
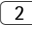



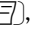
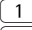
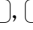
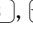
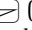
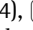
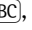
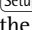



Error Codes

Error codes are displayed in the main display or text line for 2 seconds.
The program then returns automatically to the previous status.

Display	Cause	Solution
No segments appear on the display	No AC power is available The AC adapter is not plugged in Automatic shutoff configured in Setup (code $\overline{B} \ 7 \ 1$)	Check the AC power supply Plug in the AC adapter Press \overline{ON} to switch on the scale or select code $\overline{B} \ 7 \ 2$ in Setup ("no automatic shutoff")
H	The load exceeds the scale capacity	Unload the scale
L or Err 54	The weighing pan is not in place	Place the weighing pan on the scale; may have to be turned off and on again if "L" is displayed
Err 01 > Display range	Data output not compatible with output format	Change the configuration in the Setup menu
Err 02 Cal. n. possible	Calibration/adjustment condition not met, e.g., – The scale was not tared – The scale is loaded	Calibrate only when zero is displayed Press \overline{Tare} to tare Unload the scale
Err 03 Cal./adj. interrupt	Calibration/adjustment could not be completed within a certain time	Allow the scale to warm up again and repeat the adjustment process
Err 05	Built-in calibration weight does not return to initial position	Contact your local Sartorius Service Center
Err 06 Int. wt. defective	Built-in calibration weight is defective	Contact your local Sartorius Service Center
Err 07 Function blocked	Function not allowed in scales verified for use in legal metrology	Contact your local Sartorius Service Center for information on having the settings changed
Err 08* <> zero range	The load on the scale is too heavy to zero the readout	Check whether the "power-on zero range" is set
Err 09* < 0 not allowed	Taring is not possible when the gross weight is \leq zero	Zero the scale
Err 10 Tare fct. blocked	Tare key and 2nd tare memory are blocked when there is data in the tare memory for the formulation application	Press \overline{CF} to clear the formulation application; the tare key and 2nd tare memory are then accessible
Err 11 Tare2 blocked	Tare memory not allowed	Check the tare value entered
Err 12 Tare2 > Max.	Tare memory greater than weighing range or range limits	Check sample/container
Err 17 Adj.-wt. > Max.	Internal adjustment is not possible because preload is too heavy	Reduce the preload or change the configuration
Err 30 Print fct. blocked	Interface port for printer output is blocked	Contact your local Sartorius Service Center

* = occurs only via the SBI interface (ESC f3_/f4_)


Display/Problem	Cause	Solution
Err 31 Print fct. blocked	Interface handshake interrupted (XOFF, CTS)	Transmit XON, then CTS
Ref.wt. too light	Error in storing reference weight (with the counting or weighing-in-percent application)	Weight too light or there is no sample on the scale
Cannot update	Reference updating not possible (with the counting application)	See "Counting" in "Operation" for reference updating criteria
Not a number xxxxx Too low xxxxx Too high	Input wrong (for any application program), e.g., alphabetic input not allowed	Follow the instructions for the application programs
Too many char.	Input text too long	Allowable text lengths, incl. decimal point: – S ID, NUM, L ID, ID: max. 20 characters – W ID: max. 14 characters
Equation too long	Equation exceeds 28 characters	Limit equation to 28 characters
Limits different than unit	The unit for the tolerance limits entered during checkweighing is different than that required for the current application	Change the tolerance limits to fit the application
Err 10x x = 1 : x = 2 : x = 3 : x = 4 : "Checkerboard" pattern displayed continuously	Key is stuck Key pressed when switching on the scale:  (F1, F2, F5, F6),  (CF)  (F3),  (0),  (3),  (4),  (9)  (2),  (5),  (6),  (.),  (Tare)-right  (1),  (7),  (8),  (F4),  (ABC),  (Tare)-left  (Setup) key was pressed when turning on the scale, or is stuck	Release key or Contact your local Sartorius Service Center
Err 340	Operating parameter (EEPROM) is wrong	Contact your local Sartorius Service Center
No WP	Weighing platform is defective	Contact your local Sartorius Service Center
blocked	Function blocked	none
The special code  remains displayed	None of the keys has been pressed since the scale was turned on	Press a key
The weight readout changes constantly	Unstable ambient conditions Too much vibration, or the scale is exposed to a draft A foreign object is caught between the pan and the scale housing	Set up the scale in another area Change Setup configurations to adapt the scale to the ambient conditions Remove the foreign object
The weight readout is obviously wrong	The scale has not been calibrated/adjusted The scale was not tared before weighing The scale is not level The dust cover is caught under the weighing pan	Calibrate/adjust the scale Tare before weighing Level the scale See "Replacing the Dust Cover" in the chapter "Care and Maintenance"

If any other errors occur, contact your local Sartorius Service Center!

Overview

General Views of the Scales

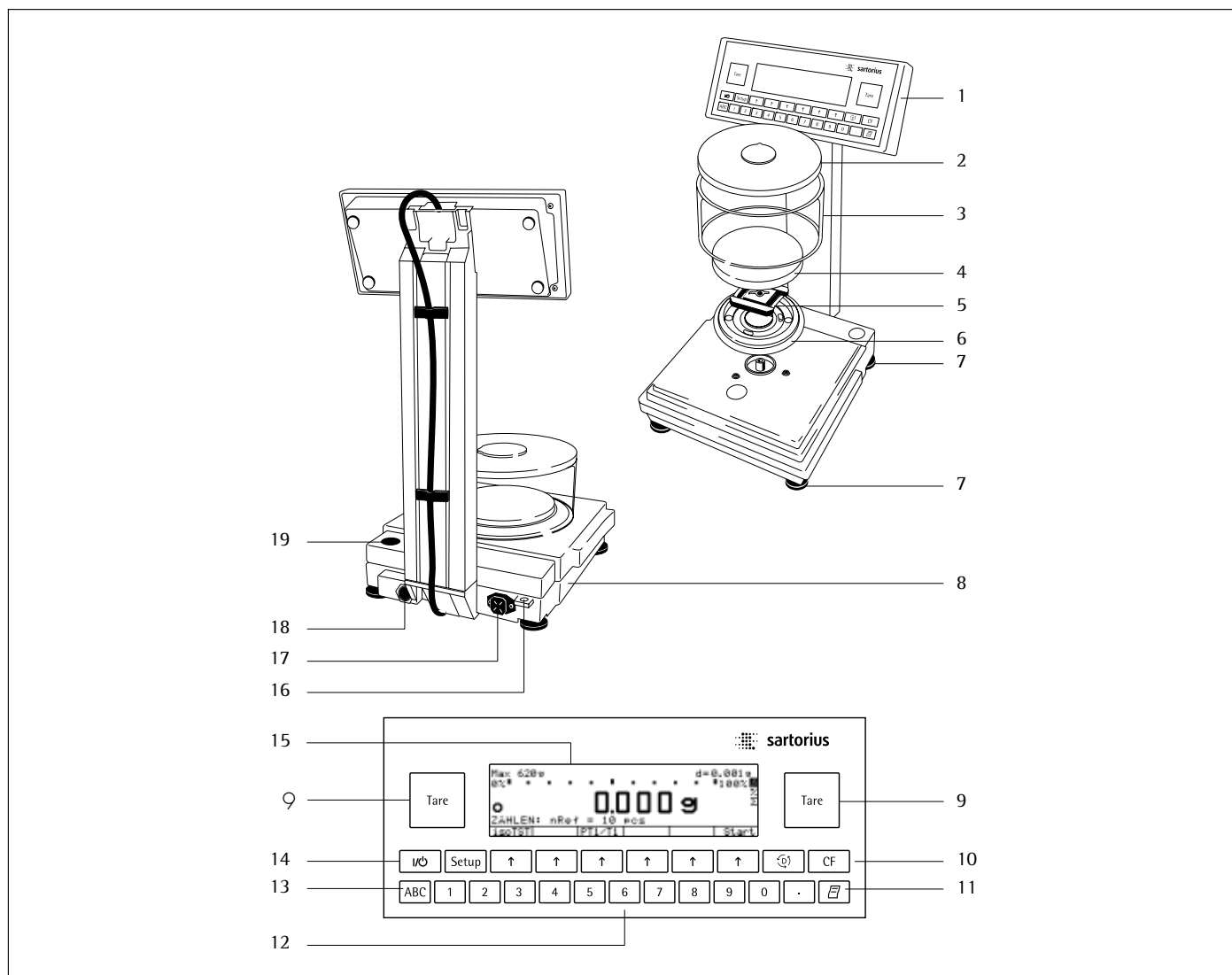FC06BBE-SX(CE)




Pos.	Designation	Spare Part Order No.	Pos.	Designation	Spare Part Order No.
1	Display and control unit		13	Numeric keys	
2	Draft shield cover	69 LP0002	14	Shift key for entering letters	
3	Glass cylinder	69 14290	15	On/off key	
4	Weighing pan	69 LP0004	16	Display	
5	Pan support	69 LP0005	17	Lug for attaching an anti-theft locking device	
6	Shield disk	69 LP0003	18	Main grounding terminal	
7	Leveling foot	69 B20005	19	AC power socket	
8	Menu access switch		20	Data interface port	
9	Metrological ID label (only on scales verified for legal metrology in the EU)		21	Level indicator	
10	Tare key			Not shown:	
11	Function keys			Dust cover for weighing platform	69 60FB01
12	Print key			Dust cover for display&t control unit	69 60LP03
				Caps and plugs for covering ports (set)	69 B20009

Overview of the Scales

FC6CCE-HX(CE), FC2CCE-SX(CE), FC12CCE-SX(CE), FC12CCE-IXCE, FC6CCE-SX(CE)

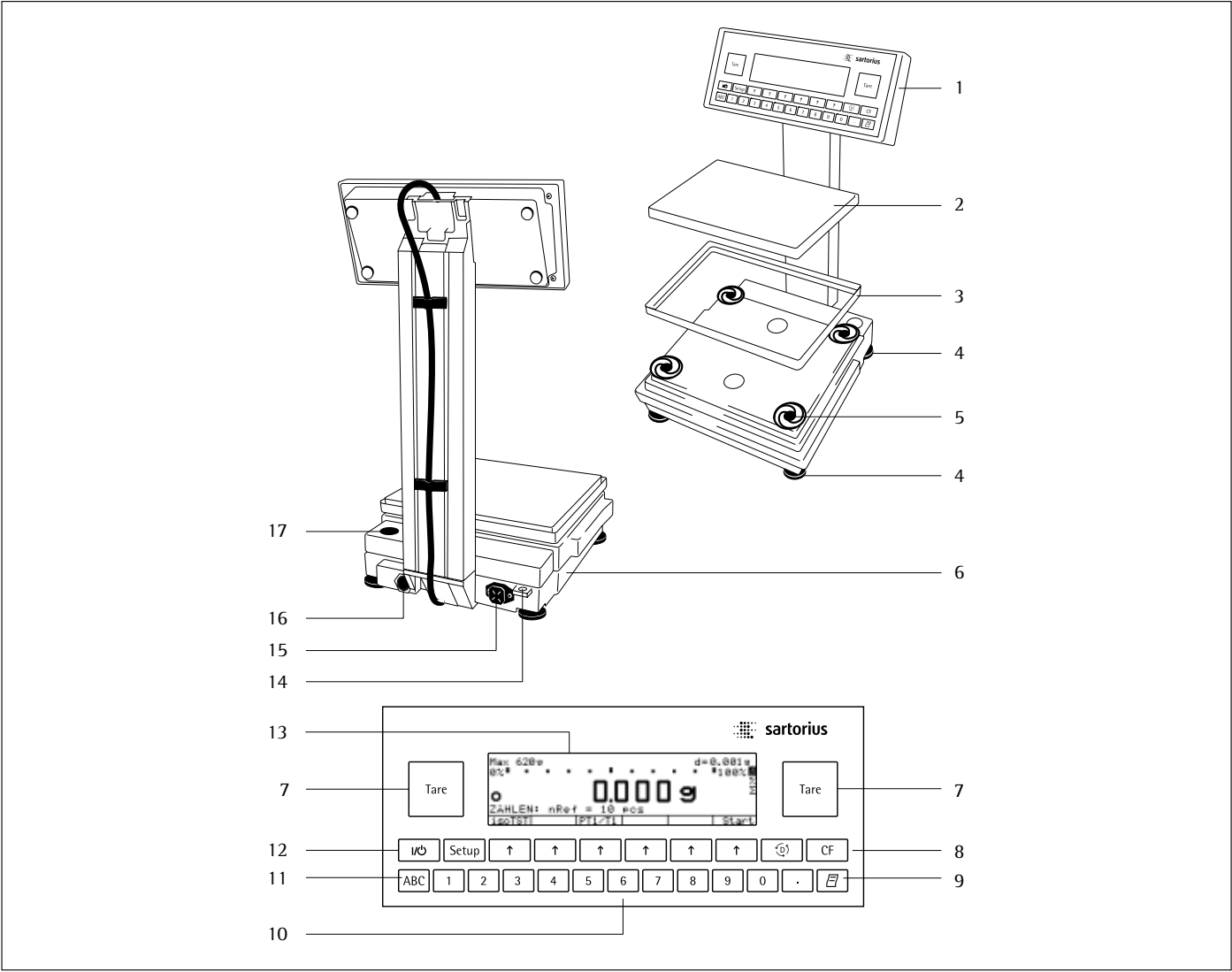“-CE” identifies the scale as verified for legal metrology in the EU*



Pos.	Designation	Spare Part Order No.	Pos.	Designation	Spare Part Order No.
1	Display and control unit		12	Shift key for entering letters	
2	Weighing pan/load plate	69 LP0007	13	On/off key	
3	Pan shield (depends on type of model)	69 LP0008	14	Display	
4	Leveling foot	69 B20005	15	Lug for attaching an antitheft locking device	
5	Menu access switch		16	Main grounding conductor	
6	Shock absorber	69 LP0010	17	AC power socket	
7	Metrological ID label (only on scales verified for legal metrology in the EU)		18	Data interface port	
8	Tare key		19	Level indicator	
9	Function keys				
10	Print key				
11	Numeric keys				

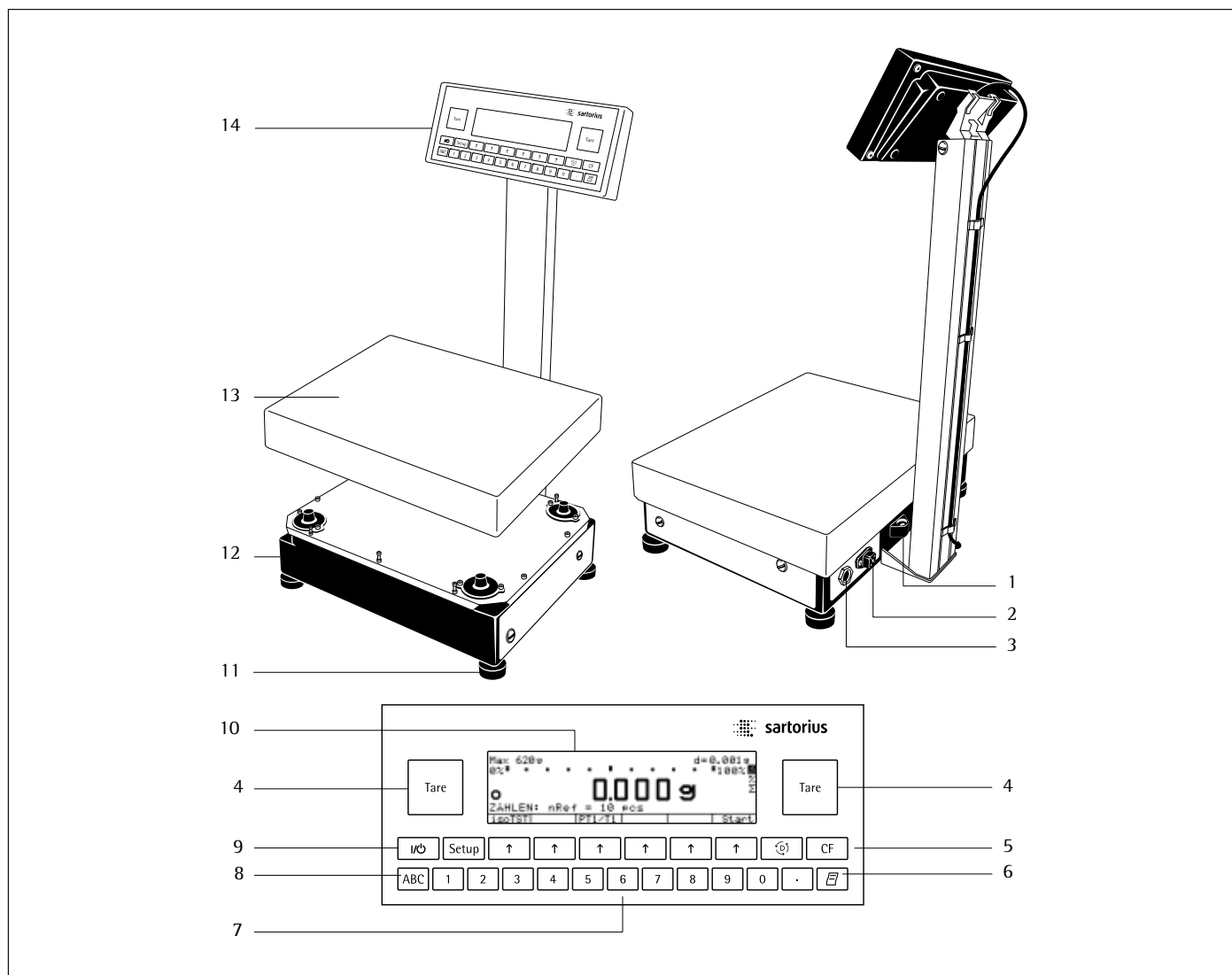
* including the Signatories of the Agreement on the European Economic Area

Not shown:
Dust cover for weighing platform 69 60FB02
Dust cover for display&control unit 69 60LP03
Caps and plugs for covering ports (set) 69 B20009

Overview of the Scales

FC34EDE-HX(CE), FC34EDE-PX(CE), FC16EDE-HX(CE), FC12EDE-PX(CE), FC64EDE-SX(CE), FC64EDE-HX

“CE” identifies the scale as verified for legal metrology in the EU*




Pos.	Designation	Spare Part Order No.	Pos.	Designation	Spare Part Order No.
1	Main grounding terminal		11	On/off key	
2	Level indicator		12	Display	
3	AC power socket	69 14290	13	Leveling foot	69 LC0093
4	Data interface port	69 LP0004	14	metrological ID label (only on scales verified for legal metrology in the EU)	
5	Manu access switch		15	Load plate	
6	Tare key		16	Display and control unit	
7	Function keys				
8	Print key				
9	Numeric keys				
10	Shift key		Not shown: Dust cover for display&control unit		69 60LP03

* including the Signatories of the Agreement on the European Economic Area

General View of the Terminal

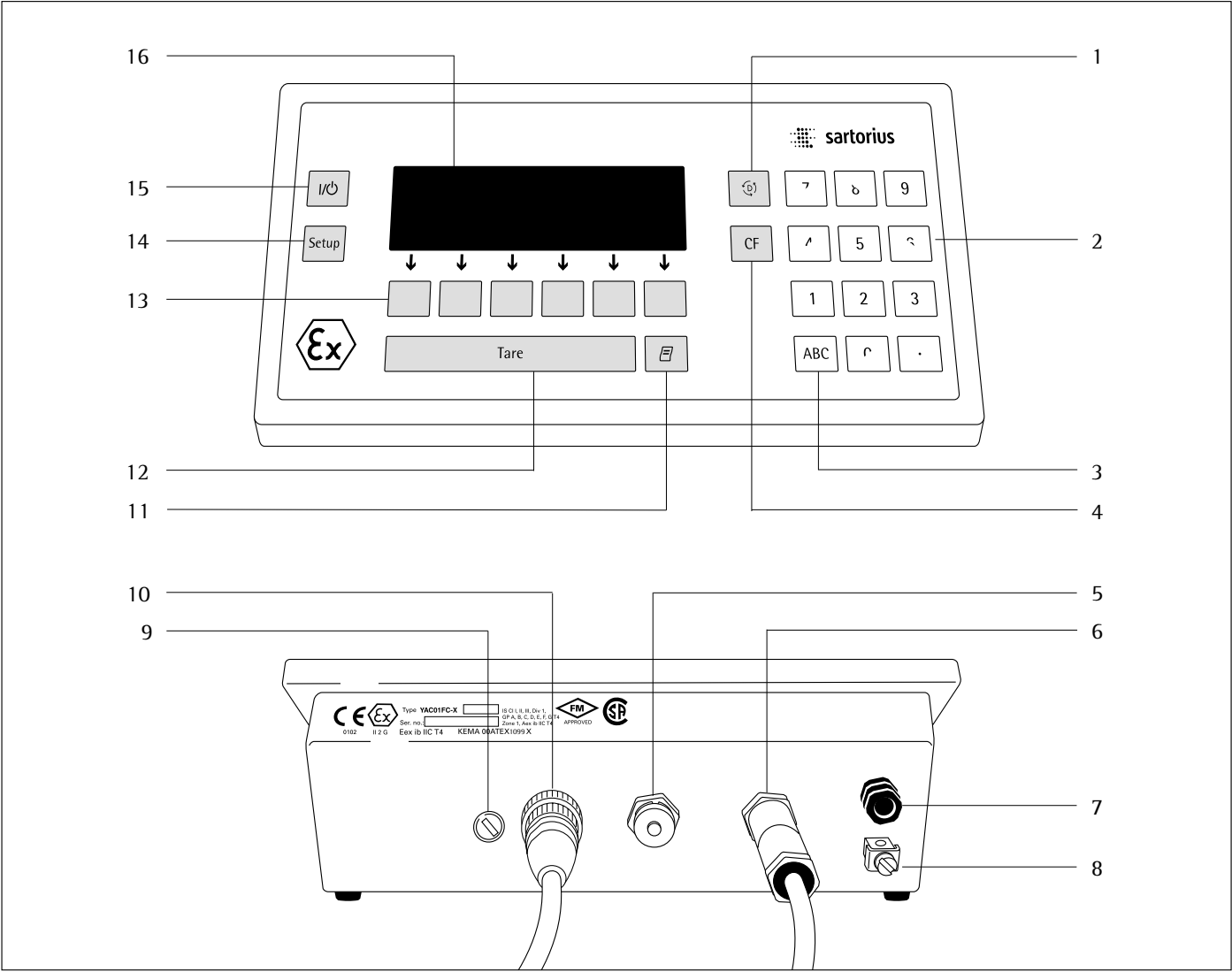FCA Models (this example: YAC01FC-X terminal)



Pos.	Designation
1	Toggle key for changing application programs
2	Numeric keys
3	Shift key for entering letters
4	CF key (clear function)
5	Serial communications port (14-pin)
6	Power jack
7	Vent valve
8	Terminal for an equipotential bonding conductor
9	Plug covering access switch

Pos.	Designation
10	Cable gland for connecting a weighing platform
11	Print key
12	Tare key
13	Function keys
14	Setup key
15	On/standby key
16	Display

Specifications

Standard Models

General Specifications:

AC power source/power requirements	AC adapter, 90 V AC (min.) to 264 V AC (max.)
Frequency	48 – 60 Hz
Allowable ambient operating temperature	0 ... +40 °C (273 ... 313 K, 32 °F ... 104 °F)
Adaptation to ambient conditions	By selection of 1 of 4 optimized filter levels
Display update (depends on the filter level selected)	0.1 – 0.4
Power consumption	25 VA
Hours of operation with fully charged YRB 06 Z external battery pack, approx.	14 h
Selectable weight units	Grams, kilograms, carats, pounds, ounces, Troy ounces, Hong Kong tael, Singapore tael, Taiwanese tael, grains, pennyweights, milligrams, parts per pound, Chinese tael, mommes, Austrian carats, tola, baht and mesghal
Selectable application programs	Mass unit conversion, counting, weighing in percent, animal weighing, calculation, density determination, over/under checkweighing, time-controlled functions, totalizing, statistics, 2nd tare memory, identifiers, product data memories
Built-in interface	RS-232C
Format:	7-bit ASCII, 1 start bit, 1 or 2 stop bits
Parity:	odd, even or space
Transmission rates:	150 to 19,200 baud
Handshake:	Software or hardware

Specifications of the Individual Models:

Model		FC06BBE-SX	FC6CCE-HX	FC2CCE-SX	FC12CCE-SX	FC6CCE-SX
Readability	g	0.001	0.01	0.01	0.1	0.1
Weighing capacity	g	620	6,200	2,200	12,000	6,200
Max. capacity	kg	3	25	10	50	50
Tare range (subtractive)	g	–620	–6,200	–2,200	–12,000	–6,200
Electronically compensated preload (without restricting weighing range)	g	93	–	110	1,200	1,240
Max. preload when starting calibration/adjustment (scale must be zeroed)	g	110	5,200	1,300	10,000	6,400
Repeatability (standard deviation)	≤±g	0.001	0.01	0.01	0.05	0.05
Linearity	≤±g	0.002	0.02	0.02	0.2	0.1
Sensitivity drift within +10 ... +30 °C	≤±/K	2 · 10 ⁻⁶	2 · 10 ⁻⁶	2 · 10 ⁻⁶	4 · 10 ⁻⁶	4 · 10 ⁻⁶
Response time (average)	s	1.5	1.5	1.5	1	1
External calibration weight (of at least accuracy class...)	g	500 (E2)	5,000 (E2)	2,000 (F1)	5,000 (F1)	5,000 (F2)
Load plate	mm	Ø 130	218×200	218×200	218×200	218×200
Dimensions (W×D×H)	mm	240×294×86	240×294×86	240×294×86	240×294×86	240×294×86
Net weight, approx.	kg	7	8.4	7.3	7.3	7.3
Dust and water protection rating according to EN 60529		IP54				

Specifications

Model		FC64EDE-H, FCA64EDE-HX	FC34EDE-H, FCA34EDE-HX	FC16EDE-H, FCA16EDE-HX
Readability	g	0.1	0.1	0.1
Weighing capacity	kg	64	34	16
Max. overload capacity	kg	130	130	130
Tare range (subtractive)	kg	–64	–34	–16
Electronically compensated preload (without restricting weighing range)	kg	13	4	4
Max. preload when starting calibration/adjustment (scale display does not have to be zeroed)	kg	approx. 45	approx. 21	approx. 19
Repeatability	≤±g	0.1	0.1	0.05
Linearity	≤±g	0.5	0.2	0.2
Sensitivity drift within +10 ... +30°C	≤±/K	3 · 10 ^{–6}	2 · 10 ^{–6}	2 · 10 ^{–6}
Response time (average)	s	1.5	1.5	1.5
External standard calibration weight (of at least accuracy class...)	kg	10 (F1)	10 (F1)	10 (F1)
Platform dimensions	mm	300×400		
Net weight, approx.	kg	16.0		
Dust and water protection rating according to EN 60529		IP44	IP65	IP65

Model		FC34EDE-PX, FCA34EDE-PX	FC12EDE-PX, FCA12EDE-HX	FC64EDE-SX, FCA64EDE-SX
Readability	g	0.1/0.2/0.5	0.1/0.2	1
Weighing capacity	kg	8/16/34	6/12	64
Max. overload capacity	kg	130	130	130
Tare range (subtractive)	kg	–34	–12	–64
Electronically compensated preload (without restricting weighing range)	kg	4	4	13
Max. preload when starting calibration/adjustment (scale display does not have to be zeroed)	kg	approx. 21	approx. 10	approx. 45
Repeatability	≤±g	0.05/0.05/0.1	0.05/0.05	0.3
Linearity	≤±g	0.2	0.2	1
Sensitivity drift within +10 ... +30°C	≤±/K	2 · 10 ^{–6}	2 · 10 ^{–6}	3 · 10 ^{–6}
Response time (average)	s	1.5	1.5	1.5
External standard calibration weight (of at least accuracy class...)	kg	10 (F2)	10 (F1)	10 (F2)
Platform dimensions	mm	300×400		
Net weight, approx.	kg	16.0		
Dust and water protection rating according to EN 60529		IP65		

Model		FCA150IGG-HX	FCA300IGG-HX
Readability	g	1	2
Weighing capacity	kg	150	300
Max. overload capacity	kg	600	600
Tare range (subtractive)	kg	150	300
Electronically compensated preload (without restricting weighing range)	kg	30	60
Max. preload when starting calibration/adjustment (scale display does not have to be zeroed)	kg	approx. 120	approx. 300
Repeatability	≤±g	1	2
Linearity	≤±g	4	8
Sensitivity drift within +10 ... +30°C	≤±/K	$2.5 \cdot 10^{-6}$	$2.5 \cdot 10^{-6}$
Response time (average)	s	1.5	1.5
External standard calibration weight (of at least accuracy class...)	kg	50 (F2)	50 (F1)
Platform dimensions	mm	800×600	
Net weight, approx.	kg	70	
Dust and water protection rating according to EN 60529		IP67	

Specifications

Models Verified by the Manufacturer, with EC Type Approval

General Specifications

AC power source/power requirements	AC adapter, 90 V AC (min.) to 264 V AC (max.)		
Frequency	48 – 60 Hz		
Allowable ambient temperature range	0...+40°C (273...313 K, 32°F...104°F)		
Adaptation to ambient conditions	By selection of 1 of 4 optimized filter levels		
Display update (depends on the filter level selected)	0.1–0.4		
Power consumption	25 VA		
Selectable weight units	Grams, kilograms		
Selectable application programs	Mass unit conversion, counting, weighing in percent, animal weighing, calculation, density determination, over/under checkweighing, time-controlled functions, totalizing, statistics, 2nd tare memory, identifiers		
Built-in interface	RS-232C Format: 7-bit ASCII, 1 start bit, 1 or 2 stop bits Parity: Space, odd, even Transmission rates: 150 to 19,200 baud Handshake: Software or hardware		

Specifications of the Individual Models:

Model		FC06BBE-SXCE	FC6CCE-HXCE	FC2CCE-S0CE
Type		isoTEST in conjunction with BD BF		
Accuracy class *		Ⓔ	Ⓔ	Ⓔ
Scale interval d*	g	0.001	0.01	0.01
Max. weighing capacity*	g	620	6,200	2,200
Verification scale interval e*	g	0.01	0.1	0.1
Min. capacity*	g	0.05	0.5	0.5
Max. overload capacity	kg	3	25	10
Tare range (subtractive)		≤100% of the max. weighing capacity		
Electronically compensated preload (without restricting weighing range)	g	93	–	110
Max. preload when starting calibration/adjustment (scale must be zeroed)	g	110	5,200	1,300
Application range according to CD*	g	0.02–620	0.5–6,200	0.5–2,200
Response time (average)	s	1.5	1.5	1.5
Allowable ambient operating temperature		+10° to +30°C (50° to 86°F)		
Pan/load plate size	mm	Ø 130	218×200	218×200
Dimensions (W×D×H)	mm	240×294×86	240×294×86	240×294×86
Net weight. approx.	kg	7	8.4	7.3
Dust and water protection rating according to EN 60529		IP 54 (protected against harmful dust deposits and splashes of water)		

* CD = Council Directive 90/384/EEC on non-automatic weighing instruments used within the European Economic Area

Model-Specific Specifications:

Model		FC12CCE-SXCE	FC6CCE-SXCE	FC12CCE-IXCE
Type		isoTEST in conjunction with BD BF		
Accuracy class *		Ⓔ	Ⓔ	Ⓔ
Scale interval d*	g	0.1	0.1	0.5
Max. weighing capacity*	g	12,000	6,200	12,000
Verification scale interval e*	g	1	1	0.5
Min. capacity*	g	5	5	25
Max. capacity	kg	50	50	50
Tare range (subtractive)		≤100% of the max. weighing capacity		
Electronically compensated preload (without restricting weighing range)	g	1,200	1,240	1,200
Max. preload when starting calibration/adjustment (scale must be zeroed)	g	8,200	2,440	8,200
Application range according to CD*	g	5–12,000	5–6,200	25–12,000
Response time (average)	s	1	1	1
Allowable ambient operating temperature		+10° to +30°C (50° to 86°F)		
Pan size	mm	218×200	218×200	218×200
Net weight, approx.	kg	7.3	7.3	7.3
Dust and water protection rating according to EN 60529		IP54		

Model		FC34EDE-HXCE, FC16EDE-HXCE, FC34EDE-PXCE, FC12EDE-PXCE, FC64EDE-SXCE, FCA34EDE-HXCE FCA16EDE-HXCE FCA34EDE-PXCE FCA12EDE-PXCE FCA64EDE-SXCE				
Type		isoTEST in conjunction with BF BF				
Accuracy class *		Ⓔ	Ⓔ	Ⓔ	Ⓔ	Ⓔ
Scale interval d*	g	0.1	0.1	0.1/0.2/0.5	0.1/0.2	1
Max. weighing capacity*	kg	34	16	8/16/34	6/12	64
Verification scale interval e*	g	1	1	1	1	10
Min. capacity*	g	5	5	5	5	50
Tare range (subtractive)		≤100% of the max. weighing capacity				
Max. capacity	kg	130				
Electronically compensated preload (without restricting weighing range)	kg	4	4	4	4	13
Max. preload when starting calibration/adjustment (scale does not have to be zeroed)	kg	ca. 21	ca. 19	ca. 21	ca. 10	ca. 45
Application range according to CD*	g	5–34,000	5–16,000	5–34,000	5–12,000	50–64,000
Response time (average)	s	1.5				
Allowable ambient operating temperature		+10° to +30°C (50° to 86°F)				
Pan size	mm	300×400				
Net weight, approx.	kg	16.0				
Dust and water protection rating according to EN 60529		IP65				

* CD = Council Directive 90/384/EEC for non-automatic weighing instruments used within the European Economic Area


Specifications

Model		FCA150IGG-HXCE	FCA300IGG-HXCE
Type		isoTEST in conjunction with HC BF	
Accuracy class *		Ⓔ	Ⓔ
Scale interval d*	g	1	20
Max. weighing capacity*	kg	150	300
Verification scale interval e*	g	10	20
Min. capacity*	g	50	1000
Tare range (subtractive)		≤100% of the max. weighing capacity	
Max. overload capacity	kg	600	600
Tare range (subtractive)	kg	150	300
Electronically compensated preload (without restricting weighing range)	kg	30	60
Max. preload when starting calibration/adjustment (scale display does not have to be zeroed)	kg	approx. 120	approx. 300
Application range according to CD*	kg	0.05-150	1-300
Response time (average)	s	1.5	1.5
Allowable ambient operating temperature		0 to +40°C (273 to 313°K; +32° to 104°F)	
External standard calibration weight (of at least accuracy class...)	kg	50 (F2)	50 (F1)
Platform dimensions	mm	800×600	
Net weight, approx.	kg	70	
Dust and water protection rating according to EN 60529		IP67	


* CD = Council Directive 90/384/EEC for non-automatic weighing instruments used within the European Economic Area

Dimensions (Scale Drawings)

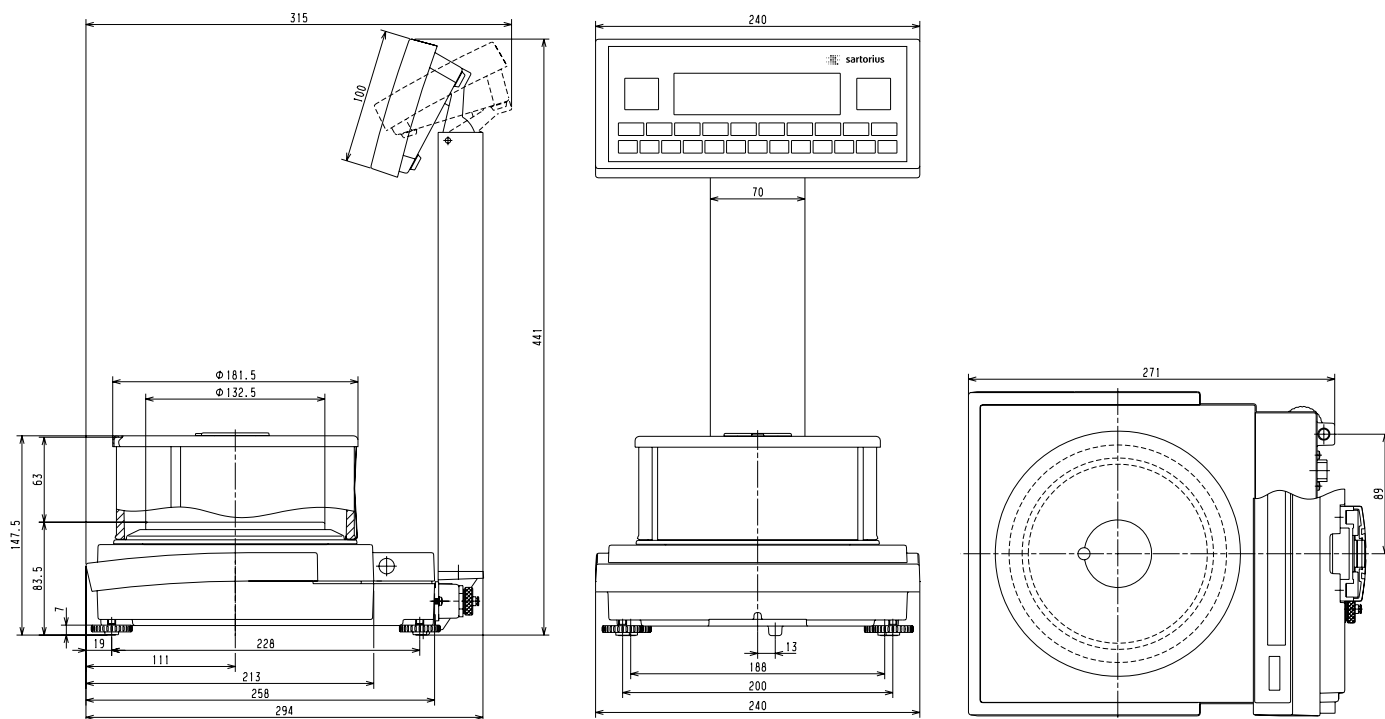FC06BBE-SX (CE)




FC6CCE-HX (CE), FC2CCE-SX (CE)

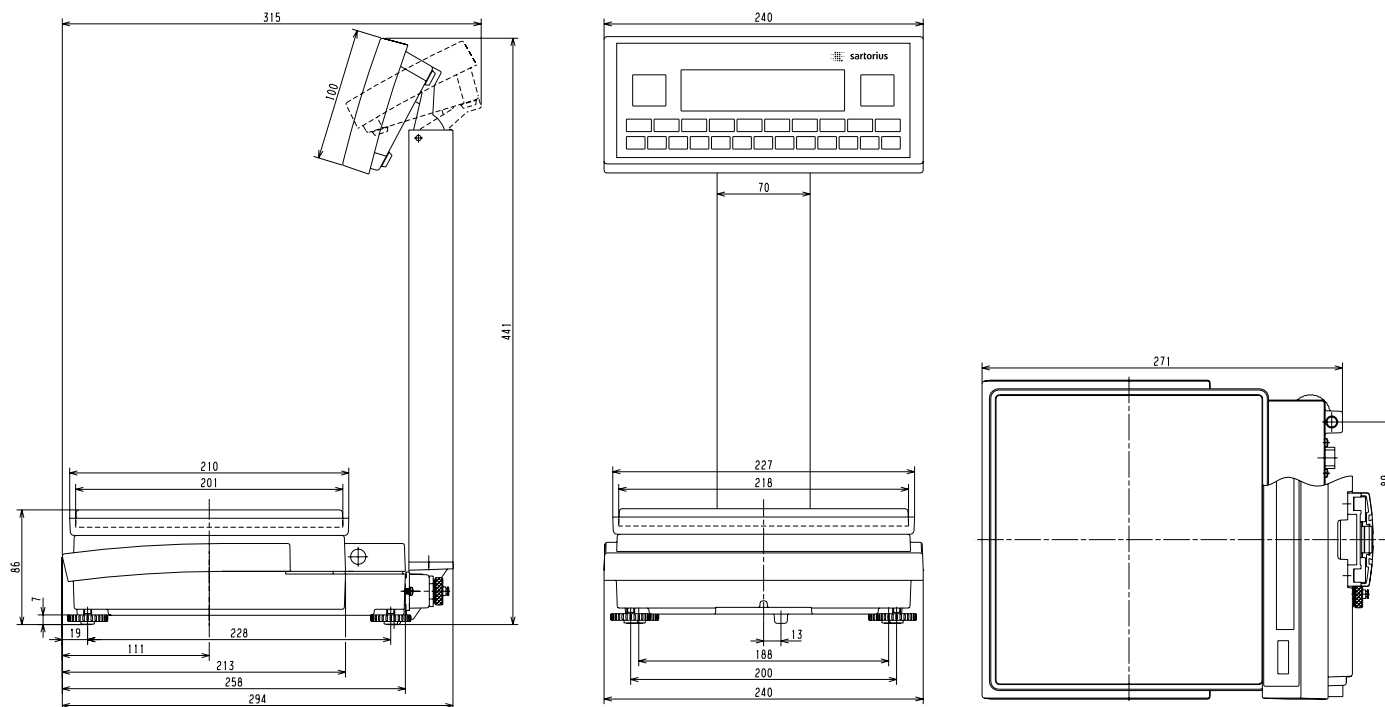


All dimensions are given in millimeters

Dimensions (Scale Drawings)

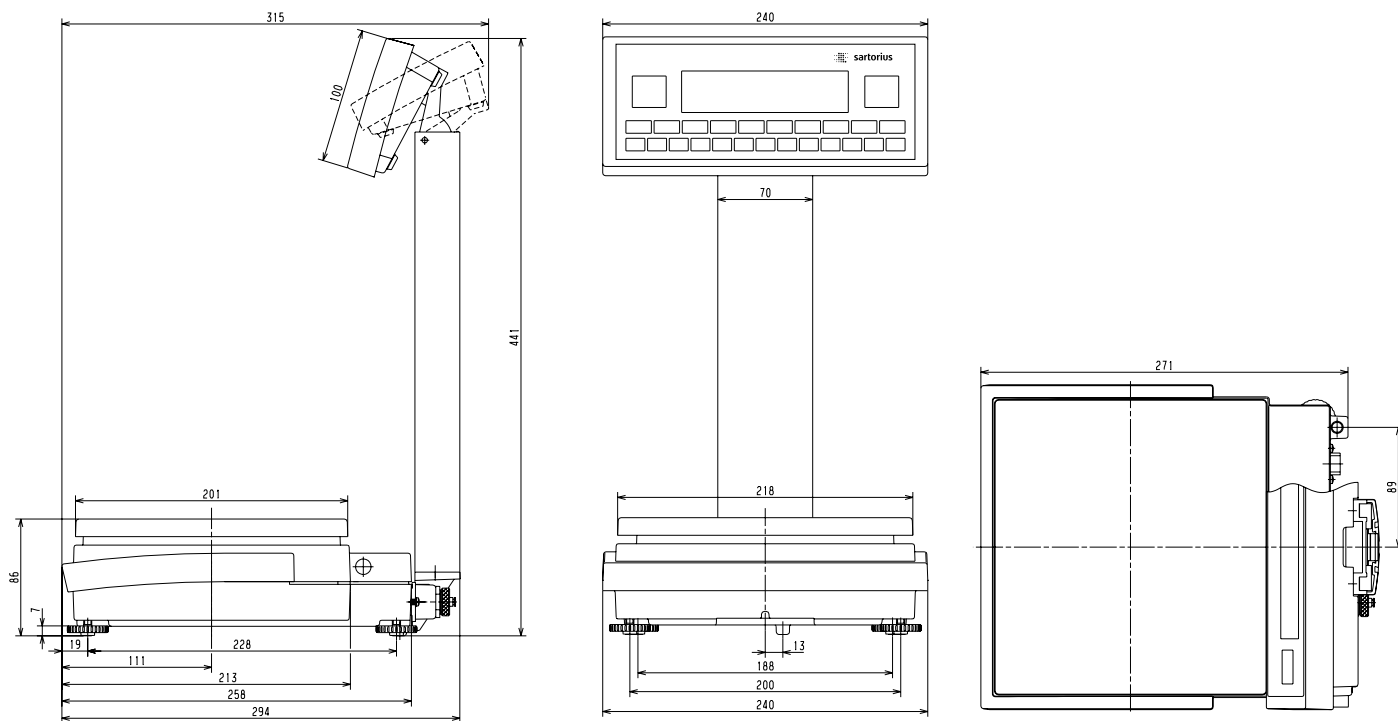FC12CCE-SX (CE), FC12CCE-IXCE, FC6CCE-SX (CE)




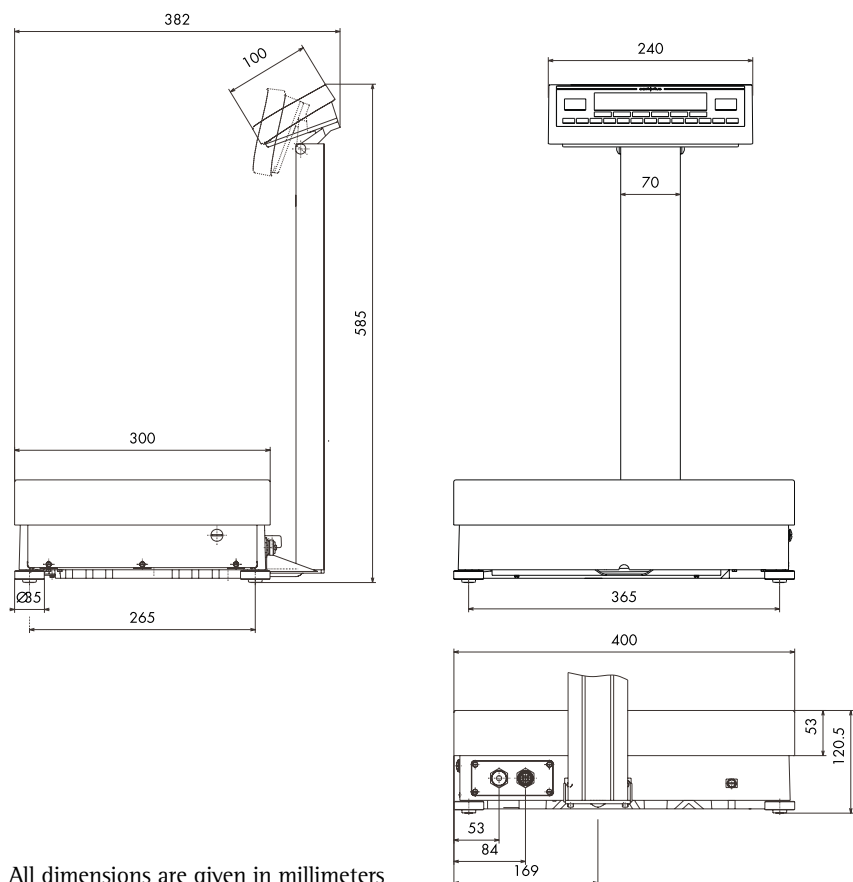FC64EDE-HX(CE), FC34EDE-HX(CE), FC16EDE-HX(CE), FC34EDE-PX(CE), FC12EDE-PX(CE),FC64EDE-SX(CE)




All dimensions are given in millimeters

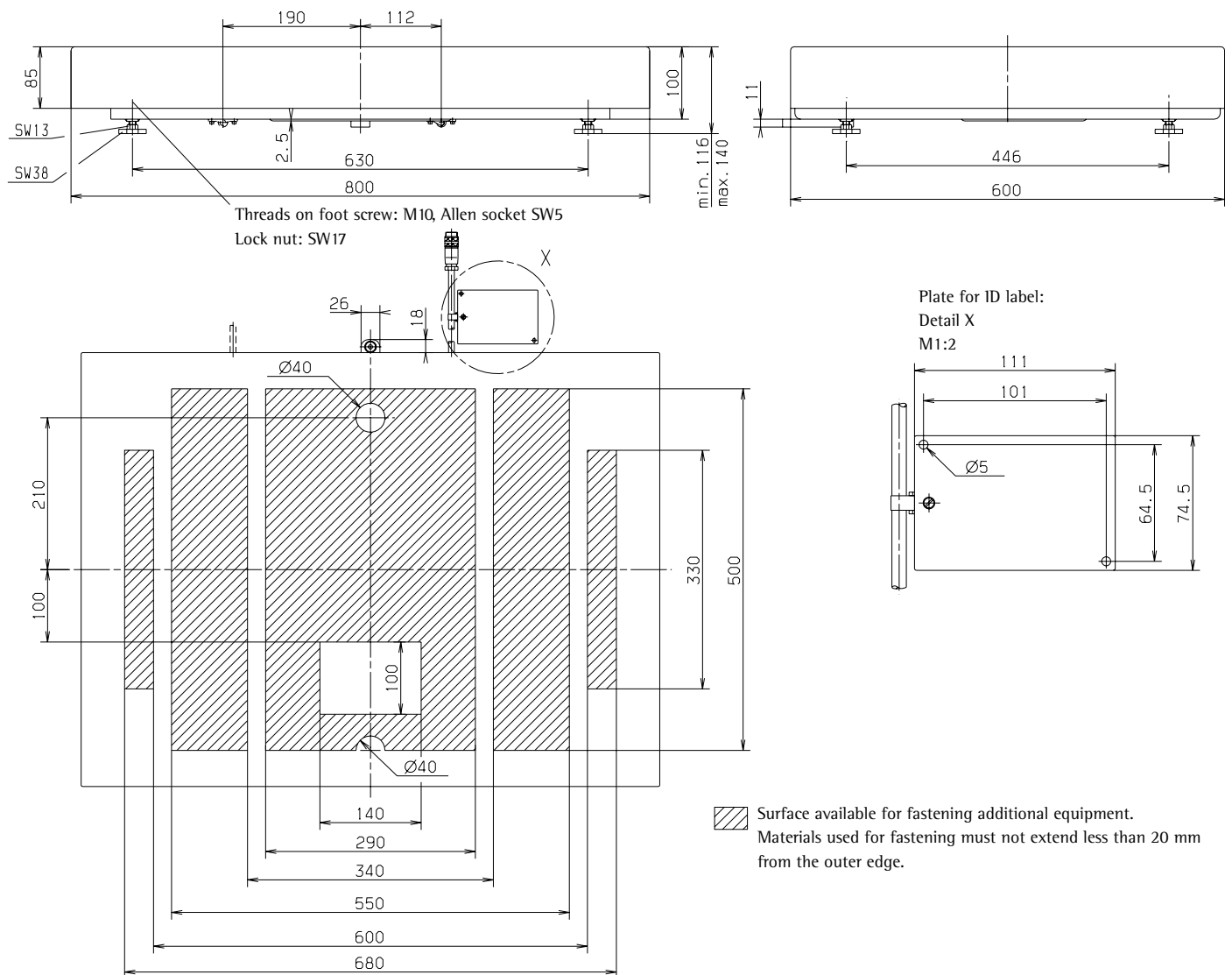Platform on FCA...IGG Models




Platform on FCA...EDE Models

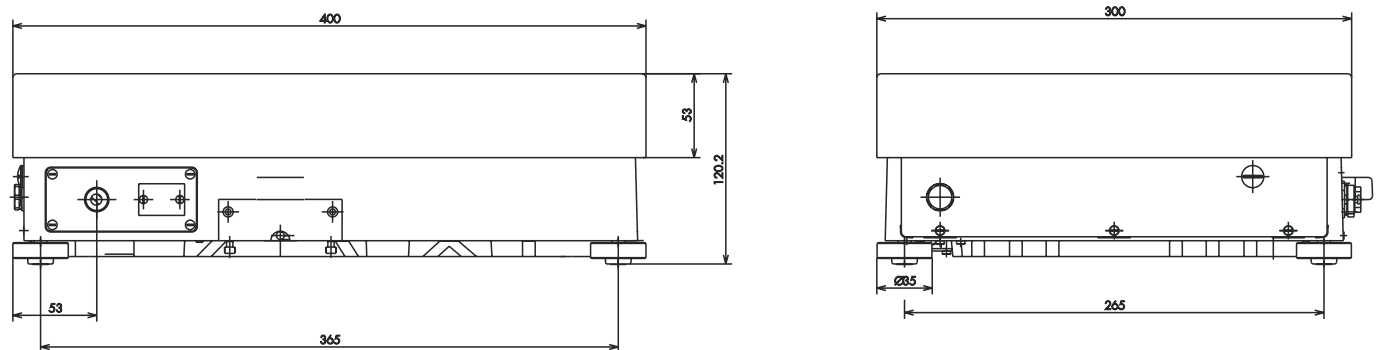

All dimensions are given in millimeters


Dimensions (Scale Drawings)

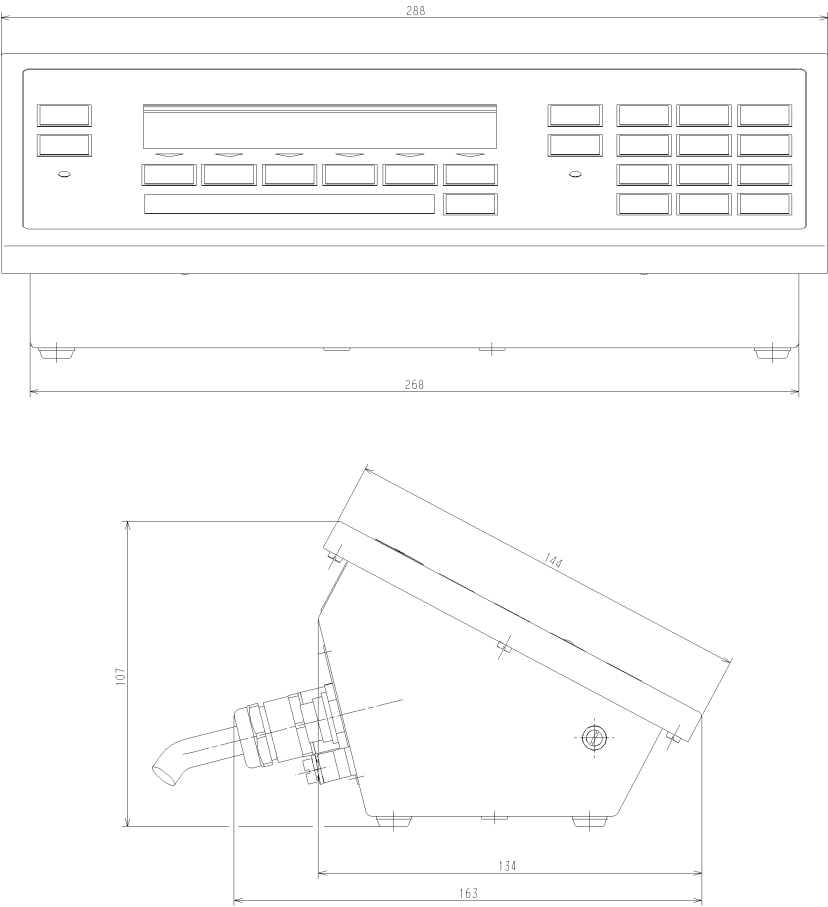FCA Display and Control Unit



Control Box Installation:
Make sure you use screws of the required length. Install additional panels or other materials as necessary to ensure that the housing is impenetrable in accordance with the IP protection rating.

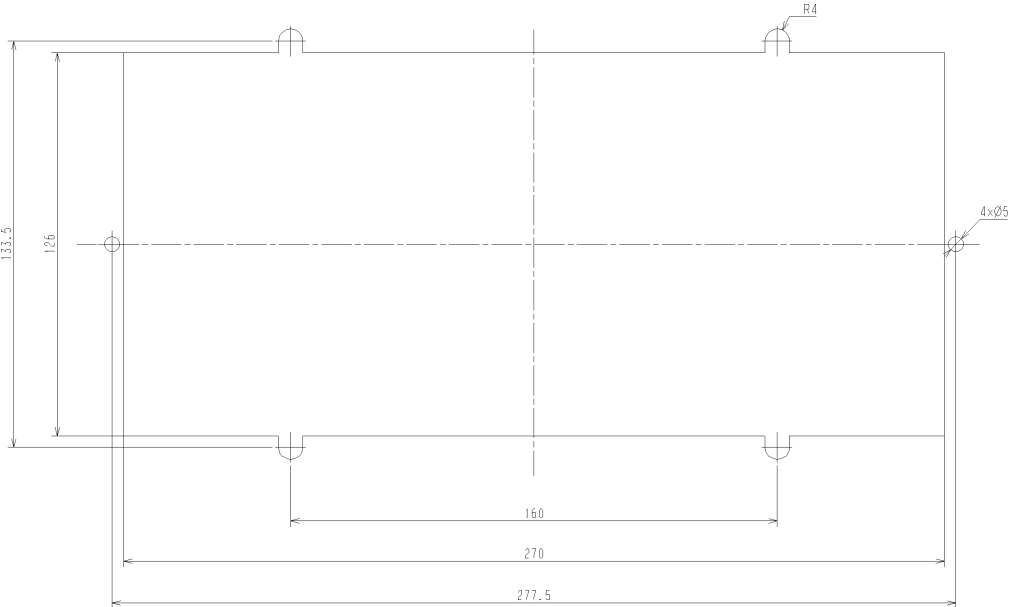Installation dimensions:



All dimensions are given in millimeters

Accessories (Options)



Product

Order No.


Power supplies

AC adapters for use in safe area (110 to 240 V)

EU excluding UK
UK
USA

YPS02-ZDR
YPS02-ZGR
YPS02-ZKR

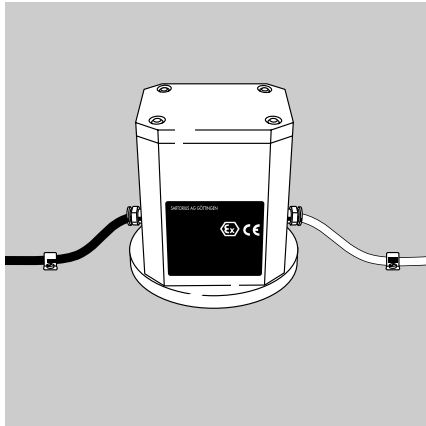- ⚠ Install the cable so that it is protected from damage.
Connect the terminal of the housing to a terminal for equipotential bonding (PE = protective earthing or grounding conductor)



Power supplies for use in hazardous area (110 to 240 V)

EU except UK
UK
USA

YPS02-XDR
YPS02-XGR
YPS02-XKR

- ⚠ Install the cable so that it is protected from damage.
Connect the terminal of the housing to a terminal for equipotential bonding (PE = protective earthing or grounding conductor)

External rechargeable battery pack, for use in Zone 1, 21 or 20 hazardous areas


YRB02-X

Battery charger for YRB02-X battery pack (for use in safe area)

YRB02LD

24V industrial power supply module for use in safe or hazardous area, switching cabinet installation: 14"

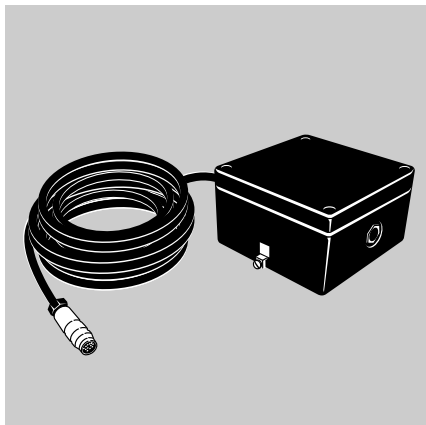YPS02-XV24



Electrical accessories

Zener barrier (2 data lines) for connecting up to 8 weighing instruments (model YD002F-X output port (RS 485) required); cabling supplied in accordance with customer requirements

YDI01-Z

Zener barrier (4 data lines) for connecting accessories installed in the safe area (e.g., printer, PC); cable for connecting Zener barrier to weighing instrument: 20 m



YDI02-Z

Zener barrier (4 data lines, 4 control lines) for connecting accessories installed in the safe area (e.g., PLC); cable length: 20 m


YDI03-Z

- ⚠ Install the cable so that it is protected from damage.
Connect the terminal of the housing to a suitable ground electrode (PE = protective earthing or grounding conductor).

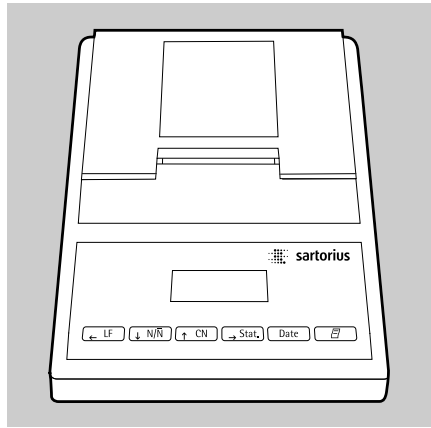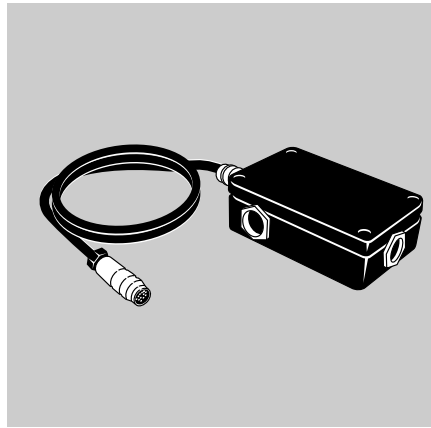Accessories (Options)



Product	Order No.
T-connector (1 input, 2 outputs), 14-pin round male connector; cable length: 1 m	YTE02-X
⚠ Install the cable so that it is protected from damage.	
Foot switch for hazardous areas/locations (Choice of functions to emulate keys on the display and control unit: [E], [Tare], [Q], [isoTEST], [CF]	YPE05-X
Printer (Data printer (not for use in hazardous areas/locations) for hard copies with the date, time, statistical evaluation, transaction counter; with an LCD (can be connected using: adapter cable YCC01-0016M3 or directly, without the printer cable included, using the interface cable YCC01-0019M3) Can be used in legal metrology. Please note: For verified scales with e ≠ d, the last digit on the printout must be specially identified as shown on the scale display. Please set this parameter in the printer menu (see instructions manual for the printer)	YDP03-OCE

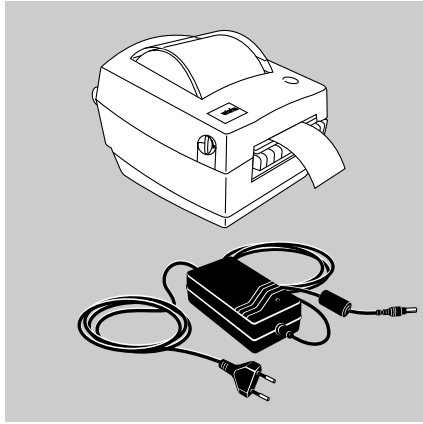


Product**Order No.**

Verifiable strip and label printer (for FC models only)
with “thermo-direct” print head; paper width: up to 108 mm;
with 100–240 V external power supply and power cord
(EU and USA); only for use with flexible printout configuration
software; adapter cable 69Y03142 required

YDP12IS-OCEUV

Verifiable strip and label printer (for FC models only)
with “thermo-transfer” print head; paper width: up to 108 mm;
with 100–240 V external power supply and power cord
(EU and USA); only for use with flexible printout configuration
software; adapter cable 69Y03142 required

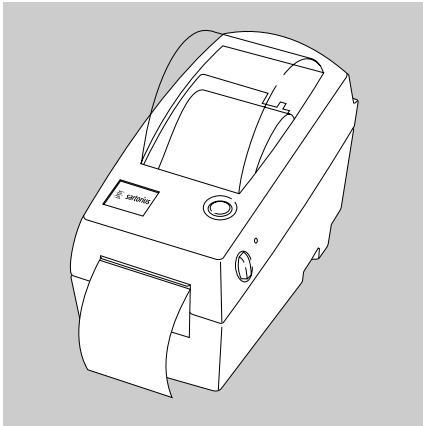YDP12IS-OCEUVTH



Verifiable strip and label printer
with “thermo” print head; paper width: up to 60 mm;
with 100–240 V external power supply;
adapter cable 69Y03142 required

YDP04IS-OCEUV

Accessories (Options)

Product	Order No.
Data interface (must be installed at the Sartorius Fast Factory or by a Sartorius service technician) for FC...BBE/CCE models: TTY/10 mA for FC...BBE/CCE models: RS-485 for FC...EDE models RS-485 for FC...EDE models TTY/10 mA for FCA models: TTY/10 mA for FCA models: RS-485	YD001F-X YD002F-X YD003FC-X YD004FC-X YD001FT-X YD002FT-X
Sartorius serial Profibus interface for RS-232 for the installation in safe areas behind Zener barrier	YSPI-232
Sartorius serial Profibus interface for RS-485 for the installation in safe areas behind Zener barrier	YSPI3-485
Cable for connecting weighing platform to separate display and control unit (length: 2.70 m ~ 9 ft.) – for FC...EDE models – for models with a weighing capacity ≤ 12 kg	Ask for information YCC01-19M3
Tiltable retainer – for installing the FC...EDE or FCA display and control unit on the front of the scale or on the wall	YDH01F
Wrap-around load plate – for models FC6CCE-SX, FC12CCE-SX, FC12CCE-I X	YLP01
Hook for under-scale weighing for FC-EDE models	69EA0040
Calibration weights for all FC/FCA scales, extensive assortment, optionally available with officially recognized DKD certificate	Ask for information
Extension cable 12-pin round male connector/12-pin round female connector (6 m; ~ 20 ft.)	YCC01-01ISM6
Interface cable for connecting a PC to the YDI03-Z, YDI01-Z Zener barrier	YCC01-03ISM5
Interface cable for directly interfacing the scale to the YDP03-OCE data printer (via Zener barrier)	YCC01-0019M3
Adapter cable From round male connector to RS232-D-Submini female interface connector; for directly connecting Sartorius accessories to the scale (via Zener barrier)	YCC01-0016M3
Adapter cable, from 25-pin D-Submini male connector to 9-contact D-Submini female connector; length 0.25 m (~ .8 ft.)	6965619
Adapter plate for fastening the FCA display and control unit to the floor stand	YAS01FCT-X
Floor stand (stainless steel), height: 1.1 m (approx. 3.5 ft) – for FCA display and control units (adapter plate YAS01FCT-X required)	YDH03IS
Retainer for mounting the YRB02-X external rechargeable battery pack on a YDH03IS/YDH02CIS column or on the wall	YBPH01

Product	Order No.
SartoConnect data transfer software for connecting your Sartorius scale to a PC running the Windows 95, 98 or NT operating system. For loading measurement data from the scale into any PC for further processing (e.g., with Excel or Access) Includes 9/12-pin cable for connecting the scale to the PC.	YSC01I
Sartorius Win Scale Scale driver software for use under Windows 95/98/2000/NT. Display scale readout on the PC monitor and provides secure memory for storing data that is subject to legal control.	YSW03
Configuration software (PC, DOS) For storing and loading scale settings	YAD01IS
Sartorius "NICE Label Express" For creating customer-specific printouts	YAD02IS

Declarations of Conformity

The CE Mark on Sartorius Equipment

In 1985, the Council of the European Community approved a resolution concerning a new approach to the technical harmonization and standardization of national regulations. The organization for monitoring compliance with the directives and standards concerning the CE marking is governed in the individual EU Member States through the implementation of the EC Directives adopted by the respective national laws. As of December 1993, the scope of validity for all EC Directives has been extended to the Member States of the European Union and the Signatories of the Agreement on the European Economic Area.

Sartorius complies with the EC Directives and European Standards in order to supply its customers with weighing instruments that feature the latest advanced technology and provide many years of trouble-free service.

The CE mark may be affixed only to weighing instruments and associated equipment that comply with the applicable Directive(s):

Council Directive 89/336/EEC "Electromagnetic compatibility (EMC)"

1. Electromagnetic compatibility:

1.1 Reference to 89/336/EEC:
Official Journal of the European Communities, Nr. 2001/C105/03

EN 61326-1
Electrical equipment for measurement control and laboratory use EMC requirements

Parts 1: General requirements
Defined immunity to interference:
Industrial areas, continuous, un-monitored operation
Limitation of emissions:
Residential areas Class B

Important Note:
The operator shall be responsible for any modifications to Sartorius equipment and for any connections of cables or equipment not supplied by Sartorius and must check and, if necessary, correct these modifications and connections. On request, Sartorius will provide information on the minimum operating specifications (in accordance with the Standards listed above for defined immunity to interference).

Council Directive 73/23/EEC "Electrical Equipment Designed for Use within Certain Voltage Limits"

Applicable European Standards:

EN 60950
Safety of information technology equipment including electrical business equipment

EN 61010
Safety requirements for electrical equipment for measurement, control and laboratory use
Part 1: General requirements

If you use electrical equipment in installations and under ambient conditions requiring higher safety standards, you must comply with the provisions as specified in the applicable regulations for installation in your country.

Weighing Instruments for Use in Legal Metrology: Directive 90/384/EEC "Non-automatic weighing instruments"

This Directive regulates the determination of mass in legal metrology.

For the respective Declaration of Conformity for weighing instruments that have been verified by Sartorius for use as legal measuring instruments and that have an EC Type-Approval Certificate, see the page after next.

This Directive also regulates the performance of the EC verification by the manufacturer, provided that an EC Type-Approval Certificate has been issued and the manufacturer has been accredited by an officer or a Notified Body registered at the Commission of the European Community for performing such verification.

The legal basis allowing Sartorius to perform EC verification is constituted by the EC Council Directive No. 90/384/EEC on non-automatic weighing instruments that has been in effect since January 1, 1993, in the Internal Market as well as by the Certificate of Accreditation of the Sartorius AG Quality Management System issued by the Metrology Department of the Regional Administration Office of Lower Saxony, Germany ("Niedersächsisches Landesverwaltungsamt -Eichwesen") on February 15, 1993.

For information on the CE mark on Sartorius equipment and legal regulations currently applicable in your country, and to obtain the names of the persons to contact, please ask your local Sartorius office, dealer or service center.

"New Installation" Service

Initial verification is covered in our "New Installation" service package. In addition to initial verification, this package provides you with a series of important services which will guarantee that you achieve optimal results with your weighing instrument:

- Installation
- Startup
- Inspection
- Training
- Initial verification

If you would like Sartorius to perform initial verification of your weighing instrument, contact an authorized service representative.

"EC Verification" – A Service Offered by Sartorius

Our service technicians are authorized to perform verification* of your weighing instruments that are acceptable for legal metrological verification and can inspect and verify the metrological specifications at the place of installation within the Member States of the European Union and the Signatories of the Agreement on the European Economic Area.

Subsequent Verifications within the European Countries

The expiration date of the verification depends on the national regulations of the country in which the weighing instrument is used. For information on verification and legal regulations currently applicable in your country, and to obtain the names of the persons to contact, please contact your local Sartorius office, dealer or service center.

* in accordance with the accreditation certificate issued to Sartorius AG


Declaration of Type Conformity to Directive No. 90/384/EEC

This declaration is valid for non-automatic electromechanical weighing instruments for use in legal metrology. These weighing instruments accepted for legal metrological verification have an EC Type-Approval Certificate. The model(s) concerned is(are) listed below along with the respective type, accuracy class, and number of the EC Type-Approval Certificate:

Model	Type	Accuracy Class	EC Type Approval No.	In Conjunction with Test Certificate	
				Type	Certificate No.
FB/FC.....-OCE	iso-TEST	II	D97-09-018	BD BF	D09-96.30
FB/FC.....-OCE	iso-TEST	II	D97-09-018	BB BD	D09-95.08
FBG/FCG.....-OCE	iso-TEST	II	D97-09-018	BF BF	D09-96.30
FC.....-XCE	iso-TEST	II	D97-09-018	MA BF	D09-96.30
FC.....-XCE	iso-TEST	II	D97-09-018	BA BF	D09-96.30
FC.....-XCE	iso-TEST	II	D97-09-018	BF BF	D09-96.30
FB/FC.....-XCE	iso-TEST	II	D97-09-018	BD BF	D09-96.30
FCA/FCB.....-XCE	iso-TEST	II	D97-09-018	BF BF	D09-96.30
FCA...-HXCE	iso-TEST	II	D97-09-018	HC BF	D09-96.30

SARTORIUS AG declares that its weighing instrument types comply with the requirements of the Council Directive on non-automatic weighing instruments, no. 90/384/EEC of 20 June 1990; the associated European Standard "Metrological aspects of non-automatic weighing instruments," No. EN 45501; the amended, currently valid versions of the national laws and decrees concerning legal metrology and verification in the Member States of the European Union, the EU, and the Signatories of the Agreement on the European Economic Area, which have adopted this Council Directive into their national laws; and with the requirements stipulated on the Type-Approval Certificate for verification. This Declaration of Type Conformity is valid only if the ID label on the weighing instrument has the CE mark of conformity and the green metrology

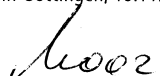Sartorius AG
37070 Goettingen, Germany
Signed in Göttingen, 19.11.2002



Dr. G. Maaz
(Head of Technical Operations)

sticker with the stamped letter "M" (the two-digit number in large print stands for the year in which the mark has been affixed):



If these marks are not on the ID label, this Declaration of Type Conformity is not valid. Validity can be obtained, for example, by submitting the weighing instrument for final action to be taken by an authorized representative of SARTORIUS AG. The period of validity of this Declaration of Type Conformity shall expire upon any tampering with, repair or modification of this weighing instrument or, in some Member States, on the date of expiration.

The operator of this weighing instrument shall be responsible for obtaining an authorized renewal of the verification, such as subsequent or periodic verification, of the weighing instrument for use as a legal measuring instrument.


J. Rehwald
(Head of Quality Management & Services)



(1) **EC-TYPE EXAMINATION CERTIFICATE**

(2) Equipment or protective system intended for use in potentially explosive atmospheres – Directive 94/9/EC

(3) EC-Type Examination Certificate Number: **KEMA 01ATEX1099 X**

(4) Equipment or protective system:

**Weighing Module series, Type IS.....-X.... and
Weighing System series, Type FC.....-X.... and
Weighing System series, Type FCA.....-X.... and
Weighing System series, Type FCB.....-X....**

(5) Manufacturer: **Sartorius A.G.**

(6) Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

(7) This equipment or protective system and any acceptable variation thereto is specified in the schedule to this certificate and the documents therein referred to.

(8) KEMA Quality B.V., notified body number 0344 in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994, certifies that this equipment or protective system has been found to comply with the Essential Health and Safety Requirements relating to the design and construction of equipment and protective systems intended for use in potentially explosive atmospheres given in Annex II to the Directive.

The examination and test results are recorded in confidential report no. 2010745.

(9) Compliance with the Essential Health and Safety Requirements has been assured by compliance with:

EN 50014 : 1997 – EN 50020 : 1994

(10) If the sign "X" is placed after the certificate number, it indicates that the equipment or protective system is subject to special conditions for safe use specified in the schedule to this certificate.

(11) This EC-Type Examination Certificate relates only to the design, examination and tests of the specified equipment or protective system in accordance with the Directive 94/9/EC. Further requirements of the Directive apply to the manufacturing process and supply of this equipment or protective system. These are not covered by this certificate.

(12) The marking of the equipment or protective system shall include the following:

II 2 G EEx ib IIC T4

Arnhem, 18 October 2001
KEMA Quality B.V.

T. Pijker
Certification Manager

* This Certificate may only be reproduced in its entirety and without any change

KEMA Quality B.V.
Utrechtseweg 310, 6812 AR Arnhem, The Netherlands
P.O. Box 5185, 6802 ED Arnhem, The Netherlands
Telephone +31 26 3 56 20 08, Telefax +31 26 3 52 58 00

ACCREDITED BY THE
DUTCH COUNCIL FOR
ACCREDITATION



(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 01ATEX1099 X

(15) Description

The Weighing Modules series, Type IS....BBE-.X...., IS....CCE-.X...., IS....EDE-.X.... and IS....IGG-.X.... serve to measure weights and have connections for supply and a remote terminal for data handling.

Depending of the type the connections to the remote terminal and supply are via plugs or an integral cable, in the latter case e.g. via a junction box Type YAS06IS-X.

The Weighing Modules differ in size and construction, depending of the weight they can measure.

The Weighing Systems series, Type FC.....-X...., FCA.....-X.... and FCB.....-X.... serve to measure weights, handle the data and provide communication with other units (e.g. a PC). The Weighing Systems consist of a Weighing Module of one of the above mentioned types and a terminal with a display and a keyboard.

Type FC.....-X.... has an integral terminal, for Type FCA.....-X.... and FCB.....-X.... the terminal is remote mounted.

All types are supplied by four independent intrinsically safe supplies.

All types have the same supply and data communication circuits.

All types have external connections for passive equipment (e.g. a foot switch).

Type FCB.....-X.... has additional I/O circuits for connection of the outputs to associated circuits via barriers and input circuits to passive circuits (e.g. switches).

Ambient temperature range -20 °C ... +40 °C.

Electrical data

All types

Unless otherwise specified, the connection-identifications are shown below per type.

Supply circuits in type of explosion protection intrinsic safety EEx ib IIC, only for connection to a certified intrinsically safe circuit, with following maximum values (for each circuit):

Circuit	U_i	I_i	P_i	C_i	L_i
V-1	12,6 V	133 mA	1,68 W	0 nF	0 mH
V-2	12,6 V	133 mA	1,68 W	0 nF	0 mH
V-3	8,6 V	187 mA	1,61 W	300 nF	0 mH
V-4	12,6 V	150 mA	1,89 W	100 nF	0 mH

Data communication circuit in type of explosion protection intrinsic safety EEx ib IIC, with the following maximum values (per circuit unless otherwise specified):

$$U_o = 8,6 \text{ V}$$

$$I_o = 23 \text{ mA}$$

$$P_o = 50 \text{ mW}$$

The maximum allowed external capacitance $C_o = 6 \mu\text{F}$,
the maximum allowed external inductance $L_o = 60 \text{ mH}$.

(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 01ATEX1099 X

Electrical data (continued)

Circuits only for connection to a certified intrinsically safe circuit, with following maximum values (for each circuit per applicable connection combination):

Circuit	U _i [V]	I _i [mA]	P _i [mW]	C _i [nF]	L _i [mH]	Connections	Note
RS232	12,6	85	270	3	0	A/J/K/N/M (BU9) or 9/10/11/12/8 (BU6)	1)
RS232	12,6	10	30	4	0	C/D/E/F/G/M (BU9) or 2/3/4/5/6/8 (BU6)	1)
RS485	12,6	85	270	110	0	J/K/L/M (BU9) or 3/4 (BU30) + 9/10 (BU6)	1)
RS485	12,6	10	30	4	0	C/D/E/F/G/M (BU9) or 6/11/8 (BU6)	1)
TTY	14,7	50	265	0	0	G/K/D/F/J (BU9) or 3/7/2/4/6 (BU30)	1)
TTY	14,7	130	100	0	0	C/E/D/F/J (BU9) or 1/5/2/4/6 (BU30)	1) 2)
I/O-out	30	120	900	0	0	1/2, 3/4, 5/6, 7/8 (ST20)	1) 3)

Notes:

- 1) The current I_i must be resistively limited
- 2) May alternatively be connected to converter Type 725324, certified per Certificate of Conformity PTB No. Ex-85.B.2089.
For this combination the following data applies:
Maximum allowed external capacitance C_o = 620 nF and maximum allowed external inductance L_o = 2 mH.
- 3) For Type FCB.....-X..... only.

Maximum values of the output data (for each circuit per applicable connection combination):

Circuit	U _i [V]	I _i [mA]	P _i [mW]	C _i [μF]	L _i [mH]	Connections	Note
RS232	12,6	28	88	1,15	50	B/O/M (BU9) or 1/7/8 (BU6)	1)
RS485	12,6	28	85	1,15	50	B/O/M (BU9) or 7/8 (BU6)	1)
TTY	12,6	28	85	1,15	50	7/8 (BU6) or 7/8 (BU6)	1) 2)
I/O-in	12,6	30	95	1,15	50	9/10, 9/11, 9/12 (ST20)	3)

Notes:

- 1) Only for use in combination with passive intrinsically safe circuits
- 2) Same connector used with or without junction box.
- 3) For Type FCB.....-X..... only.

(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 01ATEX1099 X

Electrical data (continued)

Connections per type

Weighing Modules series,

Type IS...BBE-X..., IS...CCE-X..., IS...EDE-X... and IS...IGG-X...

Connections:

- Supply circuits: ST8 pins 6/1 (V-1), 4/2 (V-2), 5/3 (V-3), 7/8 (V-4),
in combination with termination box: BU2 pins 6/1 (V-1), 4/2 (V-2), 5/3 (V-3), 7/8 (V-4)
- Data communication circuits: BU9 (pins as specified above),
in combination with termination box: BU6 + BU30 (pins as specified above)
- When termination box used:
Maximum cable length between Weighing Module and termination box Type YAS06IS-X
is 15 m.

Weighing System series Type FC.....-X....

Connections:

- Supply circuits: BU2 pins 6/1 (V-1), 4/2 (V-2), 5/3 (V-3), 7/8 (V-4)
- Data communication circuits: BU9 (pins as specified above)

Weighing System series Type FCA.....-X....

Connections:

- Supply circuits: ST8 pins 6/1 (V-1), 4/2 (V-2), 5/3 (V-3), 7/8 (V-4)
- Data and communication circuits: BU9 (pins as specified above)
- Maximum cable length between Weighing Module and Terminal is 15 m

Weighing System series Type FCB.....-X....

Connections:

- Supply circuits: ST8 6/1 (V-1), 4/2 (V-2), 5/3 (V-3), 7/8 (V-4)
- Data and communication circuits: BU9
- Maximum cable length between Weighing Module and Terminal is 15 m

(16)

Report

KEMA No. 2010745

(17)

Special conditions for safe use

1. The output circuits of the I/O of Type FCB.....-X.... may only be connected to shunt diode zener barriers with one side grounded. This ground and the external ground connection of the enclosure of the terminal of Type FCB.....-X.... must be connected to the potential equalising system within the hazardous area.
2. For ambient temperature range and electrical data, see (15).

(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 01ATEX1099 X

(18) **Essential Health and Safety Requirements**

Essential Health and Safety Requirements not covered by the standards listed at (9)	
Clause	Subject
1.0.5	Marking
1.0.6 b) and d)	Instructions

These Essential Health and Safety Requirements are examined and positively judged.
The results are laid down in the report listed at (16)

(19) **Test documentation**

1. EC-Type Examination Certificate KEMA 98ATEX0612 X
EC-Type Examination Certificate KEMA 00ATEX1012 X

dated

2. Description (34 pages) 02.02.2001
3. Drawing No. 35520-000-60-A4 (3 sheets) 20.02.2001
4. Samples

AMENDMENT 1

to EC-Type Examination Certificate KEMA 01ATEX1099 X

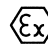

Manufacturer: **Sartorius A.G.**

Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

Description

In future the Weighing Modules series, Types IS....EDE-.X...., IS....IGG-.X...., the Weighing Systems series, Types FCA.....-X...., FCB.....-X.... and the junction box Type YAS06IS-X are alternatively also constructed in accordance with EN 50281-1-1:1998 and the documentation as listed below.

The marking for these versions becomes

 II 2 G EEx ib IIC T4 and  II 1 D T 135 °C

The maximum surface temperature T 135 °C of the enclosure is based on an ambient temperature of 40 °C.

The electrical data of the Weighing Modules series, Types IS....BBE-.X...., IS....CCE-.X...., IS....EDE-.X...., IS....IGG-.X.... and the Weighing Systems series Types FC.....-X...., FCA.....-X...., FCB.....-X...., is amended as specified below.

Electrical data

Data output circuit

for units with RS485 circuit in type of explosion protection intrinsic safety EEx ib IIC, (BU9)

with the following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 85 \text{ mA} \\P_o &= 270 \text{ mW}\end{aligned}$$

The maximum allowed external capacitance $C_o = 1,15 \mu\text{F}$, the maximum allowed external inductance $L_o = 5 \text{ mH}$.

Only for connection to a certified intrinsically safe circuit, with following maximum values.

Connections J/K:

$$\begin{aligned}U_i &= 12,6 \text{ V (each line to earth)} \\I_i &= 810 \text{ mA (resistive limited)} \\P_i &= 2,5 \text{ W}\end{aligned}$$

The effective internal capacitance $C_i = 10 \text{ nF}$, the effective internal inductance L_i is negligibly small.

Connections D/E/F/G & M (values for each circuit):

$$\begin{aligned}U_i &= 12,6 \text{ V} \\I_i &= 10 \text{ mA (resistive limited)} \\P_i &= 30 \text{ mW}\end{aligned}$$

The effective internal capacitance $C_i = 4 \text{ nF}$, the effective internal inductance is negligibly small.

Connections B/C/O/L are only for use in combination with passive intrinsically safe circuits, the effective internal capacitance $C_i = 4 \text{ nF}$, the effective internal inductance L_i is negligibly small.



AMENDMENT 1

to EC-Type Examination Certificate KEMA 01ATEX1099 X

Electrical data (continued)

All other data remain unchanged.

Special conditions for safe use

For applications of the Weighing Modules series, Types IS....EDE-.X...., IS....IGG-.X...., the Weighing Systems series, Types FCA.....-X...., FCB.....-X.... and the junction box Type YAS06IS-X in explosive atmospheres caused by air/dust mixtures, the dust layer may not exceed a thickness of 5 mm.

All other data remain unchanged

Test documentation

	<u>dated</u>
1. Description	16.11.2001
2. Drawing No.	
65638-000-39-A4	10.10.2000
65651-000-55-A3	25.09.2001
65638-000-55-A3 (2 sheets)	11.09.2001
65638-000-55-A3	11.09.2001
98100-036-30 Rev.00	04.10.2000
98100-036-31 Rev.00	04.10.2000
65638-700-90-A4 Rev.00	06.10.2000
65552-000-05-A4	02.02.2001
35520-700-30-A3	14.02.2001
35632-210-02-A4	28.05.2001
98100-031-25-A4 Rev.01	19.09.1995

Arnhem, 8 May 2002
KEMA Quality N.V.

T. Pijpker
Certification Manager

[2018412/2]

Page 2/2



(1) **EC-TYPE EXAMINATION CERTIFICATE**

(2) Equipment or protective system intended for use in potentially explosive atmospheres
- Directive 94/9/EC

(3) EC-Type Examination Certificate Number: **KEMA 98ATEX0892 X**

(4) Equipment or protective system: **Power Supply Type YPS02-X..**

(5) Manufacturer: **Sartorius AG**

(6) Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

(7) This equipment or protective system and any acceptable variation thereto is specified in the schedule to this certificate and the documents therein referred to.

(8) KEMA, notified body number 0344 in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994, certifies that this equipment or protective system has been found to comply with the Essential Health and Safety Requirements relating to the design and construction of equipment and protective systems intended for use in potentially explosive atmospheres given in Annex II to the Directive.

The examination and test results are recorded in confidential report no. 80892.

(9) Compliance with the Essential Health and Safety Requirements has been assured by compliance with:

EN 50014 : 1992 + prA1 EN 50018 : 1994 EN 50020 : 1994

(10) If the sign "X" is placed after the certificate number, it indicates that the equipment or protective system is subject to special conditions for safe use specified in the schedule to this certificate.

(11) This EC-Type Examination Certificate relates only to the design and construction of the specified equipment or protective system. If applicable, further requirements of this Directive apply to the manufacture and supply of this equipment or protective system.

(12) The marking of the equipment or protective system shall include the following:

II 2 (2) G EEx d [ib] IIC T4

Arnhem, 15 June 1998
by order of the Board of Directors of N.V. KEMA

C.M. Boschloo
Certification Manager

° This Certificate may only be reproduced in its entirety and without any change





SCHEDULE

(13)

(14)

to EC-Type Examination Certificate KEMA 98ATEX0892 X

(15) **Description**

The Power Supply Type YPS02-X.. provides four independent intrinsically safe output channels for intrinsic safe scales. The power supply is mounted in an enclosure in type of protection flameproof enclosure "d". The maximum length of the interconnection cable type LiYC-Y-CY 4x0.5 between the Power Supply and the scales is 50 m.

Ambient temperature range 0 °C ... +40 °C.

Electrical data

Main supply non intrinsically safe circuit, suitable for connection to electrical equipment with working voltages up to 264 V.

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_1, white) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 133 \text{ mA} \\P_o &= 1,68 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_2, brown) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 133 \text{ mA} \\P_o &= 1,68 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_3, green) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 8,6 \text{ V} \\I_o &= 187 \text{ mA} \\P_o &= 1,61 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 4 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_4, yellow) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 150 \text{ mA} \\P_o &= 1,89 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$



SCHEDULE

(13)

(14)

to EC-Type Examination Certificate KEMA 98ATEX0892 X

(15) **Electrical data** (continued)

Cable type LiYC-Y-CY 4x0.5

Maximum cable capacitance $C_{\text{cable}} = 28,2 \text{ nF}$

Maximum cable inductance $L_{\text{cable}} = 8,4 \mu\text{H}$

Maximum allowed cable length = 50 m

The intrinsically safe circuits are infallible galvanically isolated from the non-intrinsically safe circuits up to a sum of peak voltages of 375 V.

Installation instruction

The Power Supply provided with certified cables entries Type ADE can only be used for fixed installation. Certified cable entries Type AGRO can be used for flexible- and for fixed installation.

Routine tests

The transformer shall, before mounting into the apparatus, withstand per Clause 8.1.5 of EN 50 020 - 1994 without breakdown the application of 2500 V between the primary and secondary winding.

Routine tests according to Clause 16 of EN 50018 are not required since the type test has been made at a static pressure of four times the reference pressure.

(16) **Report**

KEMA No. 80892

(17) **Special conditions for safe use**

None

(18) **Essential Health and Safety Requirements**

Essential Health and Safety Requirements not covered by standards listed at (9)	
Clause	Subject
1.0.6	Instructions

These Essential Health and Safety Requirements are examined and positively judged. The results are laid down in the report listed at (16).



SCHEDULE

(13)

(14)

to EC-Type Examination Certificate KEMA 98ATEX0892 X

(19) Test documentation

1. EC-Type Examination Certificate PTB 98ATEX1023 U
KEMA 98ATEX0611 X
Certificate of Conformity ISSeP 92.C.103.997
LCIE 92.C6125 X

signed

2. Description 65473-000-06-A4,
Rev. 00 (11 pages) 05.05.1998

3. Drawing No. 65473-700-06-A4)
65473-000-33-A3) 05.05.1998
65473-700-05-A3)

65473-000-05-A2 02.02.1998

65473-120-90-A3)
65473-120-95-A3) 09.01.1998
65473-120-01)

4. Samples

AMENDMENT 1


to EC-Type Examination Certificate KEMA 98ATEX0892 X

Manufacturer: **Sartorius AG**
Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

Description

In future the Power Supply Type YPS02-X.. is alternatively also constructed in accordance with EN 50281-1-1:1998 and the documentation listed below.

The marking for this version becomes

 II 2 G EEx d[ib] IIC T4 and  1 D T 135 °C

The maximum surface temperature T 135 °C of the enclosure is based on an ambient temperature of 40 °C.

Electrical Data

The value of the maximum power of each supply and input circuit changes as follows:

Terminal V_1, white : P_o = 1,46 W
Terminal V_2, brown : P_o = 1,46 W
Terminal V_3, green : P_o = 1,51 W
Terminal V_4, yellow : P_o = 1,68 W

All other data remain unchanged

Installation instructions

1. For use in potentially explosive atmospheres caused by the presence of flammable gases, fluids or vapours: The cable entry devices shall be of a certified flameproof type, suitable for the conditions of use and correctly installed.
2. For use in potentially explosive atmospheres caused by the presence of combustible dust: The cable entry devices shall be of a certified flameproof type, suitable for the conditions of use and correctly installed. The minimum ingress protection requirement of IP6X according to EN 60529 must be satisfied.
3. With the use of conduit entries a sealing device shall be provided either in the flameproof enclosure or immediately on the entrance thereto. The sealing device shall be of a certified flameproof type, suitable for the conditions of use and correctly installed.

Special conditions for safe use

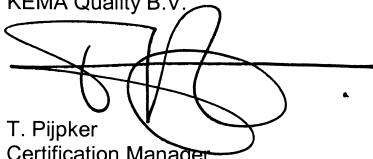
1. For applications in explosive atmospheres caused by air/dust mixtures, the dust layer shall not exceed a thickness of 5 mm.
2. For ambient temperature range and electrical data see (15) and above.

Test documentation

dated

1. Drawing No. 65473-000-36-A4 10.10.2000

Arnhem, 11 April 2002
KEMA Quality B.V.



T. Pijpker
Certification Manager

[2018412]



(1) **EC-TYPE EXAMINATION CERTIFICATE**

(2) Equipment or protective system intended for use in potentially explosive atmospheres
- Directive 94/9/EC

(3) EC-Type Examination Certificate Number: **KEMA 98ATEX0611 X**

(4) Equipment or protective system: **Power Supply Type YPS02-Z..**

(5) Manufacturer: **Sartorius AG**

(6) Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

(7) This equipment or protective system and any acceptable variation thereto is specified in the schedule to this certificate and the documents therein referred to.

(8) KEMA, notified body number 0344 in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994, certifies that this equipment or protective system has been found to comply with the Essential Health and Safety Requirements relating to the design and construction of equipment and protective systems intended for use in potentially explosive atmospheres given in Annex II to the Directive.

The examination and test results are recorded in confidential report no. 80611.

(9) Compliance with the Essential Health and Safety Requirements has been assured by compliance with:

EN 50014 : 1992 + prA1

EN 50020 : 1994

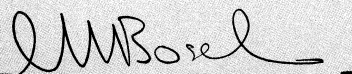
(10) If the sign "X" is placed after the certificate number, it indicates that the equipment or protective system is subject to special conditions for safe use specified in the schedule to this certificate.

(11) This EC-Type Examination Certificate relates only to the design and construction of the specified equipment or protective system. If applicable, further requirements of this Directive apply to the manufacture and supply of this equipment or protective system.

(12) The marking of the equipment or protective system shall include the following:

 II (2) G [EEEx ib] IIC

Arnhem, 18 May 1998
by order of the Board of Directors of N.V. KEMA



C.M. Boschloo
Certification Manager

• This Certificate may only be reproduced in its entirety and without any change





SCHEDULE

(13)

(14)

to EC-Type Examination Certificate KEMA 98ATEX0611 X

(15) **Description**

The Power Supply Type YPS02-Z.. provides four independent intrinsically safe output channels for intrinsically safe scales. The Power Supply Type YPS02-Z.. is mounted in an enclosure providing a degree of ingress protection of IP 54. The maximum length of the interconnection cable type LiYC-Y-CY 4 * 0.5 between the Power Supply and the scales is 50 m.

Ambient temperature range 0 °C ... +40 °C.

Electrical data

Main supply non intrinsically safe circuit, suitable for connection to electrical equipment with working voltages up to 264 V.

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_1, white) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 133 \text{ mA} \\P_o &= 1,68 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_2, brown) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 133 \text{ mA} \\P_o &= 1,68 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_3, green) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 8,6 \text{ V} \\I_o &= 187 \text{ mA} \\P_o &= 1,61 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 4 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

Supply and input circuit in type of explosion protection intrinsic safety
terminal (V_4, yellow) EEx ib IIC with following maximum values:

$$\begin{aligned}U_o &= 12,6 \text{ V} \\I_o &= 150 \text{ mA} \\P_o &= 1,89 \text{ W}\end{aligned}$$

Maximum allowed external capacitance $C_o = 1 \mu\text{F}$
Maximum allowed external inductance $L_o = 300 \mu\text{H}$

SCHEDULE

(13)

(14) to EC-Type Examination Certificate KEMA 98ATEX0611 X

(15) **Electrical data** (continued)

cable type LiYC-Y-CY 4 * 0.5 Maximum cable capacitance $C_{\text{cable}} = 28,2 \text{ nF}$
 Maximum cable inductance $L_{\text{cable}} = 8,4 \mu\text{H}$
 Maximum allowed cable length = 50 m

The intrinsically safe circuits are infallible galvanically isolated from the non-intrinsically safe circuits up to a sum of peak voltages of 375 V.

Installation instructions

The Power Supply Type YPS02-Z.. must be installed outside the hazardous area.

After installation of the Power Supply, all applicable separations shall meet the requirements per Clause 6.4 of EN 50 020 - 1994.

Inside the hazardous area, the power supply must be mounted in an enclosure which is suitable for this purpose. This combination shall be separately investigated and certified.

Routine test

The transformer shall, before mounting into the apparatus, withstand per Clause 8.1.5 of EN 50 020 - 1994 without breakdown the application of 2500 V between the primary and secondary winding.

(16) **Report**

KEMA No. 80611

(17) **Special conditions for safe use**

None

(18) **Essential Health and Safety Requirements**

Essential Health and Safety Requirements not covered by standards listed at (9)	
Clause	Subject
1.0.6 b	Instructions

These Essential Health and Safety Requirements are examined and positively judged. The results are laid down in the report listed at (16).



SCHEDULE

(13)

(14)

to EC-Type Examination Certificate KEMA 98ATEX0611 X

(19) **Test documentation**

1. Certificate of Conformity KEMA No. Ex-97.D.1279 X

signed

2. Product Compliance Report ANNEX II,
65463-700-70-A4, Rev. 00

24.03.1998

3. Drawing No. 65463-000-34-A3

03.03.1998

4. Samples



AMENDMENT 1

to EC-Type Examination Certificate KEMA 98ATEX0611 X

Manufacturer: **Sartorius AG**

Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

Description

In future the Power Supply Type YPS02-Z.. may also be constructed in accordance with the documentation listed below.

The marking becomes Ex II (2) G (1) D [EEx ib] IIC

Item (18) is amended with

Essential Health and Safety Requirements not covered by the standards listed at (9)	
Clause	Subject
2.2.2.2. and 2.2.2.4	Explosive atmospheres caused by air/dust mixtures

These Essential Health and Safety Requirements are examined and positively judged.

Electrical Data

The value of the maximum power of each supply and input circuit changes as follows:

Terminal V_1, white : $P_o = 1,46 \text{ W}$
Terminal V_2, brown : $P_o = 1,46 \text{ W}$
Terminal V_3, green : $P_o = 1,51 \text{ W}$
Terminal V_4, yellow : $P_o = 1,68 \text{ W}$

All other data remain unchanged

Test documentation

dated

1. Drawing No. 65463-000-38-A4 10.10.2000

Arnhem, 11 April 2002
KEMA Quality B.V.

T. Pijpker
Certification Manager

[2018412]



(1) **EC-TYPE EXAMINATION CERTIFICATE**

(2) Equipment or protective system intended for use in potentially explosive atmospheres – Directive 94/9/EC

(3) EC-Type Examination Certificate Number: **KEMA 03ATEX2164 X**

(4) Equipment or protective system: **Power Supply type YPS02-XV24**

(5) Manufacturer: **Sartorius AG**

(6) Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

(7) This equipment or protective system and any acceptable variation thereto is specified in the schedule to this certificate and the documents therein referred to.

(8) KEMA Quality B.V., notified body number 0344 in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994, certifies that this equipment or protective system has been found to comply with the Essential Health and Safety Requirements relating to the design and construction of equipment and protective systems intended for use in potentially explosive atmospheres given in Annex II to the Directive.

The examination and test results are recorded in confidential report no. 2029029.

(9) Compliance with the Essential Health and Safety Requirements has been assured by compliance with:

EN 50014 : 1997 EN 50020 : 2002 EN 50028 : 1987 EN 50281-1-1 : 1998

(10) If the sign "X" is placed after the certificate number, it indicates that the equipment or protective system is subject to special conditions for **safe** use specified in the schedule to this certificate.

(11) This EC-Type Examination Certificate relates only to the design, examination and tests of the specified equipment or protective system in accordance with the Directive 94/9/EC. Further requirements of the Directive apply to the manufacturing process and supply of this equipment or protective system. These are not covered by this certificate.

(12) The marking of the equipment or protective system shall include the following:

II 2G 1D EEx m[ib] IIC T4 T 80 °C

Arnhem, 2 December 2003
KEMA Quality B.V.

C.G. van Es
Certification Manager

* This Certificate may only be reproduced in its entirety and without any change

KEMA Quality B.V.
Utrechtseweg 310, 6812 AR Arnhem, The Netherlands
P.O. Box 5185, 6802 ED Arnhem, The Netherlands
Telephone +31 26 3 56 20 08, Telefax +31 26 3 52 58 00

ACCREDITED BY THE
DUTCH COUNCIL FOR
ACCREDITATION



(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 03ATEX2164 X

(15) Description

The Power Supply type YPS02-XV24 contains four separate intrinsically safe supply circuits. The main supply and the intrinsically safe circuits are provided with a permanently connected cable. The length of the cable for the intrinsically safe circuits is maximum 100 m.

Ambient temperature range -20 °C ... +40 °C.

The maximum surface temperature T 80 °C of the enclosure is based on an ambient temperature of 40 °C.

For the exclusion of dust, the unit has a degree of ingress protection of at least IP6x in accordance with EN 60529.

Electrical data

Main supply circuit 24 Vdc nominal, 30 V maximum, 6,5 W typical
Breaking capacity of fuse is 1500 A
 $U_m = 253 \text{ Vac}$

Supply circuit V1 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 1 and 6) with the following maximum values:

U_o	=	12,6	V
I_o	=	133	mA
P_o	=	1,43	W
C_o	=	0,5	μF
L_o	=	0,3	mH

Supply circuit V2 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 2 and 4) with the following maximum values:

U_o	=	12,6	V
I_o	=	133	mA
P_o	=	1,43	W
C_o	=	0,5	μF
L_o	=	0,3	mH

Supply circuit V3 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 3 and 5) with the following maximum values:

U_o	=	8,6	V
I_o	=	186	mA
P_o	=	1,49	W
C_o	=	0,8	μF
L_o	=	0,2	mH

Supply circuit V4 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 7 and 8) with the following maximum values:

U_o	=	12,6	V
I_o	=	150	mA
P_o	=	1,65	W
C_o	=	0,5	μF
L_o	=	0,3	mH

(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 03ATEX2164 X

Electrical data (continued)

The intrinsically safe supply circuits of Power Supply type YPS02-XV24 are also suitable for connection to the applicable supply circuits (V1 ... V4) of the following certified equipment:

CIXS3	approved per KEMA 03ATEX1157 X
CW.XS3-.....	approved per KEMA 03ATEX1157 X
FC.....-X....	approved per KEMA 01ATEX1099 X
FCA.....-X....	approved per KEMA 01ATEX1099 X
FCB.....-X....	approved per KEMA 01ATEX1099 X
FCT01-X	approved per KEMA 00ATEX1012 X
FCT01-XV1	approved per amendment 1 to KEMA 00ATEX1012 X
IS.....-X....	approved per KEMA 01ATEX1099 X
EB... -IX....	approved per KEMA 98ATEX1722 X

Routine test

The encapsulated parts of the unit shall be subjected to the routine tests in accordance with EN 50028, Clause 7.1 and 7.3.

(16)

Report

KEMA No. 2029029.

(17)

Special conditions for safe use

1. When the cable of the main supply is connected within the hazardous area a certified junction box suitable for the application and environment shall be used. For instance an enclosure and terminals in type of protection increased safety "e". For applications in explosive atmospheres caused by air/dust mixtures the termination box shall have an ingress protection of at least IP6x per EN 60529.
2. The enclosure shall be connected to the potential equalizing system, either direct or via the screen of the permanently connected cable. When the connection is made via the screen, it shall be assured that equalization currents stay within the current rating of the screen.
3. The connection cable without or with one outer screen shall be installed in such a way that mechanical damage is excluded.
4. For applications in explosive atmospheres caused by air/dust mixtures and where category 1D apparatus is required, the dust layer may not exceed a thickness of 5 mm.
5. For ambient temperature range and electrical data, see (15).

(18)

Essential Health and Safety Requirements

Covered by the standards listed at (9).

(19)

Test documentation

dated

Documentation list 65684-000-45-A4

09.09.2003



(1) **EC-TYPE EXAMINATION CERTIFICATE**

(2) Equipment or protective system intended for use in potentially explosive atmospheres – Directive 94/9/EC

(3) EC-Type Examination Certificate Number: **KEMA 03ATEX2137 X**

(4) Equipment or protective system: **Ex rechargeable battery pack type YRB02-X**

(5) Manufacturer: **Sartorius AG**

(6) Address: **Weender Landstraße 94-108, 37075 Göttingen, Germany**

(7) This equipment or protective system and any acceptable variation thereto is specified in the schedule to this certificate and the documents therein referred to.

(8) KEMA Quality B.V., notified body number 0344 in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994, certifies that this equipment or protective system has been found to comply with the Essential Health and Safety Requirements relating to the design and construction of equipment and protective systems intended for use in potentially explosive atmospheres given in Annex II to the Directive.

The examination and test results are recorded in confidential report no. 2021992.

(9) Compliance with the Essential Health and Safety Requirements has been assured by compliance with:

**EN 50014 : 1997 EN 50019 : 2000 EN 50020 : 2002
EN 50028 : 1987 EN 50281-1-1 : 1998**

(10) If the sign "X" is placed after the certificate number, it indicates that the equipment or protective system is subject to special conditions for safe use specified in the schedule to this certificate.

(11) This EC-Type Examination Certificate relates only to the design, examination and tests of the specified equipment or protective system in accordance with the Directive 94/9/EC. Further requirements of the Directive apply to the manufacturing process and supply of this equipment or protective system. These are not covered by this certificate.

(12) The marking of the equipment or protective system shall include the following:

II 2G 1D EEx em[ib] IIC T4 T 80 °C

Arnhem, 15 August 2003
KEMA Quality B.V.

C.G. van Es
Certification Manager

* This Certificate may only be reproduced in its entirety and without any change





(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 03ATEX2137 X

(15) Description

The Ex rechargeable battery pack type YRB02-X contains four separate intrinsically safe supply circuits, powered by two integral rechargeable batteries (each 6,5 Ah or 10 Ah). The Ex rechargeable battery pack is provided with a permanently connected cable of maximum 50 m length.

Ambient temperature range -20 °C ... +40 °C.

The maximum surface temperature T 80 °C of the enclosure is based on an ambient temperature of 40 °C.

For the exclusion of dust, the applicable constructions of the Ex rechargeable battery pack have a degree of ingress protection of at least IP6x in accordance with EN 60529.

Electrical data

Supply circuit V1 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 1 and 6) with the following maximum values:

U_o	=	12,6	V
I_o	=	133	mA
P_o	=	1,43	W
C_o	=	0,5	μF
L_o	=	0,3	mH

Supply circuit V2 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 2 and 4) with the following maximum values:

U_o	=	12,6	V
I_o	=	133	mA
P_o	=	1,43	W
C_o	=	0,5	μF
L_o	=	0,3	mH

Supply circuit V3 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 3 and 5) with the following maximum values:

U_o	=	8,6	V
I_o	=	186	mA
P_o	=	1,49	W
C_o	=	0,8	μF
L_o	=	0,2	mH

Supply circuit V4 in type of explosion protection intrinsic safety EEx ib IIC,
(connections 7 and 8) with the following maximum values:

U_o	=	12,6	V
I_o	=	150	mA
P_o	=	1,65	W
C_o	=	0,5	μF
L_o	=	0,3	mH



(13)

SCHEDULE

(14)

to EC-Type Examination Certificate KEMA 03ATEX2137 X

Electrical data (continued)

The intrinsically safe supply circuits of the Ex rechargeable battery pack may be connected to the applicable supply circuits (V1 ... V4) of the following certified equipment:

CIXS3	approved per KEMA 03ATEX1157 X
CW.XS3.....	approved per KEMA 03ATEX1157 X
FC.....-X....	approved per KEMA 01ATEX1099 X
FCA.....-X....	approved per KEMA 01ATEX1099 X
FCB.....-X....	approved per KEMA 01ATEX1099 X
FCT01-X	approved per KEMA 00ATEX1012 X
FCT01-XV1	approved per amendment 1 to KEMA 00ATEX1012 X
IS.....-X....	approved per KEMA 01ATEX1099 X
EB.....-IX....	approved per KEMA 98ATEX1722 X

Routine test

The routine tests in accordance with EN 50028, Clause 7.1 and 7.3 shall be performed on the encapsulated parts of the unit.

(16) **Report**

KEMA No. 2021992.

(17) **Special conditions for safe use**

1. The batteries shall only be charged outside the potentially explosive area. The precautions as described in the instruction manual shall be taken into account as well.
2. The Ex rechargeable battery pack may only be used as moveable or stationary equipment.
3. The enclosure shall be connected to the potential equalizing system, either direct or via the screen of the permanently connected cable.
4. For applications in explosive atmospheres caused by air/dust mixtures and where category 1D apparatus is required, the dust layer may not exceed a thickness of 5 mm.
5. For the ambient temperature range and electrical data, see (15).

(18) **Essential Health and Safety Requirements**

Covered by the standards listed at (9).

(19) **Test documentation**

dated

Documentation list 65656-000-45-A4

06.05.2003



Declaration of Conformity **to Council Directives 89/336/EEC , 73/23/EEC** **and 94/9/EEC**

The Ex rechargeable battery pack type YRB02-X

meets the requirements of the test standards listed below, in conjunction with the associated auxiliary peripheral devices and installation equipment listed in Annex A2 (see Annex A1 for a technical description).

1. Electromagnetic Compatibility
 - 1.1 Source for 89/336/EEC: EC Official Journal, No. 2001/C105/03
EN 61326-1 Electrical equipment for measurement, control and laboratory use
EMC requirements
Part 1: General requirements
Limitation of emissions: Residential areas, Class B
Defined immunity to interference: Industrial areas, continuous unmonitored operation
2. Safety of Electrical Equipment
 - 2.1 Source for 73/23/EEC: EC Official Journal, No. 2001/C106/03
EN 61010 Safety requirements for electrical equipment for measurement, control and laboratory use
Part 1: General requirements
EN 60950 Safety of information technology equipment
3. Equipment or protective systems or components intended for use in potentially explosive atmospheres and for use in presence of combustible dust
 - 3.1 Source for 94/9/EEC: EC Official Journal, No. 2001/C51/03
EN50014 General requirements
EN50019 Increased safety "e"
EN50020 Intrinsic safety "i"
EN50028 Encapsulation "m"
EN50281-1-1 Electrical apparatus protected by enclosures – Construction and testing
 - 3.2 Type Examination: KEMA/Arnheim (NL) (Notified Body, Reg. No. 0344):
KEMA 03ATEX2137 X
 - 3.3 Production Quality Assessment Notification:
Certified by PTB/Braunschweig (Notified Body, Reg. No. 0102)
Notification No.: PTB 97ATEX Q021-1

Sartorius AG
37070 Goettingen, Germany
2003

W. Obermann
(Senior Vice President, R&D
Electronics Engineering
Mechatronics Division)

Dr. K. Klein
(Senior Vice President, R&D
Mechanical Engineering
Mechatronics Division)



Certificate of Compliance

Certificate: 1153350 (LR 56628-34)

Master Contract: 167555

Project: 1268586 (Edition 5)

Date Issued: December 14, 2001

Issued to: Sartorius AG
Weender Landstraße 94-108
37075 Göttingen
Germany

The products listed below are eligible to bear the CSA Mark shown



Issued by:

Andrew Redeker

Andrew Redeker, C.E.T.
Certification Specialist

Authorized by:

John Verwey
John Verwey, P.Eng.
Operations Manager

CLASS

2258 03 - PROCESS CONTROL EQUIPMENT - Intrinsically Safe and Non-Incendive- For Hazardous Locations

PRODUCTS

Part A:

Class I, Div 1, Groups A,B,C, and D; Class II, Div 1, Groups E, F and G; Class III; Ex ib IIC T4

Precision Weighing System, Model FCabc-dXef and Precision Weighing Module, Model ISbc-dXef; Rated Input 15.5 V dc; 800 mA; Temperature Code T4; Intrinsically Safe when connected as per Control Drawing 33956-000-07-A4.

Notes:

- (1) Suffix "a" in Model FC may be "blank", "A" or "B".
- (2) Suffix "c" in the above Models may be
"BBE", "CCE" certified for Class I, Div 1, Groups A,B,C and D; Ex ib IIC T4 only
"EDE", "IGG" certified for Class I, Div 1, Groups A,B,C, and D; Class II, Div 1, Groups E, F and G; Class III;
Ex ib IIC T4
- (3) Options include RS485 communications, 10mA data communication, or any of the following simple apparatus;
Junction Box (YAS06IS-X), T-conector (TYE02-X), foot switch (YPE05-X)
- (4) The suitability of the final installation is to be determined by the local authorities having jurisdiction.



Certificate: 1153350

Master Contract: 167555

Project: 1268586

Date Issued: December 14, 2001

Part B:

Power Supplies, Model YPS02-ZKR, rated at input 100-240 vac, 50/60 Hz, 25VA, output 15.5 V dc; 800 mA. for use in Ordinary (Non-Hazardous) Locations only and Model YPS02-XKR, rated at input 100-240 vac, 50/60 Hz, 25VA, output 15.5 V dc; 800 mA for use in Class I, Groups B, C and D; Class II, Groups E, F and G; Class III Hazardous Locations; provides intrinsically safe outputs for connection to Precision Weighing System, Model FC and Precision Weighing Module, Model IS when connected as per Control Drawing 33956-000-07-A4.

Notes:

- (1) The above model numbers include suffixes designating indicator type, output signal, electrical connection, and options.
- (2) The suitability of the final installation is to be determined by the local authorities having jurisdiction.

APPLICABLE REQUIREMENTS

- | | | |
|-------------------------------|---|--|
| CSA Standard C22.2 No 0-M1991 | - | General Requirements - Canadian Electrical Code Part II. |
| 0.4-M1982 | - | Bonding and Grounding of Electrical Equipment (Protective Grounding). |
| 0.5-M1982 | - | Threaded Conduit Entries |
| 25-M1966 | - | Enclosures for Use in Class II Groups E, F and G Hazardous Locations. |
| 30-M1986 | - | Explosion-proof Enclosures for Use in Class I Hazardous Locations. |
| 142-M1987 | - | Process Control Equipment. |
| 157-M1992 | - | Intrinsically Safe and Non-Incendive Equipment for Use in Hazardous Locations. |

Factory Mutual Research

1151 Boston-Providence Turnpike
P.O. Box 9102 Norwood, MA 02062 USA
T: 781 762 4300 F: 781 762 9375 www.fmglobal.com

CERTIFICATE OF COMPLIANCE

HAZARDOUS (CLASSIFIED) LOCATION ELECTRICAL EQUIPMENT

This certificate is issued for the following equipment:

Isab-cxde Weighing Module series with standard RS232 Data Communication

IS/I,II,III/1/A,B,C,D,E,F,G – 33956-000-07-A4/ T4

I/1/ AEx ib IIC/T4–33956-000-07-A4

Maximum Entity Parameters: See control drawing 33956-000-07-A4

a = Figure for maximum mechanical load (03 thru 300)

b = Weighing unit type EDE or IGG (DIP); BBE or CCE (not DIP)

c = Mechanical sensitivity (H, I, P, S)

d = Calibration/Custody (CE)

Options:

e = Specials, color, applied software (V1, V2), RS485 Data Communication (instead of RS232), TTY Data Communication with 10 mA output (instead of RS232), YAS06IS-X Junction Box Type YAS06IS-X

Fcab-c-dXef Weighing System series with standard RS232 Data Communication

IS/I,II,III/1/A,B,C,D,E,F,G – 33956-000-07-A4/ T4

I/1/ AEx ib IIC/T4 –33956-000-07-A4

Maximum Entity Parameters: See control drawing 33956-000-07-A4

a = Terminal type: FC is with terminal YAC01LA-X00FC

FCA is with terminal YAC01FC-X

FCB is with terminal YAC02FC-X

b = Figure for maximum mechanical load (03 thru 300)

c = Weighing unit type EDE or IGG (DIP); BBE or CCE (not DIP)

d = Mechanical sensitivity (H, I, P, S)

e = Calibration/Custody (CE)

Factory Mutual Research

Options:

f = Specials, color, applied software (V1, V2), RS485 Data Communication (instead of RS232), TTY Data Communication with 10 mA output (instead of RS232), YAS06IS-X Junction Box Type YAS06IS-X

Equipment Ratings:

Intrinsically safe for Class I, Division 1, Groups A, B, C and D; suitable for Class II and III, Groups E, F and G; IS, Class I, Zone 1, AEx ib IIC, T4, hazardous locations, for indoor only use
Approved for:

SARTORIUS AG
WEENDER LANDSTRASSE 94-108
D-37070 GÖTTINGEN
GERMANY

Factory Mutual Research

This certifies that the equipment described has been found to comply with the following
Factory Mutual Research Approval Standards and other documents:

Class 3600	1998
Class 3610	1999

Original Approval Job Identification: 3012724 Approval Granted: April 23, 2002

Subsequent Revision Reports / Date Approval Amended

Factory Mutual Research Corporation



Nicholas Ludlam
Technical Team Manager
Approvals Division

May 9th, 2002
Date

Physikalisch-Technische Bundesanstalt

Braunschweig und Berlin



(1) Production Quality Assessment Notification

(Translation)



- (2) Equipment or protective systems or components intended for use in potentially explosive atmospheres - **Directive 94/9/EC**
- (3) Notification Number: **PTB 97 ATEX Q021-1**
- (4) Product group(s): Balances, load cells and power supply units in the determining types of protection "Intrinsic Safety" and "Flameproof Enclosure"

A list of the EC-Type Examination Certificates covered by this notification is held by the notified body.

- (5) Applicant: Sartorius AG
Weender Landstraße 94-108, D-37075 Göttingen
- (6) Actual manufacturer: Sartorius AG
Weender Landstraße 94-108, D-37075 Göttingen
- (7) The Physikalisch-Technische Bundesanstalt (PTB), notified body No. 0102 for Annex IV in accordance with Article 9 of the Council Directive 94/9/EC of 23 March 1994 notifies to the applicant that the actual manufacturer has a production quality system which complies to the Annex IV of the Directive.
- (8) This notification is based on the confidential audit report No. 01QS003, issued the 2001-01-29. This notification is valid until 2003-12-18 and can be withdrawn if the actual manufacturer no longer satisfies to the requirements of Annex IV.

Results of periodical reassessment of the quality system are a part of this notification.

- (9) According to Article 10 (1) of the Directive 94/9/EC the CE-Marking shall be followed by the identification number 0102 of PTB as the notified body which is involved in the production control stage.

Zertifizierungsstelle Explosionsschutz
By order

Braunschweig, January 29, 2001

Dr.-Ing. U. Johannsmeyer
Regierungsdirektor



Sheet 1/1

Notifications without signature and official stamp shall not be valid. The notification may be circulated only without alteration. Extracts or alterations are subject to approval by the Physikalisch-Technische Bundesanstalt. In case of dispute, the German text shall prevail.

Physikalisch-Technische Bundesanstalt, Bundesallee 100, D-38116 Braunschweig

Physikalisch-Technische Bundesanstalt

Braunschweig und Berlin



EG-Bauartzulassung

EC type-approval certificate

Zulassungsinhaber:

Issued to:

Sartorius AG
Weender Landstraße 94-108
37075 Göttingen
Bundesrepublik Deutschland

Rechtsbezug:

In accordance with:

§ 13 des Gesetzes über das Meß- und Eichwesen (*verification act*)
vom/dated 23. März 1992 (BGBl. I S. 711) in Verbindung mit Richtlinie
(*in connection with council directive*) 90/384/EWG, geändert durch (*amen-*
ded by) 93/68/EWG

Bauart:

In respect of:

Nichtselbsttätige elektromechanische Waage
Nonautomatic electromechanical weighing instrument
Typ/type: iso-TEST
Genauigkeitsklasse/class I, II, III, IIII Max 0,05 kg ... 300 t
Option: Mehrteilungswaage, Mehrbereichswaage
Multi-interval instrument, multiple range instrument

Zulassungsnummer:

Approval number:

D97-09-018 4. Revision

Gültig bis:

Valid until:

2007-06-26

Anzahl der Seiten:

Number of pages:

13

Geschäftszeichen:

Reference No.:

1.14 – 02000518

Benannte Stelle:

Notified Body:

0102

Im Auftrag

By order


Link



Braunschweig, 2002-03-26

Siegel
Seal

Physikalisch-Technische Bundesanstalt

Braunschweig und Berlin



Prüfschein

Test certificate

Ausgestellt für:

Issued to:

Sartorius AG
Weender Landstraße 94 – 108
37075 Göttingen
Bundesrepublik Deutschland

Prüfgrundlage:

In accordance with:

EN 45501 (1992), Nr.8.1, OIML R 76-1 (1992)

Gegenstand:

Object:

Lastaufnehmer mit Wägezelle und Auswerteelektronik mit digitalem Ausgang als Modul einer elektromechanischen Waage zum Anschluß an geeignete Anzeige- und Bedienterminals
Load receptor with load cell and electronic device with digital output as module of an electromechanical weighing instrument for connection to suitable display- and operator-terminals
Typ / type **BA BF, BC BF, BD BF, BF BF, HC BF, MA BF und MD BF**

Kennummer:

Serial number:

Prüfscheinnummer:

Test certificate number:

D09-96.30 7. Revision / *Revision 7*

Datum der Prüfung:

Date of Test:

Anzahl der Seiten:

Number of pages:

12

Geschäftszeichen:

Reference No.:

1.14 – 02001430

Benannte Stelle:

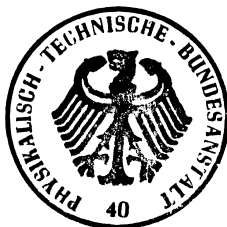
Notified Body:

0102

Im Auftrag

By order

Link



Braunschweig, 2002-11-13

Siegel

Seal

alternative

alternative type BC BF external electronics

Type BF BF only

locked

unlocked

locked

unlocked

below weighing pan type MA BF only

Legend:

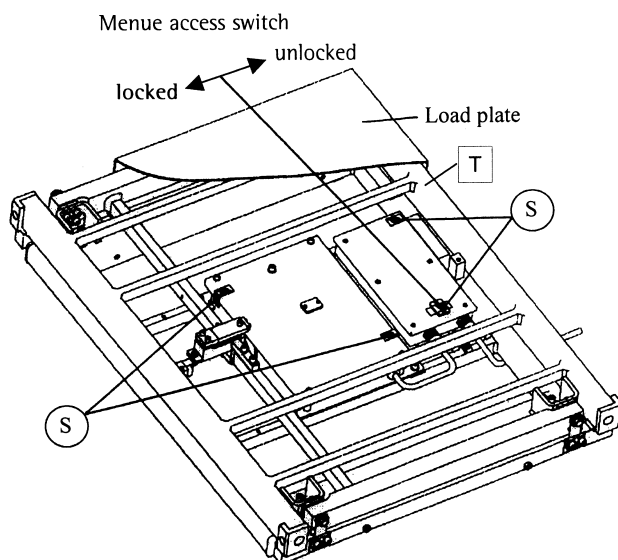
- K** Descriptive plate with CE-conformity mark
- M** Mark for EC verification (green metrology sticker)
- S** Protective seal, not for accuracy class **I**
- T** Plate with model designation
- P** Menu access switch

Indicating and operator terminals isi..., YAC01LA..., YAC01LP..., YAC01FC..., YAC02FC..., TN front-mounted, raised (post-mounted) or positioned separately.
Alternative to terminal: PC with Sartorius Win Scale YSW03 software


SARTORIUS AG GÖTTINGEN Germany **iso-TEST** 12345678 D97-09-018
CE00 0111 **M** **(II)** +10°C / +30°C
 Max 8200 g d= 0,01 g
 Min 0,5 g e= 0,1 g

Types: BA BF, BC BF, BD BF, BF BF, MA BF, MD BF
EC Type-approval D97-09-018 + EC Test certificate D09-96.30

Plates and Markings




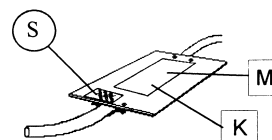- M** Mark for EC verification (green metrology sticker)
- K** Descriptive plate with CE-conformity
- S** Protective mark
- T** Plate with model designation




Affix the ID label of the weighing instrument to the delivered tag plate. Affix the ID tag plate to the data cable of the weighing module near the indicator. The verification officer or an authorized Sartorius representative must then seal the ID tag plate to the fastener.

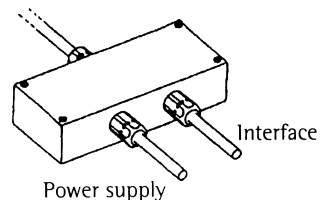Alternative: Junction box



Example of descriptive plate of the verified weighing instrument (weighing module+terminal)

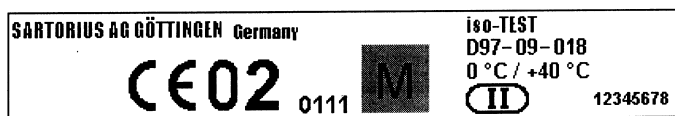K



Example of plate with model designation

T



PPHX251102e

[HC BF Ex]

Type: iso-TEST EC type-approval certificate D97-09-018
+ Type: HC BF EC test certificate D09-96.30

Index

	Page		Page		Page
Accessories (options)	131	Factory setting	11, 29	Recalculation	58
Acoustic signal	16	FlexPrint printout function	93	Reference sample updating	49
Adapt filter	11	Formulation	76	Repeatability test (reproTEST)	45
Animal weighing	54				
Application filter	11	General view of scales	117	Sample ID	33
Application parameter	17			Scale operating menu	31
Application parameter overview	18	Handshake mode	105	Scale-specific information	28
Application programs	46			Second tare memory (preset tare)	84
Application menu settings	18	Individual identification codes	86	Selecting the calibration	
Applications, combining	94	Input	3	parameter	41
Applications, configuring	17	Interface port	101	Software handshake	105
Automatic shutoff	16	Internal calibration	40	Specifications	121
Auto print functions	124	ISO/GMP printout/record	100	Statistics	80
Auto zero	11	isoCAL function	43	Symbols in display	4, 35
				Synchronization	105
Bar graph	4	Keys	3, 16, 35		
Basic weighing function	31			Time-controlled functions	70
Block printout	44, 100	Language setting	7	Toggling between weight units	46
		Line format	24, 102	Totalizing	72
Cabling diagram	107	Loading storad data	91		
Calculation	63	Lot ID	24	Under scale weighing	31
Calibration/adjustment	36			Universal remote control switch	34
Calibration/adjustment, automatic	43	Manual storage – M+	90	User data, input	7
Calibration/adjustment, external	42	Manual/automatic print mode	25	User ID	15, 33
Calibration/adjustment, internal	40	Metrological data, in display	97		
Certificates of compliance	162	MP8 interface emulation	29	Weighing in percent	51
Checkweighing	65			Weight set ID	11
Combining applications	94	Navigating in the Setup	8		
Configuring applications	17				
Configuration	7, 34	Operating design	3		
Consultation	2	Operation	31		
Contents	2	Options (accessories)	131		
Contrast	16	Overview of application menu	15		
Control lines	66				
Control port	34	Parameter settings	6		
Counting	48	Password, enter/change	13		
Counting and checkweighing		Password, general	Appendix		
with stat. evaluation	95	Pin assignment chart	106		
		Plates and markings	170		
Data output functions	97	Power-on mode	35		
Data ID codes	103	Preset tare (2nd tare memory)	84		
Data input format	104	Print application parameters	87		
Data interface port	101	Print manually/automatically	87		
Data output format	102	Printout	23		
Data output, following		Printout configuration	25		
print command	105	Printout, animal weighing	99		
Data output, activate	105	Printout, calculation	99		
Declarations of conformity	136	Printout, time-controlled			
Device parameters	15, 34	functions	70, 99		
Dimensions (scale drawings)	127	Printout, totalizing,			
Display	4	formulation, statistics	99		
Display backlighting	65	Product data memory	91		
EC type-approval certificate	168				
EC type examination certificates	138				
Enter/change password	13				
Enter time/date	9				
Error codes	115				
External calibration	42				
Extra functions					
(application menu)	84				

Entering the General Password

Enter/Change Password

- Select the Setup menu: Press **Setup**
 - > **SETUP SELECTION** is displayed
 - Select the user input function:
Press the **Input** soft key
 - > The password prompt is displayed:
- | | | | |
|--|--|--------------|---|
| SETUP | | PASSW. CHECK | |
| Enter password: XXXXXXXXXX | | | |
| << | | | ↓ |
- Enter the General Password (see below)
 - Confirm password: Press the **↓** soft key
 - > User data is displayed

- § Select the password setting function:
Press the **↵** soft key repeatedly until
- > **Enter password:**
is displayed, together with the current password setting
- Define a new password:
Enter letters/numbers for the new password (8 characters max.)
To delete the current password:
press **⌫** and confirm
- To confirm the new password:
press the **↓** soft key
- Exit the Setup menu:
Press the **<<** soft key
- > Restart your application

Detach the General Password card at the perforations and store it in a secure location

General Password:
40414243

Sartorius AG
Weender Landstrasse 94–108
37075 Goettingen, Germany

Phone +49.551.308.0
Fax +49.551.308.3289
www.sartorius.com

Copyright by Sartorius AG,
Goettingen, Germany.
All rights reserved. No part
of this publication may
be reprinted or translated in
any form or by any means
without the prior written
permission of Sartorius AG.
The status of the information,
specifications and illustrations
in this manual is indicated
by the date given below.
Sartorius AG reserves the
right to make changes to
the technology, features,
specifications and design of
the equipment without notice.

Status:
June 2004, Sartorius AG,
Goettingen, Germany